


# FIREFLY®

Firefly® Play Action™ Toothbrush, Featuring Sonic the Hedgehog™  
the Smart Interactive Sonic Toothbrush That Turns Brushing Teeth Into a Game!



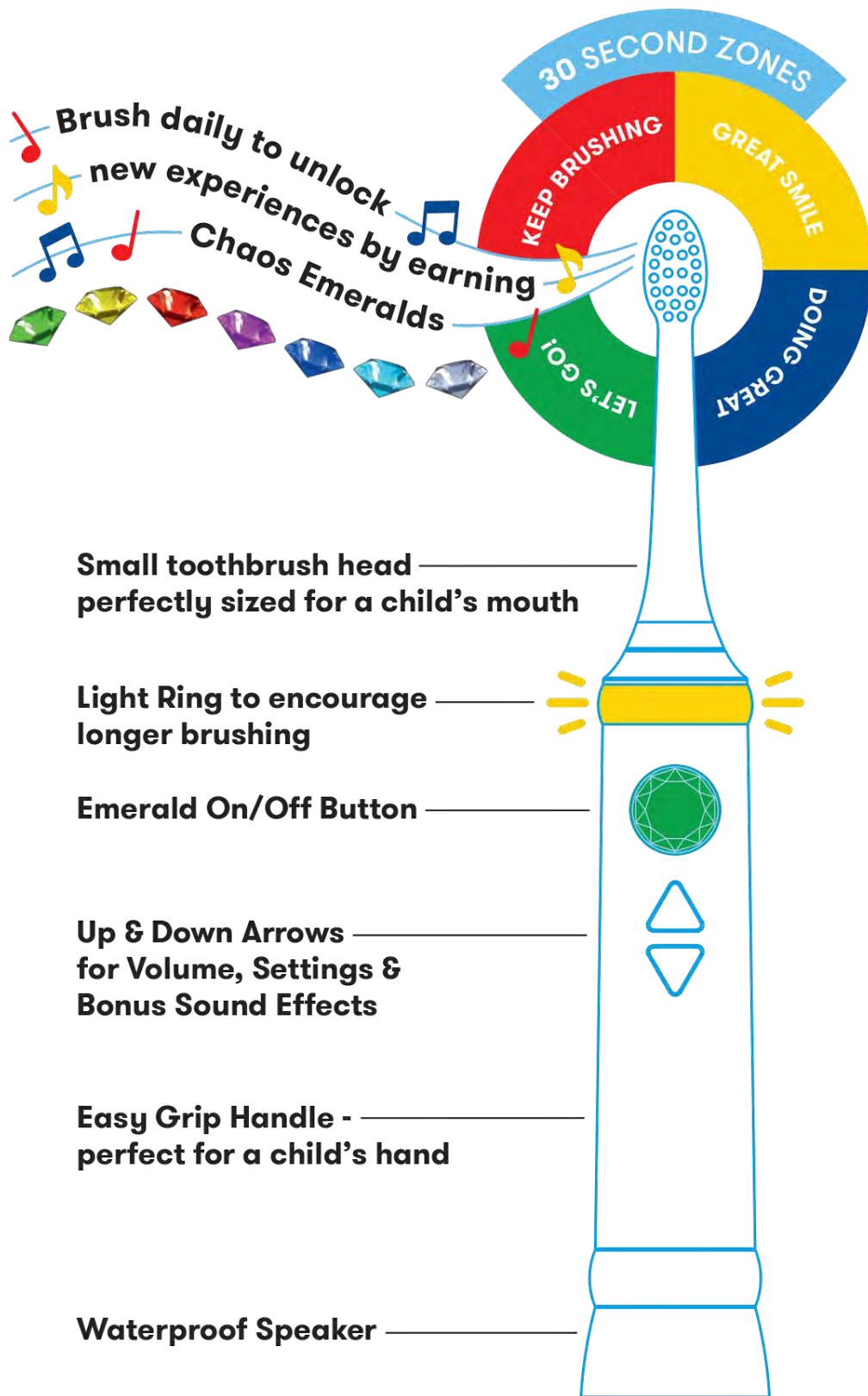

## DIRECTIONS FOR USE OF FIREFLY® PLAY ACTION™ TOOTHBRUSH

### FIREFLY® PLAY ACTION™ TOOTHBRUSH KIT INCLUDES:

1. Toothbrush Handle with:
  - Light Ring to encourage longer brushing and provide feedback
  - Emerald On/Off Button
  - Up & Down Arrows to select levels, play bonus sounds, and adjust volume (in parent mode)
  - Easy Grip Handle perfectly sized for a child's hand
  - Waterproof Speaker with high quality sound
2. (2) Replacement Brush Heads with tongue cleaner
3. 3D Antibacterial Cap\* featuring Sonic the Hedgehog

### INSTRUCTIONS FOR USE

#### HOW TO TURN FIREFLY® PLAY ACTION™ TOOTHBRUSH ON AND OFF

1. Press Emerald Button to power on
2. Long-press Emerald Button to turn off (toothbrush will automatically turn off after inactivity)

#### PARENT MODE - TO CHECK BRUSHING PROGRESS & ADJUST VOLUME

1. Long-press Emerald Button while powered off to access Parent Mode
2. The brush will automatically provide feedback on previous brushing performance including:
  - All 4 zones completed or not completed
  - Recent Progress
  - Time to change brush head (when enough brushing sessions have been completed)
3. The LED ring will illuminate the number of brushing intervals completed as represented by up to 4 of the quadrants of the ring being illuminated. Each illuminated quadrant represents up to 30 seconds of completed brushing time with a maximum of 2 minutes of use shown
4. The brush will then prompt to adjust the volume, "Use the arrows to change the volume." The brush volume can be adjusted to one of three settings:
  - Loud
  - Normal
  - Quiet
5. To exit parent mode, either press the emerald on/off button or wait and the brush will automatically turn off and save the settings



  
**INTERACTIVE  
ADVENTURES**  
Unlock game levels  
by brushing with  
Sonic the Hedgehog

\* Contains silver which naturally inhibits bacterial growth to help protect the cover



## TO USE FIREFLY® PLAY ACTION™ TOOTHBRUSH

1. Wet brush head and apply toothpaste
2. Press Emerald Button to turn on and start brushing, Sonic the Hedgehog says the greeting, "Welcome to the FIREFLY PLAY ACTION toothbrush. Let's clear all four zones. Start brushing to begin."
3. Toothbrush will gently vibrate as music plays and ring lights up
4. Sonic the Hedgehog & friends provide a guided brushing experience starting with zone 1, until all 4 zones of the mouth are complete (zones are approximately 30 seconds each)
5. During brushing session, Sonic the Hedgehog will share encouraging messages, e.g., "Keep Brushing!"
  - Flip toothbrush at any time to hear a "cool flip" line from Sonic the Hedgehog
  - Press Up/Down Arrows to play bonus sound effects - including Sonic the Hedgehog's signature spin and ring collection sound
6. Upon completing all 4 zones, Sonic the Hedgehog celebrates, saves progress, and says "Keep brushing each day to unlock more chaos emeralds. See Ya Next Time!"
7. Once each level has been unlocked, the preferred adventure can be selected with the Up/Down Arrows at the beginning of each brushing

## PRACTICE MODE – FOR FIRST TIME

1. There will be a vocal prompt to press Emerald On/Off Button to begin brushing in practice mode
2. Once complete, use Up/Down Arrows to select color of Emerald, then press Emerald On/Off Button to begin new brushing session

## TO ATTACH AND REMOVE THE TOOTHBRUSH HEAD

1. To attach brush head, align dot on back of head with the arrow on top of handle
2. Firmly place brush head onto handle and twist clockwise until dots on both handle and brush head are aligned
3. To remove brush head, twist brush head counterclockwise until dot on the back of the head aligns with the arrow on the top of the handle and pull straight

## TO REMOVE AND CHANGE THE BATTERY

1. To access the battery compartment, locate screw on the back of handle and remove using a small Phillips head screwdriver with a counterclockwise motion
2. Gently pull down on the red plastic cap to access inside of battery compartment
3. Remove each battery from inside and recycle according to local requirements
4. Replace Alkaline AAA batteries with new of the same type (do not mix old and new batteries) and note the batteries direction according to the battery symbol
5. Replace red plastic cap, ensuring parts are aligned with the screw hole on the back of the handle. Screw hold should be visible and parts firmly secured
6. Replace the screw using a clockwise motion until screw head is flush with the handle being careful to not apply excessive force

