


FIREFLY®

Firefly® Play Action™ Toothbrush featuring Trolls is a Smart Interactive Sonic Toothbrush that turns brushing teeth into a game because at Firefly®, we believe in making healthy habits fun!

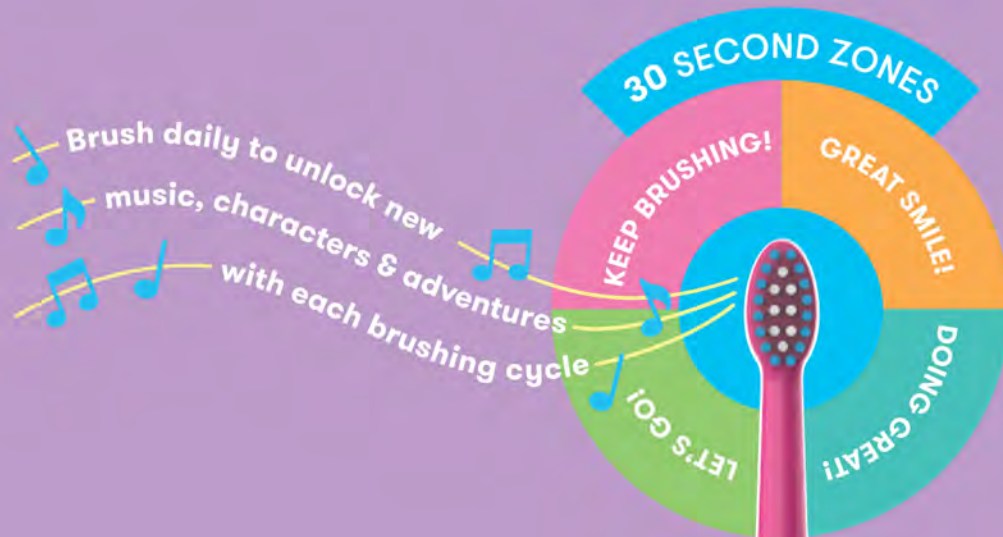

Small toothbrush head — perfectly sized for a child's mouth

Light Ring to encourage longer brushing

Gem Shaped On/Off Button

Up & Down Arrows for volume, settings & bonus sound effects

Easy Grip Handle perfect for a child's hand

Waterproof Speaker





DIRECTIONS FOR USE OF FIREFLY® PLAY ACTION™ TOOTHBRUSH

FIREFLY® PLAY ACTION™ TOOTHBRUSH KIT INCLUDES:

1. Toothbrush Handle with:
 - Light Ring to encourage longer brushing and provide feedback
 - Gem On/Off Button
 - Up & Down Arrows to select levels, play bonus sounds, and adjust volume (in parent mode)
 - Easy Grip Handle perfectly sized for a child's hand
 - Waterproof Speaker with high quality sound
2. (2) Replacement Brush Heads with tongue cleaner
3. 3D Antibacterial Cap* featuring Poppy from Trolls

INSTRUCTIONS FOR USE

HOW TO TURN FIREFLY® PLAY ACTION™ TOOTHBRUSH ON & OFF

1. Press Gem Button to power on
2. Long-press Gem Button to turn off (toothbrush will automatically turn off after inactivity)

PARENT MODE - TO CHECK BRUSHING PROGRESS & ADJUST VOLUME

1. Long-press Gem Button while powered off to access Parent Mode
2. The brush will automatically provide feedback on previous brushing performance including:
 - All 4 zones completed or not completed
 - Recent Progress
 - Time to change brush head (when enough brushing sessions have been completed)
3. The LED ring will illuminate the number of brushing intervals completed as represented by up to 4 of the quadrants of the ring being illuminated. Brushing time will build up from 10 seconds per zone to the dentist recommended 30 seconds per zone with progression through each level. Each illuminated quadrant represents up to 30 seconds of completed brushing time with a maximum of 2 minutes of use shown
4. The brush will then prompt to adjust the volume, "Use the arrows to change the volume." The brush volume can be adjusted to one of three settings:
 - Loud
 - Normal
 - Quiet
5. To exit parent mode, either press the Gem on/off button or wait and the brush will automatically turn off and save the settings

* Contains silver which naturally inhibits bacterial growth to help protect the cover



TO USE FIREFLY® PLAY ACTION™ TOOTHBRUSH

1. Wet brush head and apply toothpaste
2. Press Gem Button to turn on and start brushing
3. Toothbrush will gently vibrate as music plays and ring lights up
4. Poppy & friends provide a guided brushing experience starting with zone 1, until all 4 zones of the mouth are complete (zones are approximately 30 seconds each)
5. During guided brushing session, child receives Trolls themed messages and encouragement, e.g., "Zone 1 complete; Move your brush to the opposite side on top; Music is life!"
6. Press Up/Down Arrows to play bonus sound effects
7. Upon completing all 4- zones, a surprise is unlocked. "Press the Gem Button to unlock your surprise!"
8. Once Gem Button is pressed, it unlocks the next adventure with a Trolls themed reminder, "Play your new adventure next time you brush."
9. Once each level has been unlocked, the preferred adventure can be selected with the Up/Down Arrows



“show me your moves!”

PRACTICE MODE – FOR FIRST TIME

1. There will be a vocal prompt to press Gem On/Off Button to begin brushing in practice mode
2. Once complete, use Up/Down Arrows to select color of Gem, then press Gem On/Off Button to begin new brushing session

TO ATTACH AND REMOVE THE TOOTHBRUSH HEAD

1. To attach brush head, align dot on back of head with the arrow on top of handle
2. Firmly place brush head onto handle and twist clockwise until dots on both handle and brush head are aligned
3. To remove brush head, twist brush head counterclockwise until dot on the back of the head aligns with the arrow on the top of the handle and pull straight

TO REMOVE AND CHANGE THE BATTERY

1. To access the battery compartment, locate screw on the back of handle and remove using a small Phillips head screwdriver with a counterclockwise motion
2. Gently pull down on the red plastic cap to access inside of battery compartment
3. Remove each battery from inside and recycle according to local requirements
4. Replace Alkaline AAA batteries with new of the same type (do not mix old and new batteries) and note the batteries direction according to the battery symbol
4. Replace red plastic cap, ensuring parts are aligned with the screw hole on the back of the handle. Screw hold should be visible and parts firmly secured
5. Replace the screw using a clockwise motion until screw head is flush with the handle being careful to not apply excessive force

