Fit Plus

User Manual

EN

Smart Body Fat Scale

1- 4

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FR

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

Precautions

Read carefully before use:

- ⊘ Do not use the body fat scale if you have a pacemaker or any other internal medical device in your body. When you take a measurement, the scale passes a very low, harmless electrical current through your body. The current may cause your internal medical devices to malfunction.
- ⊘ Do not use the scale on a wet surface. This could lead to slipping.
- ⊘ Do not use a damaged charging cable or a faulty charger to charge your scale.
- ⊘ Dispose of the scale properly following your local regulations.
- ⊘ Do not jump on the scale. This scale has a weight limit of 396 lbs.


Scale Overview

GE Fit Plus Smart Body Fat Scale measures weight, BMI, body fat %, and other body compositions using bioelectrical impedance analysis (BIA) technology. This smart scale is intended for home use only.

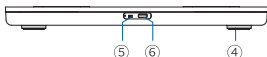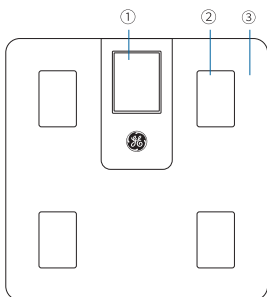Scale Layout

- ① HD LCD Display
- ② Electrodes
- ③ Tempered Glass Surface
- ④ Scale Feet
- ⑤ Reset Button
- ⑥ Charging Port

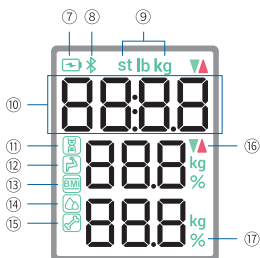LCD Display

- ⑦ Charging Indicator (Always on when the scale is charging and working)
- ⑧ Bluetooth Connection
- ⑨ Unit
- ⑩ Weight Measurement
- ⑪ Body Fat Percentage Measurement
- ⑫ Muscle Mass Measurement
- ⑬ BMI Measurement
- ⑭ Body Water Measurement
- ⑮ Bone Mass Measurement
- ⑯ Data Change
- ⑰ %

Scale Overview

Model:	Fit Plus
Dimensions:	11" × 11" × 1" / 280 × 280 × 25 mm
Weight:	2.9 lb / 1.3 kg
Weighing Range:	0.4 - 396 lb / 0.2 - 180 kg
Power:	200mAh polymer lithium battery
Rated Input:	5V $\overline{\text{=}}$ 200mA

Get Started

It needs you to create a profile in the app and open it when weighing at first use to get more data displayed on the screen.

1. Download the app and create an account

- (1) Download the app. Go to the App Store for iOS or Google Play Store for Android, and search for "Fit Profile", or scan the QR code below.
- (2) Sign up and complete your profile.

Scan QR code to download the app

Sign up

2. Setup your scale

If you use an Android phone, please open "Location". And if iPhone, please agree with the app "Fit Profile" to use your phone's Bluetooth first.

- (1) Charge your scale by using the enclosed USB-C cable with a power adapter (not included).
- (2) Place the scale on a hard flat floor surface.
- (3) Open the "Fit Profile" app and turn on the Bluetooth.
- (4) Turn on the scale by briefly stepping on and off it.



Your scale will be paired with your smartphone automatically within seconds. If not,

- (1) Tap "+ Pair the Device"; Or click the icon at the top-right-hand corner and tap "+ Add device"
- (2) Turn on the scale by briefly stepping on and off it. On the app, tap the smart scale icon, then select "Fit Plus" and click "Connect."

Your scale has now been successfully paired with your smartphone and is ready to start taking your measurements.

3. Take measurements

- (1) Open the "Fit Profile" app, and then keep the app open on the "Home" screen.
- (2) Place the scale on a hard flat floor surface, do not place it on carpet or any other soft surfaces.

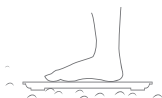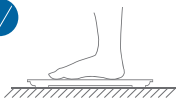- (3) Turn on the scale by briefly stepping on and off it. Then, wait for the scale to set to 0.
- (4) Step on the scale while barefoot till the weight numbers flash twice and display on the scale. Your measurement is over, and the results will also display on the app.

Multiple User Profiles

Separate Fit Profile Accounts:

Users can create or log in to their individual Fit Profile accounts and pair the smart scale with their personal accounts. Weight data for each user remains private and can be synced with third-party apps.

Shared Fit Profile Account:

Users with a Fit Profile account can create multiple user profiles within the same account. All users can log in to the shared account on their own smartphones without needing others to log out.

To add more users and share your account, click the icon at the top left corner on the Home screen in the app.

Note: Users sharing a Fit Profile account can access each other's data recorded under this account. Only the main user's data can be synced with third-party apps.

Athlete Mode

If you have athletic build, you may turn on the Athlete Mode:

Open the Fit Profile app and go to "Me" page -> your profile name -> "Athlete Mode".

Switching Units

By default, the scale uses pounds (lb) as the unit of measurement, but it also supports kilograms (kg).

To switch the unit of measurement on the scale, first connect it to the Fit Profile app, then go to the app "Me" page -> "Settings" -> "Unit".

Recommendations

To have accurate, consistent, and comparable measurement results, please follow the recommendations below.

1. DO' s

- Ensure the personal information in your profile is correct. Choose "Athlete Mode" if you work out intensively and regularly, or you are a muscular person.
- Place the scale on a hard flat floor surface when you weigh.
- Make sure the scale displays a 0 before stepping on to measure. While barefoot, stand still on the scale until the measurement is completed.
- Weigh yourself under similar conditions and around the same time every day.
- Empty your bladder before weighing.

2. DON' Ts

- Do not place the scale on a carpet or a soft surface when you weigh.
- Do not measure with socks or shoes on.
- Do not measure when your body is wet. For example, immediately after a shower, bath, sauna, or strenuous exercise.
- Do not measure after drinking a lot of water or alcohol.
- Do not measure if you have a cold or fever.

3. This smart scale is intended for home use only.

- For people ages 18 years and above, it provides full body composition measurements.
- For children under 10, it only provides weight and BMI measurements.
- For children between the ages of 10 and 17, it provides full body composition measurements. (Please note that the body mass-related items like body fat%, etc., are only for your reference.)

4. This smart scale is not intended for use by pregnant women.

Scale Display

Measurements

Lo Low Power	EEr Overloaded	CAL Calibrating	FULL Fully Charged
EEr 1 Unstable Measurement	---- Unable to calibrate	No Display Low Battery, Auto off	

Messages Displayed

Troubleshooting

Symptom	Cause	Solution
Lo	Low power	Please charge your scale.
EEr	Overload	Use the scale within 396 lbs.
EEr1	Unstable measurement	Place the scale on the flat, hard surface and measure again.
EEr2	Standing on the scale for too long	Step off the scale; your measurement is complete when the weight number flashes twice.
----	Unable to calibrate	Place the scale on the flat, hard surface and let the scale calibrate.
No Display	Auto-off, or no power	Auto-off: Step on then off to turn on the scale; No power: Please charge the scale.
Display Weight Only	Unable to determine body composition	No User Profile: Create an app profile and weigh once with Bluetooth connected. No Electrical Current: Need to stand on the scale barefoot. Users: Not suitable for children under 10. Other Issues: Contact Customer Support for assistance.
Unable to Sync with App	Lost connection, or scale turns off	Reconnect the scale: For iPhones, authorize bluetooth; for Android, enable location. Click the icon on the top-right corner and tap "My Device" -> "+ Add device" on the app.
Unable to Turn on, or Disfunction	Abnormal	Reset the scale: the reset button is next to the charging port.

Warranty

This product is warranted by the manufacturer against defects in materials and workmanship for one (1) year from the original date of purchase. Proof of purchase is required.

This warranty does not cover accidental damage, mistreatment, or unauthorized repair. Detailed warranty policy: www.fit-profile.com/warranty

If you have any issues or questions, please contact Global eFit Inc.

Contact: Global eFit Inc.

cs@fit-profile.com

(866) 786-5988 (Mon – Fri. 9:00 am – 5:00 pm, US Central Time)

FCC Compliance Statements

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, under Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Caution: Changes or modifications not expressly approved by the party responsible for compliance could void the user authority to operate the equipment.

This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment. This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter.

ISED Compliance Statements

This device contains license-exempt transmitter(s)/receiver(s) that comply with Innovation, Science, and Economic Development Canada's license-exempt RSS(s). Operation is subject to the following two conditions:

- (1) This device may not cause interference.
- (2) This device must accept any interference, including interference that may cause undesired operation of the device.

This equipment complies with IC RSS-102 radiation exposure limits set forth for an uncontrolled environment.

Global eFit Inc.

17304 Preston Road, Suite 800, Dallas, TX 75252, USA

www.fit-profile.com

cs@fit-profile.com

1-866-786-5988

(toll-free US & Canada; M-F 9:00am - 5:00pm Central Time)

Model of Independent Module Used in This Product: C21007C10

FCC ID: 2BAI6-GEX00

IC: 30163-GEX00

GE is a registered trademark of General Electric Company.

Manufactured under trademark license.

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Si tiene algún problema o pregunta, comuníquese con Global eFit Inc.