## Fit Prime

**User Manual**

EN

**Digital Body Weight Scale**

1- 4

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FR

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## Precautions

Read carefully before use:



-  Do not use the scale on a wet surface. This could lead to slipping.
-  Dispose of the scale properly following your local regulations.
-  Do not jump on the scale. This scale has a weight limit of 400 lbs.

## Scale Overview

GE Fit Prime Digital Body Weight Scale measures weight and BMI.


This digital scale is intended for home use only.

BMI(Body mass index ) is a person's weight in kilograms divided by the square of height in meters.This indicator will be calculated and shown in Fit Profile App.
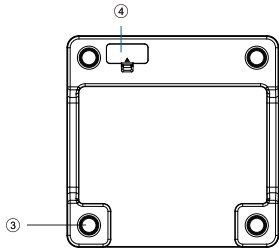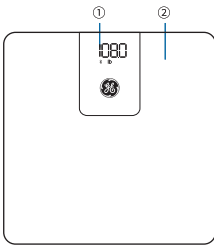
### Scale Layout

- ① LED Display
- ② Tempered Glass Surface
- ③ Scale Feet
- ④ Battery Box/Compartment
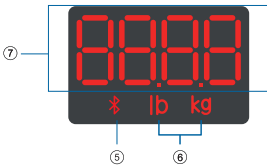
### LED Display

- ⑤ Bluetooth Connection
- ⑥ Unit
- ⑦ Weight Measurement

## Scale Overview

|                 |                             |
|-----------------|-----------------------------|
| Model:          | Fit Prime                   |
| Dimensions:     | 11 x 11 x 1 inch            |
| Weight:         | 2.6lb / 1.18kg              |
| Weighing Range: | 0.4 - 400 lb / 0.2 - 180 kg |
| Power:          | 2*1.5V AAA battery          |

## Explore more useful functions of the smart app!

### 1. Setup your scale

---

- (1) Before putting your digital scale into operation, please open the battery compartment on the underside of the device and insert the two 1.5 V batteries (type AAA) included.
  - **Do not** mix new and old batteries.
  - If the scale not in use for a long period of time, remove the batteries to extend the lifetime of the scale.
- (2) Place the scale on a hard flat floor surface.

### 2. Take measurements

---

- (1) Place the scale on a hard flat floor surface, do not place it on carpet or any other soft surfaces.

- (2) Turn on the scale by briefly stepping on and off it. Then, wait for the scale to set to 0.
- (3) Step on the scale while barefoot till the numbers display on the scale. Your measurement results will also display on the app.

### 3. Download the app and create an account

---

- (1) Download the app. Go to the App Store for iOS or Google Play Store for Android, and search for "Fit Profile", or scan the QR code below.
- (2) Sign up and complete your profile.

Scan QR code to download the app

Sign up

- (3) Open the "Fit Profile" app and turn on the Bluetooth.
- (4) Turn on the scale by briefly stepping on and off it.

Your scale will be paired with your smartphone automatically within seconds. If not,

- (1) Tap "+ Pair the Device".
- (2) Turn on the scale by briefly stepping on and off it. On the app, tap the smart scale icon, then select "Fit Prime Silver" and click "Connect."

Your scale has now been successfully paired with your smartphone. Open the "Fit Profile" app and the app open on the "Home" screen while do the weighing, it will record your data automatically.

You can track your data history freely or sync to Google fit or Apple health, or try other exciting functions with the App.

### Multiple User Profiles

#### Separate Fit Profile Accounts:

Users can create or log in to their individual Fit Profile accounts and pair the smart scale with their personal accounts. Weight data for each user remains private and can be synced with third-party apps.

### Shared Fit Profile Account:

Users with a Fit Profile account can create multiple user profiles within the same account. All users can log in to the shared account on their own smartphones without needing others to log out.

To add more users and share your account, click the icon at the top left corner on the Home screen in the app.

**Note:** Users sharing a Fit Profile account can access each other's data recorded under this account. Only the main user's data can be synced with third-party apps.

### Athlete Mode

If you have athletic build, you may turn on the Athlete Mode:

Open the Fit Profile app and go to "Me" page -> your profile name -> "Athlete Mode".

### Switching Units

By default, the scale uses pounds (lb) as the unit of measurement, but it also supports kilograms (kg).

To switch the unit of measurement on the scale, first connect it to the Fit Profile app, then go to the app "Me" page -> "Settings" -> "Unit".

## Recommendations

**To have accurate, consistent, and comparable measurement results, please follow the recommendations below.**


### 1. DO' s

- Ensure the personal information in your profile is correct.
- Place the scale on a hard flat floor surface when you weigh.
- Make sure the scale displays a 0 before stepping on to measure. While barefoot, stand still on the scale until the measurement is completed.
- Weigh yourself under similar conditions and around the same time every day.
- Empty your bladder before weighing.
- Please use 1 5V AAA batteries
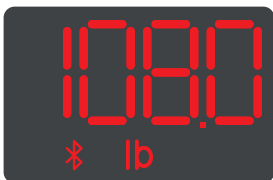
### 2. DON' Ts

- Do not place the scale on a carpet or a soft surface when you weigh.
- Do not measure with socks or shoes on.
- Do not measure when your body is wet. For example, immediately after a shower, bath, sauna, or strenuous exercise.
- Do not measure after drinking a lot of water or alcohol.
- Do not measure if you have a cold or fever.

## Scale Display

### Measurements

|  |   |   |   |
|--|---|---|---|
| <br>Normal working | No display<br>Auto-off, or no power   | <br>No wireless connection | <br>Low power                         |
| <br>Overload       | <br>Unstable measurement | <br>Unable to calibrate    | <br>The scale needs to be calibrated |

## Messages Displayed

## Troubleshooting

| Symptom  | Cause                            | Solution  |
|--|----------------------------------|---|
|    | Low power                        | Please change and use the new batteries.  |
|    | Overload                         | Use the scale within 400 lbs.   |
|    | Unstable measurement             | Place the scale on the flat, hard surface and measure again.  |
|  | Unable to calibrate              | Place the scale on the flat, hard surface and let the scale calibrate.  |
| No Display   | Auto-off, or no power            | Auto-off: Step on then off to turn on the scale;<br>No power: change and use the new batteries.   |
|  | No wireless connection           | Reconnect the scale. On the app tap the device icon located in the top-right-hand corner, and then tap "+Add device". Next, turn on the scale then tap the smart scale icon on the app screen. Then, select "Fit Prime Silver" and tap "Connect." |
|  | The scale needs to be calibrated | Step on the scale with one foot to turn it on, and step on it again after it shows "0.00"   |

## Warranty

This product is warranted by the manufacturer against defects in materials and workmanship for one (1) year from the original date of purchase. Proof of purchase is required.

This warranty does not cover accidental damage, mistreatment, or unauthorized repair. Detailed warranty policy: [www.fit-profile.com/warranty](http://www.fit-profile.com/warranty)

If you have any issues or questions, please contact Global eFit Inc.

**Contact:** Global eFit Inc.

[cs@fit-profile.com](mailto:cs@fit-profile.com)

(866) 786-5988 (Mon – Fri. 9:00 am – 5:00 pm, US Central Time)

## FCC Compliance Statements

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, under Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Caution: Changes or modifications not expressly approved by the party responsible for compliance could void the user authority to operate the equipment.

The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction

## ISED Compliance Statements

This device contains license-exempt transmitter(s)/receiver(s) that comply with Innovation, Science, and Economic Development Canada's license-exempt RSS(s). Operation is subject to the following two conditions:

- (1) This device may not cause interference.
- (2) This device must accept any interference, including interference that may cause undesired operation of the device.

This equipment complies with IC RSS-102 radiation exposure limits set forth for an uncontrolled environment.

Global eFit Inc.

17304 Preston Road, Suite 800, Dallas, TX 75252, USA

[www.fit-profile.com](http://www.fit-profile.com)

[cs@fit-profile.com](mailto:cs@fit-profile.com)

1-866-786-5988

(toll-free US & Canada; M-F 9:00am - 5:00pm Central Time)

FCC ID: 2BAI6-SC500

IC: 30163-SC500

GE is a registered trademark of General Electric Company.

Manufactured under trademark license.

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Si tiene algún problema o pregunta, comuníquese con Global eFit Inc.