



Triple Action Eye Care Therapy Wand
(User Manual)



Precautions

Thank you for purchasing this product. Please read this instruction manual carefully before using this product and keep it for future reference.


1. This product should not be used by small children, pregnant women, patients with epilepsy, tumors and cancer patients, patients with cardiac pacemakers, implanted with cardiac defibrillators or other implanted devices. If you have any medical questions, please consult your physician before use.
2. Stop using this device if you feel any discomfort. Do not use this product while wearing contact lenses. If you feel abnormalities in your eyes, eyelids or skin during use or right after use and it persists, consult your medical specialist.
3. Avoid direct contact with the eyes.
4. Do not use while charging.
5. For first time users or people with sensitive skin: Apply for a short duration to ensure there are no side-effects.

6. Do not repair, disassemble or modify this product without authorization.

Function Introduction

This Eye Care Therapy Wand is based on LED, Thermal, and high-frequency Vibration was developed specifically for the delicate skin around the eyes. It assists in promoting skin elasticity, minimizing wrinkles and fine-lines, smoothing the skin around the eyes, and reducing under-eye puffiness. For best results apply some Eye Cream before using this wand.

Product Structure

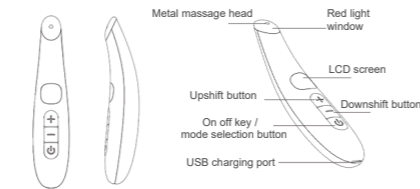

Product Specifications

Product name	Intelligent eye beauty device	Product model no	1203
Net weight	47.5g	Power	0.7W
Rated voltage	DC 3.7V	Material	Zinc alloy+ABS
Temperature	37°C-45°C	Working voltage	3.7-4.2v
Vibration frequency	10000RPM		
Charging time	1.5 hours		
Unit size	139×29×28.5mm		

Three Working Modes

Mode 1: Vibration + hot compress + red light
Triple-effect care, whitening and skin rejuvenation anti-wrinkle, reshaping the elasticity of the eye muscles.


Mode 2: hot compress + red light
The hot compress wakes up the sleeping eye muscles, and the LED red light care enhances cell activity, promotes collagen regeneration, and repairs damaged cells.

Mode 3: Vibration + hot compress
10,000 times / minute high-frequency vibration massage, 37 ° C -45 ° C constant heating, wake up the sleeping eye muscles, accelerate the absorption of nutrients.

Using Method



1. Apply your favorite Eye Cream around the eyes and massage in circular motions until it is fully absorbed. (Recommended)

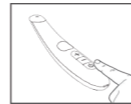

2. Long press (switch button) to turn on device and short press (switch button) to switch the mode, and press "+"/"- " to adjust the temperature as needed.



3. Holding the upper part of the device body, use the massage head to massage as shown.



4. After using each time, long pressing (on/off button) to turn off, please wipe the massage head with a dry towel or paper towel.

Charging Method

1. Charge this Beauty Wand between applications.
2. During the charging process, the display screen battery symbol flashes and the battery symbol stops flashing after charging is completed.
3. Please charge occasionally to avoid aging after not using for a long time. This device is not water-proof or water-resistant.
4. When charging is complete, please unplug the charger.

Maintenance

1. If applying Eye Cream before treatment, spread the cream evenly and let it absorb before using the device.
2. Please wipe the product clean with a dry towel after each treatment is completed.
3. Please don't pour water on the device or wash the body with water. This device is not water-proof or water-resistant.
4. Please don't use corrosive to clean the device.
5. Please store it in a cool dry place.
6. Please don't use the device if it is damaged.
7. Please don't give a strong impact on the body.

Troubleshooting

Frequent Questions	Possible Causes	Treatment
The product can't be turned on	The turning on mode is wrong	Short press the on/off button
	Power is exhausted	Please try charging, if you still can't turn on after charging, please contact the seller
The product does not work after being turned on.	The battery is exhausted	Please try charging, if it is still not available after charging, please contact the seller.
The product's vibration is not strong after touching the skin or does not work.	The resistance is big after touching the skin.	Using together with skin care products is recommended to achieve more significant effect.

Warranty Policy

GLO24K offers a limited warranty on its Beauty Devices (check website for details)



TO MAXIMIZE RESULTS SCAN HERE
For more information about the device please visit
WWW.GLO24K.COM
Dist. by Y&E Ent., LLC
Atlanta, GA 30338 USA
71-75 Shelton Street, London, UK