

## Functional Features

### 1 Different lighting effects:

Various data show that different wavelengths of light can achieve different penetration requirements. The product has both 670nm and 830nm wavelengths, which can meet the needs of different parts. 60pcs 670nm red light beads, 60pcs 830nm infrared light beads, a total of 120 light beads.

### 2 Premium-quality and light:

Low light wane and long usage lifespan. Foldable and portable, light weight, only 133g(4.7 ounce), no burden to use.

### 3 Optical Design of Convex:

Convex lamp beads can make the light more concentrated, so that the light travels in a straight line. It can reach 1-3mm hair follicle root, activate the hair follicle and let it re-enter the healthy growth phase. Pure physical therapy with no side effects.

### 4 Automatic timing function:

Once turned on, it will automatically turn off after 20 minutes, no worry about over-treatment. If you need to pause during use, you can long press the switch to turn it off. It is recommended to use 20 minutes each time, 3-4 times a week.

### 5 Function:

Red light therapy hat can promote the proliferation of brain cells, speed up cell energy metabolism, and relax the brain. It also can reduce oil production, decrease shedding, promote hair growth, restore hair growth healthy cycle, and repair damaged hair.

## Caution

- 1 Do not washed and avoid water contact with plugs and connectors. It is recommended to use professional cleaners for wiping.
- 2 Please use the matching power cord, it is recommended to use a 5VIA plug for charging.
- 3 When used with a power bank, the product will be limited by the power of the power bank device, and the power and brightness will be significantly reduced.
- 4 Avoid squeezing and over-bending the product when using it.


## Quality and Technical Certification

- 1 The product is produced by Shenzhen Yihong Lighting Co., Ltd. entrusted by Scienlodic Ltd.
- 2 FDA owner number is 10082798.
- 3 The product has passed the European patent, and the US patent is pending. Do not misappropriate and imitate without authorization.

## After Sales

If you have any problems in use, please seek help from the seller or professional guidance in time.

All Scienlodic products with quality problems can be reissued or replaced free of charge within 2 years. Please check and use the product in time, so as not to miss the post-sale period.

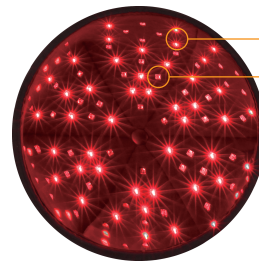

60pcs visible red light

60pcs invisible infrared light

**Tips:** Visible light for human eyes is about 400-800 nm. Half of the light on the hat is 830nm invisible infrared light. You can use a power meter to test whether the light is bright.

**Scienlodic**<sup>®</sup>  
PHOTON THERAPY PAD

## Usage Demonstration



1. Product connect to the controller




2. Plug in the power cord (plug not included)



3. Long press the switch to turn it on or off



4. Put the hat on the head

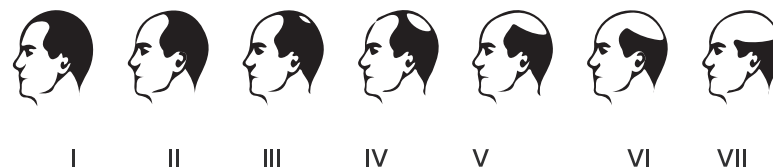

APPLICABLE

NOT APPLICABLE

Levels	3 Levels Adjustable Pulse Mode
L0:No Pulse	Suitable for short-term continuous treatment
L1: 10Hz	Helps induce brainwave entrainment
L2:40Hz	Helps restore the brainwaves to a natural order



Product name	Red light therapy hat
Wavelength	670nm&830nm
Number of led light	60pcs*670nm, 60pcs*830nm, Total 120pcs Led lights
Product perimeter	25.6in/65cm
Product diamete	7.87in/20cm
Product weight	4.7ounce/133g
Due to manual measurement, error between 0.4-0.8 inch	