



Caution

1. Products are not waterproof, pay attention to keep the product from contacting water.
2. The higher the light power generated heat will also be higher, pay attention to the use of distance and use time, to avoid burns.
3. Please use the matching charging cable, it's recommended to use the 5VIA plug.
4. Low battery power will affect the working performance. Please charge it in time.
5. If you have any questions during use, please consult the seller on the order page in time.

Quality and Technical Certification

1. The product is produced by Shenzhen Yihong Lighting Co., Ltd. entrusted by Scienlodice Ltd.
2. FDA owner number is 10082798
3. The product has passed the European patent, and the US patent is pending. Do not misappropriate and imitate without authorization.

After Sales

If you have any problems in use, please seek help from the seller or professional guidance in time.
All Scienlodice products with quality problems can replace or refund within 2 years.

Product Features

1. Different wavelengths of light:

Different wavelengths of light can achieve different penetration requirements. Red light therapy lamp in purple contains 3 wavelengths of blue light 460nm, red light 660nm and infrared 850nm. Red light therapy lamp in red contains 2 wavelengths of 660nm and 850nm.

2. Innovative pulse irradiation mode:

Red light therapy lamp in Purple has 2 pulse modes-blue light 460nm pulse mode and red light 660nm pulse mode. Red light therapy lamp in Red has 1 pulse mode - red light 660nm pulse mode.


3. Convenient, Rechargeable and Portable

4. Wide Range of Application:

It can relieve pain and relax muscles, and is suitable for use on various parts of the body such as the face, hands and joints. It helps with arthritis and inflammation. It is recommended to use it for 10-20 minutes at a time, several times a week. Keep using it for a few weeks for good results.



Usage Demonstration

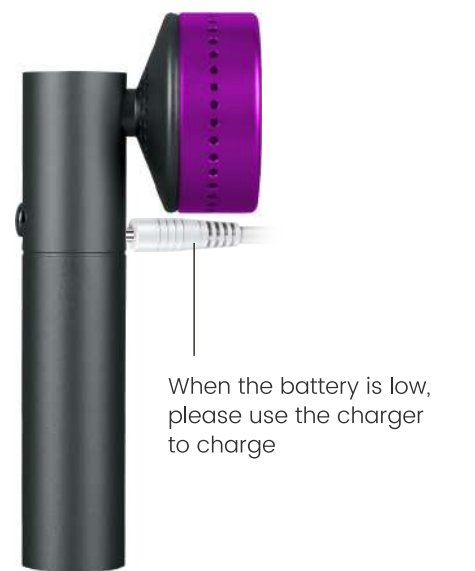

How to use red light therapy lamp in purple?

1. Press once, turn on. Simultaneously emit blue light, infrared light and red light.
2. Press again, start blue light pulse mode.
3. Press the 3rd time, start red light pulse mode.
4. Press the 4th time, turn off.



How to use red light therapy lamp in red?

1. Press once, turn on. Simultaneously emit red light and infrared light.
2. Press again, start red light pulse mode.
3. Press the 3rd time, turn off.

Product name	Multi-wavelength light therapy lamp	
Appearance color	Purple	Red
Exterior material	Aluminum alloy	
Wavelength	460nm&660nm&850nm	660nm&850nm
Number of led light	460nm*3pcs - 660nm*3pcs 850nm*3pcs - Total 9pcs	660nm*5pcs - 850nm*4pcs Total 9pcs
Light power density at high brightness mode	When measured at 1 inch: over 129mw/cm ² When measured at 3 inch: 72.5mw/cm ²	When measured at 1 inch: over 98mw/cm ² When measured at 3 inch: 56mw/cm ²
Product size	Handle: 5.1*1 inch (13*2.5*2.5cm) Lamp head: 2*2.5 inch (5*5*6.5cm) Overall: 5.9*2*2.5 inch (15*5*6.5cm) Weight: 0.29lb (133g)	
Product power	9W	
Battery information	3200mA	