MUNCHIES™

- (EN) **EPS: Quick Start Guide**
- (ES) **EPS: Guía de inicio rápido**
- (FR) **EPS: Guide de démarrage rapide**
- (IT) **EPS: Guida Rapida**
- (DE) **EPS: Schnellstartanleitung**
- (NL) **EPS: Snelstartgids**

 **EOCA MD Pty Ltd** - 16 Transvaal Ave, Double Bay, NSW
2028, Australia, info@orthomunchies.com

AF Pharma Service Europe SL
Muntaner 281, 08021, Barcelona, Spain

**AF Pharma
Service UK LTD**

Net Wt: 28g

MUNEPS2INTM-1.0.01

Align Technology GmbH,
Dürener Straße 405
50858 Köln, Germany

GTIN: 09356117000022


Customer Care:

Australia Ph: 1300 688 565,
USA/Canada Ph: 844 808 8272
www.orthomunchies.com


(01)09356117000022

MUNCHIES™ FEATURES

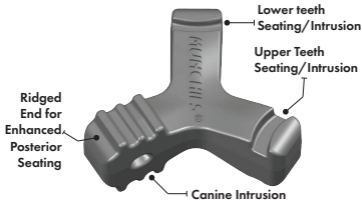EPS: Quick Start Guide

1. With your aligners in place, put your top teeth into the bigger groove on the Munchies™ EPS and bite your lower teeth against it. Bite on all your front teeth for about 20 seconds each.

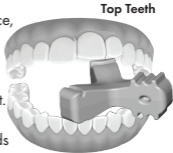Bottom Teeth

2. With your aligners in place, put your bottom teeth into the smaller groove on the Munchies™ EPS and bite your top teeth against it.

Bite on all your front teeth for about 20 seconds each.

3. With your aligners in place, put the ridged end of the Munchies™ EPS between your upper and lower back teeth on one side of your mouth and bite up and down for 30 seconds.

Back Teeth

You will feel the aligner fit far better.

Repeat on the other side of your mouth.

4. Repeat this process each time you re-insert your aligners.

MUNCHIES™ INDICATIONS FOR USE

SUPERIOR SEATING

The anatomically shaped Grooves and Ridges on the arms of each Munchies™ device have been designed to engage each tooth to deliver specific forces to ensure your aligners are seated intimately against your front and back teeth and engaging the critical attachments, without damaging them in any way.

DISCOMFORT RELIEF

Munchies™ devices have visco-elastic properties that provide relief from the discomfort felt during orthodontic tooth movement, especially in the first 24 hours after changing to the next aligner or after a wire adjustment. Patients gain relief when biting initially into a soft device, followed by a firmer device, the phenomenon is called the Bite Wafer Effect^{1,2}.

AID INTRUSION MOVEMENTS

The anatomically shaped Grooves on the arms of each Munchies™ device have been designed to engage with either your upper and lower front teeth, where biting into the device assists the aligner to deliver the light, continuous force needed to move your teeth.

RECAPTURE THE PERFECT FIT

Sometimes during treatment with clear aligners, your teeth do not move exactly as predicted, resulting in a small space appearing between the biting edge of your teeth and the aligner. If you notice a space or your Healthcare Professional identifies an issue with your aligner fit, using Munchies™, start by biting on your back teeth and then move to your front teeth, both upper and lower. Continue to gently bite down until you notice that the space has reduced.

References: 1. Saloom, H. F. Pain intensity and control with fixed orthodontic appliance therapy [A clinical comparative study on Iraqi sample]. Scientific Journal Published by the College of Dentistry-University of Baghdad, 122; 2. Farzanegan, F., Zabarjad, S. M., Alizadeh, S., & Ahrari, F. (2012). Pain reduction after initial archwire placement in orthodontic patients: a randomized clinical trial. *American Journal of Orthodontics and Dentofacial Orthopedics*, 141(2), 169-173.

Intended User: Any patient undergoing orthodontic treatment (clear aligners and braces).

Indications for Use: Munchies™ devices are used to assist in seating thermoplastic orthodontic aligners or retainers, relieve discomfort associated with orthodontic treatment (including with braces), aid intrusion movements and help recapture the perfect fit of your aligners if you get off track. Munchies™ are used to ensure that orthodontic tooth movements occur as predictably and comfortably as programmed by your orthodontic treatment plan.

Performance Characteristics: Munchies™ devices are fabricated from visco-elastic material that enables the device to compress when used and rebound afterwards to its original shape.
Precautions: For intermittent use of seating aligners or retainers and relief of discomfort caused by all forms of orthodontic treatment only.

Warnings: Use of excessive pressure to seat the aligners or retainers may cause damage to either the Munchies™ device or aligners.

Life Span: Munchies™ devices when used appropriately should last between 2 to 3 weeks, such as seating your aligners each and every time after re-insertion. Overuse will lead to premature tearing, please discard if tearing occurs.

Contraindications: Do not use if you have an allergy or hypersensitivity to silicone.

Adverse Reactions: None known - any serious adverse incident experienced should be reported to the manufacturer and the relevant competent authority.

Use in Children: This device is safe for children to use with braces or clear aligners.

Biocompatibility: ISO 10933 Biocompatibility testing has been performed on the Munchies™ with all tests passing the strict standards (tests performed cytotoxicity, skin sensitisation, irritation or subcutaneous reactivity, systemic toxicity & haemocompatibility).

Maintenance, Cleaning & Storage

Prior to use, rinse your Munchies™ devices under warm water for at least 10 seconds. The most hygienic option to sanitize your Munchies™ device is to immerse them in an ultrasonic bath for 5 minutes each day. After rinsing or an ultrasonic cleaning cycle, dry and place in your Munchies™ devices carrying case or alginer case. Munchies™ devices can also be gently scrubbed with a toothbrush and toothpaste if any debris has become attached during use. Store Munchies™ at ambient temperatures (15-25°C).

Manufacturing Quality

Munchies™ devices are manufactured using Class VI medical grade silicones, and designed in Australia and manufactured in China. Munchies™ devices are listed as a Class I medical device on the TGA (ARTG Identifier: 237778), FDA and are CE marked.

Wear & Lifespan

Chewing too aggressively or too often can cause damage to your Munchies™ device, especially in the intrusive mode. Most Munchies™ devices last about two or three weeks if the correct forces have been applied.

You will notice significant indentation or deformation on the device when you use it, but these will generally rebound within 10 minutes due to the elastic nature of the material. Premature tearing is generally due to overuse or too great a force being applied.

To avoid this, simply use less force or use less frequently.

Torn devices should be replaced immediately.

WARNING:


Excessive chewing or bite forces will lead to your Munchies™ devices tearing and should be avoided.

Premature tearing in the device generally means that you have been biting down too hard or using your device too often!

Disposal

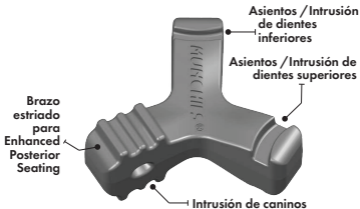Please dispose of your used Munchies™ responsibly, we recommend discarding in the General Household waste bin.

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