

 ***brio***

SmartClean

Premium Sonic Toothbrush



Contents

| | |
|------------------------------------|----|
| Safety Information | 2 |
| Components of your Brio SmartClean | 8 |
| Getting Started | 10 |
| Connecting and Charging | 10 |
| Best Practices for Li-On Batteries | 11 |
| Using the Brio SmartClean | 12 |
| Brushing Techniques | 13 |
| Brushing Modes | 15 |
| Caring for your Brio SmartClean | 18 |
| Recycling your Brio SmartClean | 19 |
| Warranty / Support | 20 |

Safety Information

IMPORTANT SAFEGUARDS: READ ALL INSTRUCTIONS BEFORE USING

The SmartClean 201 has been carefully designed to offer you and your family a new and unique brushing experience that is safe and effective. When using electrical products, basic safety precautions should always be followed.

To reduce the risk of electric shock:

- Always unplug this product after charging is complete
- Do not use while bathing
- Do not place the product where it can fall into water
- Do not place, drop, or otherwise submerge product in water or other liquid
- If charger falls into water, unplug immediately and discontinue use
- Never use charger if the cord or plug are damaged

WARNINGS To reduce the risk of burns, electrocution, fire or injury:

- Use this product only for its intended use as described in the manual.
- Do not use accessories, chargers, or attachments that are made by other manufacturers.
- Discontinue use if product is damaged in any way.


- Do not allow any portion of this product to be used as a toy.
 - Do not pull or carry by cord, use cord as a handle, or pull cord around sharp edges or corners.
 - Keep cord away from heated surfaces.
 - Do not unplug by pulling on the cord. To unplug, grasp the plug, not the cord.
 - Only connect this product to a properly grounded outlet that is GFCI protected.
 - Do not plug into an outlet with a voltage other than specified on the bottom of the charger.
 - Never force the plug into the outlet; Discontinue use if it does not easily fit into the outlet
 - Do not use the charger outdoors
 - Avoid using brush head with bent or crushed bristles. Excessively worn bristles may break off while brushing. Replace brush head every 3 months or sooner, depending on wear.
-

- This product is designed to clean your teeth, gums, and tongue only. Discontinue use of this product and contact a physician or dentist if discomfort or pain is experienced.
- This product is not intended for use by persons with reduced physical, sensory, or mental capabilities or lack of experience or knowledge of this product unless they are supervised or otherwise given instruction by a person responsible for their safety.
- Children should be supervised while using the product to ensure safe and proper use.
- Do not use attachments other than those sold by the manufacturer.
- Never place any portion of this product in the dishwasher.
- If the product is dropped, the brush head should be replaced before the next use, even if no damage is visible.

Medical Warnings

- Consult your dental professional before you use this product if you have had oral or gum surgery in the previous 2 months.
- Contact your dental professional if excessive bleeding occurs after using this product or if bleeding occurs beyond the first week of use.
- If you have a pacemaker or other implanted device, contact your physician or the device manufacturer prior to use.
- This product has been tested and is compliant with safety standards for electromagnetic devices.
- Consult your physician prior to using this product if you have medical concerns.

Components of your Brio SmartClean

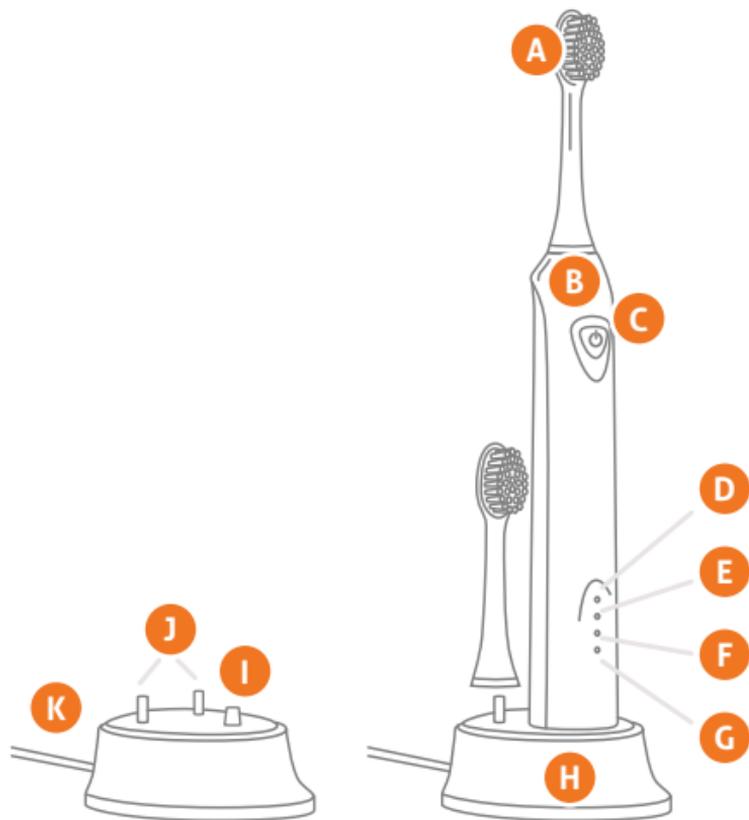

- A Brush Head
- B Toothbrush Handle
- C On/Off & Brushing Mode Button
- D Normal Working Mode Indicator
- E Soft Working Mode Indicator
- F Massage Working Mode Indicator
- G Charging Indicator
- H Charger Unit
- I Handle Stand
- J Brush Head Stand
- K Charging Cable

Getting Started

Connecting and Charging

- The SmartClean 201 has a water-resistant handle. It is designed to be rinsed lightly but not meant to be used under water or in very wet situations.
- Charge the brush before first use. The brush will come from the factory with a partial charge but should be fully charged prior to use.
- To charge the brush, plug in the charging base to a suitable outlet as described in the "Safety Information" section. Place the brush handle in the cradle. When the unit is charging, the blue "charging" indicator light on the handle will illuminate.
- A full charge takes approximately 24 hours and will sustain the brush for about six weeks of regular brushing before requiring a charge.

- The charge indicator light will switch from a solid blue light when charging to a flashing blue light when fully charged.

Best Practices for Li-ON Batteries

- The LI-ON battery cannot be overcharged, but the battery's life can be improved by more frequent charging.
- Fully charging the battery (leave on charging base for at least 24 hours) every 1-2 weeks will extend the life of the toothbrush.
- When the battery is too low, the brush will either not turn on or shut off mid-cycle. This is done to protect the battery from excessive discharge. When this happens, make sure to fully charge the battery before continuing use.
- If any diminished power or vigor is noted while using the brush, fully charge the brush before next use.
- Avoid exposing the toothbrush to excessively high or low temperatures.

Using the Brio SmartClean

Remove protective packaging and attach the brush head to the handle. Wet bristles and apply any kind of toothpaste. To avoid splashing, place the brush head in your mouth before turning on the brush.

You may notice a slight tickling or tingling sensation when first using a sonic toothbrush. As you become more accustomed to brushing with a sonic toothbrush, this sensitivity will diminish. New sonic toothbrush users may prefer to use the toothbrush on "Soft" mode for the first two weeks of use before using the brush in "Normal" mode. Replace the Brush head every 3 months, or more often if needed. Replace the head if the bristles have become much softer, bent, or misshapen.

Brushing Technique with the Brio SmartClean

It is important to note that brushing technique with a sonic toothbrush is very different than using a manual toothbrush. Please read this section carefully.

- Place the brush head in your mouth with the bristles at a 45 degree angle to your teeth.
- Start in a specific zone of your mouth – upper outside, upper inside, lower outside or lower inside.
- The 30-second timer will provide a reminder to switch zones. The SmartClean has a two-minute cycle time: After the cycle is complete, the brush will turn off.
- Gently press the button to start the brushing cycle. Hold the handle with a light grip and gently move the bristles in a small circular motion from the gums to the teeth and back. For a right-handed person, use a counter-clockwise motion, and a clockwise motion for a left-handed person.

- Spend several seconds on each tooth before moving on. Use light pressure and let the vibrations of the brush do the work.
- If you wish to whiten teeth or remove stains, extra time may be needed in those areas. Spend extra time on the areas of concern after you have completed a full brushing cycle of your entire mouth.
- If you have braces, gently turn or roll the brush head all the way around the brackets and wires in order to reach the entire tooth. Continue until you have covered your entire mouth. Fully brushing your mouth with braces may require more brushing time than a single cycle.
- Dental restorations such as implants, crowns, caps, bridges, and fillings should be brushed as you would your natural tooth surfaces unless otherwise directed by your dental professional.

Brushing Modes

- The Brio SmartClean offers three different brushing modes for differing oral care needs:
 - Normal Mode - Exceptional cleaning power for experienced sonic toothbrush users
 - Soft Mode - Gentle, yet thorough cleaning for sensitive areas or new sonic toothbrush users
 - Massage Mode - Gentle stimulation of the gums

How to Switch Between Brushing Modes

- Power On/Off and brushing mode selection are controlled by the power button. Pressing this button will cycle through the modes in order. Pressing the button again will turn off the toothbrush.
- The SmartClean 201 automatically starts in the Normal mode when pressing the power button once.

- To change from Normal mode to Soft mode, press the power button a second time.
- Pressing the power button again will switch to massage mode.
- The brushing mode indicator light indicates the current brushing mode.
- To turn off the brush, press the brushing mode button.
If the brush has been on for less than 15 seconds, it is necessary to cycle through all the modes before the brush will turn off. After the brush has been running for 15 seconds without a cycle change, it will turn off with a single push of the power button.

Caring For Your Brio SmartClean

- Thoroughly rinse brush head and bristles after each use. Towel dry and remove excess water after rinsing.
- Occasionally remove brush head to clean handle and head separately under running water. Use caution not to fully submerge the brush handle.
- Make sure to unplug the charging unit before cleaning it. The charging unit should only be cleaned with a dry or damp cloth. Never submerge the charging unit or hold it under running water.
- Do not plug in charging unit until it is completely dry.

Recycling Your Brio SmartClean

This product contains rechargeable batteries. Due to the construction and risks of dismantling this product, it is recommended that the consumer does not attempt to remove the battery from the brush handle. In the interest of protecting the environment, please recycle this unit properly with a responsible recycler in your area. For local electronics recyclers in your area visit www.e-stewards.org. If an electronics recycler is not available, please visit www.brio4life.com for instructions on how to send the brush back to us for proper recycling.

Warranty & Support

The Brio SmartClean is guaranteed to function properly from two years from the date of purchase. Defects due to faulty materials and workmanship will be repaired or replaced at Brio's expense with a qualified receipt.

To contact Brio customer service please call (248) 268-6530.

For more information on the warranty service process, please visit www.brio4life.com/service

Warranty Exclusions

Items not covered under warranty:

- Brush heads
- Damage caused by misuse, abuse, neglect, alterations, or unauthorized repair.
- Normal wear and tear, including chips, scratches, abrasions, and discoloration or fading. Slightly diminished battery capacity after first year of use.

Implied Warranties

Any implied warranties, including implied warranties of merchantability and fitness for a particular purpose, shall be limited to the duration of the express warranties set forth above. In some states, limitations on duration of implied warranties do not apply.

Limitation of Remedies

- In no event shall Brio or any of its affiliated companies be liable for any special, incidental, or consequential damages based upon breach of warranty, breach of contract, negligence, tort, or any other legal theory. Such damages include, without limitation, loss of savings or revenue, loss of profit, loss of use; the claims of third parties including, without limitation, dentists and dental hygienists, and cost of any substitute equipment or services. Some states do not allow the exclusion or limitation of incidental or consequential damages.

Troubleshooting

- If you are unable to resolve a problem of any kind based on the information within these instructions. Please contact Brio Customer Service at (248) 268-6530 to speak to a representative.

Problem

Toothbrush will not turn on

Tickling sensation when using the SmartClean

The brush head rattles

Brush seems less powerful

Solution

Recharge your brush overnight. Make sure the charging light is illuminated when the brush is placed on the charger. Remove from charger when the charging light has started blinking.

A slight tingling/tickling sensation is normal for new users of sonic toothbrushes. This sensation will diminish quickly with continued use.

Remove and re-attach brush head to handle, being sure to push down firmly.

The brush may need recharging. Place on the charger until the charge light indicates a full charge by flashing.

Problem

The brush handle got excessively wet or was submerged underwater and will no longer turn on.

Solution

Discontinue use. Dry off handle and place the brush handle in a sealed bag of uncooked rice for 7 days to assist in removing moisture from the brush. If the brush will not turn on or function properly, please properly dispose of the brush.



Be a **well** being.™