THANK YOU FOR JOINING
THE LUCENT FAMILY

We created this brand and
every product within the
range to make the world
shine a little brighter.
Starting with your smile.

STAY LUCENT

BE
LUCENT


PRISM

Sonic Electric Toothbrush

MANUAL AND WARRANTY

PRISM

INSTRUCTIONS OF USE	ENGLISH	5-6
HANDBUCH	DEUTSCH	7-8
MODE D'EMPLOI	FRANÇAIS	9-10
MANUAL DE USUARIO	ESPAÑOL	11-12
MANUALE D'USO	ITALIANO	13-14
GEBRUIKSAANWIJZING	NEDERLANDS	15-16
BRUKSANVISNING	SVENSKA	17-18
BRUGSANVISNING	NORSK	19-20
BRUGSANVISNING	DANSK	21-22
KÄYTTÖOHJEET	SUOMALAINEN	23-24
INSTRUKCJA UŻYCIA	POLSKI	25-26

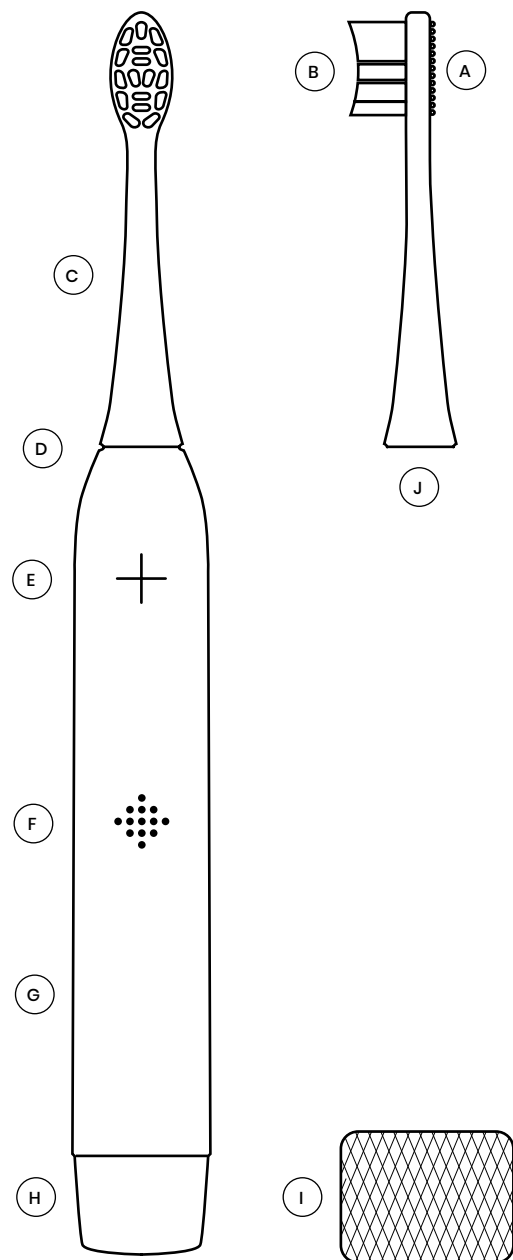INTRODUCTION

Thank you for buying the Prism Sonic Toothbrush.

Prism is developed to offer the most convenient and efficient care for your teeth. The vibrations start as soon as you put pressure on the brush and stop when you release the pressure. Just brush, smile and stay Lucent.

Of course the toothbrush is waterproof, so you don't have to worry about getting it wet. Maybe even brush in the shower.

It's important to always use a fresh brush head. Dentists recommend to change at least every 3 months. For best performance we recommend changing as soon as the bristles start bending.

Visit be-lucent.com for replaceable brush heads.

PRODUCT GUIDE/SPECS

- (A) Tongue scraper on the backside
- (B) Soft dupont bristles
- (C) Light weight and flexible head for smoother brush
- (D) Water proof (IPX8)
- (E) One-button control
- (F) Lucent LED screen
- (G) Smooth grip
- (H) Light Indicator
- (I) Wireless charging stand
- (J) Spare brush head

HOW TO USE PRISM

LOCK/UNLOCK

Lock and unlock the toothbrush
Press and hold the button for 3 seconds

ACTIVATING THE SONIC BRUSHHEAD

Start brushing

Put pressure on the brush or press the button

Stop brushing

Release pressure on the brush or press the button

THE PRESSURE SENSOR

If you brush too hard, Prism will blink and buzz, to avoid hurting your gum and enamel

CHANGE MODES

Press the button for 2 seconds when the vibrations are on to toggle between modes.

Clean ^{24 000 RPM}

High sonic frequency vibration for the daily brush

Sensitive ^{22 800/27 600 RPM}

Soft strength to protect the gum

White ^{22 800/25 800/27 600 RPM}

Pulsating vibration to remove plaque and discolouration

AUTOMATIC 2-MINUTE TIMER

VIBRATES AFTER 30 SEC

Prism has a time to guide you for the best brushing routine. Slowly move the brush towards all sides of your teeth. Focus on one quarter section of your row of teeth for 30 seconds each. Prism will do a short vibration every 30 seconds and stop automatically after 2 minutes.

CHARGING

Charging time 8 hours - lasts 60 days

Charging stand - the toothbrush stand offers wireless charging. Just plug in the cable into a USB jack and put the toothbrush inside. The bottom of the toothbrush will light up in red while charging.

INDICATIONS

The led display will indicate the battery level.

CHARGING INDICATIONS

100%

40%

Low on battery

SPECIFICATIONS

Waterproof grade:

IPX8

Rated voltage / power:

3.7V /2.5W

PRECAUTION

- Don't expose to heat above +60° degrees or cold under -20° degrees Celsius
- Please use a 5V 1A adapter for the USB connection
- If the charger is damaged, please don't use it anymore
- Do not try to detach, deform nor repair the device yourself.
- Keep the toothbrush away from children under 3 years
- For more information or customer support, please visit be-lucent.com

ONE YEAR LIMITED WARRANTY

All Be Lucent products are designed and manufactured to meet high expectations and to provide the best customer experience. However, sometimes, products may have malfunctions due to manufacturing defects. That is why we offer a full one year warranty from the date of original purchase, which covers defects in materials and workmanship.

The replacement warranty is valid only if the original receipt issued to the customer, specifying the date of purchase, is presented with the product to be replaced.

If you believe this product is defective, carefully repack the unit and return the product to your authorized seller within the warranty period, with the original receipt. Your authorized seller will replace the product, provided a defect in materials or workmanship is identified.

This limited replacement warranty does not apply in cases of abuse or misuse of the product, use contrary to the Lucent instructions, ordinary wear and tear, incorrect connection, force majeure or unauthorized repair.

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