PhysioPedal[®]

Plug-In

- EN User Manual
- ES Guía del usuario
- DE Benutzerhandbuch
- FR Mode d'emploi

by **Nobol[®]**

**Congratulations on your
new PhysioPedal**

Safety instructions	4
What's in the box	5
Specifications	6
Features	7
Remote Control / Display	8
Remote Control Buttons	9
Getting Started	10
Automated Mode	11
Manual Mode	12
Cleaning and Storage	13
Warranty	14

*Please read this user manual carefully and keep it for future reference.

The latest version of this user manual can be found on www.nobol.com

Safety Instructions

When using The PhysioPedal, basic precautions should always be followed, including the following:

Read all instructions before using The PhysioPedal.

Warning

To Reduce the risk of burns, fire, electric shock, or injury to persons:
Disconnect the charger from the cycle when not in use, and before putting on or taking off parts.

Use The PhysioPedal only for its intended use as described in this manual.
Do not use attachments not recommended by the manufacturer.

- Never operate The PhysioPedal if it's damaged, or not working properly, if it has been dropped or damaged, or dropped into water.
- Return The PhysioPedal to a service center for examination and repair.
- Never drop or insert any object into any opening.
- Do not use it in the rain.
- Do not operate where aerosol (spray) products are used or where oxygen is administered.
- For household use only.

What's in the box

A Body

B Base

C Wireless Remote

D Charging Cable

E Pedals x 2

F Handles x 2

G Base Screws x 4


H Wrench

Specifications


Product Model: P1011	Net Weight: 8.61 Kg. 19 Lbs.
Rated Power: 72W	Gross Weight: 11.25 Kg. 24.8 Lbs.
Product Size: 20.2 x 10.6 x 14.6 in.	Body: 5.4 Kg. 12 Lbs.
Output: 24V=3A	Handles: 0.3 Kg. 0.7 Lbs.
Input: 110~220V AC	Comfort Pedals: 1.6 Kg. 3.6 Lbs.

Features

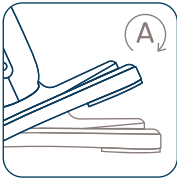
Tilt-sensitive safety feature

10° angle triggers alarm and auto-stops motor.

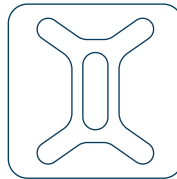Ergonomic Base

Enhanced stability, lightweight design

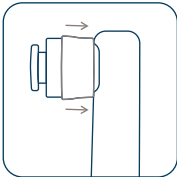
Quick Release Pedal

Pull the shaft and easily insert or release the pedal.


Remote Control

Senior friendly Intuitive Interface.

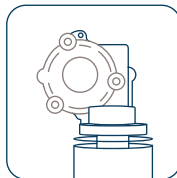Motor


Medical grade DC electric motor.

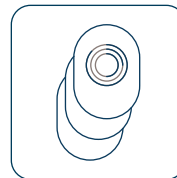Interactive Screen

Tracking and reporting of all usage, calories, distance, time used.

Remote Control / Display

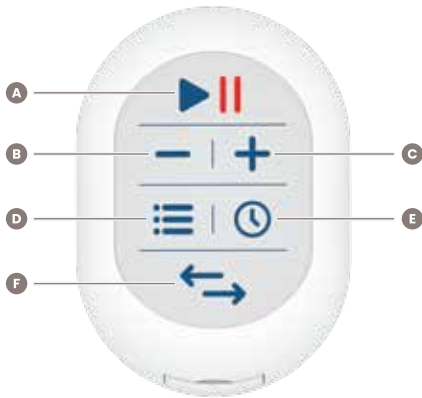- A Start/Stop
- B Decrease Speed
- C Increase Speed
- D Menu
- E Timer
- F Switch Direction

- A Rotations per minute
- B Timer
- C Timer Countdown
- D Calories
- E Steps

Remote Control Buttons

Press the **START/STOP** button once to start the Pedals or Handles to move.

Press **START/STOP** button again to stop.

The unit will start a slow speed of 30 RPM, to increase the speed press the **PLUS** button.

To decrease the speed press the **MINUS** button.

To set the timer, press the **CLOCK** button.

The PhysioPedal will start at 10 minutes.

Each press will increase the time by 10 minutes up to 90 minutes total.

If the remote control stops working, you may need to reconnect the remote control to your PhysioPedal. Turn the PhysioPedal **ON**  then, on the remote control, press and hold the **Switch Direction** button  for 3 seconds until you hear a beep.

Getting Started


Assembly

To use the PhysioPedal, use the **Tool (pic1)** to screw the **4 Screws** to attach **the base to the body** of the PhysioPedal.

Attach **the Pedals or Handles** into the Cranks. Switch from Pedals to Handles with a simple action, simply pull the knob back and slide the Pedals or Handles in or out (**pic2**).


Plug the **Power Cord (pic3)** into the Front of the Unit, screw until hand tight.

Put the **Remote Control (pic4)** over the magnetic screen.

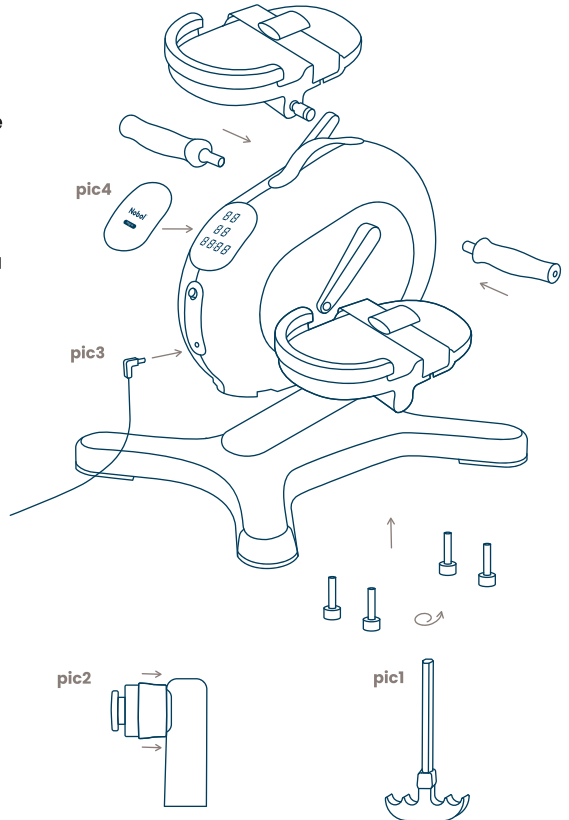
Automated Mode


To power on the PhysioPedal, make sure it is plugged into an electric outlet. Press the **Power Button** on the front of the unit to turn it **ON**.

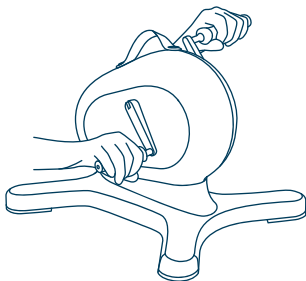To exercise the **Lower Body** sit comfortably and put both feet on the pedals.

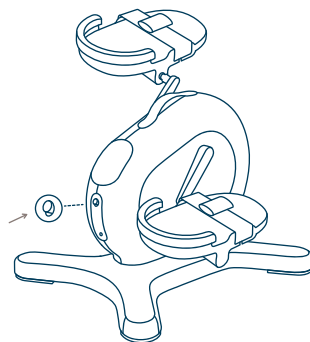
To exercise the **Upper Body** put the PhysioPedal on an elevated platform or table and put your hands on the grips.

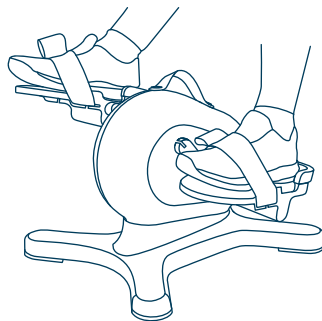Manual Mode

PhysioPedal can be used in manual mode by simply not turning on the motor, **simply Sit and Pedal**, this will NOT harm the motor.

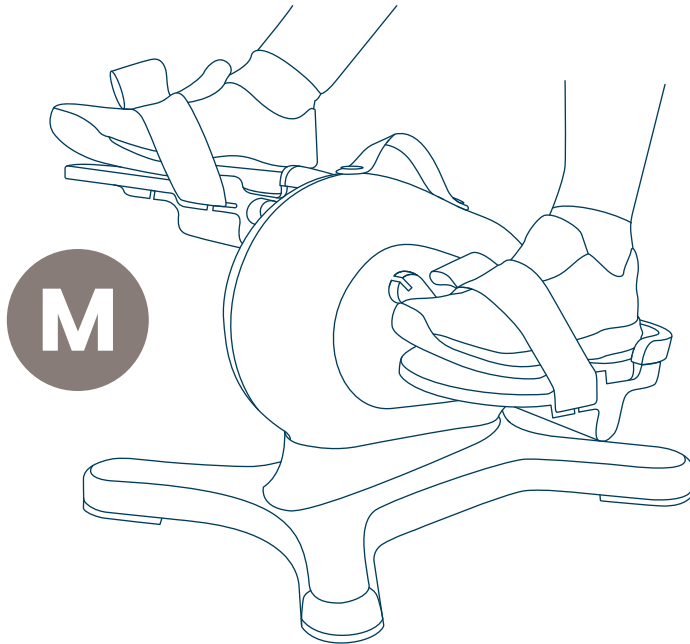Cleaning and Storing

Make sure the power is **OFF** and the battery charger is not attached. Use a damp, clean cloth and gently wipe your PhysioPedal. Store in a clean, cool, dry place, not in direct sunlight, when not in use. Any other servicing should be performed by an authorized service representative.

Warranty

Register on [Nobol.com](https://nobol.com) or contact our customer service by emailing us at cs@nobol.com to start your warranty period.

If within one year from the date of purchase, this device fails due to a defect in material or workmanship, Nobol Inc. will repair or replace the device or necessary components, free of charge.

This warranty excludes:

- A. Damage caused by accident, abuse, mishandling, or transport
- B. Units subjected to unauthorized repair
- C. Damage exceeding the cost of the product
- D. Deterioration of the delivered product resulting from abnormal storage and/or safeguarding conditions on the client's premises
- E. Products taken apart or repaired by anyone other than an authorized service representative
- F. Failure to provide the dated proof of purchase and serial number on your device
- G. Products purchased from an unauthorized retailer
- H. Please do not try and repair the unit yourself, as this will void your warranty.

*Some states and countries do not allow a limitation of damages, so the foregoing limitation may not apply to you. This warranty guarantees specific legal rights, but other rights may vary between countries and from state to state. This warranty service is available through the operating policies and procedures of Nobol Inc., and all Nobol international distributors. If you are a client, please contact the appropriate Nobol business entity dealer, or retailer, from whom you directly purchased the device for warranty and return authorization procedures. This warranty is initiated and executed by Nobol Inc.

**For more information,
please visit our website
www.Nobol.com**

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Assembled in China, Designed with ♥ in California by **Nobol, Inc.**