

# PhysioPedal<sup>®</sup>

by **Nobol<sup>®</sup>**



Effortless Exercise

# IMPORTANT SAFETY INSTRUCTIONS

When using The PhysioPedal, basic precautions should always be followed, including the following:

Read all instructions before using The PhysioPedal.

## WARNING

To Reduce the risk of burns, fire, electric shock, or injury to persons:

Disconnect the battery pack from the cycle when not in use, and before putting on or taking off parts.

Use The PhysioPedal only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.

## SAFETY INSTRUCTIONS

---

- ☐ Never operate The PhysioPedal if it's damaged, or not working properly, if it has been dropped or damaged, or dropped into water.
- ☐ Return The PhysioPedal to a service center for examination and repair.
- ☐ Never drop or insert any object into any opening.
- ☐ Do not use it in the rain.
- ☐ Do not operate where aerosol (spray) products are used or where oxygen is administered.
- ☐ For household use only.

## SAVE THIS BOOKLET

---

## SPECS



---

Power: 72W (Max)

Display: Digital

Battery Capacity: 2400mAh

Weight: 7.5Kg | 16Lbs

---

Battery Life: up to 3 hours

Output: 26V=1A

Input: 110~220V AC

Speed: 20~70 RPM

## PhysioPedal



LED Screen Display



Remote Control




Portable Battery



Oversized Comfort Pedal

## FEATURES


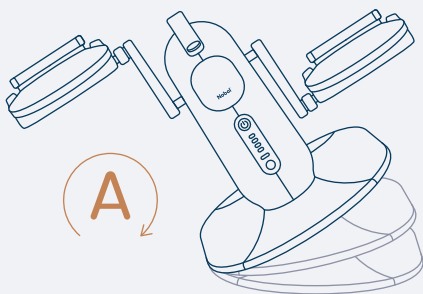


### Built in safety

If the unit tilts it automatically stops.

### Battery


High capacity rechargeable lithium-Ion battery.

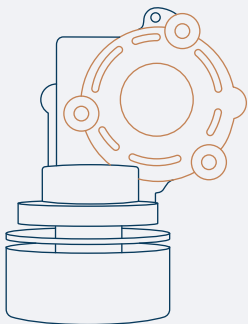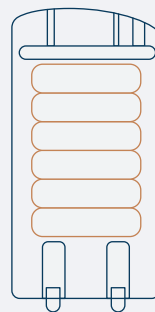

### Motor

Medical grade brushless DC electric motor.

## FEATURES


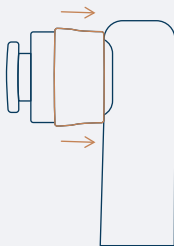


### Quick Release Pedal

Pull the shaft and easelly insert or release the pedal.

### Interactive Screen


Full tracking and reporting of all usage, calories, distance and time used.

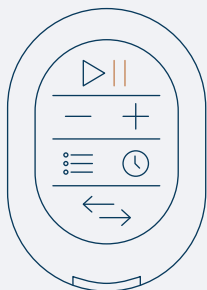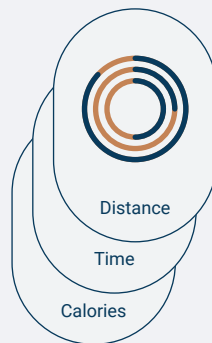

### Remote Control


Senior friendly Intuitive Interface.

DISPLAY

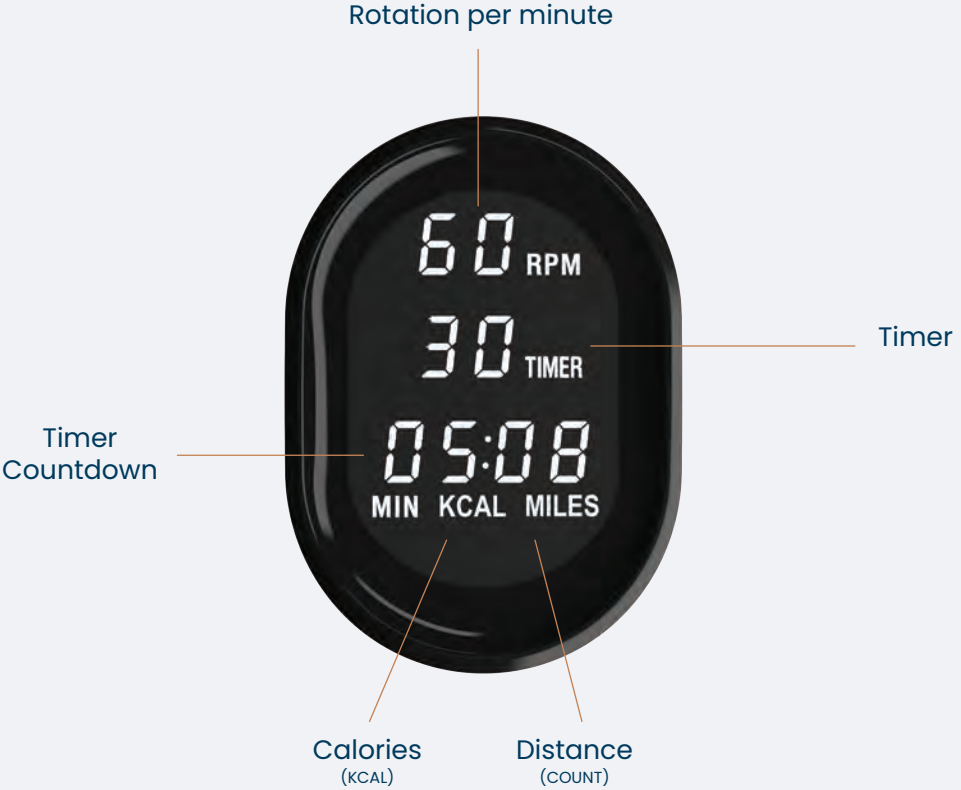



## REMOTE

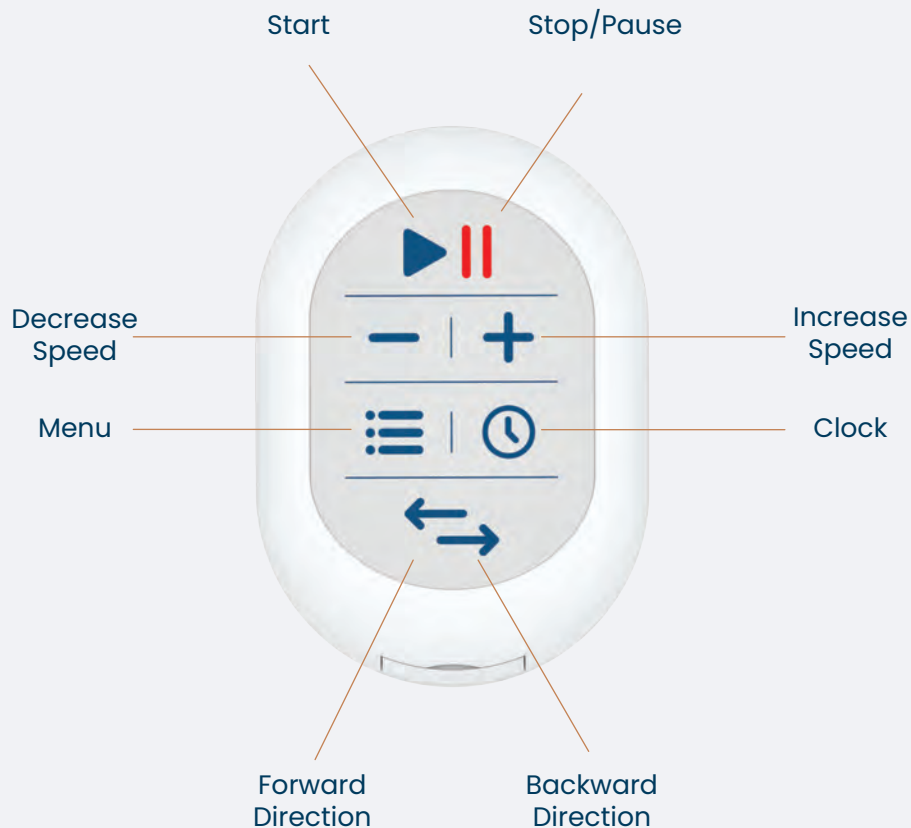

## BATTERY

lithium ion Battery Powered. Included with all PhysioPedals is the powerful rechargeable battery. Lasting for up to 8 hours on a single charge



Green lights Indicate  
how much charge is left

Red indicates low charge

Portable battery bank

Charging port



## CLEANING & STORING

---

Make sure the power is OFF and the battery charger is not attached. Use a damp, clean cloth and gently wipe your PhysioPedal. Store in a clean, cool, dry place, not in direct sunlight, when not in use. Any other servicing should be performed by an authorized service representative.

# PhysioPedal<sup>®</sup>



[www.nobol.com](http://www.nobol.com)

**Nobol<sup>®</sup>**  
*Choose Health*