

## ELECTRIC TOOTHBRUSH MANUAL

**HI THERE!  
THANK YOU FOR CHOOSING THE  
LUMINEUX® ELECTRIC TOOTHBRUSH.  
WE'RE SO HAPPY YOU'RE HERE.**

**There are several things that make this  
Electric Toothbrush so fantastically  
imperative for your oral health routine.  
Here are our top 3 reasons to use this  
Toothbrush every day:**

① The bristles of this Toothbrush are unlike those we have seen in other electric toothbrushes on the market. Five years ago, we introduced our first (manual) toothbrush. Using that experience, and after three years of research, we are excited to introduce the Lumineux® Electric Toothbrush. This electric toothbrush uses a replaceable head made with bamboo and bristles created to be extraordinarily soft for a deep clean that is super gentle on your enamel and gums. Our CEO and Founder, Dr. Maddahi, was tired of seeing unnecessary gum recession due to improper brushing or using the wrong toothbrush; he noticed even **more** gum recession with his patients that used an electric brush. This was the motivation behind his research into, and development of, a soft, flexible brushhead that cleans while helping to minimize damage to the gum. THAT is what makes the bristles of this Toothbrush so very very special.

② You may have noticed the head of the Toothbrush is made with bamboo. We chose bamboo because it is a renewable material. So, while we hope you'll keep the handle and charger of your Lumineux® Electric Toothbrush for years to come, when you replace the heads, you'll be making a better choice for our planet.

③ Electric toothbrushes in general are just better than manual ones for your oral health. When used properly, they offer a deeper clean than any manual toothbrush could dream of!

Given these facts, we're sure you're excited to use this awesome Electric Toothbrush. But before you start, please read this manual carefully and keep it somewhere safe for future reference. **If you encounter any problems during use, please contact our service center so we can help you!**

## HOW TO USE



• Carefully insert the toothbrush head onto the top of the Toothbrush body. A small gap will remain between the top of the handle and the head. This gap allows the toothbrush to vibrate correctly.

• You only need a pea sized amount of Toothpaste.



• Wet the toothbrush head before putting it into your mouth. Press the switch for 0.5 sec and choose your preferred mode (you will see each mode indicated on the handle).

• The brush bristles should be at a 45-degree angle to the gums.



Move slowly from one tooth to another tooth without going back and forth excessively. Spend about 30 seconds on each quadrant of your mouth.

## CLEANING/CARE & DISPOSAL

**Brush head:** After each use, rinse head under running water with unit switched on. Turn unit off, allow to air dry.

**Handle:** Remove brush head, run handle under water, allow to dry.

**Charging Block:** Follow instructions provided with your charging block.


**Charging Stand:** Unplug before cleaning. Wipe with a damp cloth.

**Replace brush head:** At least every 3 months. Replace sooner if brush head appears worn or bristles shed. Only use brush heads (original and replacement) from Lumineux®.

Contains batteries and/or recyclable electronic components. Do not dispose in household waste. Dispose/recycle entire unit in accordance with local environmental regulations.

## MORE ON BRUSHING

Lumineux® was founded by a renowned dentist, Dr. Kourosh Maddahi. He suggests brushing for at least 2 minutes each time you brush. During the process, your Lumineux® Electric Toothbrush will briefly stop vibrating every 30 seconds. This is to remind you to move to the next quadrant of your mouth to achieve the best cleaning results.

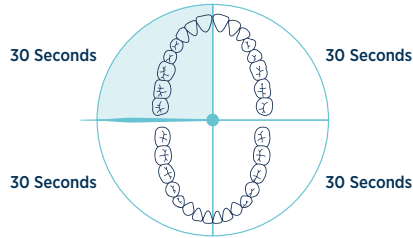

## MORE ON YOUR TOOTHBRUSH

- Brushing indicator timed for 2 minutes (30 seconds per quadrant)
- Induction charging
- 3-6 weeks of battery life
- 5 unique brushing modes for you to choose from
- Super soft bristles

**Memory function:** The toothbrush chip will remember the last brushing mode and it will start in that mode the next time it switches on.

**Cleaning mode:** Press gently and continuously on the switch button for 3 seconds to go into cleaning mode. After 3 seconds it will switch off.

**During normal operating process, press the button for more than 3 seconds and the mode intensity will switch between soft-medium-strong.**

## MODE INSTRUCTIONS

**Clean mode:** Helps remove plaque and tartar build up. Great for everyday use.

**Whiten mode:** Best paired with our Lumineux® Whitening products, this mode will help visibly whiten teeth.

**Polish mode:** This mode is designed to help smooth the surface of teeth.

**Gum care:** This mode is specifically designed to help care for your gums.

**Sensitive:** For people with sensitive gums and teeth, this mode has reduced stimulation, while still providing a dentist-approved cleaning.

## Warnings

- Place charger safe distance from water (e.g., sink). If charger falls in water, unplug immediately; do not reach into water if charger is plugged in.
- Do not use in shower or while bathing.
- Firmly insert plug into outlet; do not leave portion of plug exposed.
- When unplugging, pull the plug not the cord.
- Do not use if cord or plug is damaged. Avoid bending, jamming or squeezing the cord.
- Keep away from heated surface.
- Do not open and disassemble handle. If handle does not recharge outside warranty period, replace entire unit.
- Intended for individual use, not by multiple users.
- For ages 12 and up; not intended for use under the age of 12. Keep out of reach of children. Teenagers (age 12+) should begin use under supervision of adult to assure proper use.
- Only use with a properly attached, undamaged brush head (original or replacement) from Lumineux. Use of an improperly attached/broken brush head, or a third-party brush head, may generate small parts that can create a choking hazard for children or damage teeth.
- Use this product only for its intended use as described in this manual.
- If you are undergoing treatment for any oral care condition, consult your dental professional prior to use.
- In the first days of using any electric toothbrush, including ours, some users may experience slight gum bleeding. Typically, bleeding should stop after a few days. Should this persist, or should use cause discomfort, stop, and consult a dental health professional.
- A USB cord is provided; and is intended for use with a charging block (charger) that you provide. Since chargers do differ, consult instructions/guide provided with your charger. Use provided USB cable only with a USB-A charging block. Using a damaged cable or charger (e.g., frayed cord, damaged charger enclosure, damaged or bent plug blades), or charging when moisture is present, can cause fire, electric shock, injury or damage to your toothbrush or other property. Do not expose charger to liquids, which can cause a short circuit and overheating. Do not charge toothbrush if charger or cable gets wet. Do not attempt to dry charger or cable using an external heat source. Do not use the charger in wet locations, such as near a sink, or in or near a bathtub or shower. Never touch the charger with wet hands.
- Do not store or expose your charger to temperatures above 45°C (113° F); this may damage the charger, overheat the battery in the toothbrush or pose a risk of fire. To prevent damage to your toothbrush or charger, or their parts or internal circuits, do not use or store it in dusty, smoky, damp, or dirty environments. Keep toothbrush and charger away from heat sources and out of direct sunlight. If your toothbrush or charger becomes too hot, disconnect from power source if it is plugged in, move it to a cooler place and do not use until it has cooled.
- This toothbrush and its accessories are not toys. Children have been strangled on cords and cables. The toothbrush or charger may contain small parts or parts with sharp edges that may cause injury or create a choking hazard.
- Dispose toothbrush and charger (provided by you) in accordance with local regulations. Do not dispose in normal household waste. Improper disposal may lead to fire, explosion and/or other hazards.

## Storage and Charging

1. Before first use OR if you haven't used for a while, please fully charge.
2. For optimal performance store at room temperature (5-30°C; 40-85°F). During use or charging, temperature may increase; this is normal.
3. Charging time: 15 hours.
4. Charging Mode: Mode light will flash and turn off after 3 seconds; charging will then begin.
5. Charging Complete: The GREEN light will come on and stay on.
6. Low Battery Alert: The YELLOW light will flicker. Press the on/off button, if the yellow light goes off after 2 seconds the battery needs recharging.
7. Battery Information: 5 LED diodes: 90+%; 4 lights: 75+%; 3 lights: 50+%; 2 lights: 25+%; 1 light: 15+%; Repeating yellow light: below 10%.

## Product Specifications

**Name:** Intelligent acoustic wave electric toothbrush frequency  
**Model:** ETB-B01  
**Main power:** 2.5W  
**Battery:** Lithium-ion battery 18500/1200mAh 3.7V  
**Charging seat:** Input 5V DC 500mA

## Trouble Shooting

Condition	Possible Reason	What To Do
Vibration reduction	1. The head has been forced against the teeth while brushing. 2. Brush head is too close to the main body.	1. Place the head gently on the teeth instead. 2. Take the head out and replace on the body but leave a gap between the body and the head.
Working hours short even when the toothbrush is fully charged	1. The body was not properly placed on the charging base during charging. 2. Low battery.	1. Not enough charging time. 2. The product has not been used for over 6 months. 3. The battery has reached its usage limit.
Toothbrush doesn't work at all	The toothbrush has not been correctly charged or fully charged.	1. Charging required. 2. Contact the after-sale service to replace the battery.

## Information on Returns/Exchanges

Thank you for choosing our Toothbrush. This product should not be opened or tampered with! Refer all service and repairs to Lumineux® Customer Care. This Toothbrush has a 1-year warranty on the charger and the body. This warranty does not include the brush head. The warranty will not be valid if:

1. The Toothbrush has not been operated or maintained in accordance with Lumineux's Operating Instructions, AKA the information in this manual.
2. The Toothbrush has been used for purposes other than its intended function (uses not described in this manual).
3. The damage or malfunction of the Toothbrush is caused by any of following conditions:
  - Incorrect voltage
  - Accidents (including liquid or moisture damage)
  - Misuse or abuse of the product
  - Power supply problems including power spikes or lightning damage
  - Tampering or modification of the Toothbrush by persons other than authorized service personnel
  - Exposure to corrosive materials
  - Insertion of foreign objects into the unit
  - Use with third-party parts

Due to continuous product development, specifications are subject to change without prior notice.

## Warranty

The following limited warranty applies only to the new products. Oral Essentials Inc. "OEI" warrants the Lumineux® Electric Toothbrush for 1 year from date of purchase. Subject to the exclusions in this manual, defects due to faulty materials or workmanship will be repaired or replaced by OEI, upon receipt of proof of purchase within the applicable warranty period and a completed warranty card through [lumineuxhealth.com](http://lumineuxhealth.com)

TO THE EXTENT PERMITTED BY APPLICABLE LAW, IN NO EVENT SHALL OEI OR ITS DISTRIBUTORS BE LIABLE TO YOU OR ANY THIRD PARTY FOR ANY DIRECT, INDIRECT, INCIDENTAL, CONSEQUENTIAL, SPECIAL, EXEMPLARY OR PUNITIVE DAMAGES, OF ANY NATURE WHATSOEVER, ARISING OUT OF THE USE OF (OR INABILITY TO USE) LUMINEUX® ELECTRIC TOOTHBRUSH INCLUDING, WITHOUT LIMITATION, INJURY TO PERSONS, DAMAGE TO PROPERTY, LOSS OF VALUE OR USE OF ANY THIRD PARTY PRODUCT WITH THIS PRODUCT, EVEN IF OEI HAS BEEN ADVISED OF THE POSSIBILITY OF SUCH DAMAGES. IN NO EVENT SHALL ANY RECOVERY AGAINST OEI EXCEED THE ACTUAL PRICE PAID FOR PURCHASE OF LUMINEUX® ELECTRIC TOOTHBRUSH.

To register and/or make a claim on the warranty of your Lumineux® Electric Toothbrush, please visit [www.LumineuxHealth.com/warranty](http://www.LumineuxHealth.com/warranty)