User's Manual

CareSens® N
Blood Glucose Monitoring System

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CareSens® N

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Welcome to the CareSens® N Blood Glucose Monitoring System

Thank you for choosing the CareSens® N Blood Glucose Monitoring System. The system provides you with safe, fast, and convenient blood glucose *in vitro* (i.e., outside the body) diagnostic monitoring. You can obtain accurate results in just 5 seconds with a small (0.5 µL) blood sample. This booklet contains important information about the CareSens® N Blood Glucose Monitoring System. Please read it carefully before testing your blood glucose.

Manufacturer: i-SENS, Inc.

43, Banpo-daero 28-gil, Seocho-gu, Seoul, 06646, Korea

www.i-sens.com

Contact Customer Service:

1-800-429-5001

9 am - 9 pm EST, Monday - Saturday

At all other times or in case of emergency, please contact your healthcare professional or emergency medical response.

- No part of this document may be reproduced in any form or by any means without the prior written consent of i-SENS.
- The information in this manual is correct at the time of printing. However, i-SENS reserves the right to make any necessary changes at any time without notice as our policy is one of continuous improvement.

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Table of Contents

Information

Important Information: Read This First!	4
Important Health-Related Information	5
Specifications	6
CareSens N Blood Glucose Monitoring System	7
Inserting or Replacing the Batteries	8
CareSens N Blood Glucose Test Strip	9
CareSens N Blood Glucose Meter	11
CareSens N Blood Glucose Meter Display	12

Preparation

Setting Up Your System	13
Adjusting the Date and Time	13
Setting the Sound On/OFF	16
Setting the 'Test Result Reset'	17
Checking the System	18
Control Solution Testing	19
Comparing the Control Solution Test Results	21

Blood Glucose Testing

Using the Lancing Device	22
Preparing the Lancing Device	23
Preparing the Meter and Test Strip	24
Flagging Post-meal Test Results	24
Applying the Blood Sample	26
Discarding Used Lancets	28
Alternative Site Testing	29
HI and Lo Messages	32
Target Blood Glucose Ranges	33
Transferring Test Results	33

Additional Functions

Meter Memory	34
Viewing Test Results Stored in the Meter's Memory	34
Setting the Alarm Function	36
Setting the Post-meal Alarm (PP2 Alarm)	36
Setting the Time Alarms (Alarm 1 - 3)	37

Maintenance

Caring for Your System	39
Understanding Error and Other Messages	40
General Troubleshooting	42
Performance Characteristics	43
Warranty Information	46

Important Information: Read This First!

To receive safe and optimum system benefits, please read the entire manual contents before using the system. Please note the following instructions:

- Do not use the system for diagnosis or screening of diabetes.
- Use only fresh capillary whole blood samples for testing.
- Alternative site and fingertip test results may differ significantly due to rapid changes in the blood glucose level after eating, insulin injection or exercise.

Intended Use

The CareSens N Blood Glucose Monitoring System is intended for the quantitative measurement of glucose in fresh capillary whole blood samples drawn from the fingertips and alternative sites such as the forearm, palm, thigh, and calf. Alternative site testing should be used only during steady-state blood glucose conditions. The CareSens N Blood Glucose Monitoring System is intended for self testing outside the body (*in vitro*) by people with diabetes as an aid to monitor the effectiveness of diabetes control. It is not intended for use on neonates and is not for the diagnosis or screening of diabetes.

The CareSens N Blood Glucose Test Strips are for use with the CareSens N Blood Glucose Meter to quantitatively measure glucose in fresh capillary whole blood samples drawn from the fingertips and alternative sites.

The CareSens Control Solutions are for use with the CareSens N Meter and CareSens N Test Strips to check that the meter and the test strips are working together properly and that the test is performing correctly.

Important Health-Related Information

- The CareSens N Blood Glucose Monitoring System is intended for self-testing outside the body (*in vitro* diagnostic use).
- The CareSens N Blood Glucose Meter is designed to minimize code related errors in monitoring by using the no-coding function.
- Glucose in blood samples reacts with the chemical in the test strip to produce a small electrical current. The CareSens N Meter detects this electrical current and measures the amount of glucose in the blood sample.
- The CareSens N Blood Glucose Meter should be used only with the CareSens N Test Strips.
- An abnormally high or low red blood cell count (hematocrit level over 60% or below 20%) may produce inaccurate results.
- Inaccurate results may occur in severely hypotensive (having low blood pressure) individuals or patients in shock. Inaccurate low results may occur for individuals experiencing a hyperglycemic-hyperosmolar state, with or without ketosis. Critically ill patients should not be tested with blood glucose meters.
- Inaccurate results may occur in patients undergoing oxygen therapy.
- For questions or concerns, contact Customer Service: 1-800-429-5001 (9 am - 9 pm EST, Mon - Sat). At all other times or in case of emergency, please contact your healthcare professional or emergency medical response.

Specifications

Product specifications	
Measurement Range	20 - 600 mg/dL
Sample size	Minimum 0.5 µL
Test time	5 seconds
Sample type	Fresh capillary whole blood
Calibration	Plasma-equivalent
Assay method	Electrochemical
Battery life	5,000 tests
Power	Two 3.0 V lithium batteries (disposable, type CR2032)
Memory	250 test results
Size	93 X 47 X 15 mm (3.66 x 1.85 x 0.59 in)
Weight	51.5 g (1.82 oz) (with batteries)

Operating ranges	
Temperature	50 - 104°F
Relative humidity	10 - 90%
Hematocrit	20 - 60%

CareSens N Blood Glucose Monitoring System

CareSens N System includes the following items:

- ① CareSens N Blood Glucose Meter
- ② User's Manual
- ③ Quick Reference Guide
- ④ Batteries (2)

CareSens N System may include the following items:

- ① CareSens N Blood Glucose Test Strips (10)
- ② Lancing Device
- ③ Lancets (10)
- ④ Logbook
- ⑤ Carrying Case

- Check all the components after opening the CareSens N Blood Glucose Monitoring System package. The exact contents are listed on the main box.
- The CareSens Control Solutions (Control A & B) are necessary but not included. The CareSens Control Solutions and the cable for data management software can be ordered separately. Please contact Customer Service: 1-800-429-5001 (9 am - 9 pm EST, Monday - Saturday).


Inserting or Replacing the Batteries

The CareSens N Meter comes with two 3.0 V lithium batteries. Before using the meter, check the battery compartment and insert batteries if empty.

When the  symbol appears on the display while the meter is in use, the batteries should be replaced as soon as possible. The test results may not be saved if the batteries run out completely.



Step 1

Make sure the meter is turned off. Push the cover in the direction of the arrow to open the battery compartment.

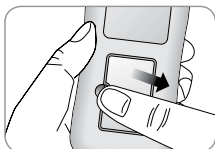
Step 2

Remove the old batteries one by one by lifting with the index finger and pulling it out with your thumb and index finger as shown in the figure on the right. Insert two new batteries with the + side facing up and make sure the batteries are inserted firmly.

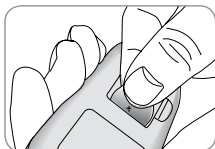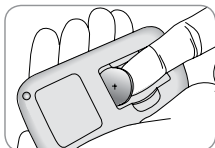Step 3


Place the cover on the battery compartment. Push down until you hear the tab click into place.

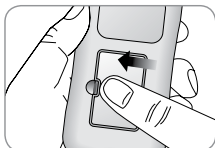Note: Removing the meter batteries will not affect your stored result. However, you may need to reset your meter settings. See pages 13-17.

CareSens N Blood Glucose Test Strip


The CareSens N Blood Glucose Monitoring System measures blood glucose quickly and accurately. It automatically absorbs the small blood sample applied to the narrow edge of the strip.

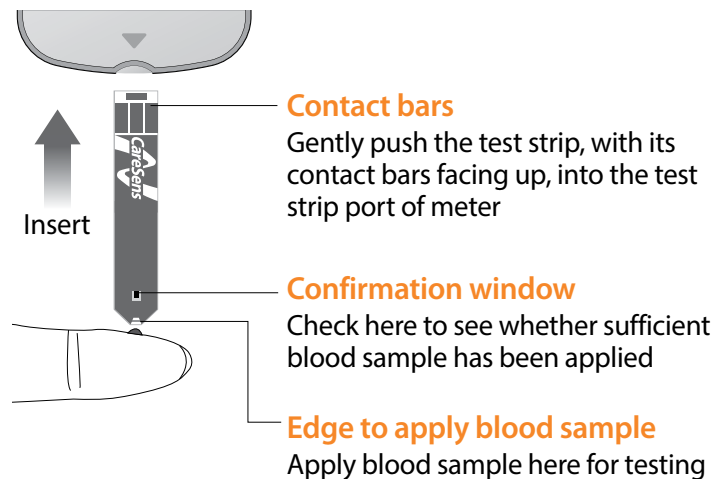CareSens N Blood Glucose Meter


Important Safety Instruction

- The CareSens N Test Strip should be used only with fresh capillary whole blood samples.
- Do not reuse test strips.
- Do not use CareSens N Test Strips beyond the expiration date. This may cause inaccurate results.
- Test strips in new, unopened vial and test strips in vial that have been opened can be used up until the expiration date printed on the test strip box and vial label if the test strips are used and stored according to its storage and handling methods.
- Store test strips in a cool and dry place at a temperature of 34-86°F (1-30°C).
- Keep test strips away from direct sunlight or heat and do not freeze.
- Store test strips only in their original vial.
- Close the vial tightly after taking out a test strip for testing and use the strip immediately.
- Handle test strips only with clean and dry hands.
- Do not bend, cut, or alter test strips in any way.
- For detailed storage and usage information, refer to the CareSens N Test Strip package insert.
- Keep the meter and testing supplies away from young children.
- The test strip vial contains drying agents that are harmful if inhaled or swallowed and may cause skin or eye irritation.

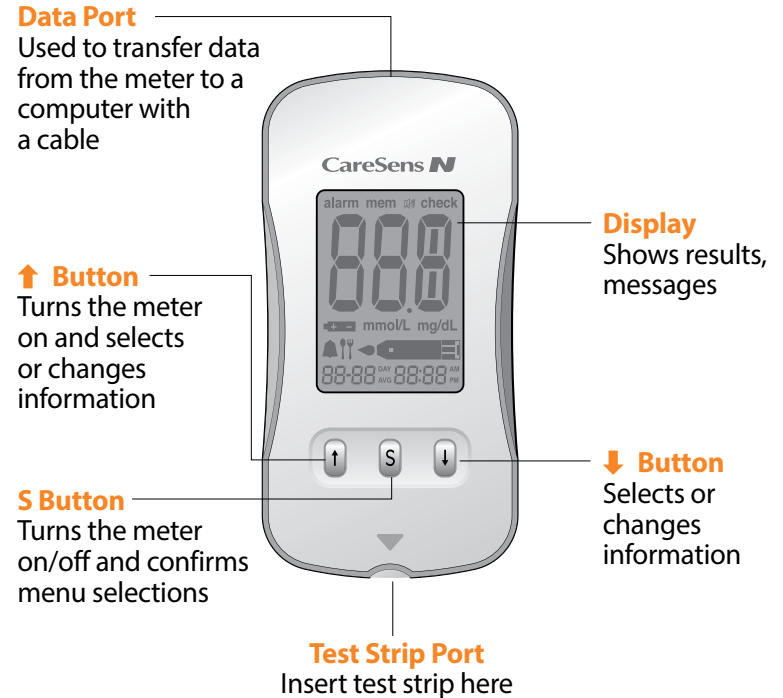Note: The cable for data management software can be ordered separately. Please contact Customer Service: 1-800-429-5001 (9 am - 9 pm EST, Monday - Saturday).

CareSens N Blood Glucose Meter Display

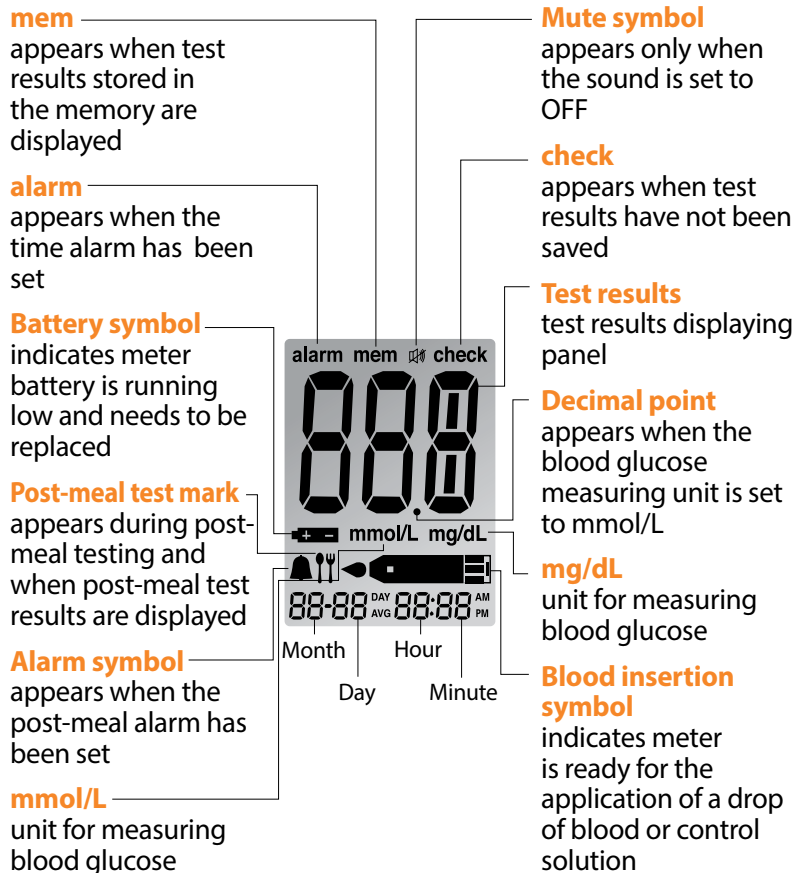Note: The unit of measurement for your CareSens N Meter is preset and fixed in mg/dL.

Setting Up Your System

You should check and update the meter-settings such as time and date before using your meter or after changing the meter batteries. Press and hold the **S** button for 3 seconds to enter the SET mode. While you are setting the meter, current settings will appear on the display. Press and hold **↓** to scroll faster. Whenever you want to complete the settings and exit from the SET mode, press and hold the **S** button for 3 seconds again.

Adjusting the Date and Time

Step 1 Entering the SET Mode

Press and hold the **S** button for 3 seconds. After all the segments flash across the screen, SET is displayed on the screen. Press the **S** button again to progress to the next step.


Step 2 Setting the Year

Press **↑** or **↓** button until correct year appears. Press and hold **↓** to scroll faster. After setting the year, press the **S** button to confirm your selection and progress to the next step.

Step 3 Setting the Month

A number indicating the month will blink on the bottom left corner of the screen. Press **↑** or **↓** button until the correct month appears. Press the **S** button to confirm your selection and progress to the next step.

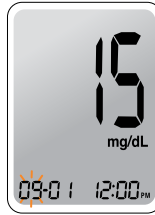
Step 4 Setting the Date

Press **↑** or **↓** button until blinking number shows the correct date. Press the **S** button to confirm the date and progress to the next step.

Step 5 Setting the Time Format

The meter can be set in the AM/PM (12-hour) or the 24-hour clock format. Press **↑** or **↓** button to select a format. The AM/PM is not displayed in the 24-hour format. After selecting the format, press the **S** button to progress to the next step.

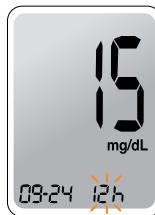Step 6 Setting the Hour

Press **↑** or **↓** button until the correct hour appears. Press the **S** button to confirm the hour and progress to the next step.

Step 7 Setting the Minute

Press **↑** or **↓** button until the correct minute appears. After setting the minute, press the **S** button to progress to the next step.


Setting the Sound On/OFF

Step 8

On pressing \uparrow or \downarrow , the screen will display the On or OFF. Press the **S** button to confirm the selection. The meter will beep in the following instances if set to On.

- When you push the **S** button or \uparrow button to turn on the meter
- When the test strip is inserted in the meter
- When the blood sample is absorbed into the test strip and the test starts
- When the test result is displayed
- When you push the \uparrow button to set the post-meal (PP2) alarm
- When it is time for a preset blood glucose test

If the sound is set to OFF, none of the sound functions will work. After setting the sound, press the **S** button to progress to the next step.

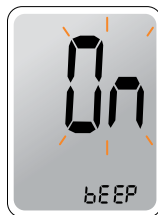Note:  symbol is displayed only when the sound is set to OFF.

Setting the 'Test Result Reset' (Deleting all the saved test results)


Step 9

In this mode all the test results stored in the meter can be deleted. Please note that if you select YES, all the stored test results will be deleted and can not be restored.

After the beeper mode is set, press the **S** button to enter the 'Test Result Reset' mode. The 'dEL' character will blink on the screen. Press \uparrow or \downarrow to alternate between 'YES' or 'no'.

To delete all the stored test results press the **S** button while the screen displays 'YES'. Then, all the test results stored in the meter will be deleted and the screen will be similar to the picture on the right.

If you do not want to delete the results press the **S** button while the screen displays 'no'. Then, the screen will return to step 2 (See page 13).

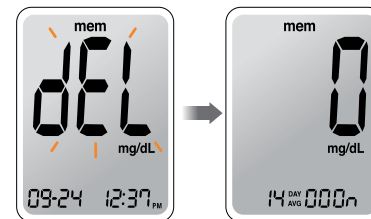Note: At any stage, if the **S** button is pressed for 3 seconds, Time and Date will finish and the meter will be turned off. Press and hold \downarrow to scroll through numbers quickly.

Checking the System

You should check your meter and test strips using CareSens Control Solution (Control A and B). The CareSens Control Solutions contain known amounts of glucose and are used to check that the meter and the test strips are working properly. The test strip vials have CareSens Control Solution ranges printed on their labels. Compare the result displayed on the meter to the control solution range printed on the test strip vial. Before using a new meter or a new vial of test strips, you should conduct a control

solution test following the procedure on page 18 with two different levels of solutions (Control A and B).

Notes:

- Use CareSens Control Solutions only.
- Check the expiration dates printed on the bottle. When you first open a control solution bottle, record the discard date (date opened plus three (3) months) in the space provided on the label.
- Make sure your meter, test strips, and control solutions are at room temperature before testing. Control solution tests must be done at room temperature (68 - 77°F, 20 - 25°C).
- Before using the control solution, shake the bottle, discard the first 1 or 2 drops and wipe the tip clean.
- Close the control solution bottle tightly and store at a temperature of 46 - 86°F (8 - 30°C).

You should do a control solution test:


- When you want to practice the test procedure using the control solution instead of blood;
- When using the meter for the first time;

- Whenever you open a new vial of test strips;
- If the meter or test strips do not function properly;
- If your symptoms are inconsistent with the blood glucose test results and you feel that the meter or test strips are not working properly;
- If you drop or damage the meter.

Control Solution Testing

Step 1


Insert a test strip into the meter's test strip port, with the contact bars facing upwards. Gently push the test strip into the port until the meter automatically turns on and the  symbol appears.

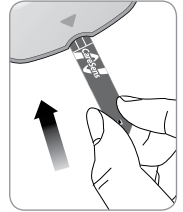Step 2

Shake the CareSens Control Solution bottle well before each test. Remove the cap and squeeze the bottle to discard the first 1 or 2 drops.

Then wipe the tip with a clean tissue or cloth. While the  symbol blinks, apply the control solution to the narrow edge of the test strip until the meter beeps. Make sure the confirmation window of test strip fills completely.

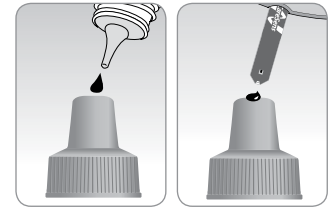Note: The meter will switch off if the control solution sample is not applied within 2 minutes of the  symbol appearing on the screen. If the meter turns off, press **S** button or **↑** button and start from step 2.


Step 3

A test result will appear after the meter counts down from 5 to 1. After your control solution result appears on the display, press **↓** for 3 seconds until the check appears on the display. When the check is displayed, the result is not stored in the meter's memory and is not included in the 14-day averages.

Step 4

Compare the result displayed on the meter to the range printed on the test strip vial. The result should fall within that range. Used strips should be discarded safely in disposable containers.

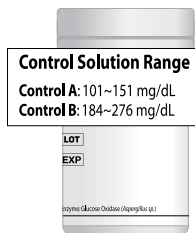Caution: The ranges printed on the test strip vial are for CareSens Control Solutions only. They do not have any connection to your blood glucose level.

Comparing the Control Solution Test Results

The test result of each control solution should be within the range printed on the label of test strip vial. Repeat the control solution test if the test result falls outside of this range. Out of range results may occur due to the following factors:


Situations	Actions
<ul style="list-style-type: none"> When the control solution bottle was not shaken well, When the meter, test strip, or the control solution were exposed to high or low temperatures, When the first drop of the control solution was not discarded or the tip of the vial was not wiped clean, When the meter is not functioning properly. 	Repeat the control solution test by referring to the "Notes" on page 18.
<ul style="list-style-type: none"> When the control solution is past the expiration date printed on the bottle, When the control solution is past its discard date (the date the bottle was opened plus three (3) months), When the control solution is contaminated. 	Discard the used control solution and repeat the test using a new bottle of control solution.

If results continue to fall outside the range printed on the test strip vial, the CareSens N Meter and CareSens N Test Strip may not be working properly. Do not use your system and contact Customer Service: 1-800-429-5001 (9 am - 9 pm EST, Monday - Saturday).

Note: The CareSens Control Solutions can be purchased separately. Contact Customer Service: 1-800-429-5001.

Using the Lancing Device

You will need a lancing device in order to collect a blood sample. You may use the lancing device contained in the CareSens N Blood Glucose Monitoring System or any other medically approved lancing device.

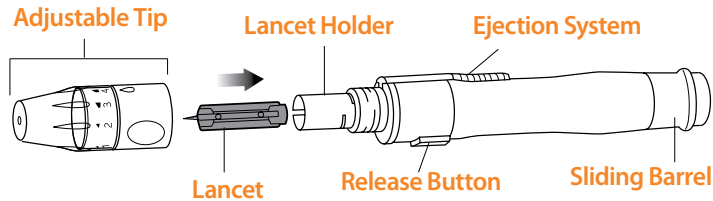- The lancing device is for use by a single user only and should not be shared with anyone including other family members.

Caution: To avoid infection when drawing a sample, do not use a lancet more than once, and:


- Do not use a lancet that has been used by others.
- Always use a new sterile lancet.
- Keep the lancing device clean.

Note: Repeated puncturing at the same sample site may cause pain or skin calluses (thick hard skin). Choose a different site each time you test.

Preparing the Lancing Device


Step 1

Wash hands and the sample site with soap and warm water. Rinse and dry thoroughly.

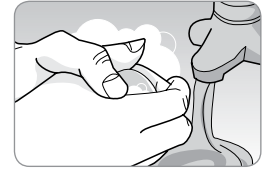
Step 2

Unscrew and remove the adjustable tip.

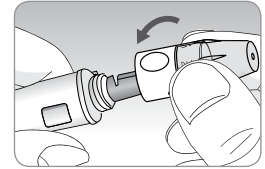
Step 3

Insert a new disposable lancet firmly into the lancet holder. Twist off the protective cover of the lancet and set it aside, then replace the adjustable tip. Keep the protective cover to replace on top of the used lancet after testing.

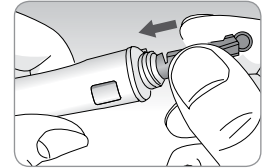
Step 4

The lancing device has five puncture depth settings, numbered 1 through 5. The smaller numbers are for a shallower puncture and the larger numbers are for a deeper puncture. Choose a depth of penetration by rotating the top portion of the adjustable tip until the setting number matches the arrow.

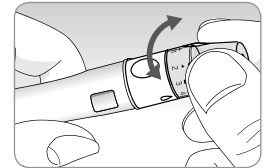Step 5

Cock the lancing device by holding the body in one hand and pulling on the sleeve with the other hand until the device clicks.

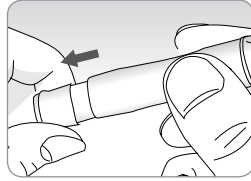
Note: The skin depth to retrieve samples will vary for various people at different sample sites. The lancing device's adjustable tip allows the best depth of skin penetration for an adequate sample size. A beginning setting of three (3) is recommended.

Preparing the Meter and Test Strip

Step 6

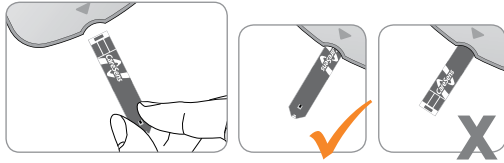Insert a test strip with the contact bars facing up into the meter's test strip port. Push the strip gently into the strip port until the meter automatically turns on and the  symbol appears.




Flagging Post-meal Test Results

CareSens N Meter allows you to flag a result of a post-meal test with  symbol. The post-meal flag () can be attached and deleted just before applying the blood sample. Once the test result is saved with the post-meal flag (), it cannot be deleted.

Step 7

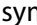If you want to attach a post-meal flag () to a test result, press and hold the  button for 3 seconds after inserting the test strip.

The post-meal flag () and the  symbol will appear on the screen. The test result will also be displayed with the post-meal flag ().

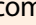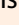
If you do not want to save the result as a post-meal test, move on to Step 8 after Step 6. If you want to delete the post-meal flag () that you attached, press and hold the  button for 3 seconds again before applying the blood sample.

Applying the Blood Sample



Step 8

Obtain a blood sample using the lancing device. Place the device against the pad of the finger. The best puncture sites are on the middle or ring fingers. Press the release button. Remove the lancing device from the finger. Wait a few seconds for a blood drop to form. A minimum volume of 0.5 microliter is needed to fill the confirmation window.

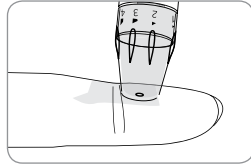Step 9

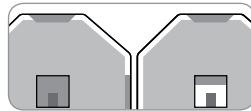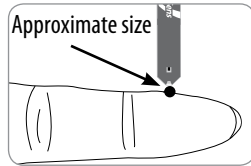After the  symbol appears on the screen, apply the blood sample to the narrow end of the test strip until the meter beeps. If the blood sample does not fill the confirmation window completely, an Er4 message may appear because of abnormal viscosity or insufficient volume.

Good Sample Insufficient Sample

Then discard the test strip, insert a new one, and apply the blood sample again after  symbol appears on the screen.

Note: The meter will switch off if the blood sample is not applied within 2 minutes of the  symbol appearing on the screen. If the meter turns off, press **S** button or **↑** button and apply the blood sample after the  symbol appears on the screen.

Step 10

The test result will appear after the meter counts down from 5 to 1. The result will be automatically stored in the meter memory.

If the test strip is removed after the test result has been displayed, the meter will automatically switch off after 3 seconds. Discard used test strip safely in a disposable container.


Alternative Site Testing

Discarding Used Lancets



Step 1

Unscrew the lancing device tip.

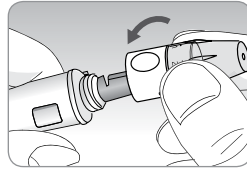Step 2

Place the protective cover on the lancet. Push the lancet ejector forward with the thumb and simultaneously pull out the sliding barrel to dispose of the used lancet in a proper biohazard container.

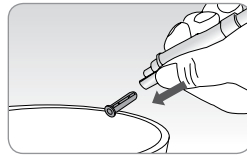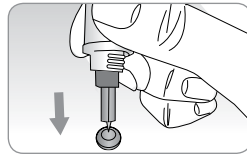Caution:


- Check for damages before using the lancet. If they have been damaged, please discard it and use other lancet.
- The lancet is very sharp. Please keep away from children.
- Keep the lancets in a cool and dry place.

Caution: The lancet is for single use only. Never share or reuse a lancet. Always dispose of lancets properly.

What is AST(Alternative Site Testing)?

When someone tests their glucose, they usually take the blood sample from the tip of the finger. However, since there are many nerve endings in the fingertip, it can be quite painful. When doing a glucose test, using different parts of the body, such as the forearms, palms, thighs, and calves, can reduce the pain during testing. This method of testing with different parts of the body is called Alternative Site Testing. While AST may reduce the pain during testing, it may not be simple for everyone and the following precautions should be observed during testing.

Alternative Sites for Testing

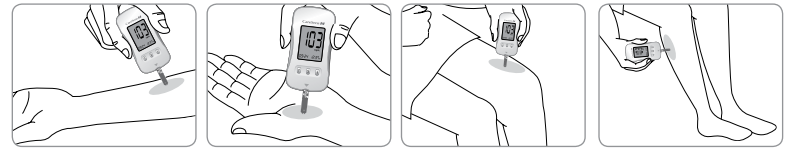Alternative Site Blood Sampling (forearm, palm, thigh, calf)

Select a clean, soft and fleshy sample site area free of visible veins and hair and away from bones. Gently massage the sample site to help blood circulation to minimize result differences between fingertip and alternative site sampling. Firmly press and hold the lancing device against site. Wait until the skin surface under the lancing device changes color. Then press the release button while continuing to apply pressure. Keep holding the lancing device against your skin until sufficient (at least 0.5 μ L) blood is drawn. Carefully lift the lancing device away from your skin.

Caution: AST results should never be used to calibrate Continuous Glucose Monitoring Systems nor entered into an insulin dosing calculator or pump for dosing recommendations.

Things to Know When Using AST

Please read the following before testing at alternative sites (forearms, palms, thighs, and calves).

The capillary whole blood of the fingertips reflects changes in glucose levels more rapidly than in alternative sites. The test results from the fingertip testing and AST may differ due to factors such as lifestyle and ingested food which affect glucose levels.

Acceptable Situations for AST

When your blood glucose levels are stable

- Fasting period
- Before a meal
- Before going to bed

Situations Requiring Fingertip Test

When your blood glucose levels are unstable

- During two (2) hours after a meal or exercise
- When sick or when glucose levels seem quite lower than test value
- When hypoglycemia is not well recognised
- When insulin has the biggest effect
- During two (2) hours after an insulin injection

AST Precautions

- Do not ignore the symptoms of hyperglycemia or hypoglycemia.
- When the results of the test do not reflect the way you feel, retest using the fingertip. If the fingertip result still does not reflect the way you feel on retesting, please consult your healthcare professional.
- Do not rely on the AST results for changing your treatment method.
- The amount of glucose in alternative sites differs from person to person.
- Before using AST, please consult your healthcare professional.

Note: Results from alternative site and fingertip samples may differ from each other, as there is a time lag for the glucose levels to reach the same value. Use a fingertip sample if you suffer from hypoglycemia or have experienced hypoglycemic shock or symptoms.

Note: If the sample drop of blood runs or spreads due to contact with hair or with lines on your palm, do not use that sample. Try puncturing again in a smoother area.

HI and Lo Messages

HI Message

The CareSens N Meter displays results between 20 - 600 mg/dL. HI appears when the blood glucose level is greater than 600 mg/dL and indicates severe hyperglycemia (much higher than normal glucose levels).

If HI is displayed again upon retesting, please contact your healthcare professional immediately.

Lo Message

Lo appears when a test result is less than 20 mg/dL and indicates severe hypoglycemia (very low glucose levels).

If Lo is displayed again upon retesting, please contact your healthcare professional immediately.

Note: If messages for hyperglycemia or hypoglycemia are displayed even if you do not have those conditions, please contact Customer Service: 1-800-429-5001 (9 am - 9 pm EST, Monday - Saturday).

Target Blood Glucose Ranges

Reminders

Time of day

Your target ranges
from your healthcare professional

Before breakfast

Before lunch or dinner

1 hour after meals

2 hours after meals

Between 2 a.m. and 4 a.m.

Expected Values: The range of a normal fasting* blood glucose level for non-diabetic adults is between 70 ~ 99 mg/dL. Two (2) hours after a meal, the range of a normal blood glucose level for non-diabetic adults is less than 140 mg/dL.

*Fasting is defined as no caloric intake for at least eight (8) hours.


Reference

American Diabetes Association. "Standards of Medical Care in Diabetes – 2015." *Diabetes Care*. January 2015; 38(1):S10, S81.

Transferring Test Results

Test results stored in CareSens N Meter can be transferred from the meter to a computer using SmartLog software and cable. The 'Pc' is displayed when the data cable connects the meter with computer.

For more information, contact Customer Service at: 1-800-429-5001.

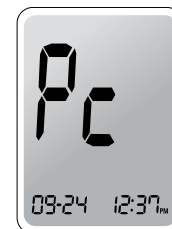Meter Memory

The CareSens N Meter can save up to 250 glucose test results with time and date. If the memory is full, the oldest test result will be deleted and the latest test result will be stored. The CareSens N meter calculates and displays the averages of total test results, pre-meal test results, and post-meal test (⏏) results from the last 14 days.

Viewing Test Results Stored in the Meter's Memory

Step 1


Press the **↑** or **S** button to turn the meter on. The current date and time will be displayed on the bottom of the screen followed by the average value and the number of the test results saved within the last 14 days.

The number of tests for the last 7 days

Step 2

Press **↓** to view the average values and the number of tests performed before eating a meal for the last 14 days. On pressing **↓** again, the average values and the number of tests performed post-meals for the same period will appear on the screen.

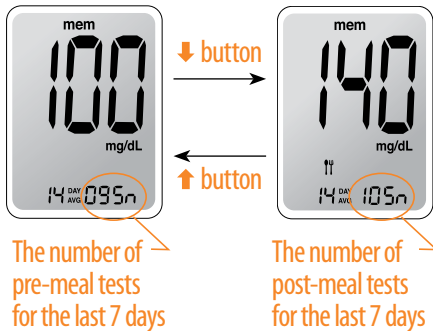Step 3

Use the **↓** button to scroll through the test results, starting from the most recent and ending with the oldest. Press **↑** to return to the result seen previously. After checking the stored test result, press the **S** button to turn off the meter.

Note: On pressing **↓**, the latest test result saved in the meter's memory will be displayed on the screen along with the date and time. Press and hold **↓** to scroll through the test results.

Setting the Alarm Function


Four types of alarms can be set in the CareSens N Meter: one post-meal alarm (PP2 alarm) and three time set alarms (alarm 1 - 3). The PP2 alarm goes off 2 hours after setting the alarm. The alarms ring for 15 seconds and can be silenced by pressing **↑**, **↓** or **S** button or by inserting a test strip.

Setting the Post-meal Alarm (PP2 Alarm)

Step 1 Setting the PP2 Alarm On


Without inserting a test strip, press and hold the **↑** button for 3 seconds to set the post-meal alarm. PP2, the **🔔** symbol and On will be displayed.

The screen will then automatically change to the memory check mode. At this time, the **🔔** symbol, indicating that the PP2 alarm has been set, will be displayed on the screen.

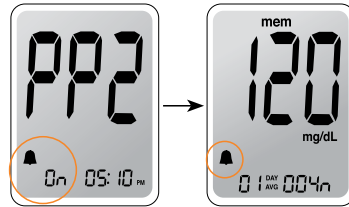Step 2 Setting the PP2 Alarm OFF

To turn off the PP2 alarm, press and hold the **↑** button for 3 seconds. PP2, the **🔔** symbol and OFF will appear on the screen. Then the screen will change automatically to the memory check mode without **🔔** symbol displayed.

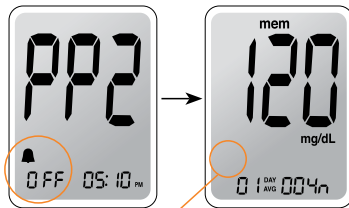🔔 symbol disappears

Setting the Time Alarms (Alarm 1 - 3)

Step 1

Without inserting a test strip, press **↑** and the **S** button simultaneously for 3 seconds to enter the time alarm mode. While OFF blinks on the screen, alarm 1 will be displayed.

Step 2

On pressing **↓**, alarm 1 is set and On is displayed on the screen. Press **↓** again to cancel alarm 1 and OFF will blink on the screen.

Step 3

Press **↑** to adjust the hour of alarm 1. A number representing the hour will blink on the screen. Press **↓** to set the hour. Press **↑** to set the minute.

Caring for Your System

Step 4

On pressing **↑**, the number indicating the minute will start blinking. Press **↓** to set the minute.

Step 5

Press the **S** button to finish and to enter the alarm 2 mode. Repeat steps 2 to 4 to set the remaining time alarms (alarm 2 and 3).

Step 6

Press and hold the **S** button for 3 seconds to finish and turn the meter off.



Use a soft cloth or tissue to wipe the meter exterior. If necessary, dip the soft cloth or tissue in a small amount of alcohol.

Do not use organic solvents such as benzene, acetone, or any household and industrial cleaners that may cause irreparable damage to the meter.

Caution: Storage and Handling

- Do not expose the meter to direct sunlight or heat for extended periods of time.
- Do not let dirt, dust, blood, or water enter into the meter's test strip port.
- Do not drop the meter or subject it to strong shocks.
- Do not try to fix or alter the meter in any way.
- Keep the meter in a cool and airy place.
- Strong electromagnetic radiation may interfere with the proper operation of this device. Keep the device away from sources of strong electromagnetic radiation.
- Store all meter components in the carrying case to prevent loss.

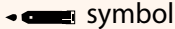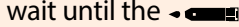Understanding Error and Other Messages

Message	What it means	What to do
	A used test strip was inserted.	Remove the strip and insert a new test strip.
	The blood or control solution sample was applied before the  symbol appeared.	Repeat the test with a new test strip and wait until the  symbol appears before applying the blood or control solution sample.
	The temperature during the test was above or below the operating range.	Move to an area where the temperature is within the operating range (50 - 104°F) and repeat the test after the meter and test strips have reached a temperature within the operating range.

Message	What it means	What to do
	The blood sample has abnormally high viscosity or insufficient volume.	Repeat the test using a new test strip.
	This error message may appear when the wrong blood glucose test strip is used instead of the CareSens N test strip.	Repeat the test with a CareSens N test strip. Insert the strip with the contact bars facing up, and push in gently until the meter turns on.
	There is a problem with the meter.	Do not use the meter. Contact Customer Service: 1-800-429-5001.

Note: If the error messages persist, contact Customer Service: 1-800-429-5001 (9 am - 9 pm EST, Monday - Saturday).

General Troubleshooting

Problem	Troubleshooting
The display is blank even after inserting a test strip.	<ul style="list-style-type: none"> • Check whether the test strip is inserted with the contact bars facing up. • Check if the strip has been inserted completely into the test strip port. • Check if the appropriate test strip was used. • Check if the batteries are inserted correctly. • Replace the batteries.
The test does not start even after applying the blood sample on the test strip.	<ul style="list-style-type: none"> • Check if the confirmation window is filled adequately. • Repeat the test after inserting a new test strip.
The test result does not match the way you feel.	<ul style="list-style-type: none"> • Repeat the test after inserting a new test strip. • Check the validity period of the test strip. • Perform control solution test.

Note: If the problem is not resolved, please contact Customer Service: 1-800-429-5001 (9 am - 9 pm EST, Monday - Saturday).

Performance Characteristics

The performance of CareSens N Blood Glucose Monitoring System has been evaluated in laboratory and clinical tests.

Accuracy: The CareSens N BGM System is calibrated to yield results equivalent to plasma glucose concentrations. The accuracy of the CareSens N System (Model GM505NA, GM505NB, GM505NC) was tested by comparing blood glucose results obtained by patients with those obtained using a YSI Model 2300 Glucose Analyzer, a lab instrument. The results below were obtained by diabetic patients at clinic centers.

	Obtained by healthcare professionals	Obtained by lay users
Slope	1.023	1.032
Y-intercept	-1.86 mg/dL	-3.63 mg/dL
Correlation coefficient(r)	0.995	0.985
Number of sample	110	105
Range tested	30 - 485 mg/dL	67 - 413 mg/dL

Accuracy results for glucose concentration < 75 mg/dL

Within ± 5 mg/dL	Within ± 10 mg/dL	Within ± 15 mg/dL
14/17 (82%)	17/17 (100%)	17/17 (100%)

Accuracy results for glucose concentration ≥ 75 mg/dL

Within ± 5%	Within ± 10%	Within ± 15%	Within ± 20%
54/93 (58%)	82/93 (88%)	93/93 (100%)	93/93 (100%)

Precision: Precision studies were performed in a laboratory using the CareSens N BGMS.

<i>Within Run Precision</i>		
Blood avg.	38 mg/dL	SD = 1.9 mg/dL
Blood avg.	86 mg/dL	SD = 3.2 mg/dL
Blood avg.	125 mg/dL	CV = 4.1%
Blood avg.	189 mg/dL	CV = 2.6%
Blood avg.	334 mg/dL	CV = 2.8%

<i>Total Precision</i>		
Control avg.	35 mg/dL	SD = 1.2 mg/dL
Control avg.	116 mg/dL	CV = 2.4%
Control avg.	324 mg/dL	CV = 2.3%

This study shows that there could be variation of up to 4.1%.

Alternative Sites Testing Evaluation

- Comparison of results using various AST sites with the results of YSI measurements.

Alternative site test results for glucose concentration < 75 mg/dL

	Within ± 5 mg/dL	Within ± 10 mg/dL	Within ± 15 mg/dL
Fingertip	1/2 (50.0%)	2/2 (100%)	2/2 (100%)
Forearm	1/2 (50.0%)	2/2 (100%)	2/2 (100%)
Palm	2/2 (100%)	2/2 (100%)	2/2 (100%)
Thigh	1/2 (50.0%)	2/2 (100%)	2/2 (100%)
Calf	1/2 (50.0%)	2/2 (100%)	2/2 (100%)

Alternative site test results for glucose concentration ≥ 75 mg/dL

	Within ± 5%	Within ± 10%	Within ± 15%	Within ± 20%
Fingertip	27/48 (56.3%)	46/48 (95.8%)	48/48 (100%)	48/48 (100%)
Forearm	19/48 (39.6%)	42/48 (87.5%)	47/48 (97.9%)	48/48 (100%)
Palm	28/48 (58.3%)	38/48 (79.2%)	46/48 (95.8%)	48/48 (100%)
Thigh	17/48 (35.4%)	36/48 (75.0%)	46/48 (95.8%)	48/48 (100%)
Calf	27/48 (56.3%)	46/48 (95.8%)	48/48 (100%)	48/48 (100%)

Warranty Information

Manufacturer's Warranty

i-SENS, Inc. warrants that the CareSens N Meter shall be free of defects in material and workmanship in normal use for a period of 5 years. The meter must have been subjected to normal use. The warranty does not cover improper handling, tampering, use, or service of the meter. Any claim must be made within the warranty period.

The i-SENS, Inc. will, at its discretion, repair or replace a defective meter or meter part that is covered by this warranty. As a matter of warranty policy, i-SENS will not reimburse the consumer's purchase price.

Obtaining Warranty Service

To obtain warranty service, you must return the defective meter or meter part along with proof of purchase to your nearest i-SENS Authorized Warranty Station.

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© 2015 i-SENS, Inc. All Rights Reserved.