

DISTRIBUTED BY AND © 2020 THREESIXTY  
SOURCING LTD. ALL RIGHTS RESERVED.  
FOR US ONLY: DISTRIBUTED BY AND © 2020  
MERCHSOURCE, LLC, IRVINE, CA, 92618. ALL RIGHTS  
RESERVED. TEL: 1-800-374-2744.  
SHARPER IMAGE® NAME AND LOGOS (INCLUDING  
THE 77 LOGO AND TOMORROW'S TOMORROW)  
ARE TRADEMARKS OR REGISTERED TRADEMARKS.  
PRINTED IN CHINA



# FIT ROLLER

## MULTIPOINT SHIATSU MASSAGER

# SHARPER IMAGE®



Please read all instructions and warnings prior to use.  
Keep this sheet as it contains important information  
for future reference.

## CONGRATULATIONS ON YOUR PURCHASE OF THE FIT ROLLER.


## RECHARGEABLE BATTERY WARNING

This product contains batteries that are non-replaceable.  
There is a risk of fire and personal injury if battery is punctured, damaged,  
or misused.  
Do not dismantle, crush, incinerate or short-circuit battery.  
Never expose battery to extreme temperatures or moisture.  
Keep away from flammable materials and direct sunlight.  
Only charge using the provided charger, cable, and/or adapters.  
Rechargeable battery is only to be charged by an adult.  
Adult supervision is required while charging the battery.  
Allow battery to cool to room temperature before charging.  
Charge on a non-flammable surface and keep away from flammable  
products.  
Discontinue charging if battery becomes hot, swells, smokes, or gives off a  
strong odor.  
Battery chargers used with the product are to be regularly examined for  
damage to the cord, plug, enclosure, and other parts. In the event of such  
damage, chargers must not be used until the damage has been repaired.  
Recycle or dispose of battery according to federal, state, and local laws.

## FCC STATEMENT

This device complies with part 15 of the FCC Rules. Operation is  
subject to the following two conditions: (1) This device may not  
cause harmful interference, and (2) this device must accept any  
interference received, including interference that may cause  
undesired operation.  
NOTE: This equipment has been tested and found to comply with the limits for  
a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are  
designed to provide reasonable protection against harmful interference in a  
residential installation. This equipment generates, uses and can radiate radio  
frequency energy and, if not installed and used in accordance with the  
instructions, may cause harmful interference to radio communications.  
However, there is no guarantee that interference will not occur in a particular  
installation. If this equipment does cause harmful interference to radio or  
television reception, which can be determined by turning the equipment off  
and on, the user is encouraged to try to correct the interference by one or  
more of the following measures:  
- Reorient or relocate the receiving antenna.  
- Increase the separation between the equipment and receiver.  
- Connect the equipment into an outlet on a circuit different from that to  
which the receiver is connected.  
- Consult the dealer or an experienced radio/TV technician for help.  
WARNING: Changes or modifications to this unit not expressly approved by the party  
responsible for compliance could void the user's authority to operate the equipment.

## GETTING STARTED

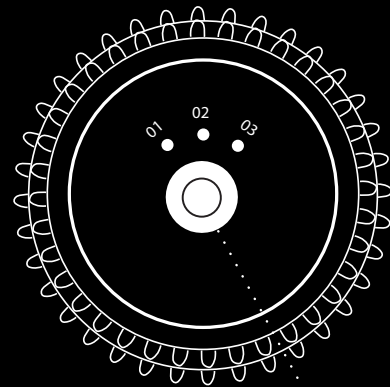


**POWER BUTTON:**  
PRESS TO TOGGLE  
THROUGH MODES

1. Fully charge the massager before first use by plugging the power adapter into the massager and then into a wall outlet.
2. Press the massager on the desired area then push the Power Button. Massage overworked muscles of the leg, arms, and feet by rolling during your warm up and cool downs at home or in the gym.
3. Choose different intensities by pressing the Power Button.
4. Use for 15 minutes at a time. Enjoy your massage!

# CONTENTS



FIT ROLLER

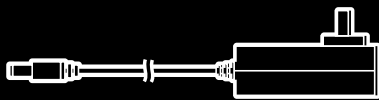


AC ADAPTER


## FEATURES

- For sore or tense muscles
- Three independent massage joints
- Four vibration levels
- Rechargeable
- Cordless


## GOOD FOR: NECK, SHOULDERS, BACK, WAIST AND THIGHS

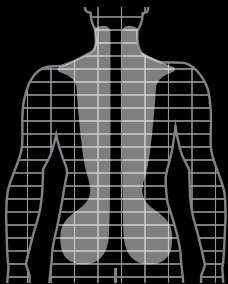


BACK

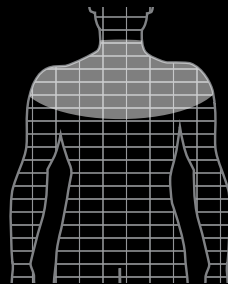

NECK + SHOULDERS

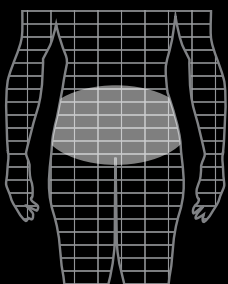

WAIST

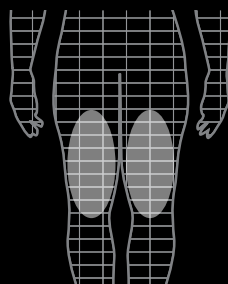

THIGHS

# CARE INSTRUCTIONS

- Spot clean only. Unplug and wipe with a damp cloth.
- Air-dry only. Never place in a dryer.
- Allow to dry fully before plugging back in.
- Do not submerge in water.
- Do not use detergent or soap when cleaning.



## RECYCLING

This product bears the selective sorting symbol for waste electrical and electronic equipment (WEEE). This means that this product must be handled pursuant to local waste electrical and electronic equipment (WEEE) directive or regulation in order to be recycled or dismantled to minimize its impact on the environment. For further information, please contact your local or regional authorities. Electronic products not included in the selective sorting process are potentially dangerous for the environment and human health due to the presence of hazardous substances.



This product comply with the applicable European directive and national regulation.



This product comply with the applicable Great Britain regulation.



This product comply with the applicable Australian/New Zealand regulation.



Batteries and accumulators must not be disposed of with household waste! Every consumer is legally obliged to hand in all batteries and accumulators, regardless of whether they contain harmful substances \* or not, at a collection point in the city or at a retail outlet so that they can be disposed of in an environmentally friendly manner. \* Cd = cadmium; Hg = mercury; Pb = lead

## DANGER



To reduce the risk of electric shock, always unplug this appliance from the electrical outlet immediately after fully charging and before cleaning.

## SAFETY WARNINGS

### PLEASE READ BEFORE USE

- Not a children's item. Keep out of the reach of children.
- Do not use this item if you are pregnant or have an injury or illness. Consult your physician if you have a medical condition that may prevent safe use of this item.
- Do not use attachments of any kind other than those that may be included.
- Do not fall asleep while using.
- Do not leave this item unattended when it is plugged in.
- Always allow item to cool if it becomes too hot.
- Never operate under blankets where air opening may be blocked. Excessive heat can occur and cause fire, electric shock, or injury to persons.
- Do not expose to extreme temperatures or moisture.
- Keep away from sources of water and heated surfaces.
- Never touch the power plug with wet hands.
- Do not carry this item by its cord.
- For indoor use only.
- Thoroughly inspect products before every use. If the product appears worn, frayed/splintered, cracked, or broken in any manner, discontinue use and discard immediately.
- The appliance is only to be used with the battery charger (OUTPUT: DC 9V 1A) provided with the appliance.