



The BLAZE.

Assembly Guide | Safety Instructions.

Questions On Assembly?

Text Us Now: (412) 423-6090

Email: Support@Recover.Red

 **RECOVERRED**





Warning



This device may cause personal injury if not adjusted or used as required by this manual. Both the operators and the technicians who use this device must be familiar with the safety requirements. The patient should also have knowledge of the basic principles. The patient should wear the protective glasses and know the precautions.

Blaze Device Introduction

The Blaze LED Red Light Therapy device is a LED light therapy device, which is suitable for skin and inflammation treatments and more. The Blaze is manufactured in our FDA registered facility. The Blaze has passed all ETL / UL Certifications for Electrical Safety Assurance. Please read these instructions carefully, and the related operators should have professional training for light therapy devices. Our LED devices are official FDA Class II medical devices. Our products are low-risk, general wellness / fitness products in accordance with the "General Wellness: Policy on Low Risk Devices" draft released January 20, 2015.

General Safety Rules



Read all instructions. Failure to follow all instructions may result in skin irritation, burns, eye damage, and/or serious injury. The term "device" in the warnings listed refers to your-operated BlazeRed Light Therapy Device. **SAVE ALL WARNINGS AND INSTRUCTIONS FOR FUTURE REFERENCE**



The safety of the patient lies on the operators, technicians and the facility this device is used at. The patient should also have knowledge of the basic principles. The patient should wear the protective glasses and know the precautions.



Do not operate device closer than 6", carefully monitor skin temperature. Maximum Permissible Exposure (M.P.E.) can be considered as dangerous if the operator is too close to the light source.



Do not stare into the Light. Staring into this device may cause burns to the eye. Some temporary color distortion can occur if eye protection is not worn.



Remove reflectors such as jewelry, watches, beauty devices or mirrors.



Clean the skin with soap and water, treatment should begin only when the skin is dry after cleansing.



Avoid any flammable items close to the device.



The Remote contains a battery. Follow all battery safety precautions. Keep batteries out of reach of children.



Do not subject the device to extreme temperatures or variations in temperature. Protect the device against moisture and direct sunlight. Device is to be used in a climate controlled environment.



Avoid Impacts or dropping of the device. Be sure when setting the device up it is not at risk to tip over. Be sure to Setup device on flat, level floor.



Device cords are only designed for device use. Do not us power cords with any other device.

Note: *There are no user serviceable parts inside. Disassembling the device will void all warranties on the product. Do not modify the device in any way. Modifying the device may result in hazardous condition.*

Disclaimer

By using this device you are agreeing to our terms of service.

A) WE ARE NOT A MEDICAL SERVICE PROVIDER, AND WE DO NOT PROVIDE ANY MEDICAL TREATMENT. We merely sell Red Light therapy devices through our website that are adjustable to any prescribed treatment level. If you require medical attention or professional advice, we would suggest you to seek such medical aid, or contact the relevant professional for the same. You must consult your physician before purchasing or using any of our products.

B) Our products are not intended to, nor do we make any claim or warranty that any of our products will diagnose, cure, or prevent any specific medical condition or disease, or will work to your expectation, or provide you any particular results. Neither any of our products, nor anything mentioned on our website shall be deemed as a substitute for professional medical advice

Limitation of Liability: To the maximum extent permitted by law, in no event shall our Company (or our directors, partners, licensors or affiliates) be liable to you or any third party for any damage, costs of procurement of substitute products, cost of medical opinions or tests or treatments or expenditures, adverse effects, illness, loss of health, medical issue, or any indirect, consequential, exemplary, incidental, special or punitive damages arising from or relating to the purchase or use of any of our products, even if our Company has been advised of the possibility of such damages. We strongly advise you to consult a registered medical practitioner before using our products. Access to, and use of, the products is at your own discretion and risk, and you will be solely responsible for any damage or loss, whatsoever.

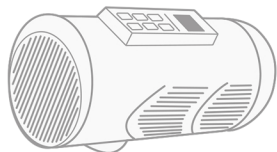

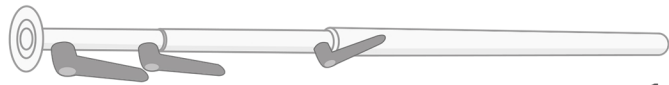
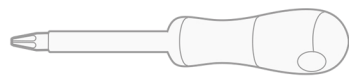



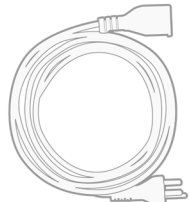
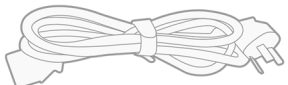

You can read more at <https://recovered.com/policies/terms-of-service>

Regular Maintenance

To ensure your device is working as best as possible we recommend inspecting the device every month.

- Check to make sure Fans are clean and not obstructed
- Check power cords and electrical connections for damage or corrosion
- Clean surface of LED lens with damp cloth or basic cleaner.







Part List

1 x LED Unit	 A cylindrical LED unit with a control panel on top and a lens on the front.	x1
Main Base	 A circular main base with a small hole in the center.	x1
Vertical Neck	 A long, thin vertical neck with a circular end and three adjustment points.	x1
Screw Driver	 A screwdriver with a pointed tip and a handle.	x1
Bracket to light screws	 Four small, cylindrical brackets used for light screws.	x4
Remote	 A rectangular remote control with several buttons.	x1
Eye protection	 Two eye protection lenses and a pair of safety glasses.	x1
Extension cord	 A coiled extension cord with a power plug and a standard outlet.	x1
Standard Power cord	 A standard power cord with a power plug and a standard outlet.	x1
Velcro Straps	 A Velcro strap used for securing equipment.	x1

INSTALL INSTRUCTIONS


STEP 1

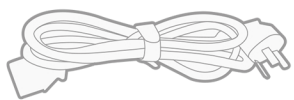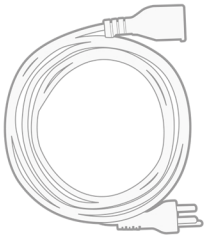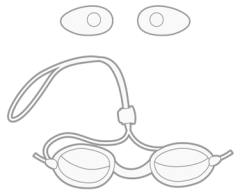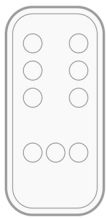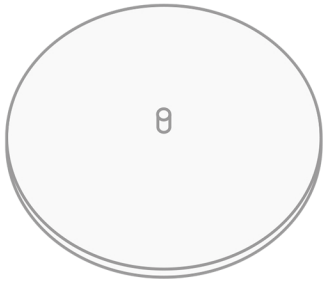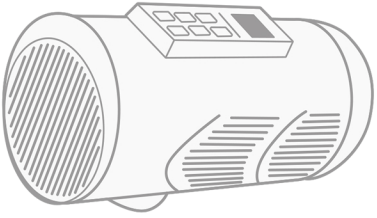Unpackage everything and locate the hardware and screws. pus the wheels into the base legs.




STEP 2

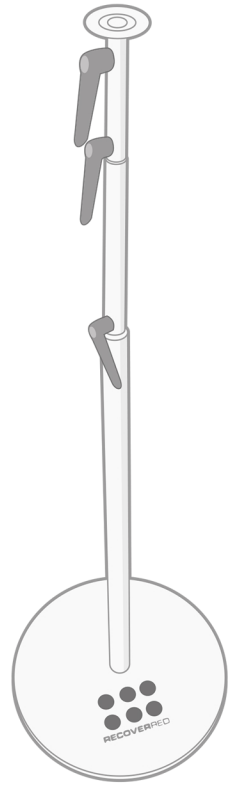Adjust the light arm straight and tighten the lever. Install the light arm onto the base by rotating onto the center bolt of the base.



STEP 3

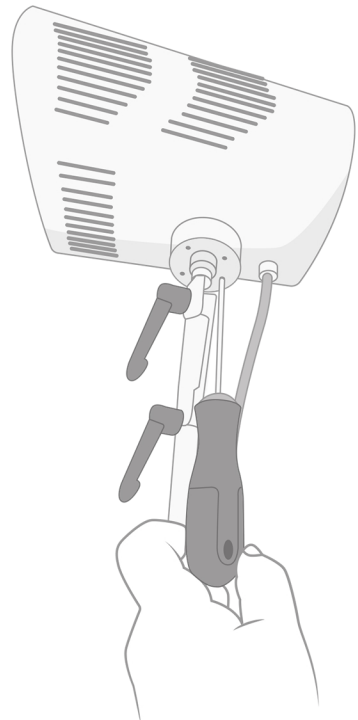
Carefully rest the light unit onto the base and attach using 4 x screws. You can rotate the light if the adjustment lever is preventing you from tightening a screw.



STEP 4

Never adjust the light while someone is underneath.

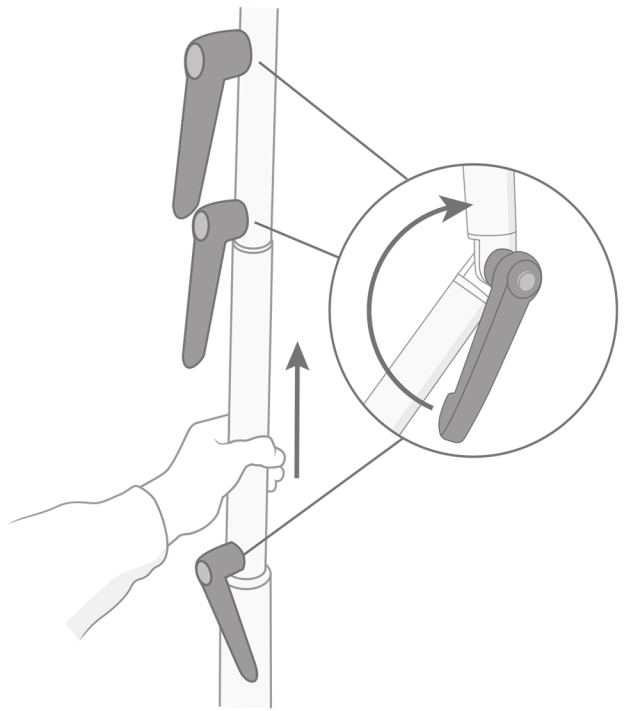Adjust the arm height as needed by lifting on the upper portion. To adjust the arm handle, slightly loosen the hand lever for that point and hold in desired position. Retighten by turning the lever clockwise, once you feel tension the hand should be turned 1 additional inch. Use the cord wrap to route the cord as needed.



More on adjustment handles:

1. Always physically check the levers and give the device a nudge after adjustment and before every use to ensure it is tight.
2. Re-position the lever for a better grip by resting thumb on the venter button and pulling the handle away from the device. You can now freely move the handle to the desired position for storage or better tightening grip.

Warning: Failure to properly secure the adjustment arms can lead to injury.

