# R1

Knee Care Device  
Instruction Manual

**KNEE JOINT CARE SOLUTION**  
Vibration Massage  
Hot Compress  
Electrical Muscle Stimulation


## R1 AT A GLANCE


- 1) Main Device
- 2) Auxiliary Device
- 3) Control Display Panel
- 4) Heating pad
- 5) Vibrators


## PREPARATION

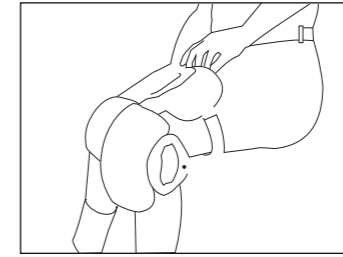1) Uses the strap to hold the main device on the knee joint, and adjusts the tightness.

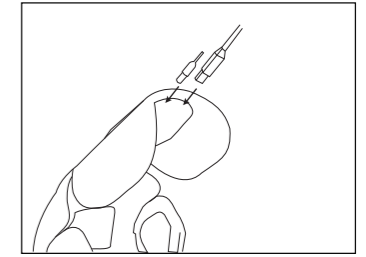2) Links the auxiliary device with the main device.

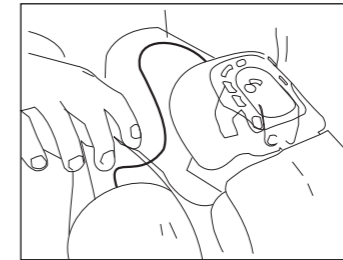3) Presses the control button to open the device.

## FEATURES

- Adjustable temperature 40°C-70°C (104°F-158°F)
- Enjoys knee massage on both legs simultaneously
- Designed to fit the knees, shoulders and elbows
- Detachable and adjustable knee straps
- EMS pads for massaging other body parts

## NOTE

- The default time for heat therapy is 30 minutes. Press the on/off button on the control panel again if you need more cycles.
- In different mode, some buttons for different functions will be locked.

## OPERATING

- 1) Press power control button to turn on the device. The default setting is Auto Mode (50 degrees Celsius, mixed strength vibration, level I EMS). Press the button again to turn off.
- 2) Press temperature control button to turn on the heating therapy. Press the temperature control buttons to adjust the temperature. The heating therapy is turned off if the temperature is lower than 40 degrees Celsius.
- 3) Press EMS button to turn on the electric impulse. Press the intensity control buttons to adjust the intensity level. Press EMS button again to toggle between different electrotherapy modes.
- 4) Press vibration massage button to toggle between low, high and mixed vibration strength. The Auto Mode will be turned off.
- 5) Press Auto Mode button to turn off the Auto mode. Press again to restart the device, and adjust temperature and EMS level..
- 6) To insert an external connector into the electrical pulse output, attach the pads to the corresponding location you would like to massage (arms, legs etc.).
- 7) Reattach pulse pads to protective paper after used.

## USE OF SILICON ELECTRODE SLICES

- Insert one end of the 2-core connector wires into the electrode port and connect the other end to the electrode slices.
- The electrode slices should be clean on their surfaces without any dust, oily or sticky substance. If their surfaces become dirty, flush with clean water and then air dry them before repeated use.
- The electrode slices can be used for around 50 times of 20 minutes session. If their viscosity cannot be recovered after several times of cleaning, please replace new ones.
- Only use the electrode slices provided with the device.
- Please dispose electrode slices according to your local environmental regulations rather than discarding them at will.

NOTE: It is recommended to use the EMS massager once or twice a day. Please give your skin a break from the electrode pads in every 20 minutes. Leaving pads on the same skin location may cause skin irritation.

## MAINTENANCE, CLEANING & STORAGE

- Store this device in a moisture-free, dry place when not in use.
- Do not wrap the power cord around the device.
- Do not use the power cord as a handle.
- Keep away from all solvents and harsh detergents.
- Unplug the device and allow it to cool before cleaning. Wipe with a soft, dry cloth as needed.
- Do not submerge the device in water or other liquids.
- Do not use abrasive cleaners, brushes, gasoline, kerosene, glass/furniture polish, or paint thinner to clean the massager.

## SAFETY INSTRUCTIONS

- 1) Only use the adapter that was supplied with the device.
- 2) Do not let the device fall from high place, and do not keep it around sharp objects.
- 3) Do not use this product when you are sleeping.
- 4) Do not use other medical electronic devices at the same time, such as pacemakers, heart-lung machines or electrocardiographs. If the patient is simultaneously using high frequency surgical equipment or other related equipment, the electrical pulse patches in the device may cause burns or damage the equipment. If the equipment is used near a microwave, the output of the device may become unstable.
- 5) Keep away from children.
- 6) When using the electric pulse function, do not move or change the position.

## DISCLAIMER

Please do not use the ELEEEELS R1 without first obtaining approval from your doctor if any of the following apply: pregnant, have skin sensitivities, heart disease, blood pressure abnormalities, malignant tumors, cerebrovascular diseases, diabetes, recent joint replacements, metal pins or plates or you have any concerns about your physical health. Do not use R1 if you feel pain, discomfort, or have a wound at or near the site of application.

## WARRANTY

IF WITHIN ONE YEAR FROM THE DATE OF PURCHASE, THIS PRODUCT FAILS DUE TO A DEFECT IN MATERIAL OR WORKMANSHIP, ELEEEELS WILL REPAIR OR REPLACE THE PRODUCT, OR NECESSARY COMPONENTS, FREE OF CHARGE

## THIS WARRANTY EXCLUDES:

- (a) damage caused by accident, abuse, mishandling, or transport;
- (b) units subjected to unauthorized repair;
- (c) units not used in accordance with ELEEEELS care instructions;
- (d) damage exceeding the cost of the product;
- (e) deterioration of the delivered product resulting from abnormal storage and/or safeguarding conditions on the client's premises, and
- (f) failure to show the dated proof of purchase

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