(BABYCARE)
WE CARE

Manufacturer: Hangzhou ChunEr Technology Co., Ltd.
Address: Room 336, Building 1, No.611 Jianghong Road, Changhe Street, Binjiang District, Hangzhou, Zhejiang, China


UKCA Name: SPEC Solutions ltd
Address: Building 1000, Lakeside North Harbour, Western Road, Portsmouth, PO6 3EZ
Email: SPEC.OFFICE02@gmail.com

EC REP: KD EUROPE
Address: Stara Debszczyzna 4, 16-424, Filipow, PL
Email: Kdeurope.office@gmail.com

T-REX TOBY ELECTRIC TOOTHBRUSH

BC2203062-1

「B22V1」

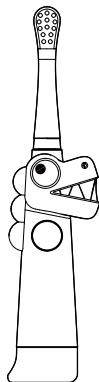Please read this manual carefully before using this product. Keep it for future reference. Failure to follow the instructions in this manual may pose a risk to the safety of your children.

PRODUCT INFORMATION

Model: BC2203062-1

Product name: T-REX TOBY ELECTRIC TOOTHBRUSH

Product size: 32*24*148.8mm

Age: 0-36M

Main materials: ABS, Silicone, PA

Rated voltage: 1.5V 

Rated power: 0.5W

Waterproof: IPX7


Vibration frequency: 12000 strokes

Battery: 1.5V "AAA LR03" (purchased separately)

PRODUCT FEATURES


1. Colorful lights and animal shapes attract babies to fall in love with teeth brushing.
2. 12000 high-frequency gentle vibrations per minute, fine and soft brush, easy to remove tooth stains.
3. IPX7 waterproof, can be washed directly for convenient usage.
4. 30 seconds reminder to change the cleaning area, 2 minutes smart timer to shut down, allow the baby to develop a scientific habit of brushing teeth.
5. Applicable to 0-year-old baby.

COMPONENT DESCRIPTION

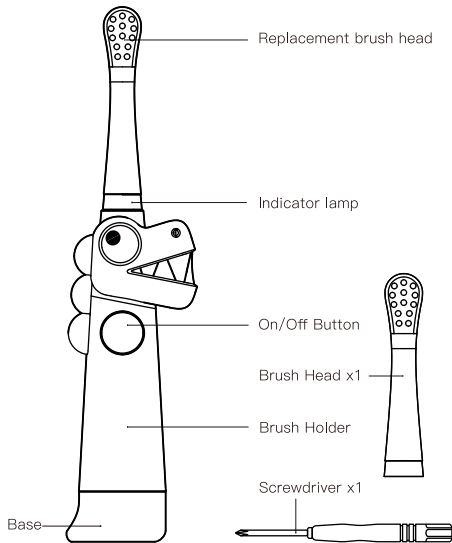HOW TO USE

1. Wet the brush and squeeze out certain amount of toothpaste.
2. Place the brush on the teeth at a small angle.
3. Press the on/off button to start the toothbrush, the colorful light will flash (about 1 second intermittent display). Then press the on/off button again to enter brushing mode, the white indicator light is on.
4. It is recommended to cover the four areas of the oral cavity (as shown in the picture), by moving the toothbrush in and out, and each take 10 seconds to brush the front, back and top of the teeth.
5. When the brushing time in one area reaches 30 seconds, the toothbrush will move intermittently 2 times to remind you to move to the next area to continue cleaning.
6. When the brushing time reaches 2 minutes, the toothbrush and indicator light will automatically turn off.
7. Rinse the brush head with flowing water, shake off any remaining water from the brush head, and wipe the handle with a towel.

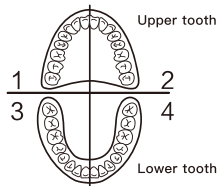(Divide the oral cavity into four areas)


HOW TO USE

KINDLY REMINDER

1. The operation sequence of the on/off button is: Colorful light mode (flashing) → brushing mode (white light is steady on) → power off.
2. When in colorful light mode, if there is no further operation within 60 seconds, the product will automatically shut down.
3. In brushing mode, press the on/off button twice (within 1 second) to turn off the product. If there is no further operation within 2 minutes, the product will automatically shut down.
4. In brushing mode, the toothbrush will move intermittently twice every 30 seconds for three times, and it automatically shuts down for the fourth time.


BRUSH HEAD REPLACEMENT METHOD

1. The brush head rotates anticlockwise to unlock, then pull to remove it.
2. Insert the positioning slot of the brush head into the bump of the host. Turn clockwise to lock and ensure installation.

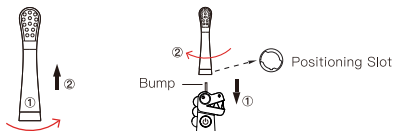
BATTERY INSTALLATION INSTRUCTIONS

Use a screwdriver to loosen the screw at the bottom of the base. Then, rotate the base clockwise and pull down.

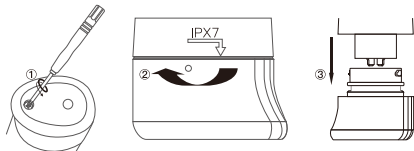
BATTERY INSTALLATION INSTRUCTIONS

With the appropriate positive and negative terminal direction, add a 1.5V "AAA LR03 " battery.

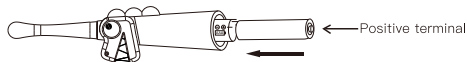Position the dot on the base in line with the arrowhead, rotate anticlockwise to reach the position. Tighten the screw in place to ensure the base cannot be pulled out.

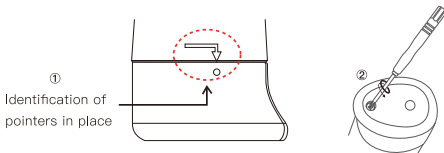NOTICE

Removing the battery is contrary to this process.
For adult installation only. For infants, please stay away from the installation site.

CLEANING AND MAINTENANCE

1. Before use, please wash the brush and ensure that the base is tightened and fixed.
2. After use, rinse the product with clean water with a temperature lower than 45°C, and place it in a ventilated and dry place.
3. Do not use steel wire, sharp instruments or corrosive detergent to clean this product.
4. Do not unscrew the base for cleaning.
5. To achieve the best cleaning effect, it is recommended to replace the brush head every three months.

WARNING

INTENDED USE: For in-home use, to promote overall oral hygiene by cleaning your teeth and gums.

1. Vigorous chewing may cause damage to the silicone cover of the toothbrush over time.
2. Babies and young children should always be supervised by a responsible adult during use.
3. Please check whether the product is damaged before use. If there is any damage, please replace the product promptly.
4. If you suspect your child has swallowed any part of the toothbrush, contact your healthcare provider immediately.

WARNING

5. This packaging is not toys, do not give them to infants to play.

6. This product is not intended for use by individuals with reduced mental, physical or sensory capabilities, or lack of sufficient experience and knowledge.

7. While battery is installed, the body shall not be immersed in water.

8. Please remove the battery when the product is not used for a long time.

9. This product is made for the functions mentioned in the instruction manual, and shall not be used for other purposes.

10. Please store it in a cool and dry place. Keep it away

WARNING

from fire and high heat source. Avoid direct sunlight.

11. Do not knock or hit the product.

12. When changing the battery, turn off the power first and dry any moisture off the body to prevent water from entering the machine. Check the positive and negative terminal of the battery, and install the battery accordingly. Make sure the base is fully tightened.

13. Do not soak the products in high-temperature water.

14. During brushing, please utilize the power of the toothbrush and do not press the teeth too hard.

15. Do not press the toothbrush on the teeth or gums too hard, or brush in the

WARNING

same position for a long time. Otherwise, it may cause damage.

16. If your oral cavity or gum has been operated on or is being treated in the past two months, please consult your dentist if you can use this toothbrush.

17. Oral bleeding may occur after the initial use of the toothbrush. If the bleeding continues for more than a week, consult the dentist's advice.

18. If your body has a heart rate regulator, or any other electronic device, consult your doctor before using this toothbrush.

19. The brush head of this toothbrush is limited to per

WARNING

sonal use and it's not recommended to share it.

20. It is recommended to change the brush head every three months.

 No modification of this equipment is allowed.

BATTERY PRECAUTIONS

1. Use only the 1.5V "AAA LRO3" battery (purchased separately).
2. The removal and installation of the battery should be correctly operated according to the positive and negative terminal and voltage mark, and the battery should not be short circuit.
3. Do not put the battery in the fire, or disassemble the battery.
4. Always follow the precautions indicated on the battery.
5. Do not use expired batteries, old batteries, or batteries that are insufficient.
6. Non-rechargeable batteries cannot be charged.
7. Only batteries of the same or equivalent to those recommended can be used.
8. The exhausted battery should be removed from the product. If not used for a long time, please remove the battery from the product.
9. Avoid use in extreme environments which would lead to reduction of battery life.
10. The battery must be removed before discarding the product.
11. When removing the battery, the appliance must be turned off.
12. The battery should be disposed of safely. Do not discard the battery as household waste.

COMMON FAULT DIAGNOSIS AND TREATMENT

During use, if you encounter a fault, in order to save time, it is recommended to check the following diagnosis before submitting for repair.

| Phenomenon | Causes | Solutions |
|-------------------------------|--|---|
| Product not working | Battery not in place | Use a 1.5V "AAA LRO3" battery |
| | Ineffective battery | Replace with new battery |
| Vibration attenuation | Wrong side of the positive and negative terminal | Install the battery according to the appropriate terminal direction |
| | Low in battery | Replace with new battery |
| Cleaning effect becomes worse | Brush head has been used for a long time | Replace the brush head |

If the product still does not work after the above inspection, please contact us. We will answer your questions as soon as possible.

Warning: Do not disassemble the appliance by yourself, improper repair may cause dangers to users!

AFTER-SALE SERVICE

Dear customer:

Thank you very much for using our products. For the product you purchased, we will offer the following replacement service, valid from the date of purchase or the date of receiving the product.

1. Replacement Policy:

1) The purchased product can be refunded or returned within 7 days without any reason (given the returned product does not affect resale).

2) If there is any quality problem after replacement, we offer to replace the product again within 15 days from the date of initial replacement.

2. The following circumstances do not fall within the scope of replacement service:

1) Fragile accessories (including but not limited to shaver blade/glass products/stewing pot/filter screen/toothbrush head/supplement solution, etc.), gifts (including but not limited to scale remover/deodorant bag/measuring cup, etc.);

2) Damage caused by irresistible natural disasters;

3) Damage caused by improper use, storage and maintenance, including but not limited to water damage and physical damage, etc.;

4) Exceeds the period of warranty defined by the replacement policy.

If you have any problems regarding the product, please contact the following email address:

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