

Hello مرحباً

Welcome to "**Lifami Family**"

You are now on your way
to better oral health

مرحباً بك في "عائلة لفمي"
أنت الآن في طريقك
لتحسين صحة الفم

SAFETY INSTRUCTIONS

DANGER

1. If you do not know the operating process and safety instructions, please do not rush to use the product to avoid hurts caused by mis-operation.
2. Do not disassemble the product or change the parts of the product to avoid damage on the product in case may caused in-valid warranty.
3. Do not use above 50°C hot water into the water tank for cleaning mouth to avoid scald.
4. If the power cable or plug is damaged, please replace it before using.
5. Do not put the product into the fire or the heating, and do not charge, use or leave it at the high temperature.

WARNING

1. Do not use this product for any other purpose other than oral cleaning.
2. Do not push the nozzle hard on the gums or teeth to avoid hurting on the gums.
3. Do not immerse this product into water for a long time to avoid incorrect operation or cause fire, electric shock and other dangers.
4. Do not damage, destroy or modify the power cable, and do not stretch, bend or twist the power cable. Do not press heavy objects on the power cable or get it caught between objects, and keep the product and cable away from hot objects.
5. Do not store the product into a hot or humid places, and avoid direct sunlight.
6. Do not charge the product in hot or humid places, such as the bathroom or toilet.
7. If not in use of the product for a long time, please unplug the power plug and charging plug to avoid fire or danger.


NOTICE

1. Do not use without a nozzle.
2. This product is not suitable for children and infants, unless they are supervised or guided by their safety responsible supervisor on how to use this product. Keep out of reach for children or infants. Do not use it as a toy.
3. If you have periodontal disease or have had oral surgery in the past 2 months, please consult your dentist before using.
4. Do not use iodine, saline (salt solution), or water insoluble concentrated essential oils in this product. Use of these can reduce product performance and will shorten the life of the product. Use of these products may void warranty.

3 PRODUCT DESCRIPTION

Parts Identification

1. Nozzle
2. Nozzle eject button
3. Charging port
4. ON/OFF Mode button
5. Pressure indicators light
6. Battery charge indicator
7. Water tank
8. Water inlet and nozzle storage

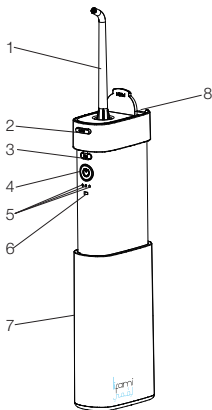

Package Inside

- . Lifami Water Flosser
- . Nozzle x 4
- . Charging cable Type-C
- . Instruction Manual & warranty
- . Bag

Nozzle Usage

A. Classic Jet Nozzle**

The everyday use

B. Perio Nozzle*

Reach crevices/problem areas
(great with PULSE mode)

C. Ortho Nozzle*

For Braces/general use

D. Tongue Cleaner Nozzle**

Fresher breath




*Replace Every 3 months

**Replace Every 6 months

GETTING STARTED


- 1 Charge the unit for **4 hours** prior first use through the charging port.

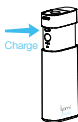

- 2 Press nozzle firmly into unit until it **CLICKS** into place. To remove nozzle, push the nozzle eject button.



- 3 Pour **warm water** into the tank and please ensure that when pouring water into the water tank, the device is switched off.

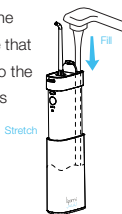


- 4 lean over sink. **Close lips** enough to prevent splashing. Keep unit upright during use to ensure continuous water flow.



- 5 Aim water at the **gumline** at a 90 degree angle. Follow the gumline and pause briefly between teeth.



- 6 Mode switch: Press the ON/OFF button to start flossing with gentle mode. **Hold for 0.5 sec** to switch to pulse mode, then hold again for 0.5 sec to switch to turbo mode.

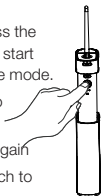


- 7 How to aim using different types of nozzles




Jet Nozzle



Perio Nozzle

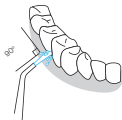

Jet Nozzle

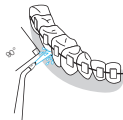

Ortho Nozzle

CLEANING

1 Nozzle

Press the nozzle button and remove the nozzle, wash it with water and then wipe dry with a cloth.

The dentist recommends to change the nozzle every 6 months.

If the top of the nozzle is deformed, replace it with a new one to avoid malfunction.



2 Product body

Remove the dirt with a cloth dampened with liquid soap, then wipe clean with a dry cloth. If the product body is dirty, wash it with water.

Do not wipe the product with thinner, gasoline, alcohol etc.

Do not wash with hot water above 50°C. Please wipe dry after cleaning.



3 Water tank

Open the tank and clean it regularly.

Wash the internal with a brush or other object and wipe with a dry cloth. When not in use of the product for a long time, please keep the water tank dry. Thank you for purchasing our product.




4 Clean

Wipe off the dirt, otherwise it will cause a leakage.



5 Nozzle Storage

After cleaning is done, store the most used nozzle in the water flosser storage.

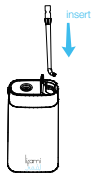

MODE INDICATION

| Mode | Mode switch | Mode indicator color | Target user | Water pressure |
|--------------------------|--|----------------------|---|----------------|
| Gentle mode | Press the switch to start flossing | Green | Gentle water pressure. Gentle cleaning. Suitable for sensitive teeth and initial time users | 40-60PSI |
| Pulse mode | Press the switch and hold for 0.5' to switch to pulse mode | Blue | Variable frequency pulse pressure alternates the strong and the weak mode | 60-110PSI |
| Turbo mode | Press the switch and hold for 0.5' to switch to strong mode | White | can deeply clean and massage the gums Strong water pressure, efficient cleaning, meeting the daily cleaning for the most of the people | 120-140PSI |
| Battery indicator | 1.The indicator turns red when run out of power. 2.The indicator flashes red when charging. 3.when fully charged, the indicator turns green. | | | |

Safety Precautions

Due to the gums are new to the product's stimulation, the gums may bleed at the initial using. (please use the Gentle mode first), the bleeding will be stopped in a week. If the situation last for more than two weeks, it may be caused by gum problem (periodontitis etc.). In this case, please stop using the product and consult your dentist.

TROUBLESHOOTING

| Problem | Possible reasons | Solution |
|--|--|--|
| Not working | Newly purchased or not used for more than three months | Please charge it for more than 4 hours |
| Low water pressure | Something inside the nozzle is blocking the nozzle Jet | Replace the nozzle |
| Battery only can be used for a short time when fully charged | The battery is at the end of its useful life | Contact your seller |
| | Insufficient charging times | It refers to recharging for more than 4 hours before using |
| The charging indicator light does not light up when charging | The power plug is not properly inserted | Please connect the power correctly |
| | Product charging circuit is faulty | Contact your seller |

Product Specifications

| | | | |
|--------------------------|-----------------|----------------------------|---------|
| Water pressure | 40---140PSI | Water tank capacity | 140ml |
| Battery capacity | 1100mAh | Noise | 60db |
| Input parameter | 5V- 1A | Rated power | 8W |
| Charging port | Type-C | Charging time | 4h |
| Waterproof level | IPX7 | Water pressure mode | 3 modes |
| Unit size (H*W*D) | 127.5*63.5*33mm | Net weight | 200g±5g |

8 WARRANTY CERTIFICATE


| | |
|---------------------------|--|
| Customer Name | |
| Country of Origin | |
| Warranty Period | |
| Invoice No. | |
| Date of Purchasing | |

Product Certification

| | | |
|---------------------------|----------------------|--|
| Product Name: | Lifami Water Flosser | |
| Model No.: | WF18 | |
| Production Time: | | |
| QC Inspector: | | |
| Inspection Result: | Qualified | |

Make water flossing a habit!

Sign up at *Lifami.com to receive exclusive discount code.

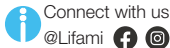

www.lifami.com