


Model B5 Quick Start

For experienced I.C.E.S. users, getting started is easy:

- 1 – Plug in the coils, up to 4 in any combination.
- 2 – Plug in the power with the supplied USB charger cable
- 3 – Switch it ON
- 4 – Let it RUN. The B5 will run automatically in 15 seconds
- 5 – The screen will look like this:
just let it run on these settings
- 6 – Use your B5 immediately without changing anything. Then over time you can learn how to completely personalize the world's most sophisticated PEMF system.

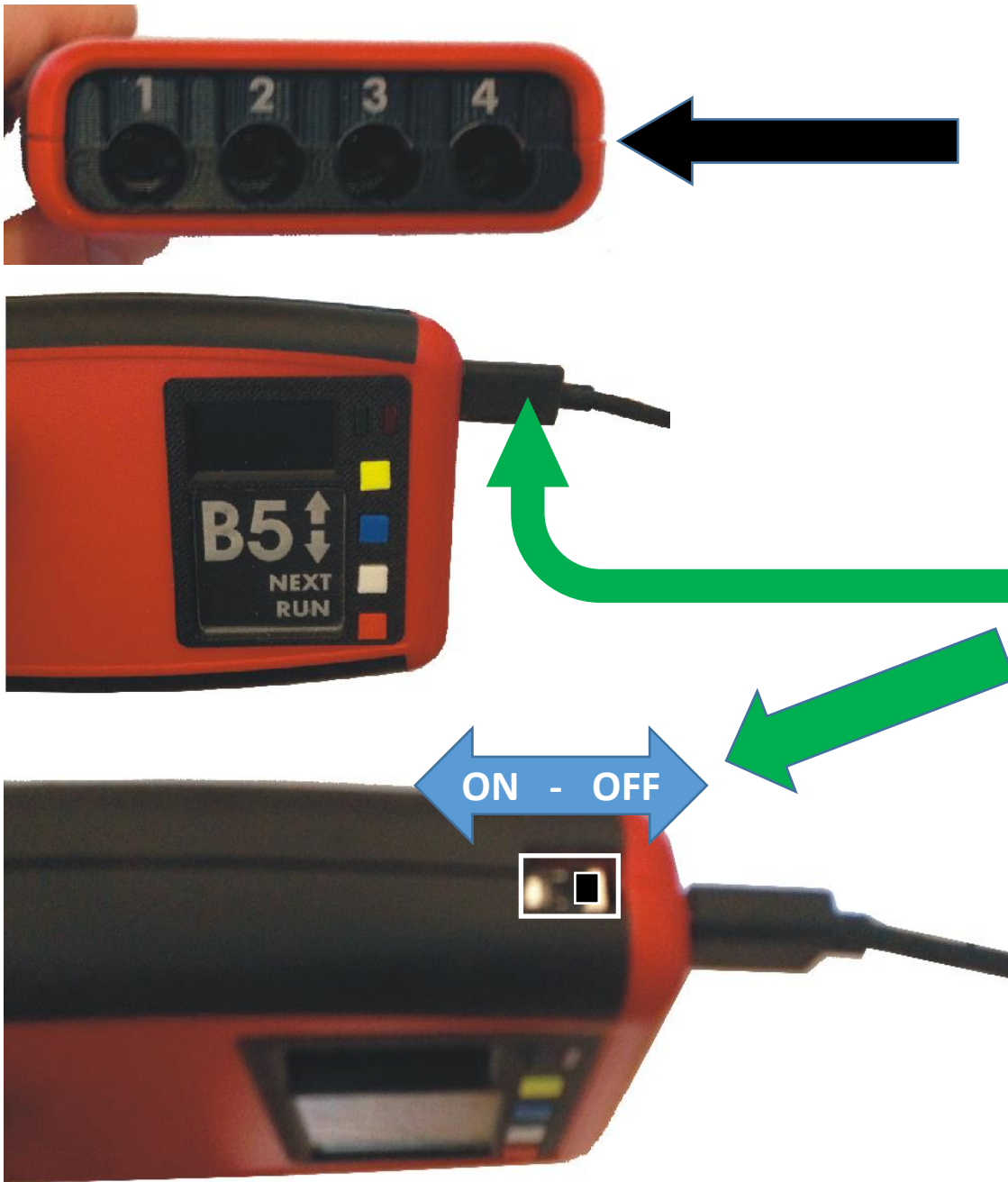

7 – Instructions at: micro-pulse.com → products → B5