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[www.thermrup.de](http://www.thermrup.de)

## Contents

1.	Important instructions.....	6
2.	Included in delivery .....	6
3.	Product features.....	7
4.	Operating.....	7
5.	Cleaning and maintenance.....	9
6.	Disposal .....	9
7.	Warranty.....	9

### 1. Important instructions

- **Non-observance of the following instructions may result in personal injury or damage to property.**
- **Do not use the device on children, disabled or sleeping persons, or heat-insensitive persons who cannot react to overheating.**
- Children shall not play with the appliance.
- **Please do not pull, twist or bend the power lines sharply.**
- The following heating pad is not intended for use in hospitals.
- **When using the heat pad for several hours, we recommend setting the lowest temperature level to avoid overheating the heated part of the body.**
- **Stop using the pad if it becomes uncomfortably hot.**

### 2. Included in delivery

Please first check the completeness of the device and inspect the product for damage. If the device is damaged, do not operate it and contact your dealer or the service department.

The scope of delivery includes:

- 1 x heating pad
- 1 x extension 18cm
- 1 x 1 x extension 60cm
- 1 x USB extension cable
- 1 x USB power adapter
- 1 x Instruction Manual


### 3. Product features

The **Thermrup Far Infrared Multifunctional Heating Pad** is suitable for heating the human body and object, the heating pad provides even and deep heat. With the included rubber bands, the heating pad can be easily attached to different body parts or objects. The product has 3 adjustable temperature levels. The special rubber heating element has a large heating surface and thus provides an overall even heat distribution.

### 4. Operation

#### 4.1 Fix to the body


Place the heating pad around your body or object. Next, close the Velcro on the rubber bands. Use the elastic strap extensions to wear the heating pad in any position on the body.
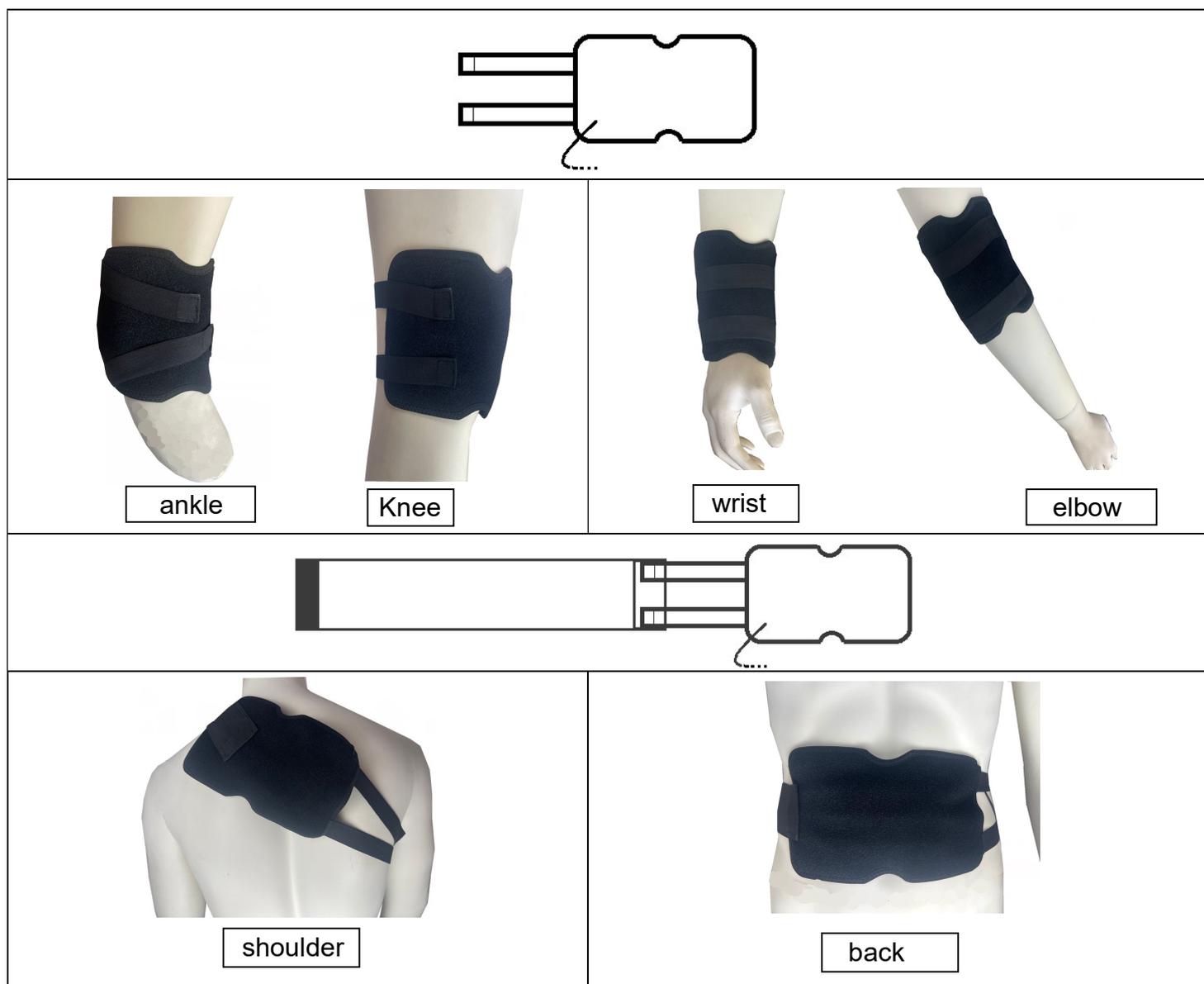
## 4.1 Operating with power bank

- Plug the USB connector of the cable into the USB socket of the power bank (not included).
- Briefly press the power button on the controller to turn on the first level.
- Next, if necessary, press the power key again to set the desired temperature level (1-2-3, red-green-blue). Each single key touch increases the temperature by one level.

The 1st level is the warmest level. (red LED indicator lights up).  
The 2nd stage is the middle stage. (green LED indicator lights up).  
The 3rd stage is the lowest stage. (blue LED indicator lights up).

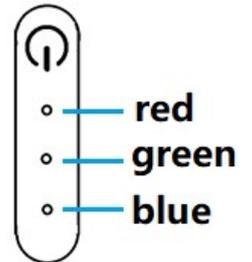- When the blue level is reached and the ON/OFF button is pressed again, the heating function is switched off and the LED indicators no longer light up. After 100 minutes of continuous operation, the heating pad will automatically turn off, restart your powerbank to continue.

**!! Discontinue use if you see reddening of the skin or if the pad becomes uncomfortably hot! Excessive heat can cause burns.**

## 4.2 mains operation

- Connect the supplied extension cable to the connection cable of the pad.
- Plug the USB connector of the extension cable into the supplied power supply.
- Press the power button of the controller short to switch on the controller. Then press the power button briefly to set the desired temperature level (1-2-3, red - green-blue).
- Press the power button until all of the LED off to switch off the controller.
- After use, remove the extension cord from the charger.

**\*\*\* Please use only the supplied charger. \*\*\***

## 5. Cleaning and maintenance

- Outside the controls on the cord, the heating pad is waterproof and hand washable.
- Never use aggressive cleaning agents or strong brushes.
- The heating pad must not be dried in a dryer or ironed. It should only be used again after the device has completely dried.

## 6. Disposal

This product must not be disposed together with the domestic waste. All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner. Consult your municipal authority or your dealer for information about disposal.

## 7. Warranty

We provide a 2-year warranty from the date of purchase for material and manufacturing defects of the product.

The warranty does not apply:

- In the case of damage caused by improper use.
- For wear parts.
- Customer's own fault.
- Once the unit has been opened by an unauthorized workshop.

For the assertion of a warranty within the guarantee period, the customer must prove the purchase.

In case of complaints please contact our service under the following contact:

E-Mail: [service@zewox.de](mailto:service@zewox.de)

**Zewox GmbH**

Oststr.159

47057 Duisburg

Germany

[www.thermrup.com](http://www.thermrup.com)