# DE Ferninfrarot Saunagürtel GB Far infrared Waist Trimmer Belt

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# **English**

#### **Contents**

1.	Important instructions	7
2.	Included in delivery	7
3.	Product features	8
4.	Operating	8
	Cleaning and maintenance	
6.	Disposal	10
	Technical data	
	Warranty	

## 1. Important instructions

- Failure to observe the following information can cause personal injury or property damage.
- Do not use the device on children, disabled people, people who are sleeping or people who are insensitive to heat and who cannot react to overheating.
- Children are not allowed to play with the device.
- This belt may only be operated in connection with the battery specified on the label.
- Do not pull, twist or sharply kink the cable.
- Only use the charger supplied to charge the battery.
- If the belt is used for several hours, we recommend that you set the temperature to the lowest possible level in order to avoid overheating of the heated body part.
- If you do not use the belt every day, pull the connection plug of the belt out of the battery.

## 2. Package Contents

- 1 x heated belt
- 1 x charger
- 1 x Battery
- 1 x adapter cable for mains operation
- 1 x instruction manual
- 1x storage bag

#### 3. Product features

The **Thermrup Far Infrared Waist trimmer belt** for indoor and outdoor use, offers deep warmth with far infrared for the waist with good freedom of movement. It's easy to put on, comfortable and safe. It has 3 adjustable temperature levels. The special rubber heating element has a large heating surface and even heat distribution. The belt is machine washable.

## 4. Operation

## 4.1 Charging the battery

- Before using it for the first time or after a long period of non-use, the battery must be fully charged for at least 4 hours.
- To charge the batteries, connect the charger cable to the battery and to a socket.
   During the charging process, the LED on the charger lights up red. When the battery is fully charged, the LED on the charger lights up green. You can then disconnect the charger from the mains and from the battery.


#### 4.2 Tightening

 Now put the sauna belt around your stomach and close the Velcro fastener. Place the heating surface of the belt around the stomach in any position, depending on where the warming effect is desired.

#### 4.3 Battery operation

- Pull the cable from the battery bag. Connect the cable to the battery.
- Place the battery in the battery bag.

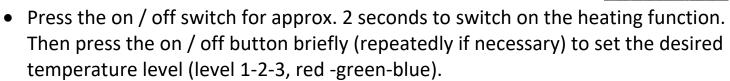- The 1st level is the warmest level. (LED display lights up red)
- The 2nd level is the middle level. (LED display lights up green)
- The 3rd level is the lower level. (LED display lights up blue)
- The four T-LED shows the remaining capacity of the battery. (If the battery is full, all four T-LEDs light up)
- After 30 seconds, the power button will go into power-saving mode and the LED indicator will darken, but the heating function will continue. If you desire a different temperature level during use, briefly press the power button again to unlock it from power-saving mode. You can then set your desired level as described above.
- Press the on / off switch for approx. 2 seconds to switch off the heating function.

\*\* Stop using it if you see redness on your skin. Excessive heat can cause burns!

## 4.4 Main Operation

- Connect the charger and the cable of the belt with the adapter cable
- Press the controller's power button for approx. 2 seconds to switch the controller on.
- Then press the power button briefly (several times if necessary) to set the desired temperature level (1-2-3, red-green-blue).
- Press the power button for approx. 2 seconds to switch off the controller.
- After use, remove the belt from the charger.

\*\*\* Please only use the supplied charger\*\*\*

## 5. Cleaning and maintenance

- · Before cleaning, disconnect the heat belts from the Battery or power supply.
- · Never use aggressive cleaning agents or strong brushes.
- · The heat belt can be washed in the washing machine up to 45 ° C in a washing bag or supplied washing bag.
- · The thermal belt must not be dried or ironed in a dryer.

## 6. Disposal

This product must not be disposed together with the domestic waste. All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner. Consult your municipal authority or your dealer for

information about disposal.

## 7. Technical data

Model: WTB620

Power supply

Heating pad :  $7.4V \sim 8.4V$ ; Max. 15W

Power supply Input :  $100-240V \sim 50 / 60Hz$ 

Output: 8.4V; 2A

Battery : Li-Po 7.4V; 3000mAh

Dimensions : approx. 115/130 x20cm

Weight : approx. 480g

## 8. Warranty

We provide a 2-year warranty from the date of purchase for material and manufacturing defects of the product.

The warranty does not apply:

- o In the case of damage caused by improper use
- o For wear parts.
- Customer's own fault.
- Once the unit has been opened by an unauthorized workshop.

For the assertion of a warranty within the guarantee period, the customer must prove the purchase.

In case of complaints please contact our service under the following

contact: E-Mail: service@zewox.de

#### **Zewox GmbH**

Oststr.159

47057 Duisburg

Germany