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The use of this device is discouraged on the

following body part or under the following

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. Bikini line

Hands and feet

Before removing hair, please cleanse off any antiperspirant, perfume, or emulsion fro

 $pigmentation. \, Level \, 1 \, is \, recommended \, before \, your \, skin \, gets \, accustomed \, to \, it.$

It is more likely to feel some pain in a bony, convex or concave skin area than others.

. Do not use this device on any body part fitting into the following

Characterized by a darker skin tone, stricken by any burn scar, chilblain, black nevus or

Within 48 hours after hair removal, protect the target skin area from exposure

▲ Be sure not to expose your eyes directly to any flashing light.

unusual to feel some mild heat on the skin exposed to flash lights.

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Danger: Erroneous operations may cause personal injuries.

Danger

o not use this product around the eyes. Protective goggles, though capable of oggles before using this product because the high-frequency flashing light might do arm to your eyesight even though such light will not be leveled at your eye

▲ Attention

oabilities or those who lack relevant experience and knowledge must use an, and maintain this product under the supervision and guidance of the

Do not try to insert any object into this device. Also, refrain from leaving this device

O Do not shed light on any clothes, hair or other objects so as to prevent burning,

ver adapter shall be kept dry. evice, please wipe off any hair or skincare product clung evice from breaking down or causing any harm to the skin riate intensity of light in line with your skin conditions. When

ng it. o ensure the surface is completely dry and further free from

ase check the room temperature before using this device, which should be used in

Do not use this product near any water-filled bathtub/basin, washbasin or

prevent any harm. ediately after shedding light into armpits.

first time, please start it up from Level 1.
r on the back, please be advised to have someone assist you ryou to operate this product single-handedly.
left unattended after it is turned on. Be sure to immediately

o not shed light on any black object so as to prevent damage to the black object an

evice.

It is device cannot be used together with any essential oil product. Please also be ure to make the target skin area clean, dry and further free from any emulsion, ibricant, or other cosmetics. ilamide) recently. son suffering from light-induced injuries (sunburn) shall not use this devir

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Do not use this device in bathrooms or other high-humidity environments. Wear protective goggles when using this device.

▲ Contraindications

number to use this device on any person suffering from severe heart disease ite or severe hypertension, malignant tumor, folliculitis, infectious skir.

being susceptible to skin rashes or allergies. This device shall not be used on any target skin area plagued by infections, herpes simplex, pathological changes, hematoma, etc. This device shall not be used on any person suffering from immunosuppressive response conditions, including HIV infection, lupus erythematosus, and porphyria. It is prohibited to use this device on any person suffering from hemorrhage disease or using external application drugs, inclusive of 13-cis-retinoic acid, alpha-hydroxy acid (AHA), beta-hydroxy acid (BHA), or azelaic acid. Use this device with caution on any person who has adnoted a chemical depilatory.

Product Profile Thank you for purchasing MEGELIN hair removal device. Please be sure to read the cautions and the instructions provided by this User Manual carefully before using this

> Parameters Product Name MEGELIN Light Sapphire Lighting Coverage 1.1*3.1=3.41cm² Wavelength 590-1200nm Power Adapter Output:24V==2.5A

tensity, high speed in emitting flashes, and a good cooling effect.The cooling surface cove ie light outlet. Whilst emitting flashes, the cooling surface simultaneously helps to reduce th

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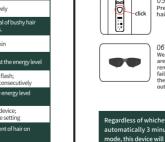
Before You Start Hair Removal 1. Cleanse the target skin area thoroughly, and make sure

Make sure the light outlet is free from any dust, hair, or dirt, has not be

Level-Intensity Checklist (Skin Tone Card)

3.Do a test on the target skin area before using this product and check the condition of the tested skin area twe hours thereafter. You can use this product safely so long as the tested skin area shows no problem. Setting Light Intensity

Skin Test

flashing lights.

to the given light intensity, and to determine intensity for each target skin area. ote: Do not use this device on a body part that is sensitive or tricky (e.g. ankles or other parts of the skeleton). Be sure to set the level at Level 1. Cleaning and Maintenance

4. When the indicator light starts blinking, press the 'Lighting button' to give out 6. Increase the intensity of light to a higher level and enable this device to giv out light for once, before sliding it to the next target skin area.

Skincare After Treatment

 Within 48 hours after the hair removal, do not expose the target skin area t amount of sunscreen cream onto the target skin area for protection purpose

After the skin test is done, wait two hours and check the skin condition If the skin area exposed to light does not show any abnormality (redness or blister), you can be reassured about using this device.

contain exfoliation particles onto the target skin area.

3. Place the device vertically onto the skin to enable its contact with the light

 The target skin area shall not get wet within 6 hours after hair removal. Thereafter, you can apply some moisturizer onto the target skin area for

. Do not expose the target skin area to any artificial tanning during the ha removal treatment. on the target skin area after hair removal.

Frequency of Use

Examine this product and make sure its body, light outlet and components ar

FAQs and Solutions This chapter summarizes those problems that may arise during the use of this device. Please contact the Customer Service enter in your country/territory should the following informa

 3. Do not use a steel wire sponge, any abrasive detergent (e.g. toothpaste) of corrosive liquid (e.g. gasoline, engine oil, and acetone) to clean this device. $\label{eq:model} \ensuremath{\mathfrak{X}} \mbox{ Note: Discontinue the use of this product when it becomes impossible to clean up its light outlet.$

Indicator is not blinking

3. Dab a few drops of water onto the cleaning cloth and then use it to clean up the cooling surface.

tion fails to help you solve the problem in question.

it cool down before cleaning it up.

2. Where necessary, use a piece of soft cleaning cloth to clean up the exterior of

ep your Payment Receipt and Warranty Card in a safe place as these documents are Special Notes

Checklist of Environmental Protection

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the normal operating condition set forth in this User Manual, the environmental fig. of this clearing will last for 10 years

The light has no effect on light golden, gr red or white hair.

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