


MEGELIN

LED Light Therapy Mask

Face & Neck Set

Please read this manual carefully before use, and keep it well, including the Product Warranty Card, and follow the warnings and precautions marked in the product manual.

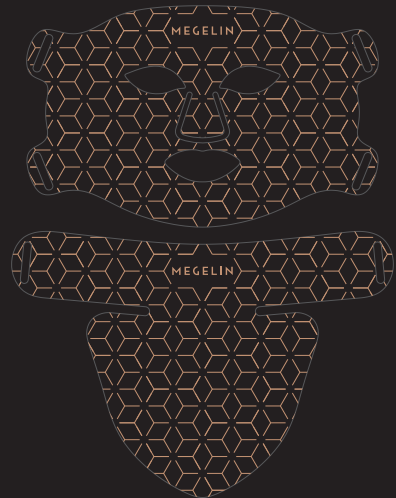

Product instruction

1 FOREWORD

Thank you for purchasing MEGELIN products. The LED light therapy mask face mask is a new product developed by our company, only to bring the theatrical care experience to consumers in a scientific, healthy and safe way.

In order to bring you a better experience and effect, please be sure to read the product manual carefully before use.

2 DEVICE LIST



Red, yellow, blue, green, purple, white, orange, seven color light options

→ Face mask interface
 → Neck mask interface
 → On-off key
 → Selection time


→ Color selection
 → Color selection




- 1 Face Mask
- 2 Neck Mask
- 3 Product controller
- 4 Face mask fixed strap
- 5 Neck mask fixed strap
- 6 Type-C charging cable

*Accessories shall be subject to the physical receipt.


3 QUICK USE INSTRUCTIONS

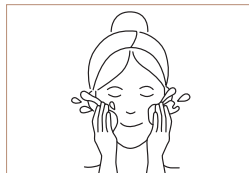

Step 1: Clean your face and neck and keep it dry.

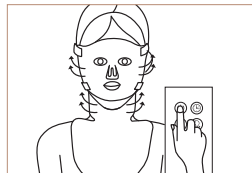

Step 2: Face and neck mask connect the controller and adjusting strap, wear the mask, long press the switch button for 1.5 seconds to turn on the machine, the default red light is turned on, short press the color light selection button to switch the color light.

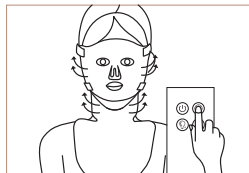

Step 3: The device contains three periods of 10, 15, and 20 minutes. The default time is 10 minutes. You can press the time button to switch the period.



Step 4: After the set time, the device will automatically shut down, remove the face and neck mask, and carry out follow-up skin care.

4 SAFETY GUIDELINES



If you are one of the following people, do not use this product

1. People who are pregnant, preparing for pregnancy and breastfeeding;
2. People who are light sensitive or who are taking any drugs that may cause light sensitivity;
3. People with genetic or skin diseases of the eye (such as lupus erythematosus, photosensitive eczema or albinism).



Charging instructions

- The controller contains lithium-ion polymer batteries. Do not use them during charging. (Note: Once used while charging, it will cause damage to the battery and device, and the mask will not light up and not work)

- Before the first use, please fully charge the controller. When charging, connect the USB charging cable to the controller, connect the plug, turn on the power, and start charging. (Note: In addition to controlling the switch, switching the color light and adjusting the time, the controller also has a "charging treasure" function, so only the Type-C charging cable and the controller can be connected when charging, without connecting to the mask, in order to reduce the damage to the mask parts)

- When the controller indicator flashes, the device is in the charging state; When the controller indicator lights on and stops blinking at the same time, the charging is complete.

5 LIGHT SENSITIVITY TEST

Please perform a light sensitivity test before using this mask for the first time

Step 1

Place the mask on your arm, turn on the switch and shine for 10 minutes. If you have red or itchy skin, which indicates that you may have photosensitive skin, it is recommended that you stop using it or consult a dermatologist before using it again.

Step 2

Before use, clean your face and keep it dry (without makeup) to prevent makeup that may contain light-sensitive ingredients that could interfere with your skin care effectiveness.

Step 3

After use, the skin may appear redness, which is a normal reaction of the skin, and will return to normal within 1-2 hours. If you feel burning or any other discomfort during the process, please stop using immediately.

6 INSTRUCTION FOR USE

How to use and steps

Step 1	Clean your face thoroughly and keep it dry.
Step 2	Connect the power cable of the face and neck mask to the controller, connect the binding band to the mask, place the mask on the face and neck, adjust and buckle the binding band, so that the mask is fixed on the face and neck, and press the on-off button for 1.5 seconds to open the mask.
Step 3	According to the different color light, the mask has a total of red, yellow, blue, green, purple, white, orange, seven kinds of color light, you can press the light selection button to switch the color light. After the device is turned on, the default color is red light. You can switch the color light by pressing the color light selection button. Press the time select button to adjust the irradiation time, the default time is 10 minutes, 10/15/20 minutes optional.
Step 4	The mask will automatically turn off when your chosen use time arrives. After closing, remove the mask, end the irradiation, and skin care can be done according to the daily process.
Remark	It is recommended to use it 3-5 times a week for 4 consecutive weeks. Each session lasts for 10-20 minutes.

7 QUALITY ASSURANCE SERVICE

Under normal use, where the product itself has quality problems, this product provides a one-year warranty for the host. Identified as damage, non-quality problems are not covered by the warranty.

If you have any questions or questions during use, you can call the phone of the purchase dealer between 9 and 18 on weekdays for consultation.

8 CARE AND MAINTENANCE

Deposit	Please store the product in a dry environment at room temperature, avoid dust, direct sunlight or damp environment.
Clean	Before and after use, use a water-based rag or a wet cleaning cloth to wipe the mask. When cleaning, ensure that the device is turned off and disconnected from the power supply.
Recycling treatment	There is a rechargeable battery in the controller, please be sure to remove the battery from the controller before the consumer decides to discard the product. The battery should be disposed of safely. When taking out the battery, the power should be cut off. The battery must be discarded in the battery recycling and processing unit stipulated by the government for decomposition.

9 CERTIFICATE

Product name: LED Face Mask

Inspection: Qualified products

Production Date: _____

Inspector: _____

10 PRODUCT PARAMETER

Product name	Megelin LED Light Therapy Mask
Product type	Face & Neck Set
Brand	MEGELIN
Product model	JML1+JML2
Operating power	9W
Electricity	2.5A
Input voltage	5V
Product size	Neck Mask 343*234*5.5mm Face Mask 309*194*5.5mm

12 SERVICE AND MAINTENANCE

1. The main engine of this product is guaranteed for one year from the date of purchase (with the bill), and accessories and consumables are not covered by the warranty; In the scope of non-warranty and beyond the warranty period, you can also be responsible for maintenance and renewal (subject to labor service fees).

2. Accessories and consumables, if you need to supplement or purchase, please contact the dealer.

3. We do not provide free warranty service for the following human factors:

- ① Failure caused by self-disassembly and modification of the instrument;
- ② Failure caused by improper use and lack of reasonable maintenance;
- ③ The failure caused by careless falls and falls.

13 WARRANTY CARD

Product name		
Order number		
Date of purchase		
Customer name and phone number		
Service time	Maintenance record	Maintainer's signature