Skin Rejuvenation Photon Mask

(Model No. PE021)

User Manual

www.projectebeauty.com

Introduction

Thank you for purchasing the Skin Rejuvenation Photon Mask (Model No. PE021). This device uses colored LED lights to provide pure light to treat a variety of skin conditions and rejuvenate skin, while reducing scar tissue. Unlike other beauty devices, this mask provides high-quality LED photon treatment that has been scientifically validated, so you can achieve efficient results faster than you ever imagined.

Photon light therapy works by using the natural healing properties present in differing wavelengths of light to treat various skin tissue issues. These state-of-the-art LEDs provide the purest light possible, especially compared to systems that use more traditional bulbs. Modern semiconductor technology has allowed dermatologists to perfect what wavelengths of light are generated, meaning more of the healing energy of the light reaches the skin tissue.

Exposing the skin's tissue surface to light activates natural photoreceptors and allows the skin cell to absorb the energy found in the light. This process accelerates biochemical reactions in the skin tissue, ultimately, leading to faster healing times, while increasing the production of collagen. This can reverse some of the more heinous effects of aging.

This technology is safe. Cosmetology experts agree these treatments have no noticeable side effects if they're administered correctly. Treatment methods provided by the LED photon mask are an important component of a complete skin care routine and should only be used for cosmetic purposes. This is not intended in any way to diagnose, treat, cure, or prevent any disease.

Benefits of Photon LED Therapy

Recent cosmetology breakthroughs have shown that LED light therapy is effective for treating several aspects of skin aging. Studies suggest this therapy can reduce the appearance of fine wrinkles, rough skin, shallowness of the flesh, and mottled hyperpigmentation. The LED light therapy increases both the proliferation and the initial production of collagen. A few studies have also suggested this therapy can cut down

Benefits of Photon LED Therapy

on the types of bacteria that cause spots and pimples. This might also help to eliminate breakouts and minimize the level of scar tissue.

Each of the colors of LED light have a different specific function. They are each selectable by way of the control box.

Red Light(~630nm): This setting allows the LED mask to penetrate the skin more deeply. Red light penetrates the skin approximately 8-10 mm, which stimulates the cells that lie at this layer. Many are fiber cells, and stimulating them can aid in the production of collagen. It can reduce wrinkles and make skin elastic and smooth. It's ideal for improving skin in the area around the eyes and forehead. It works well to cut down on fine lines. Those who need to repair damaged skin will opt for this setting.

Blue Light(~525nm): This setting allows the energy from the LED mask to penetrate the skin approximately 0.5mm. It's optimized for use with rashes and is suitable for skin that is sensitive or constantly oily. Blue light can help to reduce acne and even tighten skin that has loosened. It can quickly inhibit inflammation and possibly even destroy harmful bacteria.

Green Light(~525nm): This setting can help with mixed skin complexions, and will help to brighten dull skin that has lost its luster. Studies suggest that green light can help to clear the face's natural lymphatic system. It might even aid in providing relief for mental stress, which is of special importance, as many people use the LED photon facial mask while they relax.

Purple Light(~390nm): This setting improves blood circulation with the function of promoting the activation of skin tissue. The purple light can strengthen the protein fiber tissue and shrink pores. It can also help to get rid of yellowish complexion and remove wrinkles and acne. At the same time, it accelerates blood circulation and lymph metabolism.

Cyan Light(~490nm): This setting provides gentle acne treatment and mitigates skin inflammation. It can be used to adjust and improve the gland function of oily, black, and pink skin as well as acne. Meanwhile, the cyan light is able to promote the synthesis of protein and collagen that activate skin cells and provide the relaxation of skin contraction efficiently.

Benefits of Photon LED Therapy

Yellow Light(~590nm): This setting revitalizes the lymphatic and nervous systems, strengthens the muscular and immunological systems, whilst simultaneously improving and balancing sensitive skin. Yellow light improves microcirculation, adjusts factors, and inhibits melanin production in order to fade dark spots, freckles, and blemishes. This method of skin rejuvenation is recommended as an effective treatment for skin problems caused by age and exposure.

White Light(~Mix): This setting synergistically helps activate your skin's natural renewal system to reverse the effects of aging at the same time restoring a younger and healthier skin.

Safety Information

Precautions:

- Read all instructions before using your Skin Rejuvenation Photon Mask and make sure you understand all of the operating instructions.
- Always check the device before use to ensure it is undamaged. Use of a damaged device can be unsafe.
- Stop immediately, if you experience any adverse effect or have a concern regarding use.
- Do not look directly or shine light in the eyes while in operation.
- Do not use on eyelids; do not place the light over eyelids.
- Do not use if you are pregnant.
- Do not use if you have experienced a migraine in the last two years.
- Do not use over any suspicious or cancerous lesion. Consult your doctor.
- Do not use the LED mask while using any product that contains AHA or glycolic acid. This can cause unwanted levels of stress in the skin.
- Do not allow children to use the device unsupervised.

Safety Instructions

1. Read these instructions, hold onto them, follow them, and always heed all warnings.
2. Never use the device near water or areas of concentrated humidity.
3. Disinfect the LED surface cover with a soft cloth and a light concentration of rubbing alcohol. Carefully clean the machine of all types of skin oil.
4. Always use the device in line with the manufacturer's regulations, and never block any of the ventilation openings on the side of the mask.
5. Dispose of all products in accordance with local and national regulations and codes.
6. Do not use if you are sensitive to light. If you are unsure whether you have increased sensitivity to light, we recommend you perform the Light Sensitivity Test instructed in this manual.
7. Do not use if you are taking drugs known to make you sensitive to light. You may be sensitive to light and not know it. There are photosensitizing agents in drugs, foods, soaps, cosmetics, and dyes.
8. Drugs that are known to cause light sensitivity: Antibiotics; NSAIDS; Diures; Retinoids; Hypoglycemics, HMG-CoA reductase inhibitors (statins), Photodynamic Therapy drugs, Prophoto sensitizers, Neuroleptic drugs, or Antifungals.
9. Do not use over or near a recent steroid injection. Doing so may result in photosensitive reaction.

Safety Information

Warning

- Discontinue use if you encounter any discomfort or skin redness that lasts for more than two hours. You may have sensitivity to light or another condition. Consult your health care provider.
- Never operate this device if it has damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or submerged in water. Return the device to Project E Beauty for examination and repair.
- Do not set or keep the device where it can fall or be pulled into a tub or sink. Keep away from water (including bath and shower). If the device accidentally falls in water unplug immediately.
- Keep away from flammable aerosol products.
- Do not use other AC Power Supplies or Chargers with your LED mask.
- Turn off immediately, if the mask begins to overheat.

Contraindications

DO NOT USE if one of the following pertains to you. with these devices resulting in serious injury or death.

Pregnant: Do not use if you are pregnant or think you might be pregnant. The device has not been tested on pregnant women and the effects, if any, are unknown.

Migraine: Do not use if you have experienced a migraine in the last two years.


Seizures: Do not use if you have a history of light triggered seizures.

Lesions: Do not use this device over any suspicious or cancerous lesions. To do so may result in delaying the best medical treatment.

Failure to use and maintain the LED mask in accordance with the instructions outlined in this manual will invalidate the warranty.

What is Included

- LED Mask Device (Model No. PE021)
- Control Box
- Head Strap
- AC Mains Power Cord
- Owner's Manual (Not Depicted)

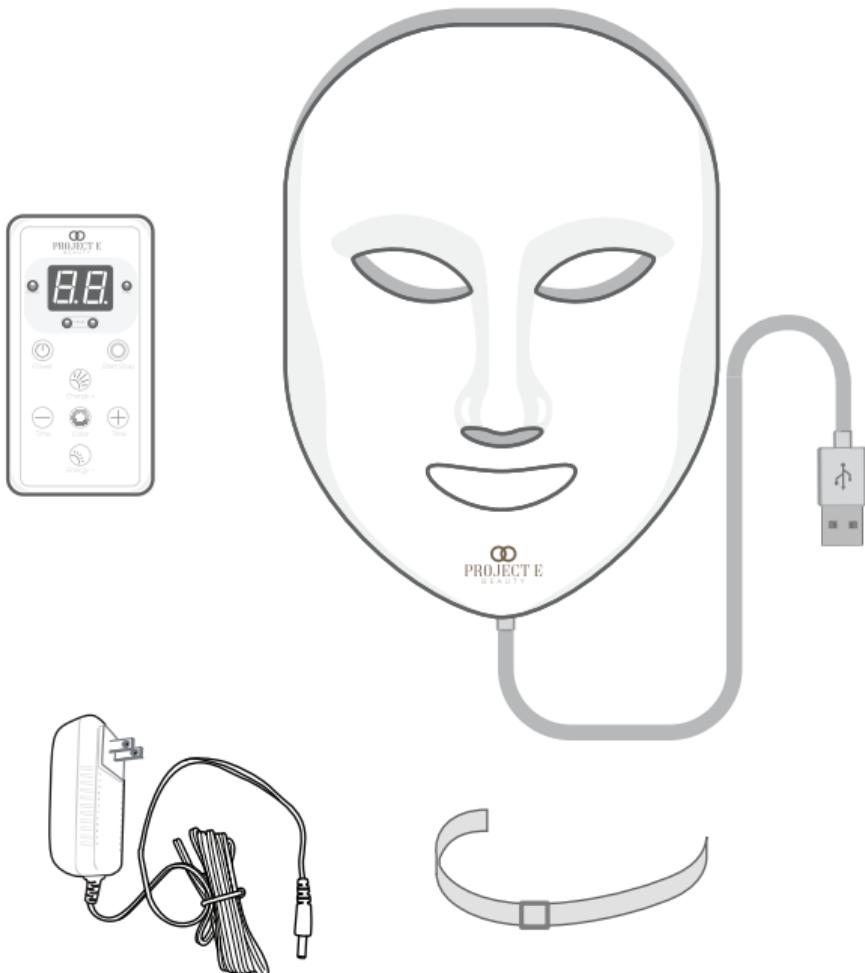Light Sensitivity Test

Recommended Before Use

In order to determine if you have a Light Sensitivity Issue please follow the instructions below. It is recommended that you test a small, less sensitive area of skin (i.e., like your forearm).

1. Press the “Power” button to turn on the mask and then press the “Start/Stop” button to activate the mask LEDs. Press “Color” button to toggle different colors.
2. Place the device on an area of the forearm. Leave in contact with your skin for a 3-minute treatment cycle.
3. Repeat by doing another 3-minute treatment on the same area of your arm.
4. Wait 24 hours. You should NOT experience any adverse reaction during the next 24 hours.

A **bad reaction** to light occurs if you have pain or discomfort in the treatment areas or any skin redness for more than two hours. It can be due to medication, supplements, or a possible skin issue.


Eye Safety

1. The design of this mask does not cover the eye area. Therefore, it is not needed to put on eye protection while performing treatments. Be sure you never look directly into the LEDs.
2. It is suggested the user selects the low energy mode at beginning and then adjusts to higher energy level once the user is adapted.
3. It is **NOT** recommended to shine the light source directly into eyes or the eyes of others.

Use of the Device

Follow these steps to ensure proper operation during a treatment:

- 1 Plug one end of the head strap into each side of the Photon LED Skin Rejuvenation Facial Mask, so the mask can be worn comfortably on the face. Adjust the strap for a good fit as desired. Make sure the Photon LED Skin Rejuvenation Facial Mask is never worn too tightly against the face.

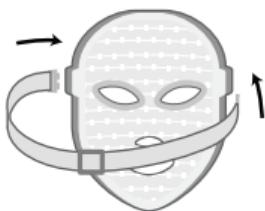- 2 Insert the USB cord into the jack on the control box. Do not force the USB plug into the jack. If it does not fit, turn the plug around and carefully try and insert in the opposite direction.

- 3 Disinfect the surface with a soft cloth. Make sure your face is clean, dry, and wiped free of any oil or cosmetics before using the mask.

- 4 Press the button marked Power to turn on the unit. A green indicator LED will switch itself on when this button is pressed.

- 5 Press the Start/Stop button to turn on the mask LEDs.


- 6** Use the color button on the control box to select the desired color setting for the mask.


- 7** Select a power level by using the Energy +/- button. The highest level is five (5).

- 8** Select a time limit for the mask, using the duration +/- buttons. The highest duration is 60 minutes.

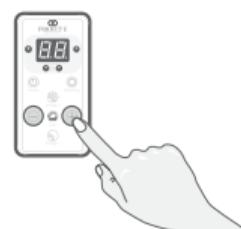- 9** Place the mask over the face, using the strap to hold the mask in place. Adjust the strap as required. You may wish to lie back to ensure good contact.

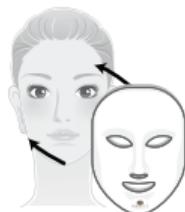- 10** Be sure you never look directly into the LEDs. These powerful semiconductor lamps are certified safe for use on skin, but can be potentially hazardous to delicate eye.

It's suggested that you do not use the mask for more than 20 minutes each time. The suggested operating time is between 10 and 15 minutes. More than likely, you will want to start with one 10- minute treatment, 4 times a week to start, unless otherwise directed by a doctor. Once you

have achieved your initial goals, you will, more than likely, want to drop down to using the device only twice a week. As the LEDs provide a pure and powerful output, these few treatments will be sufficient.

The mask will automatically shut the LEDs off, once the time limit runs out. To turn the device off before 10 minutes, use the ON/OFF switch on the control box. Other time amounts can also be specified with the control box. If you have sensitive skin, then you may wish to start with short periods of time and work your way up. Never go beyond whatever level you feel is comfortable.

Maintenance

- Be sure to unplug the device properly by holding the plug and removing it from the AC outlet.
- Disinfect the mask surface with a soft cloth after each use.
- Immediately dry the surface any time it comes into contact with moisture, chemicals, or water.
- Always keep the device clean and sterilized.
- Store the device in a clean and dry place.

Specifications

Red Light Photon Wavelength	630 nm ±5%
Blue Light Photon Wavelength.....	415 nm ±5%
Green Light Photon Wavelength.....	525 nm ±5%
Purple Light Photon Wavelength.....	390 nm ±5%
Cyan Light Photon Wavelength.....	490 nm ±5%
Yellow Light Photon Wavelength.....	590 nm ±5%
Power Output	5,1 W
Input Voltage.....	100V-240V

Warranty

Activate Your 1 Year FREE Warranty Now.

This product comes with a full one year warranty against manufacturer defects. It does not include damage caused by misuse, accident, abuse, neglect or other external causes. All we need is your name, email and order ID. That's it.

Scan the QR code or visit
www.projectebeauty.com/pages/warranty

Questions

Website: www.projectebeauty.com

Email: sales@projectebeauty.com

Tel: +1-626-759-4232

Our customer care expert team are here to answer any question you might want to ask.

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