**PROJECT E BEAUTY
ICE ROLLER
COLD THERAPY**

Operating Instructions
Model No. PE039

These instructions contain important safety and use information. Please read the instructions completely before using this product. Save the instructions for future use.

Introduction

Thank you for purchasing the Project E Beauty Ice Roller. The Roller utilizes cold therapy to rejuvenate, tighten and refresh the skin. It calms the skin especially when any discomfort caused by fatigued and irritated skin. Cold therapy has been widely utilized by estheticians and spas. The Ice Roller is a natural, durable hand-held tool that convenience your life with effortless rolling, which is easy to be added into your daily skincare routine to revive your skin.

Key Features

1. Shrink pore size, improve skin tone.
2. Eliminate puffiness under eye circles.
3. Tightens skin, reduce fine lines.
4. Soothe inflammation.
5. Reduces muscle soreness and pain.

How to Use?

- 1 Put the roller in the freezer for 4 hrs, it'll keep cool for 2 hrs.

- 2 Clean the treatment area and apply gel or cream.

- 3 Roll gently around the treatment area.

- 4 Clean the roller with cotton and 70 % alcohol. Put it back into the freezer.

Caution:
Do not heat the roller.

How to Use?

Face

Tighten Skin:


Slowly move upwards from the chin to the bottom of the ear.

Lift the device off the skin and move from the chin to the corners of the mouth.

Eye Puffiness:


Slowly move from the corner of the mouth towards under eye area and hold for 5 seconds.

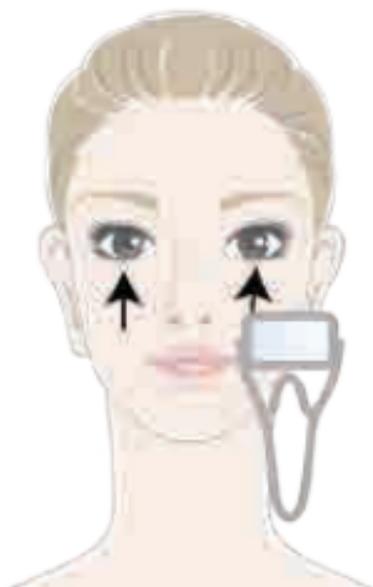Front

Slowly move upwards from above the knee to the top of the outer thigh.


Lift up and do the inner thigh.

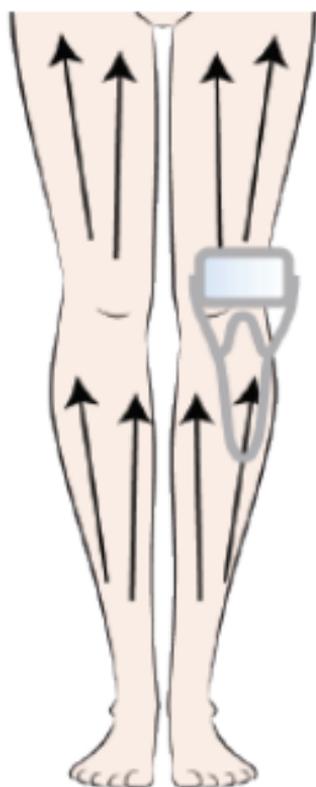Back

Slowly move upwards from above the back of the knee to the top of the outer thigh.

Lift up and do the inner thigh.

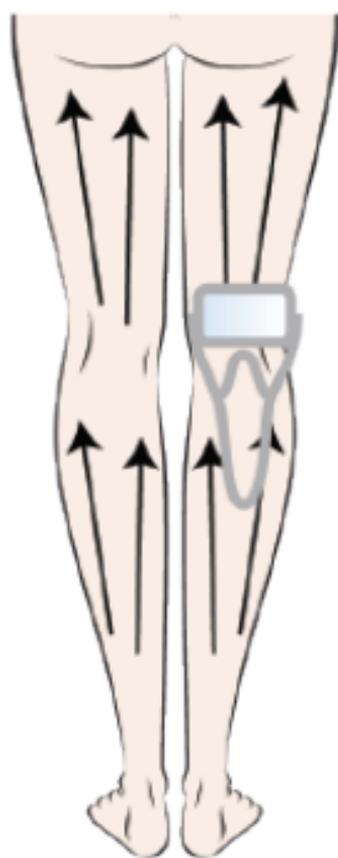
Arms

Slowly move upwards from the back if the elbow to the inside shoulder.

Lift up and do the outer shoulder.

Repeat for the front of the arms.

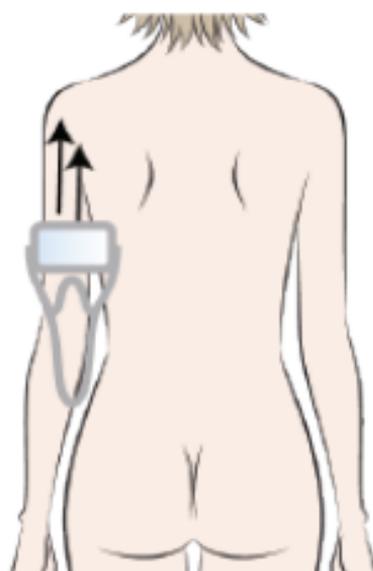Cleaning & Maintenance

Ordinary maintenance of the machine and its accessories can be carried out without opening the case.

Disinfect the roller with 70% alcohol using a cotton or soft cloth and carefully.

Warranty

The product comes with a full one year warranty against manufacturer defects. It does not include damage caused by misuse, accident, abuse, neglect or other external causes. The warranty is void if an unstable voltage or inappropriate power supply has been used.

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