Facial & Eye Photon Ultrasonic Device

E01

www.projectebeauty.com

Introduction

Thank you for purchasing the Project E Beauty Facial & Eye Photon Ultrasonic Device. This device utilizes LED Photon and Ultrasonic methods to treat a variety of skin conditions.

What separates this device from other beauty devices is the Ultrasonic feature. This function takes skin care to the next level and coupled with the LED Photon treatment, you will achieve results more effectively and faster than imaginable.

Ultrasonic therapy works using low frequency sound waves for safe mechanical exfoliation, molecular penetration and healing therapy. Photon therapy works using the natural healing power of light to treat numerous skin issues utilizing the latest in LED technology.

Exposing light to the skin surface activates the skin's photoreceptors and let the skin cells absorb its energy. It accelerates cell biochemical circulation and stimulates the production of collagen.

These safe, non-invasive technologies can be beneficial in the treatment of many skin conditions ranging from acne, wrinkles, spots removal to skin tightening. These utilized functions have been found to enhance oxygenation, stimulate dermal tissue, kill bacteria causing acne, and treat oily or combination skin problems.

Although results can vary widely by individual and skin type, skin care treatment is appropriate for all types of skin.

Cosmetologists and skin care specialists agree that these treatments, when properly performed, have no side effects. The treatment methods of this care routine and important component of your skin care routine and should only be used for cosmetic purposes. It is not intended in any way to diagnose, treat, cure, or prevent any disease.

Benefits of Ultrasonic Photon LED Therapy

This device utilizes two feature functions as follows:

LED Photon Therapy

The numerous LED photon therapy options in this device will provide different results:

BLUE (415 nm): Press the LED Select button twice: This setting will allow the device to penetrate the skin approximately 0.5 mm and is optimal to use for skin rashes. This is suitable for skin that is sensitive, oily and/ or has skin rashes. Also, effective at reducing acne and tightening loose skin.

GREEN (525nm): Press the LED Select button three times: This setting is optimal for mixed skin and will brighten dull skin.

RED (630 nm): Press the LED Select button once: This setting will allow the device to penetrate the skin approximately 8-10 mm, which will stimulate fiber cells and the production of collagen, reducing wrinkles and making skin more elastic and smooth. Optimal for surrounding eyes and forehead and works best for fine lines, reducing pore size and swelling after surgery.

Benefits of Ultrasonic Photon LED Therapy

Ultrasonic Therapy

Ultrasonic facials offer a safe and gentle alternative to dramatic plastic surgery procedures, laser resurfacing, chemical peels, collagen and Botox injections and other invasive skin rejuvenation treatments.

Results can vary by individual and skin type, are gradual, and will not just occur over night. Although Ultrasonic treatment has been shown to produce an immediate and temporary firming effect, continued application will provide more cumulative long-term, lasting results. Ultrasonic treatments can give you the following benefits:

- Stimulates Blood Circulation
- Aids in Lymphatic Drainage
- Gently Exfoliates Away Dead Skin Cells
- Increases Cell Metabolism
- Cleanses Away Toxins
- Enhances Skin Care Product Performance
- Tightens and Tones Neck and Facial Skin
- Diminishes Wrinkles and Fine Lines
- Reduces the Appearance of Congested, Puffy Eyes
- Refines Large Pores
- Improves Skin Texture and Tone
- Reduces Acne and Blemishes
- Lightens Dark Under Eye Circles
- Reduces the Appearance of Fat Deposits
- Fades Scars, Freckles, Age Spots and Redness

How Often

Perform each treatment for 10 minutes four times a week to start. Once the goal is achieved, use only twice a week.

Safety Information

Precautions:

- Read all instructions before using your Facial & Eye Photon Ultrasonic Device and make sure you understand all of the operating instructions.
- Always check the device before use to ensure it is undamaged. Use of a damaged device can be unsafe.
- Stop immediately, if you experience any adverse effect or have a concern regarding use.
- Do not look directly or shine light in the eyes while in operation.
- Do not use on eyelids; do not place the light over eyelids.
- Do not use if you are pregnant.
- Do not use if you have experienced a migraine in the last two years.
- Do not use over any suspicious or cancerous lesion. Consult your doctor.
- Do not use the device while using any product that contains AHA or glycolic acid. This can cause unwanted levels of stress in the skin.
- Do not allow children to use the device unsupervised.
- Always unplug the device correctly by holding the plug and removing it from the electrical socket.
- Remember that improper maintenance or negligence can cause damage and void the warranty. This also refers to improper storage.

Safety Instructions

1. Read these instructions, hold onto them, follow them, and always heed all warnings.
2. Never use the device near water or areas of concentrated humidity.
3. Disinfect the LED surface cover with a soft cloth and a light concentration of rubbing alcohol. Carefully clean the machine of all types of skin oil.
4. Always use the device in line with the manufacturer's regulations, and never block any of the ventilation openings on the side of the device.
5. Dispose of all products in accordance with local and national regulations and codes.
6. Do not use if you are sensitive to light. If you are unsure whether you have increased sensitivity to light, we recommend you perform the Light Sensitivity Test instructed in this manual.
7. Do no use if you are taking drugs known to make you sensitive to light. You may be sensitive to light and not know it. There are photosensitizing agents in drugs, food, soaps, cosmetics, and dyes.
8. Drugs that are known to cause light sensitivity: Antibiotics; NSAIDS; Diures; Retinoids; Hypoglycemics, HMG-CoA reductase inhibitors (statins), Photodynamic Therapy drugs, Prophoto sensitizers, Neuroleptic drugs, or Antifungals.
9. Do not use over or near a recent steroid injection. Doing so may result in photosensitive reaction.

Safety Information

Warning

- Discontinue use if you encounter any discomfort or skin redness that lasts for more than two hours. You may have sensitivity to light or another condition. Consult your health care provider.
- Never operate this device it has damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or submerged in water. Return the device to Project E Beauty for examination and repair.
- Do not set or keep the device where it can fall or be pulled into a tub or sink. Keep away from water (including bath and shower). If the device accidentally falls in water unplug immediately.
- Keep away from flammable aerosol products.
- Do not use other AC Power Supplies or Chargers with your LED device.

Contraindications

DO NOT USE if one of the following pertains to you.

Light Sensitivity: Do not use if you are sensitive to light. You are sensitive to light if you experience any discomfort or skin redness that lasts more than two hours. Discontinue using the device and consult with your healthcare provider.

Implantable Device: Do not use this device if you have a cardiac pacemaker, implanted defibrillator, shunt or other implanted electronic device. This device contains magnets that may cause interference with these devices resulting in serious injury or death.

Pregnant: Do not use if you are pregnant or think you might be pregnant. The device has not been tested on pregnant women and the effects, if any, are unknown.


Migraine: Do not use if you have experienced a migraine in the last two years.

Seizures: Do not use if you have a history of light triggered seizures.

Lesions: Do not use this device over any suspicious or cancerous lesions. To do so may result in delaying the best medical treatment.

Failure to use and maintain the device in accordance with the instructions outlined in this manual will invalidate the warranty.


What is Included

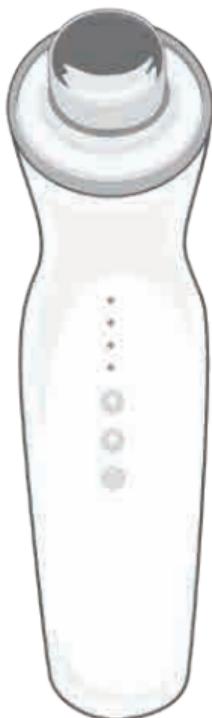Ultrasonic Photon Device


USB Cable

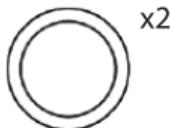x2

Two Rings (Optional)
Eye Protectors (Not Shown)
Owner's Manual (Not Shown)

Location of Controls

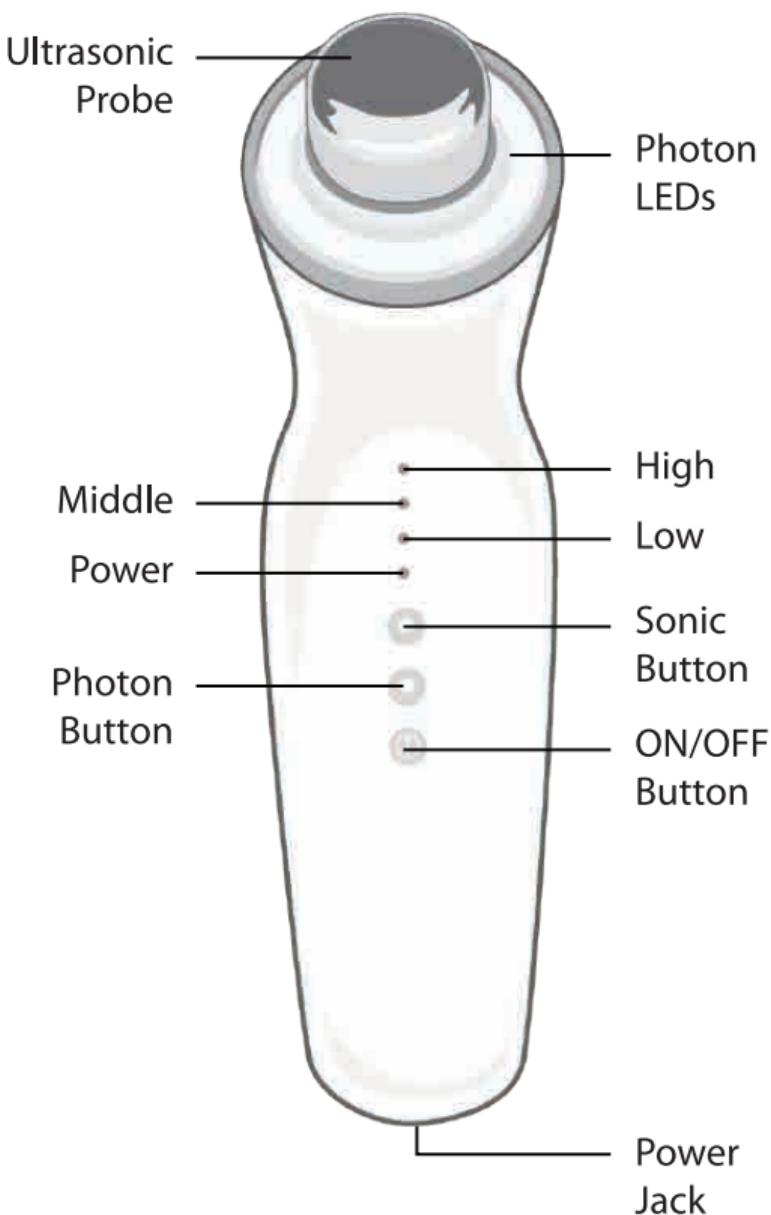Light Sensitivity Test

Recommended Before Use

In order to determine if you have a Light Sensitivity Issue please follow the instructions below. It is recommended that you test a small, less sensitive area of skin (i.e., like your forearm).

1. Press the “Photon” button to activate light.
2. Place the device on an area of the arm. Leave in contact with your skin for a 3-minute treatment cycle.
3. Repeat by doing another 3-minute treatment on the same area of your arm.
4. Wait 24 hours. You should NOT experience any adverse reaction during the next 24 hours.

A **bad reaction** to light occurs if you have pain or discomfort in the treatment areas or any skin redness for more than two hours. It can be due to medication, supplements, or a possible skin issue.

Eye Safety

1. Position the goggles over your eyes
2. Tie the elastic band so that the goggles fit comfortably.
3. It is **NOT** recommended to shine the light source directly into eyes or the eyes of others.

WARNING

Be sure you never look directly into the LEDs.

Use of the Device

Before performing a treatment you should do the following:

- Remove all jewelry and any metal objects.
- Check if absolute or relative contraindications exist (see page 8).
- Check that the device is functioning properly and the electrical cord is not damaged.
- Ensure there are not any flammable materials near the work space.


Use of the Device

Follow these steps when performing a treatment:

- 1 Before you first use the device, charge it for 3 hours.

- 2 Insert the USB cable to standard 5V charger(not included); Insert the charger into a AC outlet.

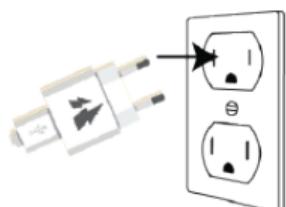Charging

This device can be used connected to AC or can be charged to be used without the AC cord.

Use of the Device

- 3** Disinfect the surface with a soft cloth.

Thoroughly clean the face with soap and water, and dry.

- 4** Use an ion gel, essence or cream on areas you will be working on to achieve the best results.

You can purchase our Hyaluronic Acid Serum at www.projectbeauty.com

It is highly recommended to put on the supplied eye protectors at this time.

- 5** Press the ON/ OFF button to turn the unit on; the Power indicator will light.


Note: The device is effective at all intensity levels including the lowest setting. It is not necessary to use the highest power level to achieve maximum results.

Use of the Device

- 6 SONIC:** Repeatedly press the Sonic button to select the desired power level (High, Middle or Low).

Optional: When using the Sonic feature, you can use a plastic ring to hold on a gel-soaked cotton cloth to the probe.

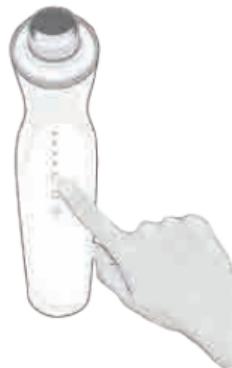Note: Ultrasound vibrations are difficult to feel on your skin, but that does not mean they are not working. To test, put a few drops of water on the ultrasound probe and you will see the water drops vibrating.

- 7 PHOTON:** Repeatedly press the LED button to select the desired LED effect. See previous page for detailed descriptions of the LED effects.

When using the Photon feature, use the eye protectors. Do not look into the LEDs directly.

- 8** Perform each treatment for 10 minutes (except nose) four times a week to start. Once the goal is achieved, use only twice a week. Place the electrode on the skin and move in the directions as shown according to the desired effects.

Use of the Device


Tighten Skin:

Slowly move from the corner of the mouth towards the area near the eyes.

Pull up Sides of Mouth:

Slowly move from the corner of the mouth slightly inwards and upwards.

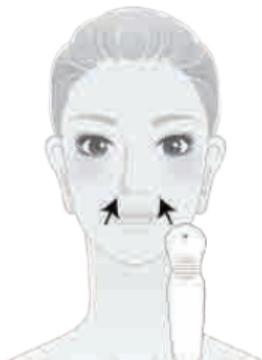Eczema and Acne:

Do not operate on the exact area of the eczema, acne, etc., but use on the surrounding skin to aid in the natural healing.

Use of the Device

Double Chin:

Slowly move upwards from the chin to the bottom of the ear, then lift the device off the skin and move from the chin to the corners of the mouth.

Note: Always move upwards; do not move the unit downwards.

- 9 After 10 minutes, the device will automatically turn off. To turn the unit off before 10 minutes, press the On/OFF button until the unit turns off.

Maintenance

- After each use, disinfect the surface with a soft cloth.
- Immediately dry off any moisture, chemicals or water on the device.
- Always keep the device clean and sterilized.
- Store the device in a dry place.

Specifications

Red Light Photon Wavelength.....	630 nm
Blue Light Photon Wavelength.....	415 nm
Green Light Photon Wavelength.....	525 nm
Ultrasonic Frequency.....	3MHz, ±5%
Power Output	2.8VA
Power Supply	100V-240V

Warranty

Activate Your 1 Year FREE Warranty Now.

The product come with a full one year warranty against manufacturer defects. It does not include deamge caused by misuse, accident, abuse, neglect or the external causes. All we need is your name, email and order ID. That's all.

Scan the QR code or visit
www.projectebeauty.com/pages/warranty

Questions

Website: www.projectebeauty.com

Email: sales@projectebeauty.com

Tel: +1-626-759-4232

Our customer care expert team are here to answer any question you might want to ask.

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