PROJECT E
BEAUTY

**Portable Rechargeable
Galvanic Device**

PE094

www.projecteauty.com

Introduction

Thank you for purchasing the Project E Beauty Portable Skin Care Galvanic Device. This device utilizes Galvanic positive and negative ionization to treat a variety of skin conditions.

Galvanic treatments work in a natural and harmless way. This treatment achieves maximum results by penetrating active substances into the tissues, where they act on the circulation to bring about an improvement in the vascular and lymphatic interchange in the area. This improves the function of the cellular membrane, and allows the trapped fluid and fat to be dispersed and eliminated.

This device utilizes a two-part treatment, first being desincrustation which will prepare the skin for the next step, iontophoresis which is used to assist the nutrients, cream, etc. into the skin for optimal results.

This safe, non-invasive technology can be beneficial in the treatment of many skin conditions such as improving skin texture, firming the skin, reducing redness, reducing circulation and more. It does this because each treatment will relax the pores (making extraction easier), deep cleanse the skin, stimulate secretions of the sweat glands as well as nerve endings, increase circulation to the area which will allow penetration of the nutrients and have an overall softening and drying effect.

Cosmetologists and skin care specialists agree that these treatments, when properly performed, have no side effects. The treatment methods of this unit are an important component of your skin care routine and should only be used for cosmetic purposes. It is not intended in any way to diagnose, treat, cure, or prevent any disease.

The device comes with a Roller Probe that is designed for the open face and neck areas as well as a Ball Probe which is the perfect size for working around the eyes and nose, giving you the best results.

Safety Information

You must read, and understand, all of the safety warnings, operating and maintenance instructions before using the Project E Beauty Portable Skin Care Galvanic Device.

Failure to follow the warnings and instructions may cause serious injury, damage to the device or void the warranty. Warnings descriptions are as follows:

WARNING!: A hazard that may cause serious injury or death.

CAUTION!: A hazard that may cause minor injury.

NOTICE: A hazard that may cause damage to the product

WARNING!

- Performing a treatment on a person with absolute contraindications may cause serious injury or death. Do not perform a treatment on a person where absolute contraindications exist.
- Use of the device near flammable materials may cause explosion, fire and serious injury or death. Always check the work space and remove any flammable materials prior to performing a treatment.
- Use of the device around water (i.e., bathtubs, wash basins, areas of high humidity) or with wet hands may result in electrical shock, serious injury or death. Never use the device around water. Always dry your hands before performing a treatment.
- Improper use of the device by unqualified people (i.e., children) may result in burns, electrical shock or other serious injury, or death. Never allow unqualified people to use this device for any purposes or play.
- Performing a treatment with a damaged device may cause serious injury or death. Always check the device for damage before performing a treatment. Never perform a treatment with a damaged device.
- Use of the device internally (i.e., mouth, nose, etc.) may cause serious injury or death. Never use the device internally.

Safety Information

CAUTION!

- Performing a treatment near areas of relative contraindications may cause minor injury. Avoid performing a treatment near areas where relative contraindications exist.
- Performing a treatment with AHA or Glycolic Acid products may cause unwanted stress to the skin. Do not use AHA or Glycolic Acid products when performing a treatment.
- Performing a treatment around broken capillaries, spider veins or rosacea may cause unwanted stress to the skin. Do not perform a treatment around broken capillaries, spider veins or rosacea when performing a treatment.
- Use of the device around jewelry and metal may cause unexpected electrical shock. Always remove all jewelry and metal prior to a treatment.

IMPORTANT NOTES:

- Improperly unplugging the plug when the device is not in use may cause damage to the cord resulting in fire or electrical shock. Always unplug the device by holding the plug and removing it from the electrical socket.
- Use of the device around water (i.e., bathtubs, wash basins, areas of high humidity) may cause damage to the device. Never use the device around water.
- Improper storage or maintenance may result in damage and void the warranty.

Safety Information

1. Read these Instructions.
2. Keep these Instructions.
3. Heed all Warnings.
4. Follow all Instructions.
5. Do not use this apparatus near water.
6. Disinfect the electrode with 70% alcohol using a soft cloth and carefully cleaning the electrodes, LED cover and Probes.
7. Do not block any of the ventilation openings. Install in accordance with the manufacturer's instructions.
8. Do not install near any heat sources such as radiators, heat registers, stoves, or other apparatus (including amplifiers) that produce heat.
9. Protect the power cord from being walked on or pinched, particularly at plugs, convenience receptacles, and the point where they exit from the apparatus.
10. Only use attachments/accessories specified by the manufacturer.
11. Unplug this apparatus during lightning storms or when unused for long periods of time.
12. Refer all servicing to qualified service personnel. Servicing is required when the apparatus has been damaged in any way, such as power-supply cord or plug is damaged, liquid has been spilled or objects have fallen into the apparatus, the apparatus has been exposed to rain or moisture, does not operate normally, or has been dropped.
13. **WARNING:** To reduce the risk of fire or electric shock, do not expose this apparatus to rain or moisture. Do not place objects filled with water, such as a vase, on the apparatus.
14. The AC mains plug is used as the power disconnect device and it should remain readily accessible and operable during intended use. In order to completely disconnect the apparatus from power, the mains plug should be completely removed from the AC outlet socket.

What is Included

Project E Beauty
Portable Skin Care
Galvanic Device

Roller Probe


Ball Probe

USB Cable

Owner's Manual
(Not Shown)

Location of Controls

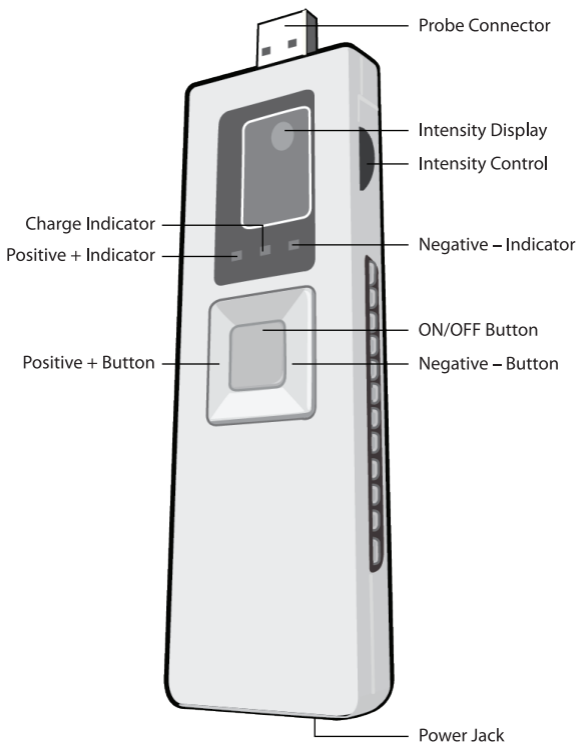Conditions for Use

WARNING!

Use of the device when the following conditions exist may cause serious damage to the skin. Do not perform a treatment if the person has any of these conditions.

Absolute Contraindications

- Cancer
- Acute infectious disease or fever
- Nausea
- Dysfunction or disorders of the nervous system
- After alcohol or drug consumption
- Inflammatory disease
- Pregnancy
- Epilepsy
- High or low blood pressure
- Heart disease
- Using a pacemaker
- Hemophilia
- Phlebitis or thrombosis

CAUTION!

Use of the device when the following conditions exist may cause potential damage to the skin. Do not perform a treatment on affected areas of the person's body.

Relative Contraindications

- Cuts or abrasions
- Vascular conditions
- Hypersensitive skin
- Excessive body hair
- Excessive metal fillings
- Heavy bridge work
- Metal pins or plates
- Sinus blockages
- Skin disorders or diseases
- Infections
- Sunburn
- Headaches or migraines

Conditions for Use

This device utilizes two feature functions as follows:

Desincrustation

Desincrustation is a deep cleansing treatment that softens sebum and is the ideal treatment for congested or oily skins. Desincrustation can be utilized as the on to prepare the skin for iontophoresis.

Iontophoresis

Iontophoresis is used to penetrate creams or solutions into the skin for a specific treatment. Iontophoresis is also used in cellulite treatments on the body to penetrate diuretic-based gels into the skin to stimulate lymphatic drainage.

This device has two probes that perform the following:

Roller Probe

The Roller Probe that is designed for the open face and neck areas.

Ball Probe

The Ball Probe is designed for working around the eyes and nose, giving you the best results.

How Often

Perform each treatment for 10 minutes four times a week to start. Once the goal is achieved, use only twice a week.

Benefits of Galvanic Therapy

Galvanic facials offer a safe and gentle alternative to dramatic plastic surgery procedures, laser resurfacing, chemical peels, Collagen and Botox injections and other invasive skin rejuvenation treatments.

Results can vary by individual and skin type, are gradual, and will not just occur over night. Although Galvanic treatment has been shown to produce an immediate and temporary firming effect, continued application will provide more cumulative long-term, lasting results. Galvanic treatments can give you the following benefits:

- Restores moisture to dry skin.
- Eliminates accumulation of oil that leads to acne.
- Hydrates skin and also increases the capacity of your skin to retain moisture.
- Wrinkles are made smooth and muscles are tightened.
- Anti-ageing properties.
- Increased Blood circulation, which nourishes the cells of the skin.
- Increases your skin's capacity to absorb nutrients.
- Aids penetration of skin care products into the deeper layers of your skin.
- Makes your skin look firm.

Use of the Device

Before performing a treatment you should do the following:

- Remove all jewelry and any metal objects.
- Check if absolute or relative contraindications exist (see page 8).
- Check that the device is functioning properly and the electrical cord is not damaged.
- Ensure there are not any flammable materials near the workspace.

Use of the Device

Charging


This device can be used connected to an AC outlet or can be charged to be used without the AC cord. To plug in and charge the unit, perform the following:

- 1 Insert the Adapter's plug into the Power jack on the bottom of the device.

- 2 Insert the USB cable to standard 5V charger(not included);
Insert the charger into a AC outlet.

The Charge LED will light to indicate the unit is charging. You can use the device while charging.

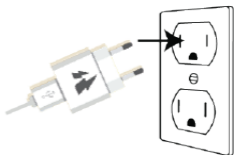Note: Recharge when the device stops working.

Use of the Device

Follow these steps when performing a treatment:

- 1 Disinfect the desired probe with a soft cloth and carefully cleaning.

- 2 To remove a probe, gently pull it off the USB-type jack.

To replace the probe, simply align the curved edge of the probe with the curved edge of the unit and gently push it onto the USB-type jack.

When using around the nose, eyes and edges of the mouth, use the Ball Probe, otherwise use the Roller Probe.


Use of the Device

3 Plug the unit in, or use a charged device (see page 11).

4 Thoroughly clean the face with soap and water, and dry.

5 Press the ON/OFF button to turn the unit on, the + and - indicators will alternatively blink.

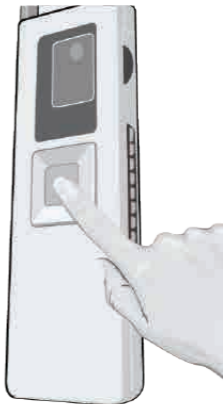6 Use an ion gel on areas you will be working on to achieve the best results.

***You can purchase our
Hyaluronic Acid Serum at
www.projectbeauty.com***

Use of the Device

- 7 Press the Negative (-) button to cleanse the skin (desincrustation) and prepare the skin for iontophoresis.

- 8 Adjust the Intensity control to a lower amount and then gradually increase as desired.

Note: The device is effective at all intensity levels including the lowest setting.

It is not necessary to use the highest power level to achieve maximum results.


- 9 To activate, hold the unit with your fingers touching the side electrodes **u** and place the whole probe **v** completely on the skin. Move in the directions as shown starting on the on the next page for best cleansing results.

When the current is passing through the skin, the + and - indicators (red and blue) will light.


Use of the Device


- Slowly move from eyebrows upwards to the top of the hairline.

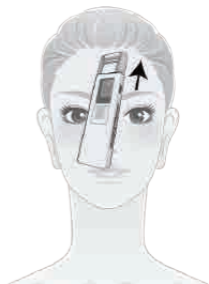- Slowly move from the bottom of the neck up the face and to the forehead.

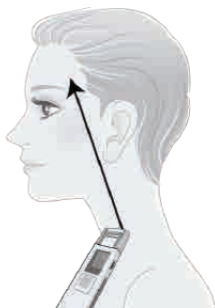- Slowly move from the nose across the cheek and towards the ear.

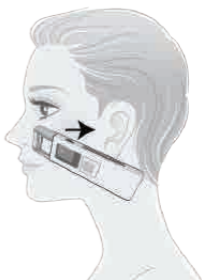
Use of the Device

- Slowly move down the chin to the neck.

For around the nose and eyes areas, make sure to use the ball probe.

Decrease the Intensity control to 0.5

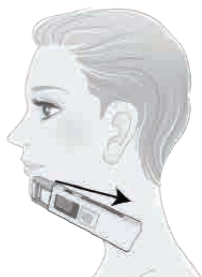- 10** For iontophoresis, apply an Ion gel on areas you will be working on to achieve the best results.

You can purchase our Hyaluronic Acid Serum at www.projectbeauty.com

- 11** Press the Positive (+) button.

Use of the Device

- 12** Adjust the Intensity control to a lower amount and then gradually increase as desired.

Note: The device is effective at all intensity levels including the lowest setting.


It is not necessary to use the highest power level to achieve maximum results.

- 13** To activate, again hold the unit with your fingers touching the side electrodes **u** and place the whole probe **v** completely on the skin. Move in the same directions as follows:

- Slowly move from eyebrows upwards to the top of the hairline.

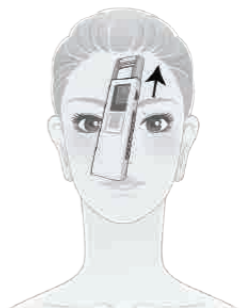
Use of the Device

- Slowly move from the bottom of the neck up the face and to the forehead.

- Slowly move from the nose across the cheek and towards the ear.

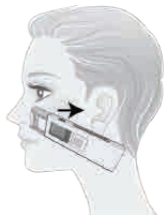- Slowly move down the chin to the neck.

Note: Always move upwards; do not move the unit downwards.

- Slowly move upwards from the chin to the bottom of the ear, then lift the device off the skin and move from the chin to the corners of the mouth.

Use of the Device

- Slowly move from the corners of the eyes to the temples to the top of the neck. Lift the device and move from the chin to the top of the neck and down to the armpits.

For around the nose and eyes areas, make sure to use the Ball Probe.

- 14** Press the ON/OFF button to turn the unit off.

Maintenance

- Properly unplug the device by holding the plug and removing it from the AC outlet.
- After each use, disinfect the electrode with a soft cloth.
- Immediately dry off any moisture, chemicals or water on the device.
- Always keep the device clean and sterilized.
- Store the device in a dry place.

Warranty

The product comes with a full one year warranty against manufacturer defects. It does not include damage caused by misuse, accident, abuse, neglect or other external causes. The warranty is void if an unstable voltage or inappropriate power supply has been used.

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