  
**PROJECT E**  
BEAUTY

**LED RF Ultrasonic  
Slimming Device  
User Manual**

PE142

[www.projecteauty.com](http://www.projecteauty.com)

## Introduction

**Thank you** for purchasing the Project E Beauty LED RF Ultrasonic Slimming Device. This device utilizes 3 different functions to treat a variety of skin and body conditions.

What separates this device from other beauty devices is the extra large ultrasound probe. Combine with Radio Frequency and Red LEDs, this takes skin care and body shaping to the next level and you will achieve results more effectively and faster than imaginable.

Features:

- 4-point focused RF
- 330K Hz Ultrasonic Technology
- USB Rechargeable
- Easy to Operate
- Noise Free

### **Ultrasonic Slimming**

Ultrasonic(Ultrasound) slimming technology is a process where you transmit waves of ultrasonic into cellulite affected areas of your body. The frequency of the ultrasonic is set to disrupt the fat cells that are causing the cellulite to appear. The continuation of the cavitation therapy causes the fat to get mopped up within the lymphatic system and disposed of via the liver. The technology is very effective and a real technological breakthrough for getting rid of cellulite

## Introduction

### RF (Radio Frequency)

For detoxicating and promoting blood circulation. With RF, the electrical energy will follow the shortest distance between two electrodes. When the electric current runs through living tissue, the high energy flow produces heat, as the living tissue will resist electricity. Tissue in both the dermal and subcutaneous layers of the skin can be heated by RF energy. As well as increasing production of collagen and elastin, RF can also be used to contract collagen fibers, resulting in skin tightening. It is also used to improve lymphatic drainage, accelerate fat metabolism and for cellulite reduction.

### RED LED

Red LED penetrates the skin approximately 8-10mm, which will stimulate fiber cells and the production of collagen, reducing wrinkles and making skin more elastic and smooth.

## Safety Information

Use of the device when the following conditions exist may cause potential damage to the skin. Do not perform a treatment on affected areas of the person's body.

### Relative Contraindications

- Cuts or abrasions
- Vascular conditions
- Hypersensitive skin
- Excessive body hair
- Sinus blockages
- Skin disorders or diseases
- Infections
- Sunburn
- Headaches or migraines

## Safety Information

You must read, and understand, all of the safety warnings, operating and maintenance instructions before using the Project E Beauty LED RF Ultrasonic Slimming Device.

**WARNING!:** A hazard that may cause serious injury.

**NOTICE:** A hazard that may cause damage to the product.

### **WARNING!**

- Performing a treatment on a person with absolute contraindications may cause serious injury or death. Do not perform a treatment on a person where absolute contraindications exist.
- Never use the device around water. Always dry your hands before performing a treatment.
- Use of the device around jewelry and metal may cause unexpected electrical shock. Always remove all jewelry and metal prior to a treatment.

Do not perform a treatment if the person has any of these conditions.

#### **Absolute Contraindications**

- Acute infectious disease or fever
- Nausea
- Dysfunction or disorders of the nervous system
- Epilepsy
- Pregnant
- Phlebitis or thrombosis
- Heart disease

**NOTE:** If you are under a doctor's care for any skin condition or if you think you may have a skin condition, please consult your physician before use.

## What is Included


Project E Beauty LED RF  
Ultrasonic Slimming  
Device

USB Cable

Owner's Manual (Not Shown)

## Location of Controls
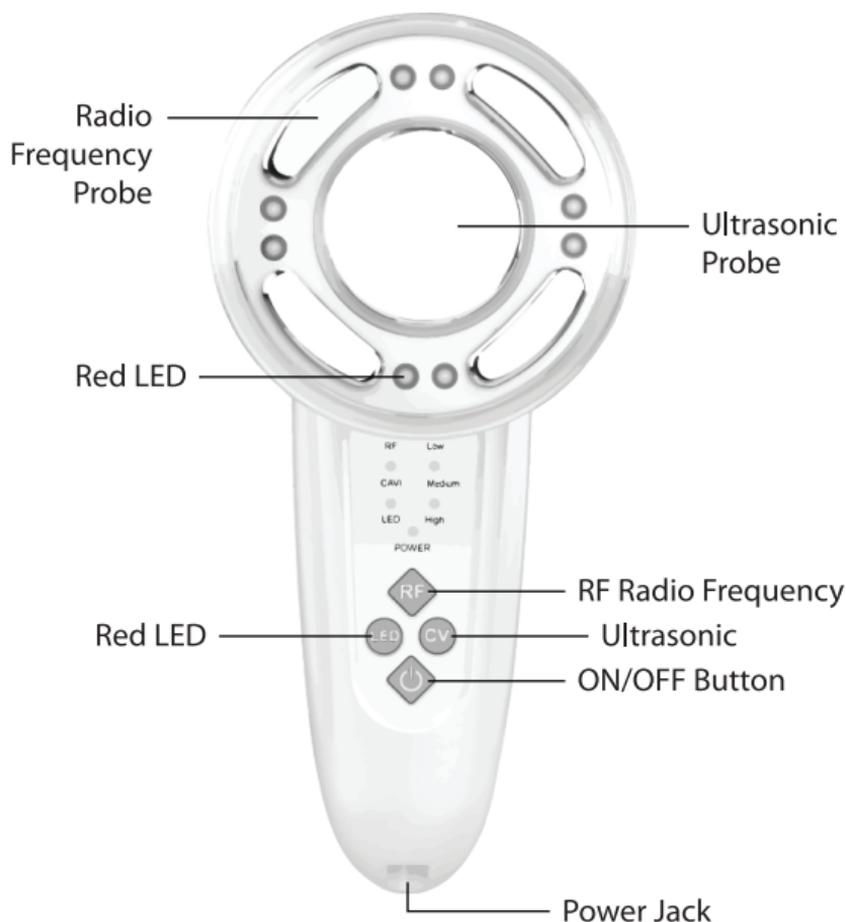
## Benefits LED RF Ultrasonic Slimming Device

**Radio Frequency Three different modes** allow the users to perform treatments on different skin types and areas, and keep them from being damaged.

**Extra large Treatment Probe** can be used to perform treatments on numerous areas, such as arms, tight or body.

# Benefits LED RF Ultrasonic Slimming Device

## How Often

Ten to fifteen minutes for each treatment.

One to two treatments a week is recommended.

## Use of the Device

Before performing a treatment you should do the following:

- Check if absolute or relative contraindications exist.
- Check that the device is functioning properly and the electrical cord is not damaged.
- Ensure there not any flammable materials near the workspace.

Follow these steps when performing a treatment:

- 1 Before your first use of the device, charge it at least 2 hours. (Note:  
Flashing Power LED: Charging  
Steady Power LED: Fully Charged

For charging:

Simply connect the USB cable to  
common 5V USB power source

Connect the cable to the charging  
jack on the massager.

### NOTICE

The massager **CANNOT** be used during charging

## Use of the Device

- 2** Disinfect the device using a soft cloth and carefully cleaning the electrodes and LED cover

- 3** Press and hold the On/Off button for 2 seconds to turn the unit on.

- 4** Press the LED button to turn on Red LED function (continuous mode).

Press again to switch to the Pulse mode. It provides short specific pulse of light energy, which enhances the treatments.

## Use of the Device

- 5 Press the CV button to turn on ultrasonic slimming function.


- 6 Press the RF button to turn on Radio Frequency function.  
Repeatedly press the RF button to select the desired RF power.

**Note:** The device is effective at all intensity levels including the lowest setting. It is not necessary to use the highest power level to achieve maximum results.

- 7 If desired, use Slimming Essence, Cream or your own products on area you will be working on to achieve the best results.

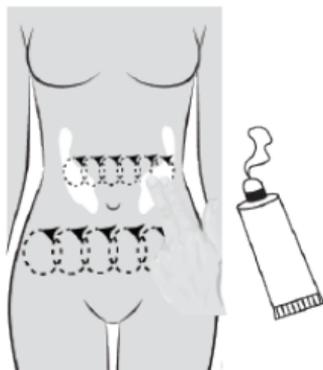
## Use of the Device


### **BODY:**

Approximate treatment time is 5-10 minutes for each area (Thighs/Legs, Buttocks, Back, Arms and Abs). If desired, use a Slimming Essence on areas you will be working on to achieve the best results.

**ABS/STOMACH:** It is recommended that the complete Ab/Stomach treatment should take 10 minutes. Slowly move in little circles, across the lower abdomen. Lift up and repeat for the upper abdomen.

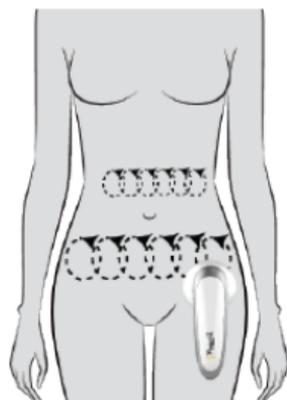
Slowly move upwards from the top of the buttocks towards the lower back. Lift up and repeat on another area with the same movement.

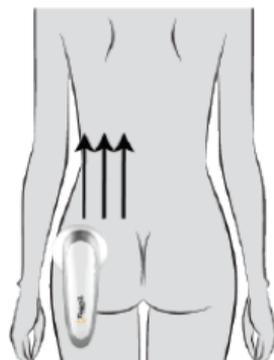**ARMS:** It is recommended that the complete arm treatment should take 5 minutes.

Slowly move in little circles, upwards from the back of the elbow to the inside shoulder. Lift up and do the outer shoulder.


Repeat for the front of the arms.
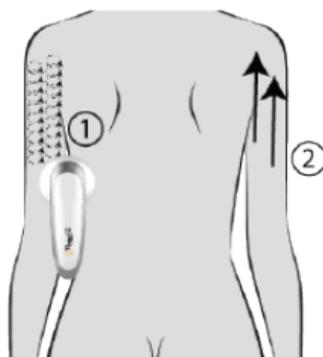
## Use of the Device


**BUTTOCKS:** It is recommended that the complete buttocks treatment should take 5 minutes.

Slowly move in little circles, upwards from the upper rear thigh to the top of the buttocks. Lift up and repeat on another area of the buttocks.

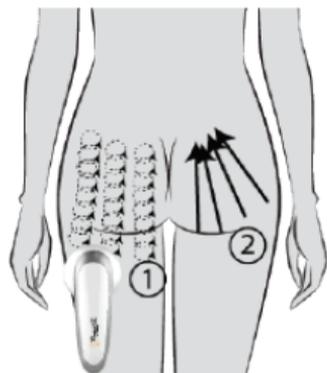Slowly move up and down in a small zig-zag motion on the buttocks.

Repeat for other cheek.

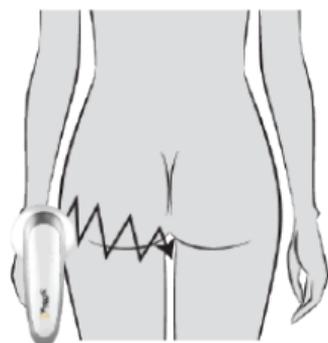**LEGS:** It is recommended that the complete leg treatment should take 10 minutes.

Slowly move in little circles, upwards from above the back of the knee to the top of the outer thigh. Lift up and do the inner thigh.

## Maintenance

- After each use, clean the surface using a soft cloth. Immediately dry off any moisture, chemicals or water on the device
- Store the device in a dry place.
- Keep out of reach of children

## Warranty

The product comes with a full one year warranty against manufacturer defects. It does not include damage caused by misuse, accident, abuse, neglect or other external causes. The warranty is void if an unstable voltage or inappropriate power supply has been used.

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