
D'arsonval High Frequency Wand

Neon & Argon

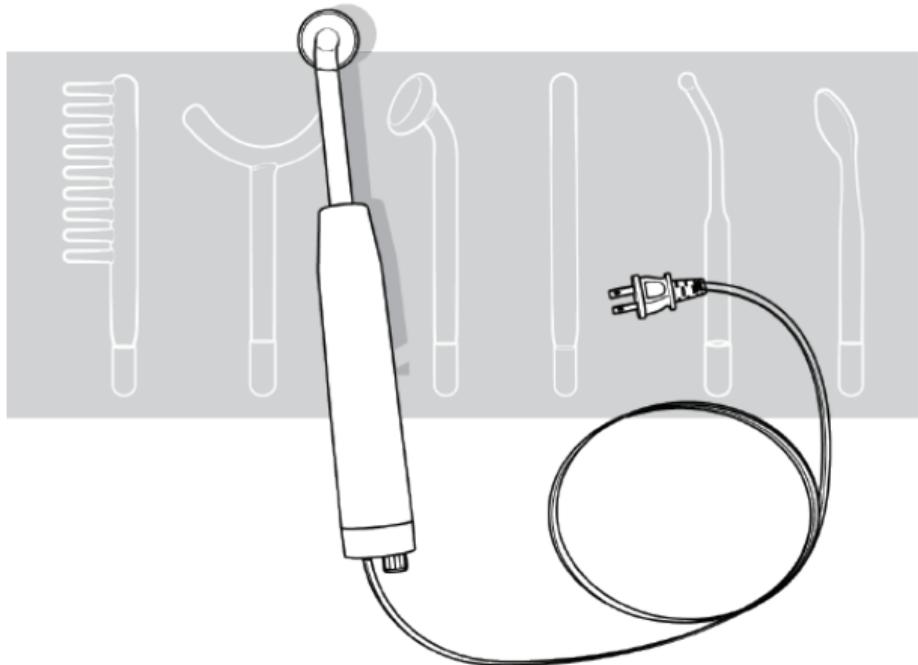OO
PROJECT E
BEAUTY

www.projectebeauty.com

Introduction

Thank you for purchasing the Project E Beauty D'arsonval High Frequency Wand Neon & Argon Device. This device utilizes high frequency electrical current to treat a variety of skin conditions.

The use of high frequency current to treat medical conditions was first developed in the late 1800's by renowned scientist Nikola Tesla. Towards the end of the century, the French biophysicist D'Arsonval introduced the use of high frequency currents to treat diseases of the skin and mucous membranes. Between 1970 and 1980 the treatment became more popular and is now widely used by skin care professionals throughout the world. This safe, non-invasive technology can be beneficial in the treatment of many skin conditions ranging from acne to wrinkles. High frequency skin care treatment has been found to enhance oxygenation, stimulate dermal tissue, kill bacteria causing acne, and treat oily or combination skin problems. It is also effective as a treatment for thinning hair and hair loss. Although results can vary widely by individual and skin type, high frequency skin care treatment is appropriate for all types of skin. Typically, the treatment has better results on oily skin.

Cosmetologists and skin care specialists agree the high frequency treatments, when properly performed, have no side effects. It is safe and beneficial to combine the high frequency treatments with other gentle cosmetic technologies such as light therapy, galvanic treatments, microcurrent, ultrasound, LED, Radio Frequency, etc. The high frequency will enhance the results of those applications.

High frequency skin care treatment is an important component of your skin care regime and should only be used for cosmetic purposes. It is not intended in any way to diagnose, treat, cure, or prevent any disease.

Safety

You must read, and understand, all of the safety warnings, operating and instructions before using the Project E Beauty Portable High Frequency Skin Care Device. Failure to follow the warnings and instructions may cause serious injury, damage to the device or void the warranty.

Warnings descriptions are as follows:

WARNING!: A hazard that may cause serious injury or death.

CAUTION!: A hazard that may cause minor injury.

NOTICE: A hazard that may cause damage to the product

⚠ WARNING!

- Performing a high frequency treatment on a person with absolute contraindications may cause serious injury or death. Do not perform a high frequency treatment on a person where absolute contraindications exist.
- Use of the device near flammable materials may cause explosion, fire and serious injury or death. Always check the work space and remove any flammable materials prior to performing a high frequency treatment.
- Use of the device around water (i.e., bathtubs, wash basins, areas of high humidity) or with wet hands may result in electrical shock, serious injury and death. Never use the device around water. Always dry your hands before performing a high frequency treatment.
- Improper use of the device by unqualified people (i.e., children) may result in burns, electrical shock or other serious injury or death. Never allow unqualified people to use this device for any purposes or play.
- Performing a high frequency treatment with a damaged device or broken electrode may cause serious injury or death. Always check the device and electrodes for damage before performing a high frequency treatment. Never perform a treatment with broken electrodes or a damaged device.
- Use of the device internally (i.e., mouth, nose etc.) may cause serious injury or death. Never use the device internally.


⚠ CAUTION!

- Performing a high frequency treatment near areas of relative contraindications may cause minor injury. Avoid performing a high frequency treatment near areas where relative contraindications exist.
- Performing a high frequency treatment with AHA or Glycolic Acid products may cause unwanted stress to the skin. Do not use AHA or Glycolic Acid products when performing a high frequency treatment.
- Performing a high frequency treatment around broken capillaries, spider veins or rosacea may cause unwanted stress to the skin. Do not perform a high frequency treatment around broken capillaries, spider veins or rosacea when performing a high frequency treatment.
- Use of the device around jewelry and metal may unexpected electrical shock. Always remove all jewelry and metal prior to a high frequency treatment.

NOTICE:

- Improperly unplugging the plug when the device is not in use may cause damage to the cord resulting in fire or electrical shock. Always unplug the device by holding the plug and removing it from the electrical socket.
- Use of the device around water (i.e., bathtubs, wash basins, areas of high humidity) may cause damage to the device. Never use the device around water.
- Improper storage or maintenance of the probe and electrodes may result in damage and void the warranty.

Parts

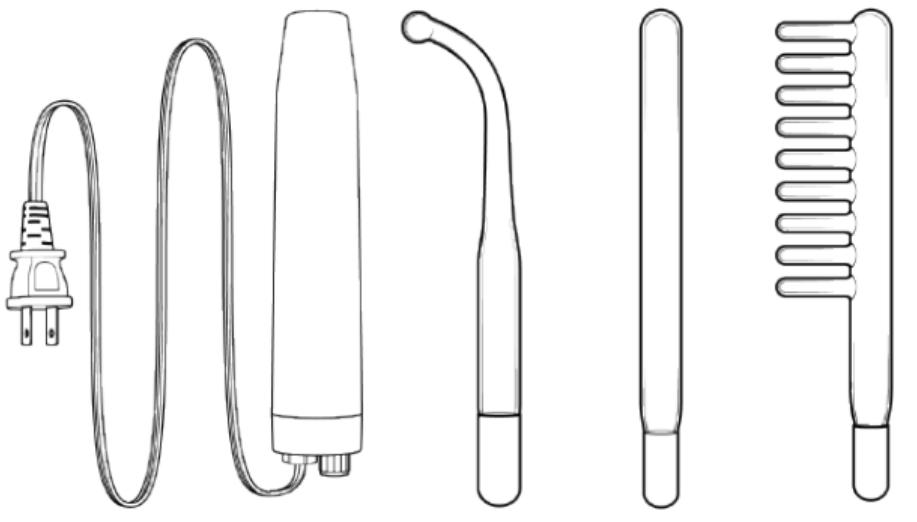
Probe

Rod (bulb)
electrode

Long electrode
(saturator)

Comb
electrode

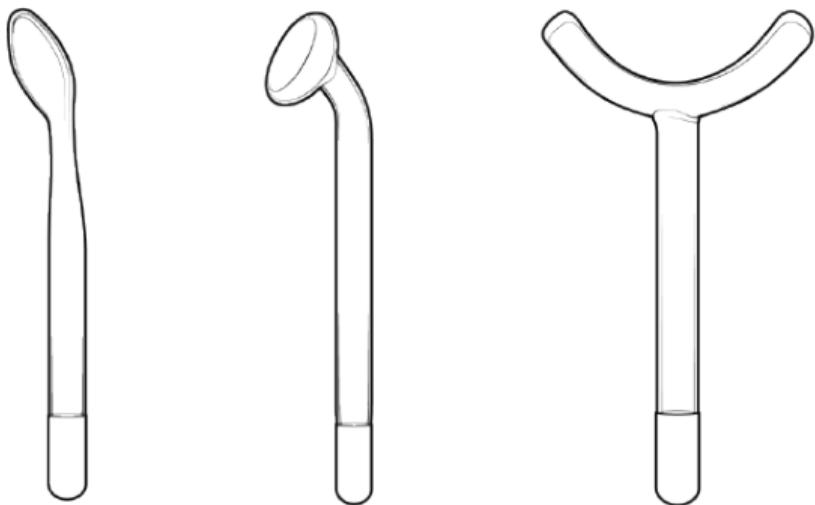Spoon
electrode

Mushroom
electrode

Wide (horseshoe)
electrode

Conditions for Use

⚠ WARNING!

Use of the device when the following conditions exist may cause serious damage to the skin. Do not perform a high frequency treatment if the person has any of these conditions.

Absolute Contraindications

- Cancer
- Acute infectious disease or fever
- Nausea
- Dysfunction or disorders of the nervous system
- After alcohol or drug consumption
- Inflammatory disease
- Pregnancy
- Epilepsy
- High or low blood pressure
- Heart disease
- Using a pacemaker
- Haemophilia
- Phlebitis and thrombosis

⚠ CAUTION!

Use of the device when the following conditions exist may cause potential damage to the skin. Do not perform a high frequency treatment on affected areas of the persons body.

Relative Contraindications

- Cuts or abrasions
- Vascular conditions
- Hypersensitive skin
- Excessive body hair
- Excessive metal fillings
- Heavy bridge work
- Metal pins or plates
- Sinus blockages
- Skin disorders or diseases
- Infections
- Sunburn
- Headaches or migraines

Only use the High Frequency Skin Care

Conditions for Use

Device for the following conditions:

Acne

The high frequency treatment can reduce, and balance, sebum (oil) production by normalizing the function of the sebum glands. The treatment loads the skin with oxygen, which cleanses the skin of unwanted toxins and prevents the growth of acne causing bacteria that causes inflammation. The high frequency treatment makes the skin more receptive to acne lotions, creams and other skin treatment products. When used with an effective acne treatment lotion, regular application of high frequency can keep acne away. Each treatment builds upon the last and the results are cumulative. Treatment may vary by individual skin type and concern.

Cellulite

The high frequency treatment can help diminish incunabula cellulite (first stages) by supplementing oxygen to the skin cells. On later stages of cellulite, the treatment can be used with anti-cellulite creams, lotions or other remedies to substantially increase penetration of the anti-cellulite components. High frequency treatment will not eliminate the problem, but with regular use, the high frequency treatment can be very effective in reducing the appearance of cellulite when used in conjunction with quality anti-cellulite products.

It should be noted that the skin care industry uses this technology to aid in lymphatic drainage, gently exfoliate

the skin, increase blood circulation, assist in the production of collagen and elastin and promote healthy cell metabolism.

Wrinkles and Aging Skin

The high frequency treatment can reduce wrinkles and improve skin tone by intensifying oxygen production, improving blood circulation and stimulate lymph drainage. The electrical impulses from the high frequency treatment produce an enriched form of oxygen, which can provide the skin with a firm, youthful, vibrant glow. They increase blood circulation, which in turn nourishes the skin's surface and renews underlying cells. The impulses also drain the congestion built up in the lymphs and reduce swelling and inflammation. Using the device with quality lotions and creams will help with the penetration of nutrients from your skin care products and improve the results of the treatments.

Hair Thinning or Loss

The lack of blood circulation is one of the main reasons for thinning hair and hair loss. High frequency treatments stimulate blood circulation which in turn provides nutrition to hair follicles. By stimulating the scalp tissue and glandular activity you can prevent further hair loss. Using the device will improve penetration and overall effectiveness of post-treatment hair growth formulas and aid in promoting healthier hair growth.

Use of the Device


Treatment Preparation

Before performing a high frequency treatment you should do the following:

- Remove all jewelry and any metal objects
- Check if absolute or relative contraindications exist
- Determine if the use of gauze is needed.

Note: Gauze is used to prevent the electrode from skipping or dragging. The face must be clean.

- Clean the electrodes
- Check that the device is functioning properly and the electrical cord is not damaged.
- Ensure there are not any flammable materials near the work space.
- Treatment Application Methods
 - There are two methods of high frequency treatment: the direct (or contact) method and the indirect (or distance) method. The direct method is the method you will use for most treatments. It is used to stimulates, sanitize and heal the skin. It generates warmth in the skin and stimulates nerve endings, which increases blood circulation and nourishes the skin. The indirect method is typically used for relaxation and sedation. Use the indirect method on scar tissue (4 to 9 months old), aging skin that is too delicate for massage and when you want to aid the penetration of creams and lotions.
 - If it is too tight to snap in the electrodes, please wrap a towel on the electrode and try to push harder to snap into the handle.
 - The electrodes may need to be activated if they are not able to light up, please try the following steps:
 1. Insert the electrode on high frequency handle.
 2. Turn on the device.
 3. Touch the electrodes to any metal object (the metal helps activate the electrodes)

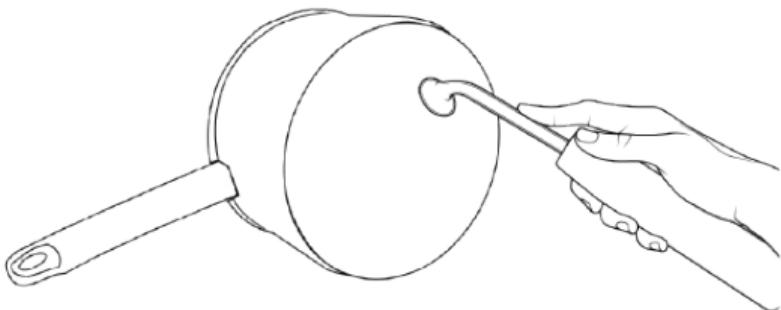Use of the Device

Using the Skin Care Device


Use this chart to choose the electrode for the area being treated.

Electrode	Area of Treatment	Condition	Timing
Rod (bulb)		For small targeted areas such as isolated acne pimples	Acne 5 - 15 minutes depending on the skin condition and reaction to the treatment
Spoon		For small areas around the eyes, nose and mouth	Wrinkles, Crow's feet, Papules, Pustules 3 - 5 minutes
Mushroom		For large areas including cheeks, forehead, chin and neck	Wrinkles, Papules, Pustules 3 - 5 minutes
Comb		Scalp area	Hair thinning, loss 3 - 5 minutes
Horseshoe (wide)		For the neck and shoulders	Varicosis treatment 3 - 5 minutes
Straight, Long Electrode (saturator)		Delicate skin, Face (use for the Indirect Method treatments)	Aging and loose skin 5 - 15 minutes depending on the skin condition and reaction to the treatment

Use of the Device

Follow these steps when performing a treatment :

1. Disinfect the electrode with a soft cloth.
2. Check that the probe is turned Off and plug the probe into the electrical outlet.
3. Insert the glass electrode into the probe; it should just snap in to place. Do not push or squeeze the electrode very hard in order to prevent the glass from breaking. Important! Insert the electrode BEFORE turning the probe On. Always turn Off the device BEFORE you remove the electrode.
4. If necessary, wrap the cord around your arm once to keep the cord from dragging across the client's skin.
5. Turn the probe ON using the On\Off knob at the end of the probe.
6. The electrode should begin to glow. Rotate the knob to the desired intensity level.
7. To begin your treatment, hold the probe and make contact with the electrode and your finger BEFORE placing it on to the client's skin. Important! Do not remove your finger until the electrode is in contact with the client's skin. Removing the finger too early creates an uncomfortable spark effect on the client's skin.

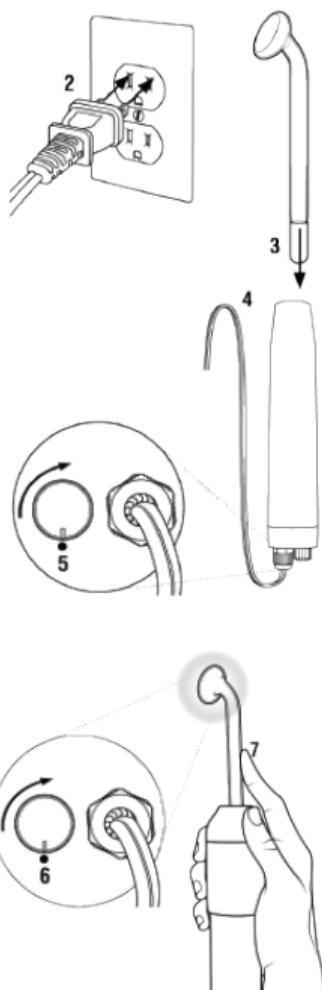Note: This device utilizes the high frequency electrical current which comes through the gas filled glass electrode. The characteristics of the current are high frequency (110-400kHz) (100,000 - 400,000 hertz), high voltage (about 20kV) and low amperage (no more than 100-200mA). The device is effective at all intensity levels including the lowest setting. It is not necessary to use the highest power level to achieve maximum results.

Use of the Device

The key to the best results is treatment consistency.


8. Place the electrode on the client's skin and remove your finger. Increase the intensity to a comfortable setting.
9. Move the electrode in small circles around the skin for a maximum of five minutes. Important! The device should be used for a maximum of 20 minutes per session with 10-15 min intervals. The overuse will void the warranty.
10. Place your finger back on the electrode and remove the electrode from the client's skin.

Note: Treat problem areas, such as blemishes, by "sparking" the spot with the edge of the electrode. Place your finger on the electrode, touch the electrode to the spot and remove your finger. Quickly replace your finger on the electrode, remove the electrode from the face and repeat in another area.

11. After the treatment is finished turn Off the On\Off knob and disconnect the power cord.
12. Carefully remove the electrode from the socket and disinfect the electrode with a soft cloth. Important! Do not immerse the electrode in water or get the metal end wet.

Dry the electrode with a clean towel and place into the storage compartment.

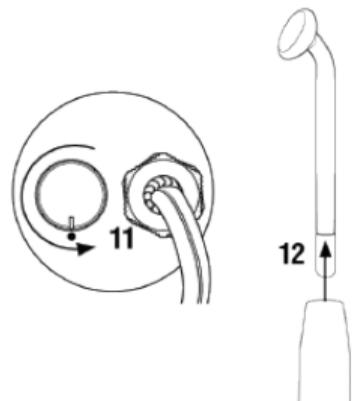Note: Argon gas Violet electrode is used for acne's skin conditions.

Neon gas Orange electrode is used for wrinkles and skin tightening.

Maintenance

- Properly unplug the device by holding the plug and removing it from the electrical socket
- Immediately dry off any moisture, chemicals or water on the device
- Always keep the electrode clean and sterilized
- To avoid damage always store the electrodes in the protective storage compartment
- Store the device in a dry place

Product Specifications

Voltage.....	110V-240V
Power	1.2 - 6.5W
Operate Frequency.....	Up to 250,000Hz
Electrode Diameter.....	11mm
Electrode Color.....	Orange X 4, Violet X 2

Warranty

Activate Your 1 Year FREE Warranty Now.

The product come with a full one year warranty against manufacturer defects. It does not include deamge caused by misuse, accident, abuse, neglect or the external causes. All we need is your name, email and order ID. That's all.

Scan the QR code or visit
www.projectebeauty.com/pages/warranty

Questions

Website: www.projectebeauty.com

Email: sales@projectebeauty.com

Tel: +1-626-759-4232

Our customer care expert team are here to answer any question you might want to ask.

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