PROJECT E
BEAUTY

**Vacuum Blackhead
Remover
User Manual**

PE711

www.projecteauty.com

Introduction

Thank you for purchasing the Project E Beauty Vacuum Blackhead Remover. This device utilizes 5 different Heads to treat a variety of skin conditions.

What separates this device from other beauty devices is the extra strong vacuum suction feature. This function takes skin care to the next level and you will achieve results more effectively and faster than imaginable.

Functions:

- 3 Different Power Suction Levels
- 5 Different Vacuum Suction Heads
- USB Rechargeable
- Easy to Operate
- Low Noise

What you need to know about blockage and inflammation of the facial sebaceous glands.

Above-average secretion by the sebaceous glands leads to oily skin. This over-activity is caused by hormonal changes, especially during puberty, inappropriate highfat or high-protein diet, excessive consumption of alcohol, etc. This over-production of sebum causes the pores to widen. Oily skin can have a shiny or dried appearance, depending on whether the flow of sebum into the pore has an oily consistency or has already dried within the follicle. Oily skin is usually "impure", because it provides an excellent feeding ground for bacteria. The oxygen supply within a healthy pore normally restricts the spread of bacteria. There is a correlation between this bacterial spread and a pore becoming visible to the eye. If the pore becomes blocked,

Introduction

the environment becomes anaerobic (i.e. oxygen-free), ideal for *Propionibacterium acne*, which spreads like wildfire. This leads to inflammation and blemishes.

Normal pore:

A normal amount of oil is secreted from the pore.

Blocked pore:

Oil trapped in the pore oxidizes and solidifies, causing clogged pores and blackheads.

Blackheads should not be removed by squeezing with fingertips or fingernails, as this can damage the skin and cause further inflammation. Cleansing two or three times a week is recommended.

Safety Information

Use of the device when the following conditions exist may cause potential damage to the skin. Do not perform a treatment on affected areas of the person's body.

Relative Contraindications

- Cuts or abrasions
- Vascular conditions
- Hypersensitive skin
- Excessive body hair
- Sinus blockages
- Skin disorders or diseases
- Infections
- Sunburn
- Headaches or migraines

Safety Information

You must read, and understand, all of the safety warnings, operating and maintenance instructions before using the Project E Beauty Vacuum Blackhead Remover.

WARNING!: A hazard that may cause serious injury or death.

NOTICE: A hazard that may cause damage to the product

WARNING!

- Performing a treatment on a person with absolute contraindications may cause serious injury or death. Do not perform a treatment on a person where absolute contraindications exist.
- Never use the device around water. Always dry your hands before performing a treatment.
- Never use the device internally.

Do not perform a treatment if the person has any of these conditions.

Absolute Contraindications

- Acute infectious disease or fever
- Nausea
- Dysfunction or disorders of the nervous system
- Epilepsy
- Pregnant
- Phlebitis or thrombosis

NOTE: If you are under a doctor's care for any skin condition or if you think you may have a skin condition, please consult your physician before use.

What is Included

Project E Beauty
Vacuum Blackhead
Remover

Small Round
Suction Tip

Medium
Round
Suction Tip


Larger
Round
Suction Tip

Oval Suction
Tip

Microcrystalline
dead skin
exfoliation Tip For
Sensitive Skin

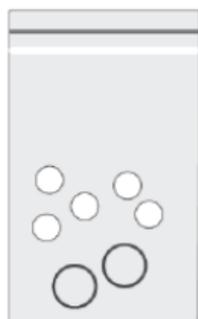Filter X 5 &
Seal Ring X 2

USB Cable

Owner's Manual (Not Shown)

Location of Controls

Benefits Vacuum Blackhead Remover

Three Different Modes allow the users to perform treatments on different skin types and areas, and keep them from being damaged.

Five Different Tips can be used to perform treatments on numerous areas, such as facial, arms or body.

Benefits Vacuum Blackhead Remover

How Often

Two to three minutes for each treatment.

One to two treatments a week is recommended.

Use of the Device

Before performing a treatment you should do the following:

- Check if absolute or relative contraindications exist.
- Check that the device is functioning properly and the electrical cord is not damaged.
- Ensure there are not any flammable materials near the work space.

Follow these steps when performing a treatment:

- 1 Before you first use the device, charge it at least 2 hours. (Note: While charging, the red light will be turned on. Once it is fully charged, the light will turn blue)

For charging:

Simply connect the USB cable to any USB power source.

Connect the cable to the Massager.

Use of the Device

2 Chose the Right Tips

For weak and gentle suction, it is used to remove small blackhead at dedicated area, best for thin, tender and sensitive skin.

Small Round Suction Tip

Medium Round Suction Tip

For strong and powerful suction, it is the best to be used to remove blackhead, grease, cosmetic residue etc.

For Body Use:

Promotes blood circulation speeds up muscle recovery. The suction allows new blood to flow into that area of tissue to enhance circulation.

Larger Round Suction Tip

Oval Suction Tip

For skin lifting, pulling it upwards to treat fine lines. This increases skin's elasticity and firmness.

This probe exfoliates dead skin that is left on the outer layer of the face leaving skin feeling instantly smooth.

Usage area: Any area on the face except around the eyes. Other areas such as knees and elbows are suitable.

Microcrystalline dead skin exfoliation Tip For Sensitive Skin

Use of the Device

- 3** Disinfect the tips with 70% alcohol using a soft cloth.


Thoroughly clean the face with soap and water, and dry. Do **NOT** use a scrub cream or exfoliator.

- 4** Use facial steamer or hot towel or blackhead remover gel/cream on your treatment area for 5 minutes to completely open your pores.

NOTICE

For best results the pores of the skin should be enlarged.

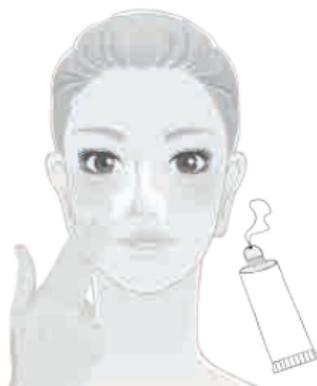- 5** Press the "On/Off" button to turn on the device.

Repeatedly press to increase the intensity

Level 1: Blue Light

Level 2: Green Light

Level 3: Red Light (Strongest suction)

Use of the Device


6 Testing:

Hold and move the suction probe around on your palm to feel the suction power and adjust strength level before using on the skin.

7 Move the appliance slowly around the tip of the nose. Hold the probe over a blackhead area to remove it.

Caution:

Always move the device around your skin slowly and frequently. Instead of moving it back and forth, move it in one direction and then reset and start again. Do not hold the machine in one place for more than 2 seconds or swelling, irritation and discomfort may occur.

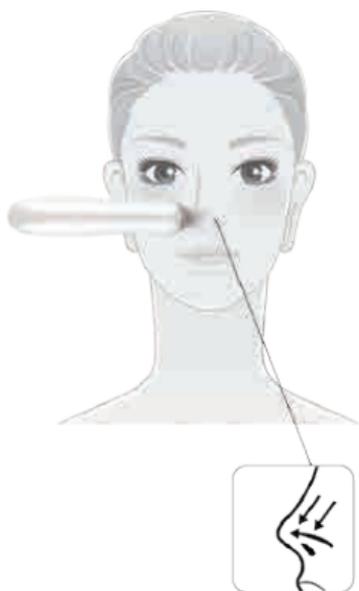8 After 5 minutes, the device will automatically turn off. To turn the unit off before 5 minutes, repeatedly press the Power button until the unit turns off.


⚠ NOTICE

Use cold steamer or cold towel together with moisturizer or other skin care products to shrink your pores after the treatment.

Use of the Device

- 9
1. Pull to remove the suction probe.
 2. Pull out the inner filter cap.
 3. Replace the filter.

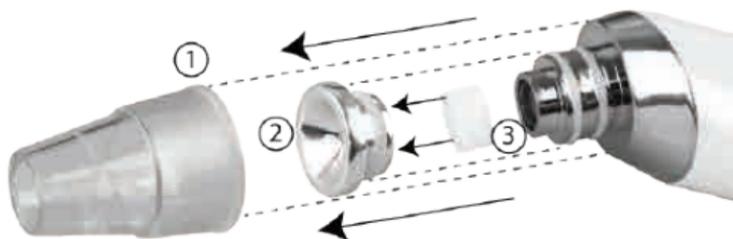⚠ NOTICE

After finish the treatment, the skin may feel tingly and flushed. This is normal because of the blood circulation and will take about 30-60 minutes to disappear.

Maintenance

- After each use, clean the surface using a soft cloth. Immediately dry off any moisture, chemicals or water on the device.
- Store the device in a dry place.
- Keep out of reach of children.

Warranty

The product comes with a full one year warranty against manufacturer defects. It does not include damage caused by misuse, accident, abuse, neglect or other external causes. The warranty is void if an unstable voltage or inappropriate power supply has been used.

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