SmoothPro+

IPL Hair Removal Device

PE720

www.projectebeauty.com

Introduction

Congratulations on purchasing Project E Beauty's SmoothPro+ (IPL Hair Removal Device). SmoothPro+ uses **Intense Pulsed Light** (IPL) technology which is known as one of the most effective methods to treat unwanted body hair with longer-lasting results. SmoothPro+ produces gentle pulses of light that penetrate the skin and are absorbed by the hair. This causes the hair follicle to go into a 'resting phase' and prevent the hair from growing any further. A few days later, the treated hair falls out naturally. By repeating the treatment according to the recommended schedule, you can continue to prevent hair regrowth.

There are so many different methods of hair removal, from traditional shaving blades to professional laser treatments. Many are time-consuming, expensive or horribly painful. SmoothPro+ is **FDA cleared** and gives you the longer-lasting results with reasonable price. You can now use this innovative light-based hair regrowth prevention technology comfortably in the privacy of your own home, whenever and wherever you want. For important information about safety and how to use the device, please read through this user manual before you use SmoothPro+.

Working Principle

The appliance works by heating up the hair and root beneath the skin. The melanin in the hair and hair root absorbs the released light pulses. The darker the colour of the hair, the more light can be absorbed. This process stimulates the hair to go into the resting phase. The hair now sheds naturally and the hair growth is inhibited.

Hairs can take one or two weeks to shed. Optimal hair removal cannot be obtained with a single session, since only hairs in the growing phase are susceptible to treatment with light. During a treatment, hairs on the treated area are usually in different stages of the hair growth cycle. To effectively treat all hairs you have to repeat the treatment every two weeks for a few times.

Benefits of SmoothPro+

This device provides the following benefits:

Intense Pulsed Light (IPL) technology for home use:

SmoothPro+ uses a light-based technology called Intense Pulsed Light (IPL). IPL is also used in the professional beauty market for hair removal and has been in use for the last fifteen years. SmoothPro+ now brings this innovative technology to a new level of home use.

IPL technology offers superior safety with low energy level:

IPL technology can achieve hair removal with energy used only at a fraction of the energy used in other light-based hair removal equipment. The low energy level used in SmoothPro+ reduces its potential to cause harm or complications, and contributes to your overall safety.

Gentle treatment, even when used on sensitive body area:

SmoothPro+ has been developed in close cooperation with leading dermatologists for effective and gentle treatment, even when used on the most sensitive body area.

Benefits of SmoothPro+

SmoothPro+ is safe for your eyes:

SmoothPro+ has a built-in safety feature for eye protection and it has been tested to international standards to ensure that its outputs are safe for your eyes.


SmoothPro+ protects your skin:

Part of the unique beauty experience with SmoothPro+ is the Detect and Set feature. This feature means you do not have to worry about choosing skin tone settings as the device automatically detects for you. It is simple and you can see skin tone setting indicated by the Power button rear side of the smart skin tone sensor window.

Is SmoothPro+ Suitable for You?

SmoothPro+ works easily, safely and effectively on a wide variety of hair colors from naturally dark blonde to brown and black; and skin tones from very white to dark brown. IPL technology cannot be used to treat white/gray, light blonder or red hair, and is not suitable for very dark skin. This is due to the high contrast required between the pigment in hair color and the pigment in skin tone.

Suitable for skin tone:

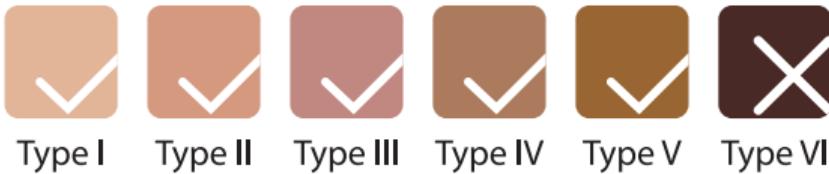Note: In the figure above, "X" means that the device is not suitable for you

Is SmoothPro+ Suitable for You?

The color charts below can help you understand whether SmoothPro+ is perfect for you!

Skin Type	Or Skin Color	Characteristics	Example
I	Pale white; fair; red or blond hair; blue eyes; freckles	Always burns, never tans	
II	White; fair; red or blond hair; blue, hazel or green eyes	Usually burns, tans with difficulty	
III	Cream white; fair with any eye or hair color; very common	Sometimes mild burn, gradually tans	
IV	Brown; typical Mediterranean caucasian skin	Rarely burns, tans with ease	
V	Dark brown; mid-eastern skin types	Very rarely burns, tans with easily	
VI	Black	Never burns, tans with easily	

Type I

Type II

Type III

Type IV

Type V

Type VI

Safety Information

Precautions:

- Read all instructions before using your Hair Removal Device. Check the device before using.
- Stop immediately, if you experience any adverse effect or have a concern regarding use.
- Use this device only for its intended use as described in this manual.
- Use with safety goggles provided with Hair Removal device.
- Do not look directly or shine light in the eyes while in operation, it may cause vision damage. Please wear provided goggles.
- Do not use if you are pregnant.
- Do not use if you have experienced a migraine in the last two years.
- Do not use over any suspicious or cancerous lesion. Consult your doctor.
- Do not allow children to use the device unsupervised.
- Do not expose the appliance to temperatures lower than 15°C or higher than 35°C during use.
- It is not recommended those people who have physical, mental disease or are under the treatment of dermatology to use this device.

Safety Information

Precautionary Instructions:

When using an electrical device, basic precautions should always be followed, including the following:

- Do not flash against any other surface than the skin. This can cause severe damage to the attachment.
- Do not open, disassemble, drop or insert any object into any openings or try to fix the device as this may cause damage to the unit, malfunction, electrical shock, fire or personal injury. There are no user-serviceable parts inside the unit.
- Dispose of all products in accordance with local and national regulations and codes.
- Do not use if you are sensitive to light. If you are unsure whether you have increased sensitivity to light, we recommend you perform the Light Sensitivity Test instructed in this manual.
- Do not use if you are taking drugs known to make you sensitive to light. You may be sensitive to light and not know it. There are photosensitizing agents in drugs, foods, soaps, cosmetics, and dyes.
- Drugs that are known to cause light sensitivity: Antibiotics; NSAIDS; Diures; Retinoids; Hypoglycemics; HMG-CoA reductase inhibitors (statins); Photodynamic Therapy drugs; Prophoto sensitizers; Neuroleptic drugs or Antifungals.
- Do not use over or near a recent steroid injection. Doing so may result in photosensitive reactions.

Safety Information

Warning

- Discontinue use if you encounter any discomfort or skin redness that lasts for more than two hours. You may have the sensitivity to light or another condition. Consult your health care provider.
- Never operate this device it has damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or submerged in water. Return the device to Project E Beauty for examination and repair.
- Do not set or keep the device where it can fall or be pulled into a tub or sink. Keep away from water (including bath and shower). If the device accidentally falls in water unplug immediately.
- Keep away from flammable aerosol products.
- Do not use other AC Power Supplies or Chargers with your hair removal device.
- Turn off immediately, if the device begins to overheat.
- You should avoid exposing the affected area to direct UV light during use.
- During use, it is strictly forbidden to block the air inlet and the heat dissipation hole to avoid affecting heat dissipation.
- If there is any stain on the filter lens of the light exit, it is forbidden to spray alcohol directly at the light exit.

Safety Information

- Continuous depilation of more than two times in the same area does not improve the course of treatment but increases the risk of adverse reactions to the skin.
- Below area is not allowed to use the device: Tanning skin; There are infections, eczema, burns, hair follicle inflammation, open lacerations, wounds or lesions, and hematoma in the hair removal area; Dark brown or black skin such as black cockroaches and freckles; Hemangioma; Cosmetic plastic surgery and the surgical site; Around the eyes and near the eyebrows; tattoo.

Contraindications

DO NOT USE if one of the following pertains to you.

Pregnant or Breast Feeding: Do not use if you are pregnant or breast feeding. The device has not been tested on pregnant women and the effects, if any, are unknown.

Unsuitable Skin or Hair Color: Do not use if your natural hair colour in the treatment areas is light blond, white, grey or red. Do not use if you are of skin type VI (you rarely or never get sunburnt and you tan very darkly; your skin colour is brownish black or darker). In this case you run a high risk of developing skin reactions, such as strong redness, burns or discolourations.

Light Sensitivity: Do not use if you are sensitive to light. You are sensitive to light if you experience any discomfort or skin redness that lasts more than two hours. Discontinue using the device and consult with your healthcare provider.

Implantable Device: Do not use this device if you have a cardiac pacemaker, implanted defibrillator, shunt or other implanted electronic device. This device contains magnets that may cause interference with these devices resulting in serious injury or death.

Contraindications

Migraine: Do not use if you have experienced a migraine in the last two years.

Seizures: Do not use if you have a history of light triggered seizures.

Lesions: Do not use this device over any suspicious or cancerous lesions. To do so may result in delaying the best medical treatment.

Failure to use and maintain the device in accordance with the instructions outlined in this manual will invalidate the warranty.

Contraindications

Never use the appliance if you suffer from any of the diseases listed below, i.e.:

- If you have a skin disease such as active skin cancer, you have a history of skin cancer or any other localised cancer in the areas to be treated.
- If you have pre-cancerous lesions or multiple atypical moles in the areas to be treated.
- If you have a history of collagen disorder, including a history of keloid scar formation or a history of poor wound healing.
- If you have a history of vascular disorder, such as the presence of varicose veins or vascular ectasia in the areas to be treated.
- If your skin is sensitive to light and easily develops a rash or an allergic reaction.
- If you have infections, eczema, burns, inflammation of hair follicles, open lacerations, abrasions, herpes simplex, wounds or lesions and hematomas in the areas to be treated.
- If you have had surgery in the areas to be treated.
- If you have epilepsy with flashlight sensitivity.
- If you have diabetes, lupus erythematosus, porphyria or congestive heart disease.
- If you have any bleeding disorder.
- If you have a history of immunosuppressive disease (including HIV infection or AIDS).

Contraindications

Never use the appliance if you take any of the medications listed below, i.e.:


- If your skin is currently being treated with or has recently been treated with Alpha-Hydroxy Acids (AHAs), Beta- Hydroxy Acids (BHAs), topical isotretinoin and azelaic acid.
- If you have taken any form of isotretinoin Accutane or Roaccutane in the last six months. This treatment thins the skin and can, therefore, make it more susceptible to tears, wounds, and irritations.
- If you are on painkillers, which reduce the skin's sensitivity to heat.
- If you are taking photosensitizing agents or medications, check the package insert of your medicine and never use the appliance if it is stated that it can cause photo-allergic reactions, photo-toxic reactions, or if you have to avoid the sun when taking this medicine.
- If you take anticoagulation medications, including heavy use of aspirin, in a manner which does not allow for a minimum 1-week washout period prior to each treatment.
- If you take immunosuppressive medications.

Contraindications


Never use the appliance on the following areas:

- On nipples, areolas, labia minora, vagina, anus and the inside of the nostrils and ears.
- Men must not use it on the scrotum and face.
- Over or near anything artificial like silicone implants, pacemakers, subcutaneous injection ports (insulin dispenser) or piercings.
- On moles, freckles, large veins, darker pigmented areas, scars, skin anomalies without consulting your doctor. This can result in a burn and a change in skin colour, which makes it potentially harder to identify skin-related diseases (e.g. skin cancer).
- On warts, tattoos or permanent make-up. This can result in a burn and a change in skin colour.
- On areas, where you use long-lasting deodorants. This can result in skin reactions.
- On sunburnt, recently tanned (last 48 hours) or fake-tanned skin. All types of tanning can influence the safety of the treatment.

What is Included

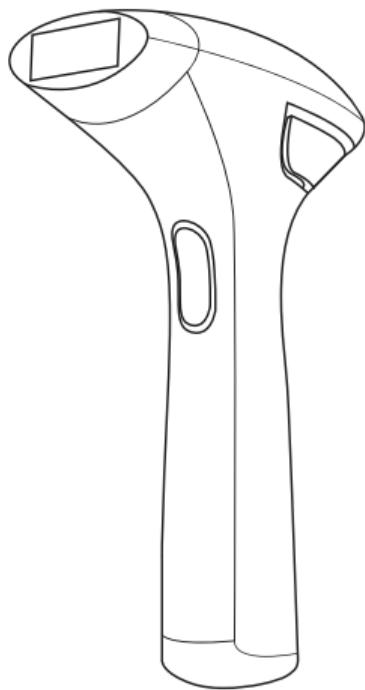SmoothPro+
IPL Hair Removal
Device

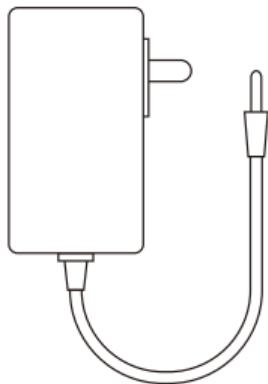
AC Adapter

Eye Protectors (Not Shown)

User Manual (Not Shown)

Location of Controls

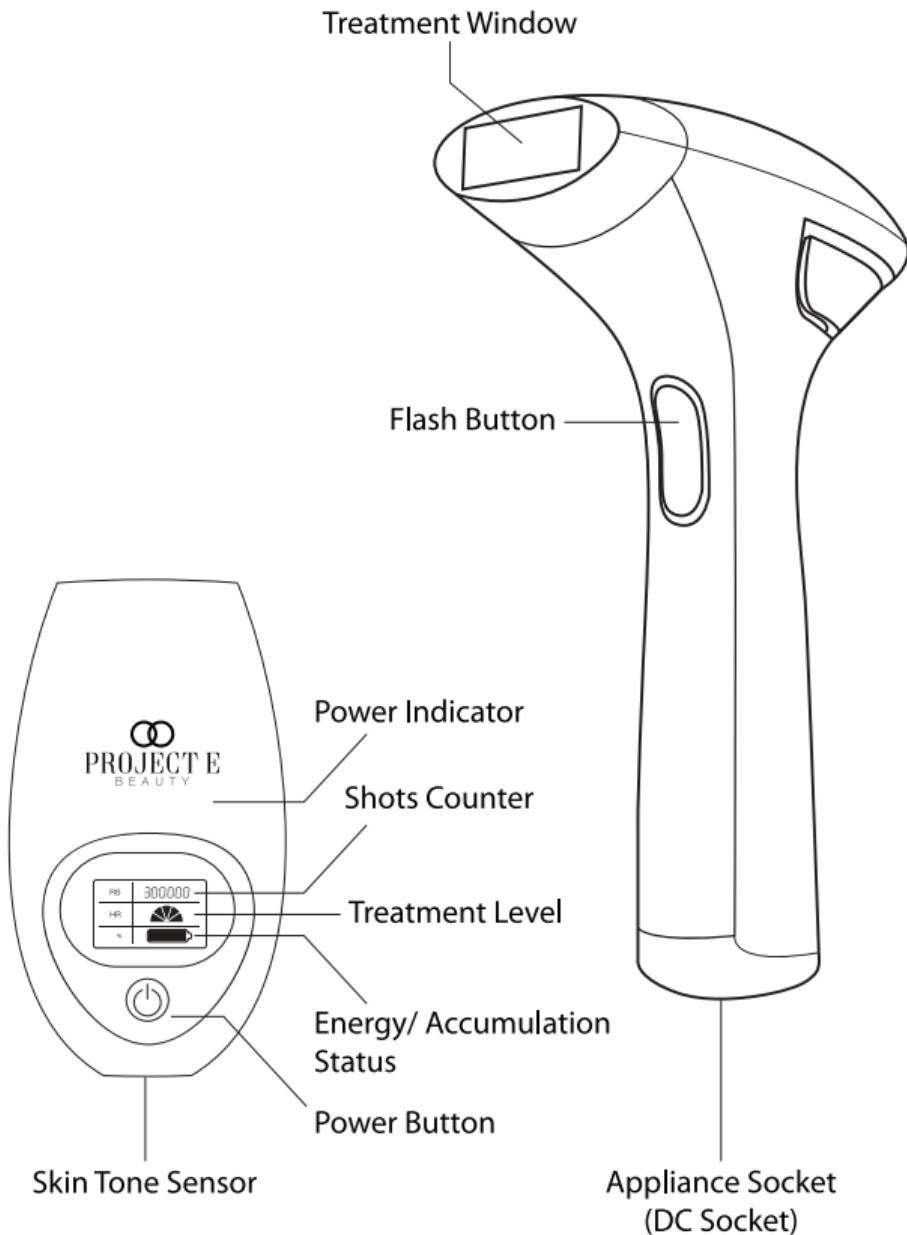Explanation of Controls

SmoothPro+ Activation Buttons and Status Indications:

1. Power button used for device activation, deactivation, and energy selection.
 - Pressing the button (after the device is connected to the power supply) will turn on the device at energy Level 1. Press again to change the energy, from the lower level to the higher, for a total of 5 different energy levels.
 - The device will switch to standby mode when left unattended or by pressing the Power button for 3 seconds.
2. Flash Button (flash pulse activation) marked on about 1/3 up the unit.
 - The flash button will activate the flash lamp if the treatment window is in full contact with the skin and the skin tone is within the allowed range. Pressing the button once will activate a single pulse, moving the device about half the distance of the window each time.

Explanation of Controls

3. The Status of Power Indicator

Indicator	Mode	Meaning
No Light	OFF	Device not working (Not connected to the power supply)
"Blinking" Red Light	Activation / Standby	Device is in standby mode (Device is activated after connected to the power supply)
"Blinking alternately" Red and Green	Flash Charging	Device is charging for the flash
"Solid" Red Light	Skin Testing	Skin test is completed, it means your skin tone is suitable for using the device
"Solid" Green Light	Hair Removal	Device can deliver a pulse (You are conducting 'Hair Removal' due to only the full skin contact and also allowed skin tone (refer to p.6-7))

Explanation of Controls

4. Other Components Explanation:

- **LCD screen:**
LCD display function, includes Pulse counter, Treatment level and Remaining electric quantity three parameters.
- **Pulse counter:**
To indicate the use of remaining number; and there are 300,000 times in total.
- **Treatment level:**
IPL's energy level selection; and it includes 1 to 5 energy levels.
- **Skin tone sensor:**
To detect your skin tone and identify whether is suitable for.
- **Treatment window:**
Place close to your skin and emit IPL, and you can start hair removal.
- **Appliance socket (DC socket):**
To connect to the AC adapter for power supply.

Light Sensitivity Test

Recommended Before Use

In order to determine if you have a Light Sensitivity Issue please follow the instructions below. It is recommended that you test a small, less sensitive area of skin (i.e., like your forearm).

1. Press the power button to turn on the device. (refer to p. 26-29)
2. Place the skin tone sensor window closely against the skin. The power light should change from green and red alternating lights into a solid red light, meaning that the skin type is appropriate for hair removal.
3. Place the device on an area of the forearm. Select the lowest level and deliver a pulse on the selected test skin area.
4. If you feel comfortable with the treatment level, increase it to Level 2 and perform an additional pulse on the next spot.
5. Only if you feel comfortable, you may continue to reach Level 3, Level 4, or the maximum treatment level of 5.
6. Wait 24 hours. You should NOT experience any adverse reaction during the next 24 hours. A bad reaction to light occurs if you have pain or discomfort in the treatment areas or any skin redness for more than two hours. It can be due to medication, supplements, or a possible skin issue.

How to Use

Preparations Before Use

1. Make sure the device is clean, dry and intact before the treatment. Check the appliance, especially the light exit window and attachment, for trapped hair, dust and fibers.

2. Shave off the surface hair on the treatment area. Clean your skin and make sure it is hair free, entirely dry and free from oily substances (e.g. present in cosmetic products like deodorants, lotions, perfumes, tanning creams and sunscreens).


How to Use

3. Plug the power adapter into power source and connect the hair removal device; the power indicator will be **red and flash**.

4. Put on the goggles over your eyes before the hair removal treatment. Never shine the light source directly into the eyes.

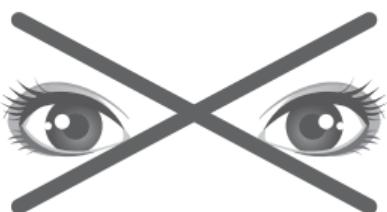How to Use

Please note that treatment on unshaved or improperly shaved areas has a number of undesired effects:

- The light exit window and attachment may become contaminated with stray hairs and dirt. A contaminated light exit window and attachment can reduce the life and efficacy of the appliance.
- Burns or black spots that appear on the light exit window and attachment due to contamination can cause the treatment to become painful or cause skin reactions like redness and discolouration.
- Hairs on your skin can burn, which can result in an unpleasant smell.

Note:

- If your last depilation method involved removing hairs by the roots (e.g epilating, waxing, etc.), wait until you notice significant hair regrowth before you use the device. Do not use depilatory creams instead of shaving.
- If shaving causes wounds or lesions, your skin may be more sensitive and the light treatment may sometimes feel uncomfortable.
- Once the hairs stop growing back, which usually happens after 4-5 treatments, you no longer have to shave before you use the appliance.

How to Use

Getting Started Steps

1. Press the power button to turn on the screen, it will show pulse counter, treatment level and remaining electric quantity. (As illustrated)

2. Place the skin tone sensor window closely against skin. The power light should change from flashing red light lights into a **solid red light**, meaning that the skin type is appropriate for hair removal. It takes up to 5 seconds to test your skin tone.

How to Use

- After testing skin tone, the device will be set to level 1. Short press power button to switch treatment level. Select the appropriate intensity level for your skin, staying within the recommended settings by using the skin chart provided. (see p.7)

- Place the device at a 90 ° angle on the skin so that the attachment is in contact with the skin.

How to Use

- For easier use, stretch the area of skin you want to treat. On areas with softer skin, you may need to press the head onto the skin with slight pressure.


- When the power light is **solid green** that means it is ready to flash.

How to Use


7. Press the flash button to release a flash. You will see a bright flash and hear a popping noise; then you may feel a mild stinging or heat sensation.
8. Let go of the flash button after the appliance has flashed.

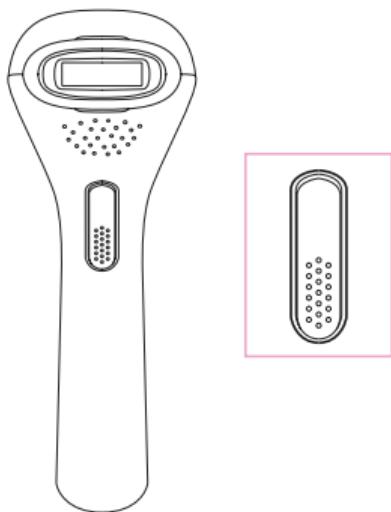9. For the next flash, repeat steps 6 to 9. After each flash, it takes up to 3 seconds until the appliance is ready to flash again.
10. Slowly move the area for hair removal, avoiding going over the same spot more than once. This does not improve the effectiveness of the treatment, but increases the risk of skin reactions.

How to Use

11. To avoid untreated areas, always make sure there is some overlap with the previously treated area when you place the appliance on the skin. The effective light only comes out of the light exit window. Make sure the flashes are produced close to each other.

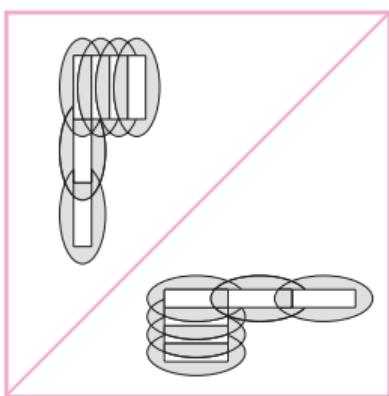12. After the treatment, press and hold the power button for 3 seconds to turn off the device; unplug the unit after each treatment session.

How to Use

13. After hair removal, it is recommended to increase moisturizing and avoid exposing the affected area to direct UV light.

14. Check the light exit window and the attachment for hairs and dirt and clean the appliance after use.

Operation Tips:

According to different treatment levels, there are some of our suggestions for operation:

- Level 1 and Level 2: It is recommended to be used on the small areas such as underarm and facial hair below the chin line.
- Level 3 to Level 5: It is recommended to be used on the large areas such as legs.
- In addition, please note the actual operation should be combined with the real situation of individuals, such as it is not recommended for dark skin.

Important NOTES:

- SmoothPro+ will not trigger a pulse if there is no full contact of the light window with the skin or your skin tone is too dark or tanned.
- Though the SmoothPro+ has been tested to international standards to ensure that its outputs are safe for your eye, it is recommended to wear goggles as added precaution. Never look directly at the light coming from the device, as it may be harmful to the eyes.
- If there is no any operation in 3 minutes, SmoothPro+ will automatically turn off.

How to Use

Skincare after treatment (Optional)


- After treatment, you can apply toner, moisturizer to the treated skin area to moisturize it. If treated area appears irritation, please use coolant or cold towel to cool down the skin.
- Please do not apply perfume, mask or makeup that contains exfoliator to the treated skin area.
- Within 24 hours after treatment, avoid UV light; during skin care, remember to apply sun cream.
- If skin irritation lasts for over 24 hours without getting better, please consult doctors.

Treatment Schedule

- Treatment Area: lips, armpit, shank, legs, arms, bikini line, etc
- Lamp life time: 300,000 flashes
- Treatment Cycle: 10 times treatment is a complete cycle
- Please see the hair reduce to decide the treatment times, it depends on individual's situation, normally it is suggested to operate again when the hair has grown 30%, it helps to achieve the best result.

First Stage (first 4 times):

- Once every two weeks

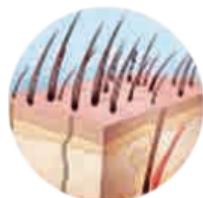Second Stage (the 5th and 7th time):

- Once a month

Third Stage (the 8th and 10th time):

- Once every 3 months

Treatment Schedule

Average Treatment Time

Area	Approx. treatment time
Underarms	1.5 minutes per underarm
Arms	3 minutes per arm
Bikini line	1.5 minutes per side
Entire bikini area	5 minutes
One lower leg	12 minutes
One full leg	20 minutes
Back	5 minutes
Lips	1 minute

Note: The average treatment times above have been observed during research and tests. Individual treatment times may vary from person to person.

Using on Different Areas

Using on the legs

- Check the recommended light intensities table on p.7 to see which settings are suitable for your legs.
- You can move the appliance up and down or across the leg, whichever is more convenient to you.
- Bony areas like shinbone and ankles are more sensitive than other areas of the body. When you press the attachment hard onto these areas, the skin can be more susceptible to skin reactions, such as skin irritation and redness.
- The appliance may leave imprints on the skin in the treated areas. These imprints may show you where you can make the treated areas overlap.

Using on Different Areas

Using in the bikini area

- Check the recommended light intensities table on p.7 to see which settings are suitable for your bikini area.
- In principle, SmoothPro+ is safe and gentle enough for treatment of the entire bikini area (e.g. for 'Brazilian' or 'Hollywood' style). Use a lower setting on darker pigmented, more sensitive skin in the bikini area.
- Do not use the appliance on inner labia, vagina and anus. Men who wish to treat their pubic area must not use the appliance on the scrotum.
- For hygienic reasons, always clean the attachment before you start treating another body area.

Using in the bikini area

- Check the recommended light intensities table on p.7 to see which settings are suitable for your underarms.
- Make sure you remove any deodorant, perfume, cream or other skin care product from your underarms before treatment.
- Raise your arm and place your hand in your neck to stretch the skin in the underarm area which makes the treatment easier to perform in this way.

Possible Skin Reactions

Common Skin Reactions

Slight Redness and Irritation: Your skin may show slight redness and/or may prickle, tingle or feel warm. This reaction is absolutely harmless and disappears quickly.

A Skin Reaction Similar to Sunburn: This may occur on the treated area for a short while. If this does not disappear within 3 days, we advise you to consult a doctor.

Dry Skin and Itching: This may occur because of the combination of shaving and light treatment. This is harmless and disappears within a few days. You can cool the area with an ice pack or a wet cloth. If dryness persists, you can apply a non-scented moisturizer on the treated area 24 hours after the treatment.

Rare side effects

Burns, Excessive Redness and Swelling: They are the result of using a light intensity that is too high for your skin colour. If these reactions do not disappear within 3 days, we advise you to consult a doctor. Wait with the next treatment until the skin has healed completely and make sure you use a lower light intensity.

Possible Skin Reactions

Skin Discolouration: Skin discolouration manifests itself as either a darker or lighter patch than the surrounding area. This is the result of using a light intensity that is too high for your skin colour. If the discolouration does not disappear within 2 weeks, we advise you to consult a doctor. Do not treat discoloured areas until the discolouration has disappeared and your skin has regained its normal skin colour.

Skin Infection and Inflammation: This occurs very rarely and may be caused by using the appliance on wounds or cuts that resulted from shaving, on already existing wounds or on ingrown hairs.

Excessive Pain: This can occur during or after treatment if you have used the appliance on unshaved skin, if you use the appliance at a light intensity that is too high for your skin colour, if you flash the same area more than once and if you use the appliance on open wounds, inflammations, infections, tattoos, burns, etc (see 'Contraindications' on p.12-16 for more information).

Cleaning

- Make sure to turn off and unplug unit before cleaning.
- It is recommended to clean the treatment window after every session by using a clean cloth with a few drops of water to gently wipe the device.
- If it is hard to clean the edges of the attachment with the cleaning cloth, use a cotton bud instead. Make sure that no fluff or fibres stay behind on the attachment or on the light exit window.
- Never use scouring pads, abrasive cleaning agents or aggressive liquids such as petrol or acetone to clean the appliance.
- Keep SmoothPro+ dry. Never submerge the device or any of its parts in water as may lead to an electrocution risk.

Storage

- To avoid breakage or missing device and accessories, please store the product in the box after each use.
- Clean the device before you store it.
- Keep it in a cool and dry place, away from sunlight. Do not store it in wet or humid environment, otherwise it may affect the product function.

FAQs on SmoothPro+

Q: When I put the device against my skin, nothing happens. Why?

A: 1) You haven't finished the skin tone testing step. Please place the skin tone sensor window against skin closely to finish the testing.
2) Treatment window is not pressed closely against skin.

Q: Why I feel heat and hurt during operation?

A: 1) The treatment level is too high for you. Please press the power button to choose a lower level.
2) There is hair left on skin. Please shave and clean the hair before use.

Q: Are there any long-term effects of SmoothPro+ usage? Is it safe for my skin?

A: Professional IPL treatments already exist for over 25 years and researches haven't shown appearance of skin cancer or fertility issues. No serious side effects or damage from long-term use have been reported. SmoothPro+ safety system ensure that light flashing is only possible if the window of the devices is in full contact with the skin. If your skin is too dark, the skin color sensor prevents the emission of flashes. Side effects and complication, while possible, are very unlikely as long as you use SmoothPro+ in accordance with this user manual.

FAQs on SmoothPro+

Q: Is SmoothPro+ safe to use on the face?

A: By checking your skin type, you can choose the specific setting that will ensure effective yet gentle treatment on your face. As the area around your eyes is very sensitive, we recommend to use SmoothPro+ only below the cheekbones in order to avoid the risk of eye damage. To achieve optimal results on the cheeks, upper lip and chin stand in front of the mirror to get a good view of the area to be treated.

Q: Is the flashing light on SmoothPro+ safe for my eyes?

A: SmoothPro+ does not hurt your eyes. The safety system prevents flashing when the device is not fully in contact with the skin. Be sure to make good skin contact to avoid unnecessary scattered light and never use SmoothPro+ in the area around your eyes or for treating your eyebrows. You may wear goggles during use. We recommend to use SmoothPro+ in a well-lit room, to reduce the perceived brightness of the flash.

Specifications

Item Name: SmoothPro+ IPL Hair Removal Device

Item ID: PE720

Type of protection: BF type

Power 36W (12V/3A)

Input voltage: 100-240V

Rated frequency range: 50/60Hz

Lamp life time: 300,000 flashes

Spot size (window size of IPL): 4.5cm²

Net weight: 386g

Package Dimension: 257.5 x 207 x 86.5 mm

Warranty

Activate Your 1 Year FREE Warranty Now.

The product come with a full one year warranty against manufacturer defects. It does not include damage caused by misuse, accident, abuse, neglect or other external causes. All we need is your name, email and order ID. That's all.

Scan the QR code or visit
www.projectebeauty.com/pages/warranty

Questions

Website: www.projectebeauty.com

Email: sales@projectebeauty.com

Tel: +1-626-759-4232

Our customer care expert team are here to answer any question you might want to ask.

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