VORA MICROCURRENT FACIAL ROLLER

PE723

INSTRUCTION MANUAL

Safety Instructions:

To ensure proper operation of this appliance, please read this instruction manual carefully before using the appliance and save it for further reference.

WARNING


- This device is not a toy, keep out of reach of children.
- Please verify the device before you start using it. Do not use if there is any damage.
- If used in a dry environment, handling this device may generate static electricity. This is a normal phenomenon. It does not affect the use of the device.

CAUTION

- The following persons should not use this device: pregnant women or while nursing or during menstruation; persons with allergies; persons with facial nerve disorders; persons who have fractures or implants with metal, silicone, or plastic; persons with other severe skin diseases.
- This device is not intended for use by persons (including children) with reduced physical, or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning the use of the device by a person responsible for their safety.
- Immediately stop using the device and consult a doctor if you feel any discomfort.
- Keep the device out of children.
- Regularly clean the device to extend its lifetime.
- Do not share the same device with others.
- Do not attempt to modify, disassemble, or repair the device by yourself to avoid accidents.
- Remove any jewelry/accessories and tie your hair when using the device, to prevent hair and jewelry from interfering with the device.
- Make sure the solar panel is exposed to the light source when using.

Device Description:

- ① Massage Rollers
- ② Solar Panel
- ③ Handle
- ④ Cleaning Cloth
- ⑤ Storage Bag

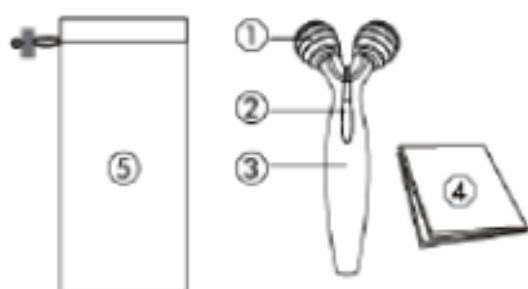Picture 1

Preparation Instructions:

- ※ This device uses light to generate microcurrent, and it does not require to be charged. The device features a solar panel on the handle, which uses light to generate microcurrent and output on the rollers.
- ※ Make sure you do not cover the solar panel while using the device.

How to Use:

- ※ Clean the device with a clean slightly wet cloth.
- ※ Use the device on the skin in different angles to increase the message intensity and efficiency.
(Picture 2)
 - Face: Gently move the rollers from the jawline to your cheekbone, or from the corner of the mouth to the ear; move the rollers in upward movements and from the center to the sides of the face. (Picture 3)
 - Neck: Gently move the rollers along the neck to the jawline. (Picture 4)
 - Arm: Gently move the rollers along the arm.
(Picture 5)
- ※ Do not use this device excessively. It is recommended to massage each area 15-20 times.

Picture 2

Picture 3

Picture 4

Picture 5

Cleaning & Maintenance:

- ※ Before and after each single use, clean the device promptly to ensure the best performance.
 - 1) Do not immerse the device in water. Clean the surface of the device with a clean slightly wet cloth.
 - 2) Never use steel wool, abrasive cleaners or corrosive liquids to clean the device, may result in its malfunction or discoloration.
 - 3) Never use lubricants, detergents, polishing agents or air fresheners on this device.
 - 4) After cleaning, please use a towel to wipe the device.

Storage:

Do not place or store the device for long period under a hot or humid environment.

Recycling:

When disposing of this device, do not put in the general household waste. Please dispose of this device in an official recycling point to contribute to environmental protection.

Device Specifications:

- Name: VORA Microcurrent Facial Roller
- Model No: PE723
- Solar Panel Rated Voltage: 200lux light $\geq 1.35V$
- Microcurrent: $\leq 10\mu A$
- Executive Standard: Q/XMMJ1002-2015

Warranty & Support:

Activate Your 1 Year FREE Warranty Now.

The product come with a full one year warranty against manufacturer defects. It does not include damage caused by misuse, accident, abuse, neglect or the external causes. All we need is your name, email and order ID. That's all.

Scan the QR code or visit
www.projectebeauty.com/pages/warranty

Customer Service

Website: www.projectebeauty.com

Email: sales@projectebeauty.com

Tel: +1-626-759-4232

Our customer care expert team are here to answer any question you might want to ask.

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