
Lumamask LED Light Therapy

(Model No. PE724)

User Manual

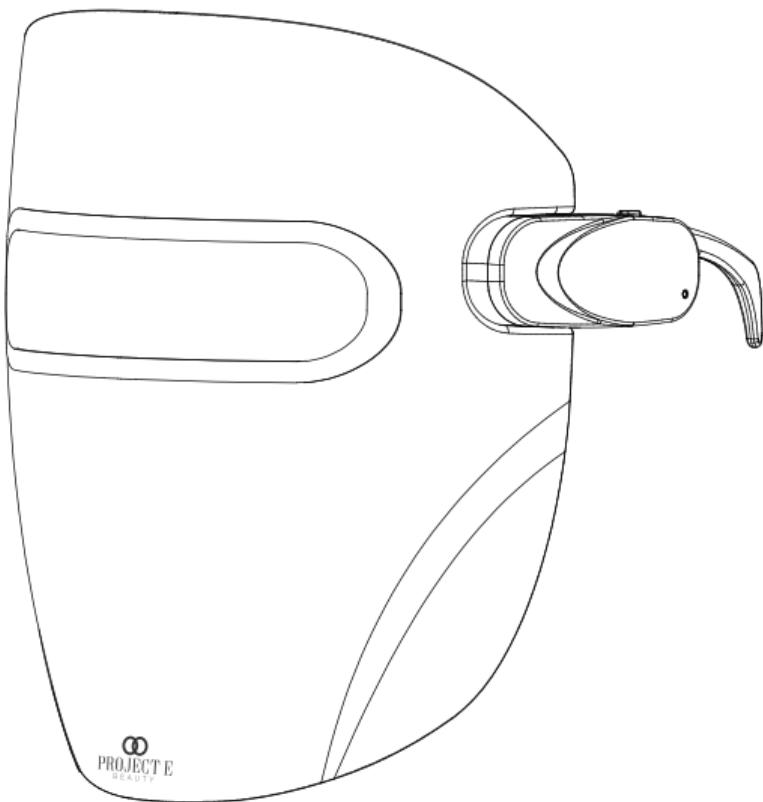www.projectebeauty.com

Introduction

Thank you for purchasing the Project E Beauty Lumamask (Model No. PE724). This device uses colored LED lights to provide pure light to treat a variety of skin conditions and rejuvenate skin, while reducing scar tissue. Unlike other beauty devices, this mask provides high-quality LED photon treatment that has been scientifically validated, so you can achieve efficient results faster than you ever imagined.

Photon light therapy works by using the natural healing properties present in differing wavelengths of light to treat various skin tissue issues. These state-of-the-art LEDs provide the purest light possible, especially compared to systems that use more traditional bulbs. Modern semiconductor technology has allowed dermatologists to perfect what wavelengths of light are generated, meaning more of the healing energy of the light reaches the skin tissue.

Exposing the skin's tissue surface to light activates natural photoreceptors and allows the skin cells to absorb the energy found in the light. This process accelerates biochemical reactions in the skin tissue, ultimately, leading to faster healing times, while increasing the production of collagen. This can reverse some of the more heinous effects of aging.

This technology is safe. Cosmetology experts agree these treatments have no noticeable side effects if they're administered correctly. Treatment methods provided by the LED photon mask are an important component of a complete skin care routine and should only be used for cosmetic purposes. This is not intended in any way to diagnose, treat, cure, or prevent any disease or diseases.

This 140g ultralight LED mask is designed to ensure a comfortable experience with a hollowed-out viewing angle. It's made of 36 LEDs, an optimal amount to maximize skin improvements, with 7 different light settings to provide treatments for various skin concerns. This mask is wireless and you can conveniently use it anywhere. Our light therapy face mask uses the latest in light irradiation technology and quality materials making it safe, comfortable, natural, and reusable.

Benefits of Photon LED Therapy

Recent cosmetology breakthroughs have shown that LED light therapy is effective for treating several aspects of skin ageing. Studies suggest this therapy can reduce the appearance of fine wrinkles, rough skin, shallowness of the flesh, and mottled hyperpigmentation. The LED light therapy increases both the proliferation and the initial production of collagen. A few studies have also suggested this therapy can cut down on the types of bacteria that cause spots and pimples. This might also help to eliminate breakouts and minimize the level of scar tissue.

Each of the colors of LED light have a different specific function. They are each selectable by way of the control box.

Red Light(~630nm): This setting allows the LED mask to penetrate the skin more deeply. Red light penetrates the skin approximately 8-10 mm, which stimulates the cells that lie at this layer. Many are fiber cells, and stimulating them can aid in the production of collagen. It can reduce wrinkles and make skin elastic and smooth. It's ideal for improving skin in the area around the eyes and forehead. It works well to cut down on fine lines. Those who need to repair damaged skin will opt for this setting.

Blue Light(~415nm): This setting allows the energy from the LED mask to penetrate the skin approximately 0.5mm. It's optimized for use with rashes and is suitable for skin that is sensitive or constantly oily. Blue light can help to reduce acne and even tighten skin that has loosened. It can quickly inhibit inflammation and possibly even destroy harmful bacteria.

Green Light(~525nm): This setting can help with mixed skin complexions, and will help to brighten dulled skin that has lost its luster. Studies suggest that green light can help to clear the face's natural lymphatic system. It might even aid in providing relief for mental stress, which is of special importance, as many people use the LED photon facial mask while they relax.

Purple Light(~390nm): This setting improves blood circulation with the function of promoting the activation of skin tissue. The purple light can strengthen the protein fiber tissue and shrink pores. It can also help to get rid of yellowish complexion and remove wrinkles and acne. At the same time, it accelerates blood circulation and lymph metabolism.

Benefits of Photon LED Therapy

Cyan Light(~490nm): This setting provides gentle acne treatment and mitigates skin inflammation. It can be used to adjust and improve the gland function of oily, black, and pink skin as well as acne. Meanwhile, the cyan light is able to promote the synthesis of protein and collagen that activate skin cells and provide the relaxation of skin contraction efficiently.

Yellow Light(~590nm): This setting revitalizes the lymphatic and nervous systems, strengthens the muscular and immunological systems, whilst simultaneously improving and balancing sensitive skin. Yellow light improves microcirculation, adjusts factors, and inhibits melanin production in order to fade dark spots, freckles, and blemishes. This method of skin rejuvenation is recommended as an effective treatment for skin problems caused by age and sun exposure.

White Light(~Mix): This setting synergistically helps activate your skin's natural renewal system to reverse the effects of aging at the same time restoring a younger and healthier skin.

Safety Instructions

Precautions:

- Read all instructions before using your Lumamask and make sure you understand all of the operating instructions.
- Always check the device before use to ensure it is undamaged. Use of a damaged device can be unsafe.
- Stop immediately, if you experience any adverse effect or have a concern regarding use.
- Do not use if you are pregnant.
- Do not use if you have experienced a migraine in the last two years.
- Do not use over any suspicious or cancerous lesion. Consult your doctor.
- Do not use the LED mask while using any product that contains AHA or glycolic acid. This can cause unwanted levels of stress in the skin.
- Do not allow children to use the device unsupervised.
- Remember that improper maintenance or negligence can cause damage and void the warranty. This also refers to improper storage.
- Dispose of all products in accordance with local and national regulations and codes.
- Do not use if you are sensitive to light. If you are unsure whether you have increased sensitivity to light, we recommend you perform the Light Sensitivity Test instructed in this manual.
- Do not use if you are taking drugs known to make you sensitive to light. You may be sensitive to light and not know it. There are photosensitizing agents in drugs, foods, soaps, cosmetics, and dyes.
- Drugs that are known to cause light sensitivity: Antibiotics; NSAIDS; Diures; Retinoids; Hypoglycemics, HMG-CoA reductase inhibitors (statins), Photodynamic Therapy drugs, Prophoto sensitizers, Neuroleptic drugs, or Antifungals.
- Do not use over or near a recent steroid injection. Doing so may result in photosensitive reactions.

Safety Information

Warning

- Discontinue use if you encounter any discomfort or skin redness that lasts for more than two hours. You may have sensitivity to light or another condition. Consult your health care provider.
- Never operate this device if it is not working properly, if it has been dropped or damaged, or submerged in water. Return the device to Project E Beauty for examination and repair.
- Do not set or keep the device where it can fall or be pulled into a tub or sink. Keep away from water (including bath and shower). If the device accidentally falls in water unplug immediately.
- Keep away from flammable aerosol products.
- Do not use other AC Power Supplies or Chargers to charge your LED mask.
- Turn off immediately, if the mask begins to overheat.

Contraindications

DO NOT USE if one of the following pertains to you.

Pregnant: Do not use if you are pregnant or think you might be pregnant. The device has not been tested on pregnant women and the effects, if any, are unknown.

Migraine: Do not use if you have experienced a migraine in the last two years.



Seizures: Do not use if you have a history of light triggered seizures.

Lesions: Do not use this device over any suspicious or cancerous lesions. To do so may result in delaying the best medical treatment.

Failure to use and maintain the LED mask in accordance with the instructions outlined in this manual will invalidate the warranty.

What is Included

- LED Mask Device (Model No. PE724)
- USB Cable
- Owner's Manual

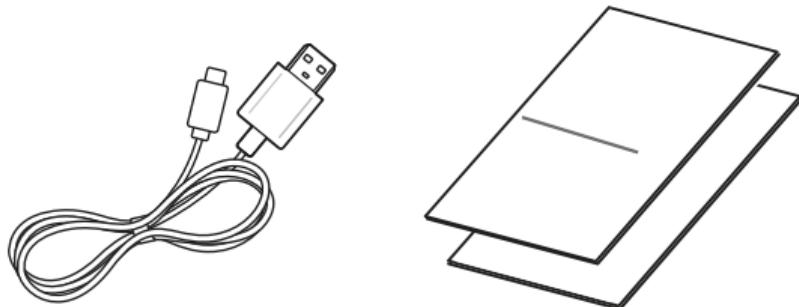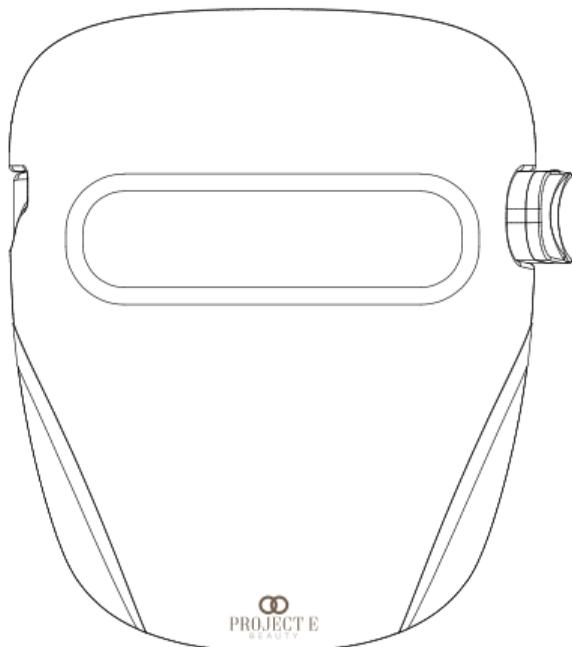Light Sensitivity Test

Recommended Before Use

In order to determine if you have a Light Sensitivity Issue please follow the instructions below. It is recommended that you test a small, less sensitive area of skin (i.e., like your forearm).

1. Press the “Power” button to turn on the mask.
2. Place the device on an area of the forearm. Leave in contact with your skin for a 3-minute treatment cycle.
3. Repeat by doing another 3-minute treatment on the same area of your arm.
4. Wait 24 hours. You should NOT experience any adverse reaction during the next 24 hours.

A **bad reaction** to light occurs if you have pain or discomfort in the treatment areas or any skin redness for more than two hours. It can be due to medication, supplements, or a possible skin issue.

Eye Safety

1. The design of this mask does not cover the eye area. Therefore, it is not needed to put on extra eye protection while performing treatments. Be sure you never look directly into the LEDs.
2. It is **NOT** recommended to shine the light source directly into eyes or the eyes of others.

How to Use

Follow these steps to ensure proper operation during a treatment:

- 1 Clean the surface with a soft wet cloth. Make sure your face is clean, dry, and wiped free of any oil or cosmetics before using the mask.

- 2 Use a serum, ion gel, essence or cream on the face and neck to achieve the best results.


You can purchase our Project E Beauty Serum at www.projectebeauty.com

- 3 Place the mask over the face and adjust the position for a comfortable fit.


- 4 Press and hold the Power button on the left side of the temples to turn on the mask.

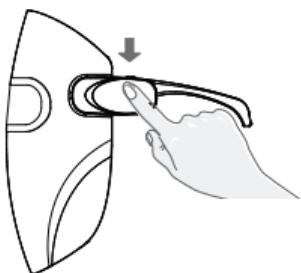
How to Use

- 5 Press again the Power button to select the desired color setting for the mask. Pressing once will select the red color, twice will select blue, three times will select green, four times will select purple, five times will select yellow, six times will select cyan, and seven times will select white.

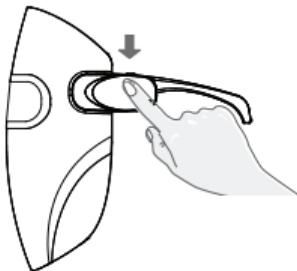- 6 The LED light will automatically be turned off after 10 minutes. Also, when the mask is out of battery, it will be turned off after blinking red light.

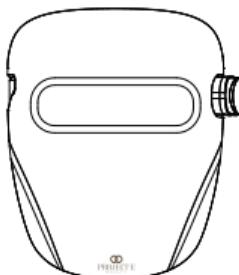- 7 Do not make direct eye contact with LED light while the device is in operation. Excessive exposure to LED light can cause eye damage.


Charging

Charging Your Device

To charge the mask please insert the USB cord into the jack on the mask and then plug it into a standard 5V USB socket. If it does not fit, turn the cord around and carefully try and insert it in the opposite direction. It takes around 2 hours to fully charge the LED mask.

| Charging Indicator Light | Status |
|--------------------------|--------------------------|
| Red | Battery is charging |
| Green | Battery is fully charged |

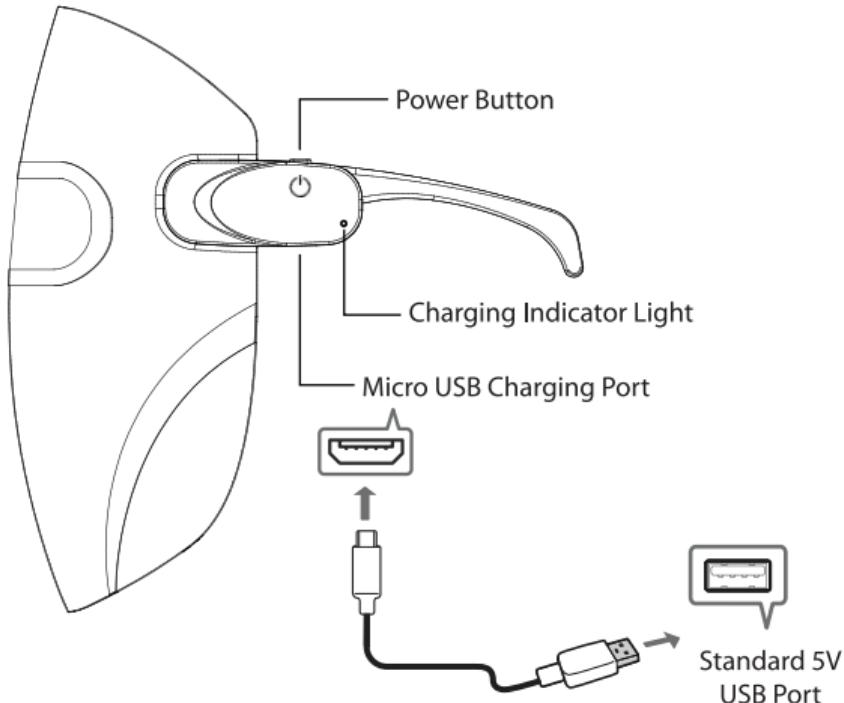Frequently Asked Questions

Q: What kind of cosmetics should I use?

A: After cleansing, please apply the foundation cosmetics.

Please pay attention to using cosmetics that contain the following ingredients as it may reduce the effect: denatured alcohol, sunscreen, retinol, oil.

Q: Can I open my eyes when I use this mask?

A: Although the eye shield is applied to prevent the LED light from being used, it is recommended that you close your eyes as much as possible because some LED light may be generated depending on the user.

Q: What should I do if I feel dry?

A: 1) Please reduce the frequency of use in a week.

2) It is recommended to use moisturizing cream or mask pack for moisturizing immediately after use.

3) If you are taking medication, consult your health care professional before using this medicine to induce skin dryness.

Q: How often should I use this mask?

A: It's suggested that you do not use the mask for more than 20 minutes each time. The suggested operating time is between 10 and 15 minutes. More than likely, you will want to start with one 10-minute treatment, 4 times a week to start, unless otherwise directed by a doctor.

Once you have achieved your initial goals, you will, more than likely, want to drop down to using the device only twice a week. As the LEDs provide a pure and powerful output, these few treatments will be sufficient.

Frequently Asked Questions

Q: How does light therapy actually work?

A: Light therapy gently penetrates the skin to clear acne. The blue light in the mask targets acne-causing bacteria, and the red light in the mask reduces acne inflammation. This is a clinically proven acne technology that is used in dermatologists' offices.

Q: Is light therapy safe? Is the light UV free? Can I get burned by it? Will the light damage my eyes?

A: Light therapy is proven to be safe on the skin and is UV free.

This product does not emit harmful, hot UV wavelengths of light. The 7 color lights that make up the core of the mask's technology are considered safe and the Light Therapy Mask delivers controlled doses. In fact, light therapy has been used by dermatologists for over a decade and is a clinically proven technology.

Q: Are there side effects?

A: Light therapy technology is safe and clinically proven to be gentle to skin.

However, discontinue use if you experience any discomfort or skin reddening or discoloration lasting for more than 24 hours.

Frequently Asked Questions

Q: Does the mask come in different sizes?

A: No. The mask is designed to fit most teens and adults.

Q: Can I use if I am pregnant?

A: Do not use the Light Therapy LED Mask if you are pregnant, may be pregnant or are nursing as the risks are unknown and have not been established.

Q: What happens if I need to end my 10-minute treatment early before the mask automatically turns off?

A: You can turn the mask off at any time by pressing the power button for a FULL second.

Q: How do I know when my treatment session is over?

A: The product will automatically turn off after 10 minutes treatment.

Q: How do I clean the device?

A: Clean the outside of the device whenever it appears dirty. Moisten a soft cloth and clean the surface of the device after each use. Do not get water or liquid inside any part of the device.

Maintenance

- Be sure to unplug the device properly by holding the plug and removing it from the USB outlet.
- Clean the mask surface with a soft wet cloth after each use.
- Immediately dry the surface any time it comes into contact with moisture, chemicals, or water.
- Always keep the device clean and sterilized. Store the device in a clean and dry place.

Specifications

| | |
|-------------------------------------|------------|
| Red Light Photon Wavelength | 630 nm ±5% |
| Blue Light Photon Wavelength..... | 415 nm ±5% |
| Green Light Photon Wavelength..... | 525 nm ±5% |
| Purple Light Photon Wavelength..... | 390 nm ±5% |
| Cyan Light Photon Wavelength..... | 490 nm ±5% |
| Yellow Light Photon Wavelength..... | 590 nm ±5% |
| Rated Power | 2 W |
| Working Voltage..... | 3.7 V |
| Charging Input Voltage..... | 5V, DC 1 A |

Warranty

Activate Your 1 Year FREE Warranty Now.

This product comes with a full one year warranty against manufacturer defects. It does not include damage caused by misuse, accident, abuse, neglect or other external causes.

All we need is your name, email and order ID. That's it.

Scan the QR code or visit
www.projectebeauty.com/pages/warranty

Questions

Website: www.projectebeauty.com

Email: sales@projectebeauty.com

Tel: +1-626-759-4232

Our customer care expert team are here to answer any question you might want to ask.

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