



REINVO

FACIAL MICRODERMABRASION
WAND

USER GUIDE

Model: PE741



MULTILINGUAL INSTRUCTIONS

MULTIPLE LANGUAGES USER MANUAL

For other languages (French, Spanish, Italian, German),
scan the following QR code



SCAN QR
CODE

TO GET USER
MANUAL

NUMÉRISER LE
CODE QR

POUR OBTENIR LE
MANUEL DE L'UTILISA-
TEUR EN FRANÇAIS

SCANSIONA
IL CODICE QR

PER OTTENERE IL
MANUALE D'USO IN
ITALIANO

ESCANEAR
CÓDIGO QR

PARA OBTENER EL
MANUAL DE USUARIO
EN ESPAÑOL

QR-CODE
SCANNEN

DEUTSCHES BENU-
TZERHANDBUCH
ERHALTEN

FUNCTION & FEATURES

FUNCTIONS

Introducing Reinvo - the ultimate at-home microdermabrasion wand. It utilizes medical-grade, diamond-encrusted tips to remove dead skin cells and impurities that clog pores, leading to blackheads, spots, and acne breakouts. Reinvo's advanced suction technology promotes circulation, skin tissue development, and collagen production, resulting in smoother and more youthful-looking skin.

Experience professional-grade exfoliation in the comfort of your own home with Reinvo - the perfect addition to your skincare routine. For safety, read the instructions carefully and consult a physician if you have sensitive skin, take medication, are undergoing treatment, pregnant, or have pre-existing medical conditions.


FEATURES

- Uses medical grade diamond encrusted tips
- Customizable suction strength
- Smart screen shows auto mode time, suction strength, and battery life
- Suitable for sensitive skin types

WHAT'S IN THE BOX?


- Reinvio Facial Microdermabrasion Wand x 1
- Charging base x 1
- USB cord x 1
- Massage head x 1
- Regular head x 1
- Fine head x 1
- Body head x 1
- Replacement filters x 30
- O-shaped rubber gaskets x 5
- Cleaning brush x 1

REINVO FACIAL MICRODERMABRASION WAND

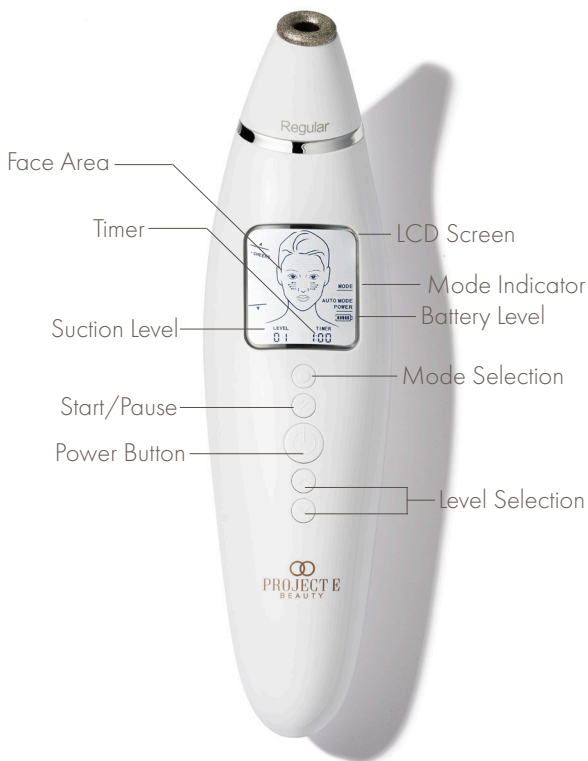

HOW TO USE REINVO?

TO CHARGE


 **Note: The device will NOT light up when charging.**

There are two ways to charge the Microdermabrasion Wand.

OPTION 1: Place the wand in the groove of the charging base, fitting the device onto the pin. Then, insert the USB cord.

OPTION 2: Insert the USB cord directly into the bottom of the wand.

It takes approximately three hours to fully charge the wand. At this point, the wand will automatically stop charging and the smart screen will show the fully charged icon.

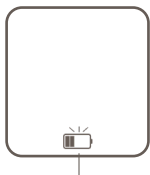

Charging Indicator




Fully Charged

When the battery status is low, the charging indicator on the screen will flash. At the same time, the device will alert with five beeps every 30 seconds until it goes into standby mode and automatically turns off after 10 seconds.

PREPARING FOR A TREATMENT

1. For best results, we suggest introducing the Micro-dermabrasion Wand into your skincare routine at night. This is because the body undergoes the natural repair process during the hours of sleep.
2. Before the treatment, remove all the makeup and wash your face. Dry thoroughly. A clean face will further enhance the exfoliation process.

HOW TO USE

 **To check your skin's sensitivity, be sure to start with the lowest suction level (on your arm).**

1. Press the power button to turn the wand on. The smart screen will light up and the device will enter standby mode.
2. Press the Mode Selection button to pick the mode you would like to use: AUTO MODE or MANUAL MODE.
3. Press the START/PAUSE button. Allow the head to gently suction your skin, making sure that it's perpendicular to the surface while moving the device as directed by the smart screen. Do not hover over a single spot for too long, and if possible, use your other hand to hold your skin taut and firm.
4. Adjust the suction level by pressing the up and down arrow while the wand is on. Level 1 is the weakest while level 3 provides the strongest suctioning power.
5. You may pause the wand any time during use via the START/PAUSE button.

AUTO MODE

This is a mode guides you through a quick four-minute session to help you achieve a cleaner, clearer, and healthier complexion. The wand will beep when it's time to move on to the next treatment area. Simply use the up and down arrow to adjust suction strength.



CHEEKS
Gently move from the
nose outward,
1 minute



NOSE
Gently move up and
down from the sides to
the middle of nose, 30
seconds



CHIN
Gently move up and
down the chin, 30
seconds



FOREHEAD
Gently move from the
middle of forehead
outwards to the
temples,
30 seconds



EYE AREA
Gently move from
side to side,
30 seconds



NECK
Gently move up and
down the neck,
1 minute

Want to skip a specific area?

- Press the START/PAUSE button.
- Press the up and down arrow to pick your preferred area.
- Press the START/PAUSE button again.

MANUAL MODE

This mode can work for a maximum of 9 minutes and is suitable for areas such as the neck, collarbone, arms, stomach, hands, legs, and feet. We do not recommend using it elsewhere.

We also do not recommend suctioning on one area for more than 90 seconds.

ENHANCE YOUR ROUTINE



REGULAR HEAD

The standard size tip that's great for use on cheeks, forehead, jaw, nose, neck; use once a week



MASSAGE HEAD

This tip enhances blood circulation and stimulates collagen production; use anytime



FINE HEAD

For hard to reach areas such as around the nose and eye areas; safe for sensitive skin, use once a week



BODY HEAD

Great for other parts of the body such as arms, legs, stomach, and back; use once a week

HOW TO CLEAN YOUR MICRODERMABRASION WAND

1. Pull the head off the device.
2. Use the cleaning brush and some rubbing alcohol or warm water to clean the surface and surrounding area. Make sure to rinse the head with warm water if you cleaned it with alcohol.
3. When the head has dried, you can now replace the used filter with a new one. Reattach the head to the wand.

FREQUENTLY ASKED QUESTIONS

Q: How does the Reinfo facial microdermabrasion wand work?

The Reinfo facial microdermabrasion wand uses a combination of diamond-tipped exfoliation and suction to remove dead skin cells and impurities from the skin's surface. This process stimulates circulation and promotes cell renewal, resulting in smoother, brighter, and more youthful-looking skin.

Q: Can I use the Microdermabrasion Wand with other skincare products in my routine?

Get the most out of your skincare routine - include the use of this wand before applying a hydrating serum or your preferred moisturizer. Avoid using this with medications and abrasive products.

Q: When is it not recommended to use the Microderm-

abrasion Wand?

If you have excessive acne, overly sensitive skin, cold sores, blisters, open wounds, fungal infections, or skin conditions such as rosacea, seborrheic dermatitis, eczema, psoriasis, convex moles, herpes, inflammation, or redness, please do not use Reinvu. If you are uncertain whether Reinvu is suitable for you, please consult your physician before use.

Q: How often should I use Reinvu ?

We suggest you use it once a week. However, if you're using the Massage Head, you can use the wand as often as you like.

Q: How long should each treatment session last?

We suggest up to four minutes or less, depending on your comfort level.

Q: When should I replace the heads, gaskets, and filters?

Replace the heads every six months. Loose or damaged gaskets should be replaced right away. The filters should be replaced after every treatment.

Q: What should be done when the suctioning power weakens?

Check the diamond head: It may become clogged with dead skin cells and impurities, reducing the suction power of the wand.

Check the filter: The filter in the wand's suction chamber may become clogged with debris, reducing the suction power of the wand.

TECHNICAL SPECIFICATIONS

Model Number	PE741
Rated Voltage	5V
Battery Capacity	850mAH
Weight	134g
Rated Power	2.5W
Dimensions	41 x 48 x 185 mm

1 YEAR WARRANTY

ACTIVATE YOUR 1 YEAR FREE WARRANTY NOW.

Reinvo Facial Microdermabrasion Wand comes with a full one year warranty against manufacturer defects. It does not include damage caused by misuse, accident, abuse, neglect or other external causes.

All we need is your name, email and order ID. That's it.



Scan the QR code or visit
www.projectebeauty.com/pages/warranty

QUESTIONS?

Website: www.projectebeauty.com

Email: sales@projectebeauty.com

Tel: +1-626-759-4232

Our customer care expert team is available to answer any questions you may have.

