



# BodySculpt

LED RF Beauty System

USER GUIDE

Model Number: PE142





# MULTILINGUAL INSTRUCTIONS

## MULTIPLE LANGUAGES USER MANUAL

For other languages (French, Spanish, Italian, German)  
scan the QR code



SCAN QR  
CODE

TO GET USER  
MANUAL

NUMÉRISER LE  
CODE QR

POUR OBTENIR LE  
MANUEL DE L'UTILISA-  
TEUR EN FRANÇAIS

SCANSIONA  
IL CODICE QR

PER OTTENERE IL  
MANUALE D'USO IN  
ITALIANO

ESCANEAR  
CODIGO QR

PARA OBTENER EL  
MANUAL DE USUARIO  
EN ESPAÑOL

QR-CODE  
SCANNEN

DEUTSCHES BENU-  
TZERHANDBUCH  
ERHALTEN

Or visit [projectebeauty.com/pages/user-manuals](https://projectebeauty.com/pages/user-manuals)



- Absolute contraindications include
  - Acute infectious disease or fever
  - Nausea
  - Dysfunction or disorders of the nervous system
  - Epilepsy
  - Pregnancy
  - Phlebitis or thrombosis
  - Heart disease
- Relative contraindications include
  - Cuts or abrasions
  - Vascular conditions
  - Skin sensitivity
  - Excessive body hair
  - Sinus blockages
  - Skin disorders
  - Infections
  - Sunburn
  - Headaches or migraines
- Never use or store the device around water
- Ensure hands and the treatment area are dry when performing a treatment
- Always remove all jewellery and metal to mitigate the risk of electric shock
- Device Handling and Care
  - Disinfection: Clean the device after every use with a soft cloth, ensuring the electrodes and LED cover are properly maintained.
  - Storage Tips: Store the device in a cool, dry place when not in use, keeping it away from moisture

# GETTING STARTED

- What's in the Box?
  - Project E Beauty BodySculpt | LED RF Beauty System x 1
  - USB Cable x 1
  - User Manual x 1


# BodySculpt

LED RF Beauty System



# CHARGING INSTRUCTIONS - HOW TO CHARGE BODY SCULPT

- Connect the USB cable to a common 5V USB power source
- Connect the USB cable to the power jack on the bottom handle of the device (see image)
- Charge for a minimum of two hours before use
- Once the device starts to charge, the power button will flash red
- Once the device is fully charged, the power button will display a steady red light
- The device cannot be used while charging. Wait for a steady light to resume before use

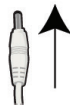



# FUNCTIONS

- Benefits
  - **Red LED:** Stimulates the production of new collagen to prevent skin sagging and enhance overall skin elasticity and firmness.
  - **Radio Frequency:** Gentle heat energy boosts the production of new collagen fibres to tighten the skin and disperse stubborn cellulite
  - **Ultrasonic Cavitation:** Ultrasonic waves directly target the fat cells to dissolve unwanted adipose tissue from areas such as the hips, torso, legs, and arms and give the appearance of a slimmer, smoother silhouette.
- How the technology works
  - **Red LED:** specific red light wavelengths penetrate the dermal skin layer to trigger the production of new collagen, elastin, and hyaluronic acid
  - **Radio Frequency:** Electrical currents create heat in the skin tissue and trigger the collagen production process to tighten loose skin and enhance skin firmness
  - **Ultrasonic Cavitation:** Through a process known as ultrasonic lipolysis, ultrasonic vibrations break down stubborn fat and cellulite by creating pressure that releases free radicals to reduce localized fat deposits

- Features
  - 10-minute auto-shut off timer
  - 3-in-1 body sculpting and firming device
  - Premium, high-performance conductor probes for advanced conductivity
  - 3 adjustable intensity levels

## HOW TO USE YOUR BODYSCULPT | LED RF BEAUTY SYSTEM

- Step-by-Step Usage Instructions
  - Preparing Your Skin
  - Clean the treatment area using a gentle, nourishing body cleanser
  - Pat dry
  - Ensure the area is free of any scented lotions
  - Apply a generous layer of Project E Beauty HydraFirm Body Gel
- Selecting Settings (Intensity, Mode, etc.) + Turning device on and off
  - Disinfect the device using a soft cloth or wipe



- **Press the Power Button:** Hold for 2 seconds to power on and off the device



- **Press the LED Button:** To turn on the red light therapy function. Press again to switch between continuous and pulsed red light.



- 'Pulsed mode' emits red light wavelengths in short bursts to enhance stimulation, allowing for deeper light penetration

- **Press the CV Button:** To activate the ultrasonic cavitation function



- **Press the RF Button:** To activate the radio frequency skin tightening function




- To navigate between radio frequency intensity levels, **press the RF button** and observe which level lights up. It is recommended to start by using the lower intensity mode and building up as you adjust




- All three modes can be used simultaneously
- Recommended Duration of Use
  - Use the device for approximately 5-10 minutes per treatment area. 5 minutes will suffice for smaller areas such as the arms, while 10 minutes is recommended for larger areas such as the torso or upper thigh.
- Recommended Frequency of Use
  - Treat each area 3-4 times per week
- Autoshut off
  - The device has an integrated shut-off timer of 10 minutes. Simply power on the device again to continue treating additional areas
- Position of the device
  - Ensure the device is positioned against the skin to allow the energy to penetrate effectively

- How to use on each body part


- **Buttocks:** Use for 5-7 minutes, slowly move in circular motions, starting at the bottom of the glute and moving upwards

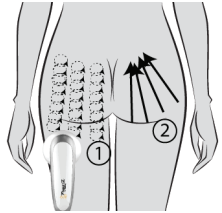


- **Legs:** Use for 10 minutes per leg, starting from the knee and working upward towards the upper thigh and buttocks. Treat the inner, outer, lower, and upper thigh

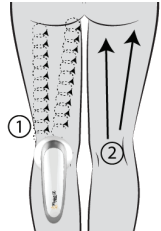

- **Arms:** Use for 5 minutes per arm, move the device in gentle circular motions starting at the elbow and moving up the back of the arm. Treat the front and back of the arm and inner and outer shoulder for a holistic treatment

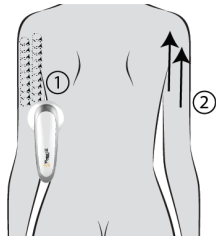

- **Stomach:** Treat for 10 minutes, starting at the lower abdomen and working in slow, gentle, circular movements up the torso until you reach the chest area

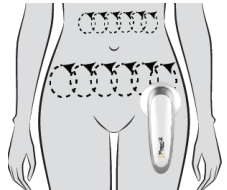

# MAINTENANCE & CARE

- Cleaning Instructions
  - Before and after each use, clean the device with a soft cloth to remove any residue from the electrodes.
- Storage Tips
  - Keep the device in a cool, dry place and avoid exposure to moisture.



# 1 YEAR WARRANTY

ACTIVATE YOUR 1 YEAR FREE WARRANTY NOW.

The BodySculpt | LED RF Beauty System comes with a full one year warranty against manufacturer defects. It does not include damage caused by misuse, accident, abuse, neglect or other external causes.

Register your warranty, all we need is your name, email and order ID. That's it.



Scan the QR code or visit  
[www.projectebeauty.com/pages/warranty](http://www.projectebeauty.com/pages/warranty)

## QUESTIONS?

Website: [www.projectebeauty.com](http://www.projectebeauty.com)

Email: [hello@projectebeauty.com](mailto:hello@projectebeauty.com)

Tel: +1-626-759-4232

Our customer care expert team is available to answer any questions you may have.

# TECHNICAL SPECIFICATIONS

Input:	DC 5V 1A
Model No/SKU:	PE142
Power Max:	11W
Ultrasound Frequency:	330KHz
Radio Frequency:	500KHz
Battery Capacity:	1400mAh
Auto Shut-Off Timer:	10min

## CERTIFICATIONS AND COMPLIANCE

CE Safety Certification, deemed safe, healthy and environmentally sound for sale in EEA.

### Legal Disclaimer

The BodySculpt | LED RF Beauty System is designed for personal cosmetic use only and is not intended to diagnose, treat, cure, or prevent any medical condition. Always follow the instructions for use and consult a healthcare professional before using this device, especially if you have a pre-existing skin condition, are pregnant, or are taking medication that increases light sensitivity.

The manufacturer and distributors are not liable for any injury, damage, or adverse effects resulting from misuse, failure to follow instructions, or unauthorised modifications of this product. Discontinue use immediately if you experience discomfort, irritation, or unusual skin reactions.

This device is not waterproof—avoid exposure to water or excessive moisture. Keep out of reach of children and store in a cool, dry place when not in use

