
BODY FAT SCALES

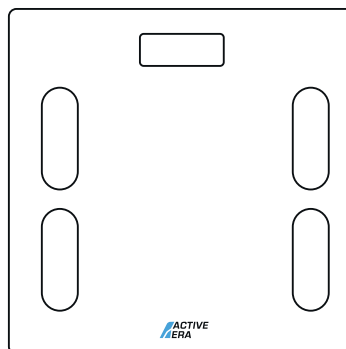MODEL - BS-02

INSTRUCTION MANUAL

BEDIENUNGSANLEITUNG
MODE D'EMPLOI
MANUALE DI ISTRUZIONI
MANUAL DE INSTRUCCIONES
HANDLEIDING
BRUKSANVISNING

IMPORTANT INSTRUCTIONS - RETAIN FOR FUTURE USE.

Body Fat Scales - BS-02

Thank you for choosing to purchase a product from ActiveEra. Please read the entire manual carefully prior to first use and keep in a safe place for future reference.

SAFETY INSTRUCTIONS

- Use this device only as described in this manual. Any other use not recommended by the manufacturer may cause fire, electric shock, or injury to persons.
- The appliance is not to be used by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction.
- Children must be supervised not to play with the appliance.
- Cleaning and user maintenance shall not be made by children unless continuously supervised.
- Children of less than 3 years should be kept away unless continuously supervised.
- Place the appliance on a flat, stable surface. Operating the product in any other position could cause a hazard. Avoid carpets and soft surfaces.
- Do not jump or stand on the outer edges of the scale.
- This device should not be used by individuals who are diagnosed with a condition or taking medication that can affect water levels in the body. The accuracy of the readings for these individuals has not been verified. Seek advice from a medical expert.
- The data provided by this device should be used as a reference only and is not a substitute for medical advice from a licenced healthcare professional
- This device should not be used to diagnose or treat any medical condition. Always consult your licenced healthcare professional prior to making changes to your diet, exercise plans, or physical activities.
- This device is slippery when wet. Make sure that the device and your feet are dry prior to use.
- Dry off the surface of the scale before use.
- Avoid placing the device in areas where it can become a tripping hazard.
- This device should not be used when it is stacked on top of another device. Doing so may cause the device to operate improperly.
- Do not place any objects on top of the scale when it is not in use.

- This device is intended for indoor household use and similar applications. Do not use outdoors.
- Do not wash the device with water or allow water or other liquids to enter into the device or contaminate electrical parts.
- Do not wipe or wash the device with chemical solvents. Prior to cleaning ensure the device is off and that the batteries have been removed.
- Do not operate the device if there are any signs of damage to the device or accessories supplied.
- Do not expose the scale to extreme temperature fluctuations.
- Do not place the scale near high heat sources e.g. ovens or radiators.
- Do not store the device in a vertical position. Doing so may damage the high precision sensors. Always store the device in a horizontal position.
- Do not attempt to repair or adjust any electrical or mechanical elements in the device. Doing so may be unsafe and will void your warranty.

BATTERIES

This product is supplied with Button Cell - CR2032 Li-ion disposable batteries.

PLEASE READ AND SAVE THESE IMPORTANT SAFETY INSTRUCTIONS

BATTERY WARNING: KEEP OUT OF REACH OF CHILDREN

Choking hazard. Do not swallow. Harmful if swallowed. Swallowing can lead to chemical burns, perforation of soft tissue, and death. Severe burns can occur within 2 hours of ingestion. If you think batteries might have been swallowed or placed inside any part of the body, seek immediate medical attention.

- This product only uses a Button Cell - CR2032 Li-ion battery. It is not compatible with other battery types.
- Do not allow children to handle batteries.
- If the battery compartment does not close securely, stop using the product and keep it away from children.
- Ensure that the battery is installed with the correct polarity (+ and -).
- Remove exhausted battery promptly.
- Always inspect new batteries before use. Damaged or weak batteries may leak acid. Discontinue use if batteries show signs of leakage, damage, or weakness at any time.
- Ensure the battery compartment is correctly and securely closed after installing the batteries.
- Wash hands after handling batteries.
- Disposable batteries are not rechargeable and must be recycled or disposed of separately from household waste. Always dispose of batteries at an official collection point and according to your local environmental laws and guidelines.

PRODUCT SPECIFICATION


Product Dimensions: 300x300x15mm
11.8"x11.8"x0.6"

Power: CR2032 3V Battery
Display: LCD

Max Weight: 180kg/400lbs/28st
Minimum effective weight: 6kg/13lbs

PRODUCT FEATURES

1. LCD display
2. Tempered glass surface
3. Battery Cover
4. Unit conversion: KG/LB/ST
5. High precision sensors
6. Metal sensor plates

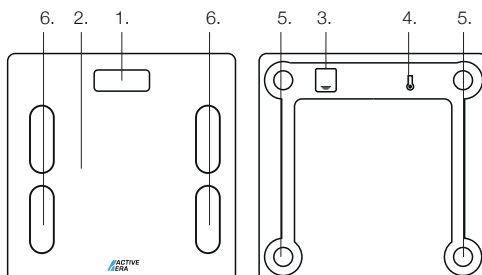SETTING UP THE PARAMETERS

The scales can store data on 10 separate users. To set up a new user or reset an existing user follow the steps below:

1. When the LCD display shows '0.0', press the 'SET' button for 3 seconds to enter the parameter input mode.
2. When the user code is flashing, press the up or down buttons to set the user number. Once selected press 'SET' to confirm your choice.
3. When the gender sign is flashing, press the up or down buttons to select the required gender. Once selected press 'SET' to confirm your choice.
4. Follow the same process for the height and age.

Remember your user code (1-10) to avoid setting the parameters again. If the scales are idle for more than 6 seconds during the parameter setting, the parameters already set will be saved and the LCD display will show "0.0" meaning that body fat measuring has now started.

WEIGHING

1. Place the scales on a solid and level surface
2. Stepping onto the scales will automatically turn them on. Once a steady measurement is taken, the final weight measurement will flash on the screen
3. Scales will automatically switch off when not in use

MEASURING BODY FAT

For an accurate result, please ensure that shoes and socks are removed before stepping onto the unit. Then align your feet so that they are both touching the metal sensor plates.

1. Press the 'SET' button to switch the unit on
2. Press the up or down button to select the user or reset the parameters
3. When the correct user number is flashing, step onto the scales (ensuring both feet are touching the sensor plates)
4. Wait until the weight reading is locked and flashing. Remain on the scales while it calculates the other measurements.
5. Once the results have been calculated they will be displayed and cycled through. The results will repeat three times and saved for next time.

Body Fat %

Visceral Fat

Body Water %

Recommended Calorie Intake

Muscle Mass

Body Mass Index

Bone Mass

To display a user's last test results, follow the steps below:

1. Press 'SET' button to turn the scales on
2. Select the correct user by pressing the up and down arrows
3. Wait for 8 seconds until the user stops flashing or press firmly on the scales platform
4. Press up and down arrows to cycle through last results

AUTO IDENTIFY USER FUNCTION

After the user has set their parameters and weighed themselves, the scales will remember the user's weight and automatically revert to their user profile, if their weight is within +/- 2kg / 4.4lbs.

BODY MEASUREMENTS EXPLAINED

These bathroom scales use bioelectrical impedance to measure body fat. A very small and safe electrical current is passed through the lower half of the body. The electrical current flows slower through bone and fat than it does through muscle and water. The scales can measure this difference in speed and estimate body fat using a multi-step, mathematical formula.

Please note, these scales have been calibrated for those individuals who have a median body type and therefore the measurements recorded by individuals with unique body types e.g. significant upper body muscle, athletes or bodybuilders may not be accurate.

Body Fat %: This measurement refers to the subcutaneous fat (outer fat) closest to the skin. Men and women have a different distribution of body fat so it is normal for body fat % to vary based on gender.

AGE	MEN	WOMEN
20-39	8%-20%	22%-33%
40-59	11%-22%	24%-34%
60-79	13%-24%	25%-36%

Body Water %: Indicates the percentage of water in the body. It is important to stay well hydrated for overall health and wellbeing.

MEN	WOMEN
50%-65%	45%-60%

Muscle Mass: Indicates the total weight of muscle in your body. The normal ranges for muscle mass are as follows:

AGE	MEN	WOMEN
20-39	34%-40%	25%-30%
40-59	33%-39%	24%-30%
60-79	33%-38%	24%-30%

Bone Mass: Indicates the total weight of bone in your body. The normal ranges for bone mass are as follows:

MEN	WOMEN
3%-5%	2.5%-4%

Visceral Fat: Indicates the amount of visceral fat in your body. Visceral fat is found in the abdomen and surrounding vital organs.

Category	Visceral Fat
Normal	1-9
High	≥10

Recommended Calorie Intake: This calculates how many calories you should consume per day based on your data and body composition.

Body Mass Index: Body Mass Index (BMI) is a number calculated from a person's weight and height. Refer to the following chart for classifications:

Category	BMI
Underweight	< 18.5
Healthy weight	18.5-24.9
Overweight	25-29.9
Obese	30-34.9
Severely Obese	35-39.9
Morbidly Obese	> 40

TO INSERT OR REPLACE BATTERY

- Carefully remove the battery cover on the underside of the scales
- Insert 1 x CR2032 3V battery into the compartment as per the directions indicated
- After inserting the batteries, replace the battery cover

UNIT CONVERSION

These scales measure weight in Kilograms (kg), Pounds (lb) and Stone (st). Press the button on the underside of the scales to switch between them.

ERROR CODES

- 'Lo' – Batteries are low and needs replacing
- 'Err' – Scales are overloaded and have exceeded maximum weight. Step off the scales to avoid damaging the sensors
- 'C' – A measuring error has occurred. Step off scales and back on again

RECYCLING & DISPOSAL

Waste electrical and electronic products (WEEE), batteries, accumulators, and packaging should not be disposed of with general household waste. This is to prevent possible harm to the environment and human well-being, and to help conserve our natural resources. Most electrical products with a plug, battery, or cable can be recycled. Contact your local council, household waste disposal service, or the store where you purchased the product for information regarding the collection schemes in your area.

Waste Electrical and Electronic Equipment (WEEE)

This symbol on the product and/or packaging indicates that Waste Electrical products (WEEE) under Directive 2012/19/EU which governs used electrical and electronic appliances, should not be disposed of with general household waste.

Waste Batteries and Rechargeable Batteries

This symbol on the batteries supplied with your product indicates that Waste Batteries under Directive 2006/66/EC which governs used batteries and accumulators should not be disposed of with general household waste. If the battery contains more than 0.004% lead, the battery symbol will be accompanied by the letters (PB).

The Mobius Loop

This symbol indicates that the product or packaging is capable of being recycled.

Take-Back Scheme

Waste electrical products (WEEE) can be sent to One Retail Group for recycling and disposal when you purchase a new product from One Retail Group, on a like-for-like basis. For more information on how to recycle this product please visit www.activeera.com

The end-user of electrical and electronic products (WEEE) is responsible for separating old batteries and lamp bulbs before they are handed over to a designated disposal and recycling service. Waste batteries can also be disposed of in collection bins, which can be found in most local supermarkets and household waste recycling centres.

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© Copyright 2022
EU Authorised Representative: Brandrep Limited