
---

**SMART SCALES**

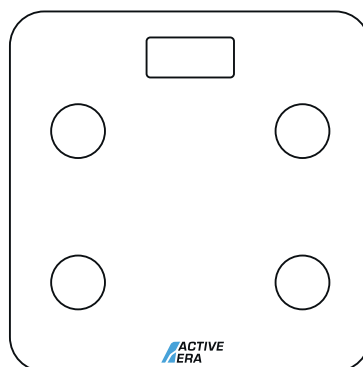**MODEL - BS-05**

---

**INSTRUCTION MANUAL**

BEDIENUNGSANLEITUNG  
MODE D'EMPLOI  
MANUALE DI ISTRUZIONI  
MANUAL DE INSTRUCCIONES  
HANDLEIDING  
BRUKSANVISNING

**IMPORTANT INSTRUCTIONS - RETAIN FOR FUTURE USE.**

# Smart Scales - BS-05

Thank you for choosing to purchase a product from Active Era. Please read the entire manual carefully prior to first use and keep in a safe place for future reference.

---

## IMPORTANT SAFETY INSTRUCTIONS

---

### **PLEASE READ AND SAVE THESE IMPORTANT SAFETY INSTRUCTIONS**

When using electrical appliances, especially when children are present, basic safety instructions should always be followed, including the following:

- Use this device only as described in this manual. Any other use not recommended by the manufacturer may cause fire, electric shock, or injury to persons.
- The appliance is not to be used by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction.
- Children must be supervised not to play with the appliance.
- Cleaning and user maintenance shall not be made by children unless continuously supervised.
- Children of less than 3 years should be kept away unless continuously supervised.
- Place the appliance on a flat, stable surface. Operating the product in any other position could cause a hazard. Avoid carpets and soft surfaces.
- This device should not be used if you have a pacemaker or any other internal medical device.
- This device should not be used if you are pregnant, have symptoms of edema, and/or are undergoing dialysis treatment.
- Do not jump or stand on the outer edges of the scale.
- This device should not be used by individuals who are diagnosed with a condition or taking medication that can affect water levels in the body. The accuracy of the readings for these individuals has not been verified. Seek advice from a medical expert.
- The data provided by this device should be used as a reference only and is not a substitute for medical advice from a licenced healthcare professional
- This device should not be used to diagnose or treat any medical condition. Always consult your licenced healthcare professional prior to making changes to your diet, exercise plans, or physical activities.
- This device is slippery when wet. Make sure that the device and your feet are dry prior to use.
- Dry off the surface of the scale before use.
- Avoid placing the device in areas where it can become a tripping hazard.
- This device should not be used when it is stacked on top of another device. Doing so may cause the device to operate improperly.

- Do not place any objects on top of the scale when it is not in use.
- This device is intended for indoor household use and similar applications. Do not use outdoors.
- Do not wash the device with water or allow water or other liquids to enter into the device or contaminate electrical parts.
- Do not wipe or wash the device with chemical solvents. Prior to cleaning ensure the device is off and that the batteries have been removed.
- Do not operate the device if there are any signs of damage to the device or accessories supplied.
- Do not expose the scale to extreme temperature fluctuations.
- Do not place the scale near high heat sources e.g. ovens or radiators.
- Do not store the device in a vertical position. Doing so may damage the high precision sensors. Always store the device in a horizontal position.
- Do not attempt to repair or adjust any electrical or mechanical elements in the device. Doing so may be unsafe and will void your warranty.

## **BATTERIES**

---

This product is supplied with AAA disposable batteries.

**PLEASE READ AND SAVE THESE IMPORTANT SAFETY INSTRUCTIONS**

### **BATTERY WARNING: KEEP OUT OF REACH OF CHILDREN**

Choking hazard. Do not swallow. Harmful if swallowed. Swallowing can lead to chemical burns, perforation of soft tissue, and death. Severe burns can occur within 2 hours of ingestion. If you think batteries might have been swallowed or placed inside any part of the body, seek immediate medical attention.

- This product only uses a AAA battery. It is not compatible with other battery types.
- Do not allow children to handle batteries.
- If the battery compartment does not close securely, stop using the product and keep it away from children.
- Ensure that the battery is installed with the correct polarity (+ and -).
- Remove exhausted battery promptly.
- Always inspect new batteries before use. Damaged or weak batteries may leak acid. Discontinue use if batteries show signs of leakage, damage, or weakness at any time.
- Ensure the battery compartment is correctly and securely closed after installing the batteries.

- Wash hands after handling batteries.
- Disposable batteries are not rechargeable and must be recycled or disposed of separately from household waste. Always dispose of batteries at an official collection point and according to your local environmental laws and guidelines.

## PRODUCT SPECIFICATION

---

Dimensions: 28cm x 28cm x 2.4cm

Net Weight: 1.53kg

Power: 3 x AAA batteries

Display: LCD

Max Weight Capacity: 180kg / 400lb / 28st


Minimum Effective Weight: 6kg / 13lb / 0.9st

Unit Increments: 0.1kg / 0.2lb

## PRODUCT FEATURES

---


1. LCD display
2. Tempered glass surface
3. Battery cover/compartament
4. Measuring unit (kg/lb/st)
5. High precision sensors
6. Anti-slip safety pads
7. Metal sensor plates
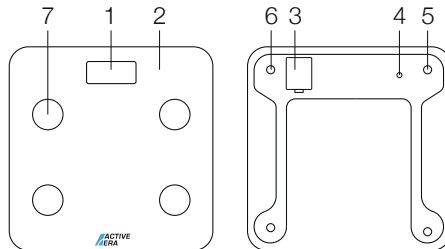
## BATTERIES

---

1. Turn the smart scales over and place on a dry and flat surface.
2. Pull the plastic tab from out of the battery compartment.
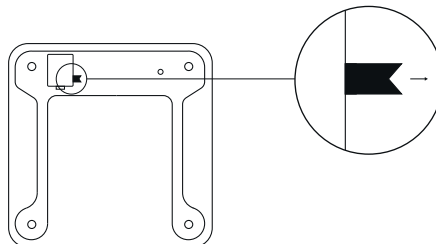
## REPLACING BATTERIES


---

1. Turn the smart scales over and place on a dry and flat surface.
2. Carefully remove the battery cover to reveal the battery compartment.
3. Remove the 3 old AAA batteries and insert 3 new AAA batteries. Ensure that the positive and negative ends of the batteries are facing the correct polarity directions as marked within the battery compartment.
4. Replace the battery cover.

## MEASURING UNITS

---

These smart scales can measure in kilograms (kg), pounds (lb), and stone (st). Press the UNIT button located on the bottom of the smart scales to cycle through and select the desired measuring unit.
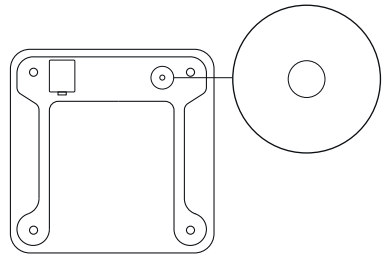
## ACTIVE ERA APP

---

### ACCOUNT CREATION

1. Open the Active Era app.
2. Press the Sign Up button.
3. Follow the on-screen instructions.

Note, entering the incorrect gender, age, or height during account creation may result in inaccurate measurements.

Note, to improve the accuracy of the measurements for individuals with unique body types (e.g. significant upper body muscle, athletes or bodybuilders) ensure that athlete mode is selected when creating an account.

Note, each main account can individually track the data of 24 users.

## PAIRING

Only pair the smart scales using the Active Era app. Do not pair the smart scales directly from your smartphone or tablet's Bluetooth menu.

To pair the smart scales with the Active Era app:

1. Enable Bluetooth on your smartphone or tablet.
2. Open the Active Era app and sign in.
3. Step onto the smart scales and stand still.
4. The smart scales will turn on and automatically pair.
5. You will know if the smart scales have successfully paired as your weight will show both on the LCD display of the smart scales and in the Active Era app.
6. Once paired the app and smart scales will automatically connect the next time the app is opened.

Note, if the smart scales do not successfully pair to the Active Era app please review the troubleshooting section of this instruction manual.

## MEASUREMENTS

---

Prior to taking measurements ensure that you have downloaded the Active Era app, created a user account, and paired the smart scales with your smartphone or tablet.

1. Place the smart scales on a dry, solid, and flat surface.
2. Open the Active Era app and sign in.
3. Remove shoes and socks and step barefoot onto the smart scales.
4. Stand still and ensure your feet are in full contact with the metal sensor plates.
5. The smart scales will turn on and an analysis of your body's composition will start automatically.
6. Once the analysis has completed a final weight measurement will flash on the LCD display of the smart scales in conjunction with the data being shown within the Active Era app.

Note, always remove shoes and socks prior to using the smart scales.

Note, to ensure accuracy of results, it is recommended to weigh yourself at the same time of day as measurements may vary after exercise, dieting, or dehydration.

## **BABY MODE**

---

Baby mode allows you to find out the weight of your baby by weighing yourself on the smart scales whilst holding your baby. To use baby mode:

1. Open the Active Era app and sign in.
2. Create a new user profile for your baby and ensure this new baby specific user profile is selected.
3. Place the smart scales on a dry, solid, and flat surface.
4. Press the menu button in the top right hand side of the screen and press Baby Mode in the drop down menu.
5. Step onto the smart scales by yourself.
6. The smart scales will turn on and an analysis of your body's weight will start automatically.
7. Once a final weight measurement has been obtained it will flash on the LCD display.
8. Immediately weigh yourself with your baby.
9. Once a final measurement has been obtained your weight plus your baby's weight will flash on the LCD display followed by your baby's weight being shown within the Active Era app.

Note, if using baby mode when your own user profile is selected within the Active Era app you will be prompted to save the data to the user profile specific to your baby or to create a new user profile for your baby.

Note, for best results, we recommend you use baby mode when your baby is calm.

## **WEIGHT ONLY MEASUREMENT**

---

These smart scales can also be used without the Active Era app to provide just a weight measurement.

1. Place the smart scales on a dry, solid, and flat surface.
2. Step onto the smart scales.
3. The smart scales will turn on and an analysis of your body's weight will start automatically.
4. Once a final weight measurement has been obtained it will flash on the LCD display.

Note, to ensure accuracy of results, it is recommended to weigh yourself at the same time of day as measurements may vary after exercise, dieting, or dehydration.

## **BODY MEASUREMENTS EXPLAINED**

---

These smart scales use bioelectrical impedance to conduct an analysis of your body's composition. A very small and safe electrical current is passed through the lower half of the body. The electrical current flows slower through bone and fat than it does through muscle and water. The smart scales can measure this difference in electrical current speed and using a multi-step mathematical formula estimate 13 individual measurements.

Please note, these smart scales have been calibrated for those individuals who have a median body type and therefore the measurements recorded by individuals with unique body types e.g. significant upper body muscle, athletes or bodybuilders may not be accurate. To improve the accuracy of the measurements for individuals with unique body types ensure that athlete mode is selected when creating an account.

### **WEIGHT**

Weight refers to the total weight of the body.

## BMI

Body Mass Index (BMI) refers to a measure that uses your height and weight to work out if your weight is healthy. Refer to the table below for reference.

| <b>CATEGORY</b> | <b>BMI</b>  |
|-----------------|-------------|
| Underweight     | < 18.5      |
| Healthy Weight  | 18.5 - 24.9 |
| Overweight      | 25 - 29.9   |
| Obese           | 30 - 34.9   |
| Severely Obese  | 35 - 39.9   |
| Morbidly Obese  | > 40        |

## BODY FAT PERCENTAGE

Body fat percentage refers to the percentage of subcutaneous fat (outer fat) closest to the skin in the body. Men and women have a different distribution of body fat so it is normal for body fat percentage to vary based on gender. Refer to the table below for the normal ranges for body fat percentage as a percentage of total body weight for reference.

| <b>AGE</b> | <b>MEN</b> | <b>WOMEN</b> |
|------------|------------|--------------|
| 20 - 39    | 8% - 20%   | 22% - 33%    |
| 40 - 59    | 11% - 22%  | 24% - 34%    |
| 60 - 79    | 13% - 24%  | 25% - 36%    |

## **VISCERAL FAT**

Visceral fat refers to a measure of visceral fat in the body. Visceral fat is found in the abdomen and surrounding the vital organs. Refer to the table below for reference.

| <b>CATEGORY</b> | <b>VISCERAL FAT</b> |
|-----------------|---------------------|
| Normal          | 1 - 9               |
| High            | > 10                |

## **BODY WATER PERCENTAGE**

Body water percentage refers to the percentage of water in the body. It is important to stay well hydrated for overall health and wellbeing. Refer to the table below for the normal ranges for body water percentage as a percentage of total body weight for reference.

| <b>MEN</b> | <b>WOMEN</b> |
|------------|--------------|
| 50% - 65%  | 45% - 60%    |

## **SKELETAL MUSCLE PERCENTAGE**

Skeletal muscle percentage refers to the percentage of skeletal muscle (muscle connected to the skeleton) in the body. Refer to the table below for the normal ranges for skeletal muscle percentage as a percentage of total body weight for reference.

| <b>MEN</b> | <b>WOMEN</b> |
|------------|--------------|
| 40% - 60%  | 30% - 50%    |

## MUSCLE MASS

Muscle mass refers to the total weight of muscle in the body. Refer to the table below for the normal ranges for muscle mass as a percentage of total body weight for reference.

| AGE     | MEN       | WOMEN     |
|---------|-----------|-----------|
| 20 - 39 | 34% - 40% | 25% - 30% |
| 40 - 59 | 33% - 39% | 24% - 30% |
| 60 - 79 | 33% - 38% | 24% - 30% |

## BONE MASS

Bone mass refers to the total weight of bone in the body. Refer to the table below for the normal ranges for bone mass as a percentage of total body weight for reference.

| MEN     | WOMEN     |
|---------|-----------|
| 3% - 5% | 2.5% - 4% |

## PROTEIN PERCENTAGE

Protein percentage refers to the percentage of protein in the body. Refer to the table below for the normal ranges for protein percentage as a percentage of total body weight for reference.

| MEN       | WOMEN     |
|-----------|-----------|
| 16% - 18% | 14% - 16% |

## **BMR**

Basal Metabolic Rate (BMR) refers to the number of calories the body requires when at rest to function effectively. BMR is relative to gender, body mass, age, weight, height, and physical activity level.

Within a healthy, balanced diet, a man needs approximately 2,500kcal a day to maintain their weight and function effectively. For a woman, that figure is approximately 2,000kcal a day.

## **METABOLIC AGE**

Metabolic age refers to a measure that compares your BMR to the average BMR of your age group. If your metabolic age is higher than your actual age, it's a sign that you need to improve your metabolic rate.

## **CLEANING**

---

Use a soft, dry or dampened cloth to clean away any dust, or debris from the exterior surface of the smart scales.

Do not wipe or clean the smart scales with chemical solvents.

## **MAINTENANCE**

---

The smart scales require little maintenance. Do not try to fix them by yourself. Refer them to a qualified service personnel if assistance is required.

## STORAGE

---

For long-term storage:

1. Remove the batteries from the battery compartment.
2. Cover the smart scales with a waterproof covering and store horizontally in a safe, dry place, out of the reach of children.

## ERROR CODES

---

| ERROR CODE | CAUSE                       | SOLUTION   |
|------------|-----------------------------|--|
| Lo         | Battery power is low        | Replace batteries  |
| Err        | Smart scales are overloaded | Step off smart scales to avoid damaging the high precision sensors |
| C          | Measuring error             | Step off smart scales and back on again                            |

## TROUBLESHOOTING

---

| PROBLEM                                       | CAUSE  | SOLUTION  |
|---|--|---|
| Smart scales not turning on                   | Batteries not installed correctly                    | Ensure batteries are properly installed   |
|   | Battery power is low                                 | Replace batteries   |
| Smart scales not connecting to Active Era app | Bluetooth not on                                     | Turn smartphone or tablet's Bluetooth on  |
|   | Smartphone or tablet operating system not compatible | Ensure your smartphone's or tablet's operating system is running iOS 8.0 and higher or Android 5.0 and higher |
| Smart scales not analysing body composition   | Wearing shoes and/or socks                           | Remove shoes and/or socks and ensure your feet are in full contact with the metal sensor plates               |

## **DECLARATION OF CONFORMITY**

---

Hereby, One Retail Group Limited declares that the radio equipment type BS-05 is in compliance with Directive 2014/53/EU. The full text of the EU declaration of conformity is available at the following internet address:  
<https://activeera.com/ec-declaration-of-conformity-pdf/>

## **FCC STATEMENT**

---

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation. Warning: Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

## RECYCLING & DISPOSAL

---

Waste electrical and electronic products (WEEE), batteries, accumulators, and packaging should not be disposed of with general household waste. This is to prevent possible harm to the environment and human well-being, and to help conserve our natural resources. Most electrical products with a plug, battery, or cable can be recycled. Contact your local council, household waste disposal service, or the store where you purchased the product for information regarding the collection schemes in your area.

### Waste Electrical and Electronic Equipment (WEEE)

This symbol on the product and/or packaging indicates that Waste Electrical products (WEEE) under Directive 2012/19/EU which governs used electrical and electronic appliances, should not be disposed of with general household waste.

### Waste Batteries and Rechargeable Batteries

This symbol on the batteries supplied with your product indicates that Waste Batteries under Directive 2006/66/EC which governs used batteries and accumulators should not be disposed of with general household waste. If the battery contains more than 0.004% lead, the battery symbol will be accompanied by the letters (PB).

### The Mobius Loop

This symbol indicates that the product or packaging is capable of being recycled.

## **Take-Back Scheme**

Waste electrical products (WEEE) can be sent to One Retail Group for recycling and disposal when you purchase a new product from One Retail Group, on a like-for-like basis. For more information on how to recycle this product please visit [www.activeera.com](http://www.activeera.com)

The end-user of electrical and electronic products (WEEE) is responsible for separating old batteries and lamp bulbs before they are handed over to a designated disposal and recycling service. Waste batteries can also be disposed of in collection bins, which can be found in most local supermarkets and household waste recycling centres.

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© Copyright 2022
EU Authorised Representative: Brandrep Limited