

The logo for SIMUSI, featuring the word "SIMUSI" in a stylized, lowercase, sans-serif font. The letters are dark gray and have a modern, rounded appearance. The logo is centered within a large, light gray circular ring that frames the text.

simusi


**Remote Control
Instructions**

Directions for use

Step One


1. Separate the device and power is on with blue indicators.
2. Press the "ON/OFF" button of the remote control for 2 seconds to turn on it, it will connect to the tens unit automatically.

Note: First time to connect the remote to the devices, you need to press the "ON/OFF" button for about 5 seconds.

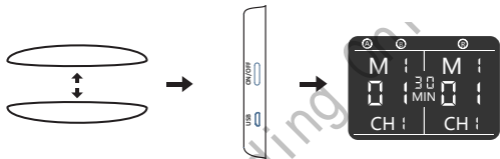


Step Two

1. Attach a gel pad to each piece of the device and remove the protective films.
2. Attach the device where you need a stunning massage.
3. Select the mode and intensity, and press "▶" to start massage.

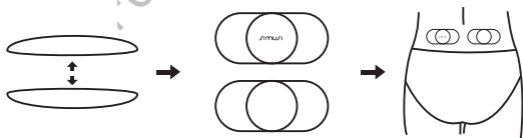

Instructions for Remote Control

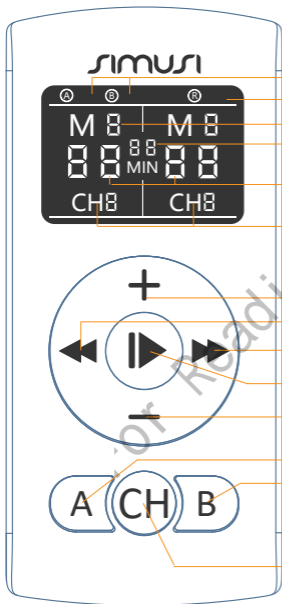

Indicators

- 1 Tens unit battery
- 2 Remote battery
- 3 Modes (18 modes)
- 4 Timing (30mins per treatment)
- 5 Intensity (15 levels)
- 6 Channels (3 channels)

Key Instructions

Intensity level up

Previous Mode

Next Mode

Play/Stop

Intensity level down

Press to select Device A

Press to select Device B

Could set different modes and intensities for each devices

Press to switch channel

CH1: 6 classic modes

CH2: 4 traditional modes

CH3: 8 scene modes

Instructions for Main Functions

Channels and Modes

CH1 - Classic Mode (6 modes)

CH1 M1: knead CH1 M2: press CH1 M3: squeeze CH1 M4: stroke
CH1 M5: elbow CH1 M6: beat

CH2 - Traditional Mode(4 modes)

CH2 M1: massage CH2 M2: acupuncture CH2 M3: cupping jar CH2 M4: scrap

CH3 - Scene Mode(8 modes)

CH3 M1: zen CH3 M2: home CH3 M3: office CH3 M4: seaside
CH3 M5: valley CH3 M6: vacation CH3 M7: rest CH3 M8: travel

Tips: Through pressing button A/B to choose A/B device
You can set different modes and intensities for each devices.



Device A
Mode: squeeze
Intensity: 3

Device B
Mode: home
Intensity: 5



Device A
Mode: scrap
Intensity: 15

Device B
Mode: beat
Intensity: 7

Directions for App Download

Method One

Search App [Relax Place] in AppStore or GooglePlay



This application is **ONLY** compatible with



iOS 7.0, Android 4.3 or above,



Bluetooth 4.0 or above.

Method Two

Scan the QR Code to download APP [Relax Place]



QR code for AppStore



QR code for GooglePlay

Note: Please refer to the main instructions book for more informations.