ULTRASONIC ORAL CLEANER

User Instruction

Thanks for choosing our product, to safely and effectively use this product, please read the instruction carefully.

Point to Note :

- ① The head of this product can be washed, it is not recommended to wash the whole body,
- ② It is normal that the temperature of the product may rise slightly when it is working or charging.
- ③ Please keep this product out of the reach of children
- ④ This product can not be placed in the sun exposure or place with too much humidity.
- ⑤ Do not disassemble this product without authorization to prevent electric shock or battery short circuit.
- ⑥ If don't use it for a long time, please make sure to charge it at least once every two months.
- ⑦ Please do not use in the following situations:
 - a. Patients with periodontitis, gingivitis and other dental diseases;
 - b. People during gums bleeding period;
 - c. People during dental fillings, veneers, and other dental procedures period;
 - d. Patients with oral diseases, pregnant women, heart disease patients, and people who use artificial heart pacemakers.

Components:

| | |
|----------------------|---|
| 1 X Host Machine | 4 X Replacement Head |
| 1 X Dental Mirror | 1 X Charging Cable |
| 1 X User Instruction | 1 X Protection Cover |
| 1 X Mouth Opener | 1 X Replacement Head Disassembly Wrench |

Specification:

| | |
|--|-----------------------|
| Charging Current: $\leq 1A$ | Charging Time: 2-3.5h |
| Operate Frequency: 30-40Hz | Power: $< 3W$ |
| Product Weight: 137g | Package Weight: 337g |
| Size: 27*27*190mm(Not including replacement head length) | |

Tips: How to distinguish between genuine and fake ultrasonic dental cleaners?

A genuine ultrasonic dental cleaner uses high-frequency vibration generated by a high-frequency ultrasonic transducer, you will not feel the vibration when you hold it in your hand. And you will only feel the vibration when you touch a hard object with your hand. When eccentric wheel motors produce false ultrasound, you will feel the vibration of your hand.

FQA:

- Q** Is it normal for my gums to bleed when I use it?
- A** If you suffer from periodontitis or gingivitis, it can easily cause your gums to bleed, even if you brush your teeth or chew hard food. If your gums bleed while using this product, stop using it and rinse your mouth with cold water.
- Q** Why the cleaner switches off sometimes when using it?
- A** If the battery indicator flashes during use, or it turns off automatically as soon as the unit is turn on, the battery is generally too low. Please use it after charging.
- Q** Is it normal the dental tip and cleaner temperature get slightly higher when using it?
- A** As the internal parts are in working condition, it is normal for the machine temperature to be slightly higher.
- Q** Is it normal if the cleaner have noises sometimes when using it?
- A** It is normal to have some noise, its because of the replacement head loose or the resonance. If the replacement head looses, please tighten it.

Charging Steps:

When charging, please use the charging cable come with the product and insert into the host machine, and insert the other end into the 5V charger for charging. The power light flashes when charging and stays on after fully charged, which takes about 3 hours. 2A charger adapter is recommended.

Operation Steps:

1. Due to safety considerations, the host machine only contains a small amount of electricity and not fully charged. Please charge it before use.
2. Know the various parts of product:

Protective Cover

Host Machine

Replacement Head

Charging Cable


Dental Mirror

Disassembly Wrench

3. Install the replacement head and tighten it with the disassembly wrench.

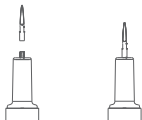4. Press and hold the button for 2 seconds to turn on the device. Enter the gear mode, you can press and hold the button for 2 seconds to turn off in any gears; to switch four gears (soft, normal, middle, strong),
5. With memory function, when you turn on the device next time, the gear will be the last time you used,
6. This product will automatically turn off after working for 10 minutes.
7. If the battery indicator flashes during use, or it turns off automatically as soon as the unit is switched on, the battery is generally too low. Please use it after charging.

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