

# ACTIVATE YOUR RESULTS

**ESSENTIAL FIRST STEP:  
MICROCURRENT  
ACTIVATOR**

Always prime your skin with NuFACE® Microcurrent Activators before treatment. Each ionized formula is specially formulated to deliver microcurrent to the skin at every level for optimal results and is supercharged with active hydrating or anti-aging ingredients.

**⚡ FIX® Line Smoothing Serum**

Diminishes the look of fine lines and wrinkles, while tightening and brightening for a radiant appearance.




**N Ū F A C E®**



Follow @MyNuFACE for Aesthetician Tips, Demos, and Live Q+A.  
Share your routine or results with #MyNuFACE

75180

**N Ū F A C E®  
QUICK START**

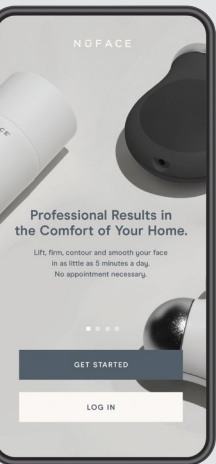

# BEFORE YOU POWER UP

## CHARGE DEVICE FOR 8 HOURS

Before your first treatment, charge your new NuFACE FIX® Line Smoothing Device and read the enclosed user manual.

## SEE HOW EASY TO USE!

Get more out of your new device by visiting [useNuFACE.com](http://useNuFACE.com) to watch step-by-step tutorials and learn more.

# GETTING STARTED

## 4 SIMPLE STEPS

### 1] PREP

Cleanse skin using an oil-free cleanser. This will wash away any excess oils that could otherwise block or prevent an effective microcurrent treatment.

### 2] ACTIVATE

Apply a mask-like layer of FIX® Serum. Use enough formula to ensure the device glides easily along the skin.

### 3] SMOOTH

Push the ON/OFF button on your FIX® Device; the white LEDs will illuminate, and the energizing hum will begin. Select a preferred treatment technique. Begin treatment using a quick erasing motion while applying light to medium pressure.

### 4] FINISH

When complete, massage in remaining FIX® Serum for added skincare benefits.

Visit [useNuFACE.com](http://useNuFACE.com) to see how easy to use!

# 3-MINUTE TREATMENTS



## AREAS

TIGHTEN FOREHEAD

SMOOTH FINE LINES AROUND EYE AREA

REDUCE SMILE LINES

ENHANCE LIP-AREA FULLNESS

## TECHNIQUES

**HORIZONTAL**  
GLIDING MOTION



**VERTICAL**  
GLIDING MOTION



**CIRCULAR**  
MOTION



Optimal results are achieved with regular use.