

EN INSTRUCTION FOR USE

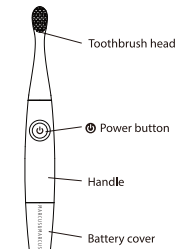Specification
Operation voltage: 1.5V DC
Battery: AAA Alkaline battery IPX7
Sonic Vibration:
Frequency 7000 oscillations / min


Warranty
1 year warranty applies to this product from the date of purchase. According to the condition of the product, defected caused by the device itself or during production will be replaced free of charge for 1 year.
Please register online at www.marcusmarcus.com

How To Place Battery To Kids 2-min Timer Sonic Toothbrush
1. Turn battery cover as direction (1) and pull down as (2).
2. Insert a new AAA alkaline battery with (+) and (-) in correct orientation.
3. Install battery cover as (3) and turn it firmly as (4).

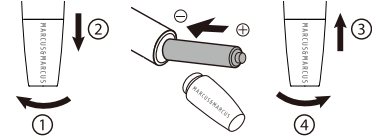How To Insert Toothbrush Head Into Kids 2-min Timer Sonic Toothbrush

1. Plug the new toothbrush head into the notch and twist clockwise until the toothbrush head is secured. Wash toothbrush head thoroughly with warm-soapy water before initial use.
2. To replace the old toothbrush, twist the old toothbrush head anti-clockwise and pull it off.

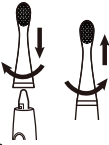How To Use Kids 2-min Timer Sonic Toothbrush

1. Apply toothpaste to the toothbrush head before pressing .
2. When the vibration of the toothbrush is gradually weaken indicating that the battery is almost drained, replace it by a new AAA alkaline battery. Before replacing the battery, ensure the handle is thoroughly dry to prevent water drops entering it.


CLEANING & CARE:

- Wash toothbrush head thoroughly with warm-soapy water before initial use. After use, rinse toothbrush head thoroughly, shake off excess water and store upright. Allow air dry.
- Do not immerse the toothbrush under water.
- Do not stroke or knock the toothbrush with too much strength.
- Do not use any other toothbrush heads other than that recommended by the manufacturer.
- Make sure the battery cover is fixed firmly when washing the toothbrush.
- Take out battery of toothbrush if no intention to use for an extended period of time.
- Store products in a dry place. Do not leave products in direct sunlight or near a source of heat.

FOR YOUR CHILD'S SAFETY AND HEALTH WARNING:

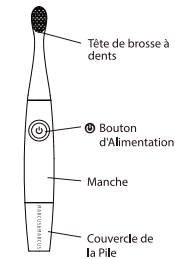- This is not a toy. This product is for tooth brushing only.
- Not suitable for children under 3. Always use with adult supervision.
- Do not leave child unattended while using the product.
- Keep away from children when not in use.
- Inspect the product before each use.
- Do not chew or bite on toothbrush. Battery choke hazard.
- Do not install batteries backwards, charge or put in fire. May explode or leak causing injury.
- Always use AAA alkaline batteries. Use of non-alkaline batteries may result in a drastically reduced operating time.
- If battery cover is not secured, water leakage may cause malfunction of the toothbrush.
- Dentists recommend to replace silicone toothbrush head every 6 months or sooner if it looks worn.
- If your child is undergoing treatment for any oral condition or experiences any discomfort, consult your dental professional prior to use.
- Discard and replace at first signs of damage or weakness.

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