



# USER GUIDE

---

TOM5 STREETER

**Synetik**  
ErgoCare

# Table of Contents

<b>Intended Use</b> .....	3
<b>Warranty</b> .....	4
<b>Safety and Warning Notes</b> .....	5
<b>Transportation</b> .....	6
<b>Overview of the Tom5 Streeer</b> .....	7
<b>Technical Specifications</b> .....	8
<b>Configuration</b> .....	8
Box Contents .....	8
Unboxing the Tom5 Streeer .....	9
<b>Assembling the Tom5 Streeer</b> .....	9
Unfolding the Stroller-style Wheelchair Frame (without the seat) .....	9
Fitting the Seat in the Forward-facing Position .....	9
Fitting the Seat in the Rear-facing Position: .....	10
Lateral Protector Installation .....	10
5-Point Belt Installation .....	11
Padded Frame Cover Installation .....	11
<b>Adjustment</b> .....	12
1. Seat-unit Tilt-in-space .....	12
2. Seat-depth Adjustment .....	12
3. Seat-base Width Adjustment .....	12
4. Backrest Angle Adjustment .....	12
5. Backrest Height Adjustment .....	13
6. Footrest Angle Adjustment .....	13
7. Footrest Height Adjustment .....	13
8. Footrest Angle Adjustment .....	13
9. Footrest Velcro Lock/S Adjustment .....	13
<b>Daily Usage</b> .....	14
Parking Brake .....	14
Swivel Lock .....	14
Handle Height Adjustment .....	14
Unfolding the Stroller (with the Seat) .....	15
Folding the Stroller (with seat attached) .....	15

Getting into the Stroller .....	16
Getting out of the Stroller.....	16
Negotiating Stairs .....	16
Negotiating Curbs or Single Steps .....	17
<b>Accessories</b> .....	<b>18</b>
Basket - Fabric.....	18
Knee Block .....	18
Calf Rest .....	19
Basket-HD (Steel) .....	19
Tilt-lock for Spastic .....	19
Grab-Bar Bracket Cover.....	19
Extended Side Protection .....	19
Footmuff - Summer .....	20
Footmuff - Winter.....	20
Work Board .....	20
Bag.....	21
Canopy.....	21
Abduction Wedge.....	21
Headrest (pair).....	22
Headwings .....	22
Neckrest.....	23
Side Trunkrest (pair).....	23
Lumbar Support (w.25cm, 35cm) .....	24
Backrest Wedge.....	24
2-point Belt.....	25
4-point Belt.....	25
Chest-belt Addition to 4-point Belt .....	26
Chest Harness .....	26
Foot Harnesses .....	27
H-style Harness .....	27
Immobilisation Waistcoat.....	28
Abduction Belt .....	29
Rain Cover.....	29

Transit Option: 010-186 .....	30
<b>Cleaning Instructions</b> .....	<b>31</b>
Removing the Upholstery .....	31
<b>Maintenance</b> .....	<b>31</b>

**IMPORTANT:**

The information in this document is subject to change without notice. For more information, please contact our Customer Service at:  
Toll free: 1-888-752-8818  
Phone: 1-450-759-9449  
Fax: 1-450-752-1644  
Email: [info@synetikergocare.com](mailto:info@synetikergocare.com)  
Website: [www.synetikergocare.com](http://www.synetikergocare.com)  
Business Hours: 8 am to 4 pm Eastern Time

Only use parts and accessories supplied by Synetik ErgoCare. Any substitution of parts with those of another company will result in the loss of warranty and may cause damage.

Read this manual carefully before using the stroller for the first time and keep it for future reference. The safety of a child could be compromised if you do not follow the instructions and safety notices mentioned in the User Guide.

In case of apparent defects, which limit the safe use of the product (for example, the chassis is damaged, there are cracks or a visible defect), it is imperative to immediately stop using the product and contact your distributor or local repair centre. Synetik ErgoCare assumes no liability for damages and accidents caused by improper use of its products.

**Intended Use**

The stroller-style wheelchair TOM5 Streeter comes in 3 sizes (mini, standard and maxi). TOM5 is designed for both short and long term use of children with special needs This stroller-style wheelchair incorporates the ALUX production method for high durability and low chassis weight. The TOM 5 is one of the lightest stroller-type wheelchairs in its category. The seat unit is detachable, reversible and entirely size and position adjustable. It also incorporates the innovative MONOSLIDE design as well as the tilt-in-space telescope-operated system, and flexible system for attaching seating restraints and positioning accessories. Fully foldable frame even with seat attached. Soft shell is covered by fully detachable, soft padded upholstery in fashionable design. Solid steel for a stable ride. You will find all the details in the assembly section and on our website: [www.synetikergocare.com](http://www.synetikergocare.com).

## Warranty

1. A 24-month warranty is applicable on this product from the date of purchase. The warranty applies to material and manufacturing defects.
2. The warranty is applicable at the location where the product was purchased, or in an authorized service repair shop.
3. A properly completed Certificate of Warranty and a copy of the proof of purchase must be presented when applying for the warranty. Do not remove the identification label. It must be legible and undamaged throughout the entire warranty period.
5. The product must be used entirely for the purpose for which it was designed. The user must follow the directions mentioned in the Instructions for Use. The manufacturer/distributor is not responsible for defects caused by any handling that is contrary to the Instructions for Use.
6. The product must be used, treated, maintained and stored properly according to directions mentioned in the Instructions for Use.
7. Warranty proceedings can be initiated, if the product under claim has been returned completely in the original package, carefully packed and protected against damage. The owner is obliged to ensure the product under claim is properly cleaned and conforms to general hygienic principles. The owner cannot ask the seller to do basic maintenance or cleaning of the product throughout the claim proceedings.
8. All defects must be claimed immediately.
9. The warranty is forfeited if the defects were caused by overloading, accident, violence or inappropriate tampering or mishandling of the product.

If you have any questions or comments, you can contact our Customer Service.

## IDENTIFICATION SHEET

Identification tag / Product type (write up):



Final inspection:

Name:

Signature:

## Safety and Warning Notes

**ATTENTION!** Read this manual carefully before first using the stroller, and keep it for further reference. The safety of a child might be in jeopardy if the instructions and safety information listed in the User Guide are not followed.

Safety points to consider when using the stroller-style wheelchair TOM5 Streeter.

- BAGS ATTACHED TO THE PUSH HANDLE AFFECT THE STABILITY OF THE STROLLER
- NEVER LEAVE A CHILD UNATTENDED
- MAXIMUM LOAD CAPACITY OF THE STROLLER:
  - Mini: 35 Kg, STD: 40 Kg, Maxi: 50 kg
- **WARNING !** Keep fingers and covers away from any moving parts when using the stroller.
- **WARNING !** Never leave a child unattended! Even if the belts and brakes are in use.
- **WARNING !** Check before use that all locking devices are properly secured, and all safety equipment is properly used.
- **WARNING !** The seat unit is not suitable for a child younger than 6 months.
- **WARNING !** This product must not be used while jogging or rollerblading.
- The attendant must be familiar with the instructions of this manual.
- Professional medical advice regarding the suitability of this product must be sought before use.
- Always make sure, before heading out, that the child is comfortable, and that clothes and blankets, etc., are tucked in, so they do not get caught in the wheels.
- Ensure the child is clear of moving parts before making any adjustments, to avoid any injury.
- When angle positioning and adjusting the seat and backrest, act accordingly to avoid sudden uncontrolled movements.

- Do not allow the child to climb unassisted into and out of the stroller.
- Always apply the parking brakes when putting the child in and taking them out of the stroller.
- Make sure that the belts and other securing accessories, safety bar and other equipment are fitted and adjusted properly.
- Do not stand on the footplate.
- The brakes must always be engaged when the stroller is not moving.
- Getting over the curb: Practice first using an empty stroller.
- The stroller is not designed to be used in the stairs.
- The attendant must take precautions against fire risk. In buildings, check for the location of fire exits and ensure all procedures are understood.
- Do not let children play with or hang onto the stroller.
- Do not overload the stroller. Overloading the chassis, seat, bag or basket causes a reduction in stability of the stroller.
- Pay full attention, when handling the stroller, so as not to injure the child or yourself.
- If it is possible, always use ramps or elevators to avoid any access restrictions.
- If possible, two people should lift the stroller.
- When you use the stroller on the road, follow the road traffic regulations.
- Do not use the stroller for any purpose other than that for which it was designed.
- The stroller is designed to carry only one child at a time.
- Do not use anything other than original PATRON accessories designed for this product.
- When the stroller is folded and unfolded, the User Guide must be followed carefully, to avoid any damage to the construction that may affect product stability.
- In case of apparent defects, which limits the product's safe use (e.g. damaged chassis, part comes off or other visible defects, cracks, etc.) it is necessary to stop using the product immediately and contact your local official distributor or service centre.
- Synetik ErgoCare does not bear any responsibility for any damages and/or accidents caused by inappropriate use of this product.

## **Transportation**

The stroller-style wheelchair TOM5 Streeter has successfully passed impact tests. The seat meets the ISO 7176-19 Standard. The tests were carried out in all possible configurations of the product.

To be compliant, the stroller must have the transit option which includes 4 tie-downs.

# Overview of the Tom5 Streeter

## SCHEMA



### SCHEMA DESCRIPTION

1. Push handle	16. Rear wheel
2. Push handle height-adjustment hinge	17. Parking brake
3. Push handle folding release	18. Basket brackets
4. Backrest	19. Seat width-adjustment plastic nuts
5. Backrest height-adjustment hex bolts	20. Footrest angle-adjustment lever



6. Backrest height-adjustment plastic nuts	21. Seat depth-adjustment hex bolts
7. Backrest recline-adjustment lever	22. Footrest height-adjustment hex bolts
8. Canopy brackets	23. Footrest angle-adjustment hex bolts
9. Interface lock safety pin	24. Footrest angle-adjustment slide plate
10. Interface locking release	25. Footrest Velcro band
11. Tilt-in-space adjustment lever	26. Footrest
12. Safety bar	27. Padded frame cover
13. Footrest	28. Lateral protectors
14. Front swivel castor wheel	29. Transport safety catch
15. Swivel lock	

## Technical Specifications

Measurements	Mini	Standard	Maxi
Backrest Height *	17.3 – 22.4	20.1 – 26.8	21.7 – 27.6
Seat Depth *	7.3 – 10.6	9.1 – 13.0	11.4 – 15.8
Seat Width *	8.7 – 11.0	9.8 – 13.8	10.6 – 15.4
Seat to Footrest *	3.5 – 10.2	4.7 – 14.2	4.7 – 15.4
Backrest Angle (°)	0 – 170	0 – 170	0 - 170
Seat Angle (°)	+38 / -28	+38 / -28	+38 / -28
Total Length *	30.3	30.3	30.3
Total Width *	23.4	25.2	26.6
Total Height/with Handle *	37.4 / 43.3	37.4 / 43.3	38.2 / 46.5
Total Weight *	18 kg	19 kg	20 kg
Maximum Load *	40 kg	40 kg	50 kg
Folded Stroller Length	33.8	34.6	35.0
Folded Stroller Width	23.4	25.2	26.6
Folded Stroller Height	14.6	15.4	16.9

\*inches

## Configuration

### Box Contents

Upon delivery, you will receive one (1) box containing the stroller and the accessories ordered. The only tool required for assembly is a hex key that is found in a small pocket on the back of the backrest. Always put the key back in its place.

### Items List:

- Stroller-style wheelchair detachable and reversible seat unit
- Stroller-style wheelchair frame
- Grab bracket covers
- Padded frame covers
- Lateral protectors
- Accessories
- 5-point Belt

### Unboxing the Tom5 Streeter

- 1- Remove the contents from the boxes and place on a flat surface.
- 2- Make sure that the material corresponds with what you ordered and that no damage was caused during transportation.
- 3- Follow the instructions in the following section (Assembling the Tom5 Streeter).

## **Assembling the Tom5 Streeter**

The instructions in the User Guide will direct you through the assembly process of the stroller-style wheelchair. However, it is recommended that the initial adjustment for the user be done in the presence of a Synetik ErgoCare representative.

**NOTE! The numbers in parenthesis refer to the Schema presented on Page 6.**

### Unfolding the Stroller-style Wheelchair Frame (without the seat)

- 1- Place folded stroller frame on its wheels.
- 2- Unsnap and release the transport safety-catch (29) on frame.
- 3- To unfold the stroller frame, stay between the rear wheels and hold the handle (1) tightly. Move the handle up and back until the handle-folding release (3) is fully locked into its end position.
- 4- To check the locking position, move the handle upward to ensure that both mechanisms are fully locked (the stroller frame must not start folding inward).

### Fitting the Seat in the Forward-facing Position

- 1- Make sure that the stroller frame is unfold and locked into the unfolded position. (see previous section).
- 2- Hold the seat in a way that it faces towards the front of the stroller.

- 3- Slide the seat interface tubular adaptors directly into the seat-frame interface sockets. Start by inserting just one side of the seat interface tubular adaptors into the OPEN interface socket (without locking pins).
- 4- Push the second side down so that the seat locks automatically (you will hear the locking sound).
- 5- Finally check if the seat is properly locked by slightly lifting the seat unit. Seat must not separate from the interface frame.

#### Fitting the Seat in the Rear-facing Position:

- 1- Make sure that the stroller frame is unfold and locked into the unfolded position. (see previous section).
- 2- Hold the seat in a way that it faces towards the back of the stroller.
- 3- Slide the seat interface tubular adaptors directly into the seat-frame interface sockets. Start by inserting just one side of the seat interface tubular adaptors into the OPEN interface socket (without locking pins).
- 4- Push the second side down so that the seat locks automatically (you will hear the locking sound).
- 5- Finally check if the seat is properly locked by slightly lifting the seat unit. Seat must not separate from the interface frame.

#### Lateral Protector Installation

- 1- Attach the two (2) plastic clips on the protector to the vertical tubing on the side of the backrest.
- 2- Repeat for the other protector.



### 5-Point Belt Installation

- a) Run the lower straps between backrest (4) and the side cushions and through the lower cam buckles and secure them by doubling the strap. In a situation where there are two accessories secured by the cam buckles, run the straps through the cam buckles only once.
- b) Lift the head cushion.
- c) Run upper (thin) straps through the holes in the backrest panel, and run them through upper cam buckles, then secure them by doubling the strap. In a situation where there are two accessories secured by the cam buckles, run straps through cam buckles only once.



### Padded Frame Cover Installation

- a) Attach the padded cover around the handle tube and over the handle-folding release (3).




# Adjustment

Before placing the child in the stroller, adjust the height of the backrest, the depth and the width of the seat and install the desired accessories.

## 1. Seat-unit Tilt-in-space

- a) Unzip the right padded frame cover (27) partially from bottom upwards to be able to touch the reclining lever (pic.5).
- b) Hold the seat backrest firmly with one hand.
- c) Pull the seat tilt-in-space lever (11) with your other hand and set the desired position of the seat (pic.5). Release the lever and rezip the padded frame cover.

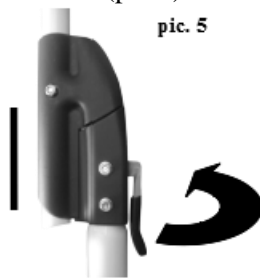

## 2. Seat-depth Adjustment

- a) Remove the seat base upholstery from the seat base (attached by Velcro fasteners).
- b) Loosen the two hex bolts (21) located under the mid-frontal part of seat using the hex key.
- c) Loosen, by hand, the four plastic adjustment nuts located under the seat (19).
- d) Set the desired seat depth by pulling or pushing the front end of the seat.
- e) Properly tighten all hex bolts and plastic nuts and reinstall the seat-base upholstery.

## 3. Seat-base Width Adjustment

- a) Remove the seat base upholstery from the side plates.
- b) Loosen, by hand, the four plastic adjustment nuts located under the seat (19).
- c) Set the seat width by pulling or pushing the side plates.
- d) Properly tighten all plastic nuts and reinstall the side panel upholstery.

## 4. Backrest Angle Adjustment

- a) Hold the seat backrest firmly.
- b) Loose backrest recline-adjustment lever (7).
- c) Set the backrest to the desired angle.
- d) Tighten the lever properly.

- e) If needed, partially unzip the rear part of the backrest cover to make handling easier.
- f) Make sure the backrest is fully locked by gently pushing and pulling on it. Backrest must stay in pre-set position.

#### 5. Backrest Height Adjustment

- a) Remove the backrest upholstery from the rear and front side of the backrest (attached by Velcro fasteners).
- b) Loosen the two hex bolts located on top-rear part of backrest (5).
- c) Loosen the two plastic nuts on the rear side of backrest (6).
- d) Set the height of the backrest by lifting or pressing down on it.
- e) Properly tighten all hex bolts and plastic nuts and reinstall the backrest upholstery.

#### 6. Footrest Angle Adjustment

- a) Hold the footrest (13) firmly.
- b) Release the footrest angle-adjustment lever (20).
- c) Set the footrest to the desired angle.
- d) Tighten the lever properly to secure it into position.

#### 7. Footrest Height Adjustment

- a) Loosen the two hex bolts located on rear side of footrest slider bar (22).
- b) Set the footrest height to the desired position.
- c) Properly retighten both hex bolts.

#### 8. Footrest Angle Adjustment

- a) Loosen the two hex bolts located under the footrest (23).
- b) Set the angle of the footrest by pulling or pushing on the footrest angle- adjustment slide plate (24).
- c) Properly retighten both hex bolts.

#### 9. Footrest Velcro Lock/S Adjustment

The Velcro security band is a safety device that prevents the footrest from folding inward.

- a) Tighten and fasten the Velcro security band located under the footrest (25) to ensure the footrest position is secured.

## Daily Usage

Synetik ErgoCare recommends performing a daily inspection before seating a child in the TOM 5 Stretcher.

### Parking Brake

- 1- To set the brake, press brake lever pad down on both sides fully to the end position (pic.9)
- 2- To disengage the parking brake, lift both brake lever pads.



### Swivel Lock

- 1- If required, the front wheel swivel can be locked.
- 2- To lock the caster forks, lift and turn the swivel lock pin by 90° on both sides and let it drop into the designated hole on the caster fork (pic.10).




### Handle Height Adjustment

- 1- To change the handle height position, press the black plastic buttons on the handle hinges on both sides simultaneously (2).
- 2- Set the handle to desired position, release knobs.
- 3- Handle will lock into new position automatically; you will hear the locking sound.

### Unfolding the Stroller (with the Seat)

- 1- Place folded stroller on its wheels.
- 2- Unsnap and release the transport safety-catch (29) on frame.
- 3- To unfold the stroller frame, stay between the rear wheels and hold the handle (1) tightly.
- 4- Move the handle up and back until the handle-folding release (3) is fully locked into its end position.
- 5- To check the locking position, move the handle upward to ensure that both mechanisms are fully locked (the stroller frame must not start folding inward).
- 6- Pull the seat tilt-in-space adjustment lever (pic.5) and move the seat to the required position.
- 7- Pull the padded frame covers (27) down over the locking mechanism and zip both covers up.
- 8- Release the backrest recline-adjustment lever (7) and set the position of the backrest.
- 9- Tighten the lever properly when the desired backrest angle is reached.
- 10- Release the footrest angle adjustment lever (20), set footrest (13) to the desired position and retighten lever.

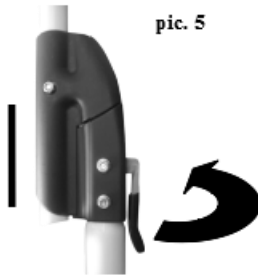


### Folding the Stroller (with seat attached)

- 1- Apply the parking brake.
- 2- Release the footrest angle-adjustment lever.
- 3- Set the position of the footrest just bit less than horizontally (pic.4).
- 4- Release the footrest Velcro security band and fold up the footrest (pic.4).
- 5- Release the backrest recline adjustment lever and push backrest as far as possible forward, towards the seat base. (pic.6).
- 6- Re-tighten the lever.
- 7- Unzip the padded frame covers partially and pull them over the locking mechanism.
- 8- Pull up the seat tilt-in-space adjustment lever (pic.5) and set the position of the seat to tilt forward to maximum (pic.6).



- 9- Lift the push handle-folding release movable parts (pic.6) on both sides and fold the push handle forward towards the front wheels.
- 10- Attach the transport safety latch (pic.14).
- 11- The stroller is folded (pic.7).

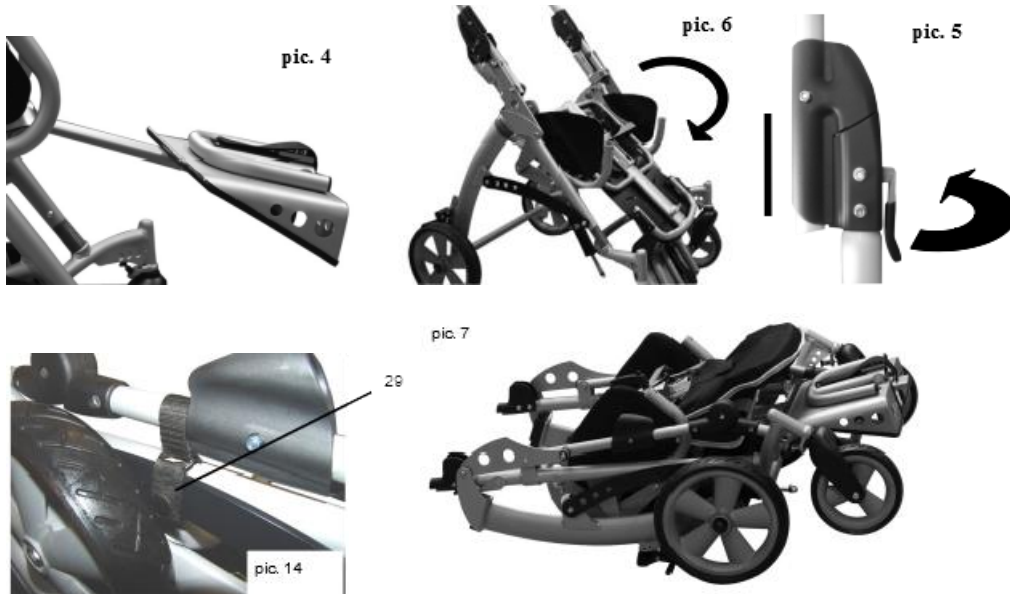

### Getting into the Stroller

- 1- Set the parking brake. If necessary, ask a capable person for help to hold the stroller firmly.
- 2- Place the child into the stroller carefully.
- 3- Fasten the safety belt.
- 4- Adjust the footrest position.
- 5- Adjust all fasteners and positioning accessories the stroller is equipped with.
- 6- Always make sure that the child is seated comfortably in the stroller.

### Getting out of the Stroller

- 1- Engage parking brake. If necessary, ask a capable person for help to hold the stroller firmly.
- 2- Remove all applied fastening and positioning accessories.
- 3- Unfasten the safety belt.
- 4- Take the child carefully out of the stroller.

### Negotiating Stairs

To Facilitate Travel:

1. Make use of ramps and elevators when available.
2. Where possible stairs should be avoided, if there is an alternative to using stairs.
3. If not, ensure at least one helper is available to provide assistance.
4. Do not attempt to take the stroller up or down a flight of stairs unless you and your helper are sure you have control of the stroller and can support it and the weight of its occupant.
5. Do not lift the stroller by the footrest.

If no alternative is possible, follow these instructions:

Going Down Stairs:

- 1- Position the stroller squarely with the front wheels on the top step.
- 2- Tilt the push handle backward slowly and carefully push the stroller forward until the rear wheels are at the top step.
- 3- The helper stands in front, steadily lowering the stroller by firmly holding on to the front frame.
- 4- Support the weight and carefully lower the rear wheels down the stairs.
- 5- Continue step by step to the bottom. (Pic.19)



Pic.19

Going Up Stairs:

- 1- Position the stroller squarely with the rear wheels against the bottom step.
- 2- Tilt the stroller backward slowly and raise toward the next step.
- 3- The helper stands in front, steadies and lifts the stroller by firmly holding on to the front frame.
- 4- Continue step by step to the top (Pic.20)



Pic.20

Negotiating Curbs or Single Steps

Down:

1. Position the rear wheels at the edge of the curb.

2. Hold the push handles firmly and tilt the stroller slightly backwards.
3. Balance on the rear wheels, carefully lower the stroller down off the curb on its rear wheels.
4. Never tip the stroller forward at a curb.

Up:

1. Position the rear wheels at the curb.
2. With a firm grip on the push handles, tilt the stroller backwards and place the rear wheels on top of the curb.
3. Position the front wheels at the curb bottom and lift to the upper level.

## Caution

1. During this operation, the attendant and helper must be well-positioned and have a secure grip on the stroller.
2. The attendant's hands must be inside the push handle straps at all time.
3. Rest at regular intervals.
4. Do not lift the stroller by its
5. footrest.

## Accessories

### Basket - Fabric

- 1- Clip the basket onto chassis-frame tubing. **Maximum load capacity is 3 kg.**



### Knee Block

- 1- Loosen the two hex bolts positioned on the rear side of footrest slider bar and remove the footrest.
- 2- Install the knee block.
- 3- Secure it into position and insert the bolts into the assembly holes.
- 4- Reinstall the footrest into the correct position and tighten both hex bolts.

### Calf Rest

- 1- Loosen the two hex bolts located on the rear side of footrest slider bar and remove the footrest.
- 2- Install the calf rest, secure it into position and insert the bolts into the assembly holes.
- 3- Reinstall the footrest into the correct position and tighten both hex bolts.

### Basket-HD (Steel)

- 1- Screw the attachment adaptors for the Basket-HD onto the rear axle.
- 2- Attach the basket to these adaptors with the screws.
- 3- Secure the basket with the belts and with the Velcro in the front part of the stroller.

### Tilt-lock for Spastic

- 1- Unscrew the backrest recline-adjustment lever.
- 2- Put tilt-lock into the groove where the lever is moving and put the second part of tilt-lock onto the screw of the lever.
- 3- Screw the lever back on.

### Grab-Bar Bracket Cover

- 1- Put the plastic plug onto the grab bar bracket.
- 2- Install the elastic cover and fasten the rubber band onto the plastic nut located under the seat
- 3- Repeat the same procedure on the other side.

### Extended Side Protection

#### Mounting –

- 1- Slip the extended side protection over the backrest and onto the tubes of the push handle from both sides with the help of the four plastic loops.
- 2- The upper loops go between plastic sleeve on the push handle and the handle knob.
- 3- The flaps of the extended side protection are secured with Velcro fastener.
- 4- The lower part of the extended side protection is secured with strips.

#### Seated Position –

- 1- In a seated position, it is important that the extended side protection is criss-crossed at the rear of the backrest and secured with a strap and Velcro lock.
- 2- This strap runs through plastic triangles near the upper loop and is secured with a Velcro fastener.

#### Semi-seated Position –

- 1- The extended side protection is secured with a strap and with a Velcro fastener that runs through a plastic triangle.

#### Flat Position –

- 1- Leave the extended side protection loose.



### Footmuff - Summer

- 1- Slip the muff over the footrest and attach it with the zippers around the push handle tubing.



### Footmuff - Winter

- 1- Attach the muff with the rubber over the backrest (take the canopy off if need be) and strap its middle part to the stroller frame.
- 2- The precut leather band sewed into the muff allows for the use of an abductor block, 4-point belt, 2-point belt, abductor belt and H-style harness.

### Work Board

- 1- Remove the push handle cover (if need be).
- 2- Place the work board attachments into the holes.

### Depth Adjustment –

- 1- Release the screws underneath the board, adjust the depth and secure the screws properly.

### Angle Adjustment –

- 1- Release the levers on each side of the work board, adjust the angle and tighten the screws.

### Bag

- 1- Attach the bag to the push handle by its straps.
- 2- **Maximum load capacity 1 kg. Bags and backpacks attached to the push handle affect the stability of the stroller.**



### Canopy

#### Fit

- 1- Remove the side paddings in the upper part of the backrest and locate the canopy brackets
- 2- Insert the canopy plastic clips into both canopy brackets.
- 3- The canopy can be adjusted by pushing onto the canopy frame.

#### Remove

- 1- Press the canopy fitting releases and pull the canopy out of the hood attachments.



### Abduction Wedge

- 1- Seat the child in the stroller and locate the optimal position for the abduction wedge.
- 2- Take the child out of the stroller.
- 3- Find the relevant fitting hole on the seat base.
- 4- Place the wedge on the seat (the narrower part should point toward the backrest) and secure the wedge screw with a nut.



### Headrest (pair)

- 1- Seat the child in the stroller and determine the appropriate height for the headrest.
- 2- Take the child out of the stroller.
- 3- Adjust the backrest height.
- 4- Find the relevant fitting holes on the rear side of the backrest base.
- 5- Insert screws through the slots in the upholstery and the backrest base.
- 6- Secure the headrest screws with nuts.



### Headwings

- 1- Seat the child in the stroller to determine the appropriate height of the headwings.
- 2- Take the child out of the stroller.
- 3- Install the headwings onto the backrest at the determined height, and secure them using Velcro straps.



### Neckrest

- 1- Slightly spread open the headrest cushion.
- 2- Attach the neckrest using the Velcro strap.
- 3- Close the cushion.



### Side Trunkrest (pair)

- 1- Seat the child in the stroller and determine the appropriate height for the trunkrests. (The trunkrests should be just a little bit lower than the occupant's arm pits.)
- 2- Take the child out of the stroller.
- 3- Unzip the backrest zippers.
- 4- Locate the relevant fitting holes on the rear side of the backrest base according to required location.
- 5- Cut slots through the upholstery as long as the metal plates of the trunkrests.
- 6- Put the metal plates of trunkrests through the slots where you just cut the upholstery.
- 7- Put the screws through the holes in the backrest base and secure them in place with nuts.
- 8- You can adjust the desired width of side trunkrests with the help of the horizontal slots in the metal plates of the side trunkrests.





### Lumbar Support (w.25cm, 35cm)

- 1- Unzip the backrest zippers.
- 2- Determine the appropriate height of the lumbar support.
- 3- Cut 5 cm slots through the upholstery.
- 4- Run the lumbar support Velcro through the backrest slots.
- 5- Lift the sides of the backrest upholstery and stick the second part of the Velcro strap onto the metal plate.



### Backrest Wedge

- 1- Peel off the backrest upholstery (with Velcro)
- 2- Stick the wedge of backrest to the desired height using Velcro
- 3- Put the upholstery back into place.



## 2-point Belt

### Fit

- 1- Secure the belt in the lowest plastic cam buckles at the rear of backrest.

### Adjust

- 1- Loosen the fastener.
- 2- Seat the child in the stroller and adjust the belt length.
- 3- To adjust the belt length, use the plastic cam buckles at the rear of backrest.
- 4- The child must always feel comfortable when using the belt.



## 4-point Belt

### Mounting

- 1- Run the lower straps between the backrest and the side cushions and secure it using the plastic cam buckles at the rear side of the backrest.
- 2- Run the upper straps through the appropriate slots in the backrest and secure them using the upper plastic cam buckles at the rear of backrest. In a situation where there are two accessories secured by the cam buckles, run straps through the cam buckles only once, otherwise run them through the buckles in reverse to secure them.

### Settings and Adjustments.

- 1- Release all three belt fasteners, seat the child in the stroller.
- 2- Set desired strap length using plastic cam buckles.
- 3- Other adjustments with the belt are to be made using the plastic cam buckle.



### Chest-belt Addition to 4-point Belt

- 1- The chest-belt addition is installed on the upper (thinner) straps of the 4-point belt.



### Chest Harness

#### Fit

- 1- Seat the child in the stroller and determine the appropriate belt height.
- 2- Take the child out of the stroller.
- 3- Attach the belt onto the backrest with the belts, fasteners and Velcro straps.

#### Usage

- 1- Undo the front Velcro strap.
- 2- Seat the child in the stroller and fasten the Velcro strap to position the child.
- 3- The child must always feel comfortable when using the belt.



### Foot Harnesses

- 1- Attach the straps of the foot harnesses by using the zipper.
- 2- The child's feet are immobilized by a Velcro strap.
- 3- The child must always feel comfortable when using the foot harnesses.



### H-style Harness

#### Mounting

- 1- Run the lower straps between backrest and side cushions and then run them through the lower plastic cam buckles placed on the backplate of the backrest. In the situation where there are two accessories secured by the cam buckles, run the straps through the cam buckles only once, otherwise run them through the buckles in reverse to secure them.
- 2- Run the upper straps through the appropriate slots in the backrest and run them through upper plastic cam buckles located on the backplate of the backrest. In a situation where there are two accessories secured by the cam buckles, run the straps through the cam buckles only once, otherwise run them through the buckles in reverse to secure them.

#### Settings and Adjustments

- 1- Open the plastic cam buckles and fold up the harness.
- 2- Seat the child in the stroller, put the harness on and secure the child using the buckles.
- 3- Adjust the length of the straps with the plastic clip.
- 4- Other adjustments are to be made using the plastic cam buckles.



### Immobilisation Waistcoat

#### Mounting

- 1- Run the lower straps between backrest and side cushions and then run them through the lower plastic cam buckles located on the backplate of the backrest. In the situation where there are two accessories secured by the cam buckles, run the straps through cam buckles only once, otherwise run them through the buckles in reverse to secure them.
- 2- Run the upper straps through the appropriate slots in the backrest and run them through upper plastic cam buckles located on the backplate of the backrest. In the situation where there are two accessories secured by the cam buckles, run the straps through the cam buckles only once, otherwise run them through the buckles in reverse to secure them.

#### Setting and Adjustments

- 1- Unzip the waistcoat.
- 2- Unlock the buckles on the higher straps.
- 3- Seat the child in the stroller,
- 4- Zip up the waistcoat and lock the buckles.
- 5- Adjust the length of the straps with plastic clip.
- 6- Other adjustments are to be made using buckles.



### Abduction Belt

#### Mounting

- 1- Run the Velcro straps through the gap between the backrest upholstery and the seat.
- 2- Place the straps behind the back tubing of the seat and secure the strap using the Velcro.
- 3- Run the thinner upper straps between the backrest and the side cushions and secure them with the upper plastic cam buckles located on the backplate of the backrest. In a situation where there are two accessories secured by the cam buckles, run the straps through cam buckles only once, otherwise run them through the buckles in reverse to secure them.

#### Settings and Adjustments

- 1- Unfasten the buckles and seat the child in the stroller.
- 2- Fasten the buckles.
- 3- Adjust the length of the straps with the plastic cam buckles.
- 4- Other adjustments are to be made using only the plastic buckles.




### Rain Cover

- 1- The rain cover should be pulled over the stroller, so that the canopy and the footrest are completely covered.
- 2- Store only when cover is clean and dry.




Transit Option: 010-186

1. During transit, the chair stroller must be forward facing.
2. The stroller tie-downs must be securely attached to the four red anchor points on the stroller. See above photo for anchor points.
3. All floor tie-down straps must be drawn tight in the front and the back to eliminate any forward/aft movement of the stroller. See pictures A1 to A6.

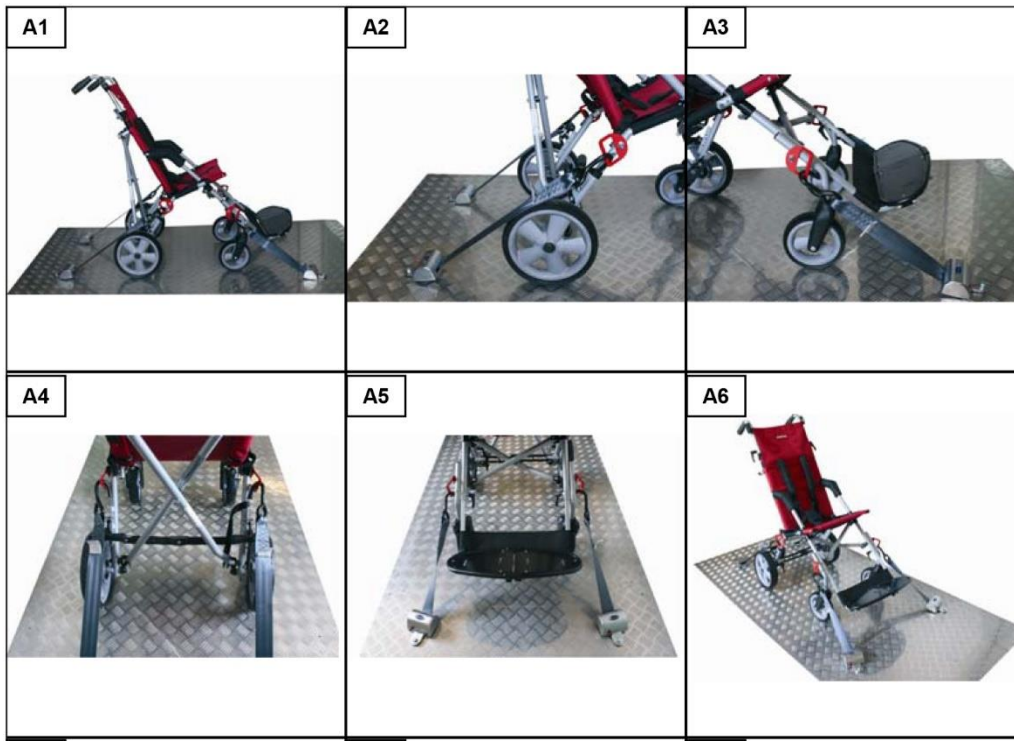

4. Install the 3-point belt of the anchor system as shown in Figures B1 to B5.

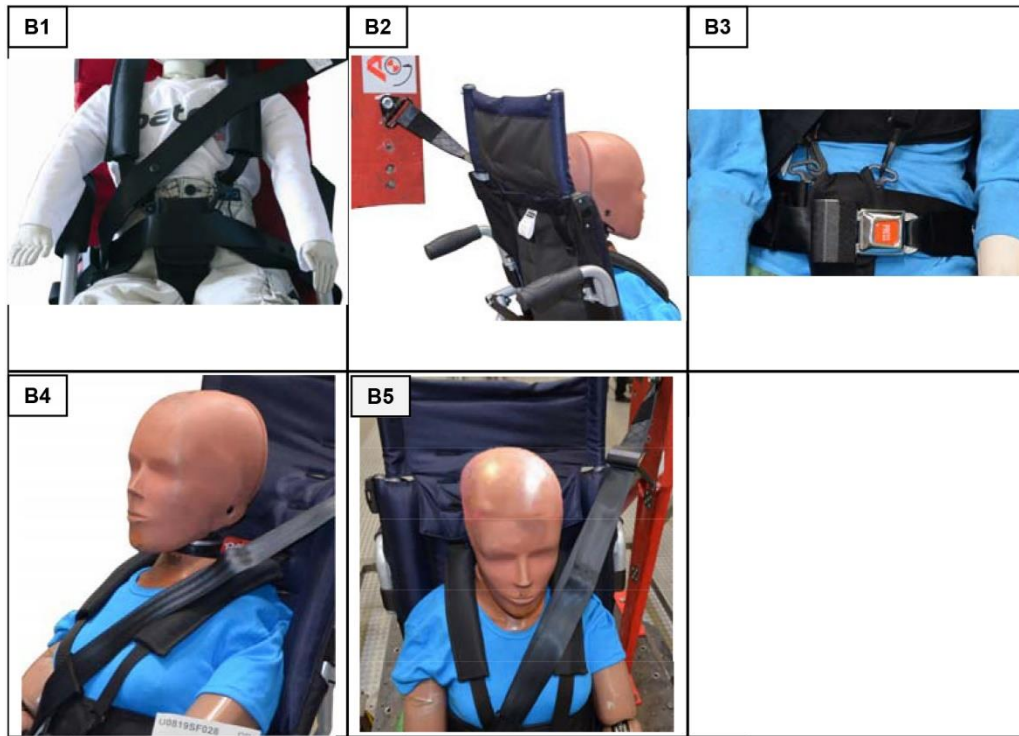

5. Stroller restraint manufacturer instructions must be followed precisely to ensure intended performance.

## Cleaning Instructions

The upholstery is detachable. It can be cleaned using warm water (Max. 40 °C) and a mild detergent. If needed, a replacement upholstery can be purchased.

### Removing the Upholstery

- 1- From the top of the frame, undo the metal snaps located on each side of the stroller.
- 2- Loosen the Velcro fasteners from the 5-point belt and remove all the straps.
- 3- Gently pull up on the backrest upholstery until it comes off the frame.
- 4- Gently pull down the seat and calf upholstery until it comes off the frame.
- 5- To re-attach the upholstery, follow steps in reverse order.

## Maintenance

Synetik ErgoCare recommends that a daily visual inspection of equipment be performed by a therapist, parent or caregiver to ensure the safe use of the product.

- Your stroller has been designed to a high standard and with correct use and maintenance will provide many years of trouble-free service.
- Regularly inspect locking devices, catches, footrest, joint strapping to ensure they are secure and in full working order.



- Brakes, wheels, tires are subject to wear and should be replaced when necessary, by an authorized repair person approved by the stroller supplier.
- We recommend that the stroller be checked by an authorized repair person every 12 months.
- We recommend that after 24 months, visit an authorized service department.
- If the stroller is damaged or in need of repair, the authorized repair person will be able to confirm what is required.
- Alternatively, Synetik ErgoCare Customer Service Department is available for product information.
- Only PATRON Bohemia a.s. replacement parts should be used on this stroller.
- Should parts of your stroller become stiff or difficult to operate (no more than every 6 months), ensure that dirt is not the cause, apply a light coating of a spray lubricant e.g. WD 40 (do not use oil or grease). Ensure that the lubricant penetrates between the components but do not overspray as this may encourage the accumulation of dirt. Wipe off any excess spray. Regular lubrication of moving parts will extend the life of your stroller and make opening and folding the stroller easier.
- Make sure that you are familiar with the operation of the stroller . It should always open and fold easily. If it does not, do not force it. Stop and read the instructions.
- Wheels, plastic and metal parts may be wiped clean with warm water and mild detergent.
- The seat fabric can be washed by hand (Max. 30 °C).
- Never clean with abrasive, ammonia-based, bleach-based or spirit-type cleaners.
- After exposure to damp conditions, dry off with a soft cloth and leave fully open in warm environment.
- Store your stroller in a safe, dry place.
- Check and tighten the screw connections regularly.
- If there are any used control cables (brakes, telescopes, tilt), check them regularly and, if necessary, oil them.
- Synetik ErgoCare reserves the right to make design changes to any of their products without prior notice as a part of their continuous development program.