
cosihome®

INSTRUCTION MANUAL
BEDIENUNGSANLEITUNG
MODE D'EMPLOI
MANUALE DI ISTRUZIONI
MANUAL DE INSTRUCCIONES
HANDLEIDING

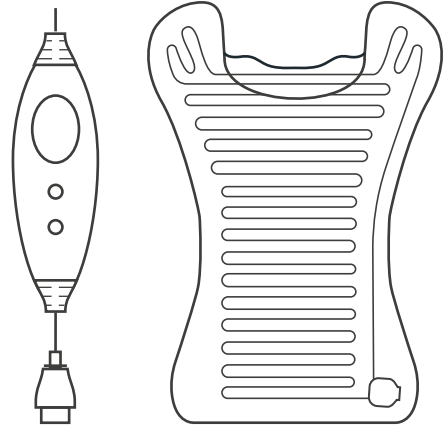MODEL CH-HP03G

Electric Heating Pad for Neck, Shoulders and Back

IMPORTANT INSTRUCTIONS – RETAIN FOR FUTURE USE

CONTENTS

EN	Safety Instructions.....	Page 3
	Specification.....	Page 7
	Symbols Explained	Page 7
	Parts	Page 7
	Operating Instructions.....	Page 8
	Error Codes.....	Page 9
	Washing.....	Page 9
	Drying	Page 10
	Storage	Page 10
	Recycling and Disposal.....	Page 11
DE	Deutsch.....	Page 12
FR	Français	Page 22
IT	Italiano	Page 32
ES	Español	Page 42
NL	Nederland	Page 52

cosihome[®]

ELECTRIC HEATING PAD FOR NECK,
SHOULDERS AND BACK
MODEL CH-HP03G

Thank you for choosing to purchase a product from Cosi Home. Please read the entire manual carefully prior to first use and keep in a safe place for future reference.

SAFETY INSTRUCTIONS

PLEASE READ AND SAVE THESE IMPORTANT SAFETY INSTRUCTIONS

When using electrical appliances, basic safety precautions should always be followed to reduce the risk of fire, electric shock, and injury to persons, including the following:

- ⚠ WARNING: KEEP AWAY FROM FIRE OR FLAMMABLE MATERIALS.**
- ⚠ WARNING: TO AVOID THE RISK OF ELECTRICAL SHOCK, FIRE OR INJURY TO PEOPLE, ALWAYS MAKE SURE THE HEATED PAD IS UNPLUGGED FROM THE ELECTRICAL OUTLET BEFORE RELOCATING, SERVICING, OR CLEANING IT.**
- ⚠ CAUTION: PROLONGED USE OF THE HEATED PAD AT A HIGH TEMPERATURE SETTING COULD LEAD TO SKIN BURNS. DO NOT USE FOR EXTENDED PERIODS OF TIME.**
- ⚠ CAUTION: IF YOU FEEL UNWELL WHILE USING**

THIS HEATED PAD, TURN OFF THE HEATED PAD IMMEDIATELY.

- Do not use by persons insensitive to heat and other very vulnerable persons.
- Use this product only as described in this manual. Any other use not recommended by the manufacturer may cause injury to persons.
- Do not use the heated pad when folded or creased.
- Do not cover or trap the controller under pillow, mattress or blanket.
- Do not place the controller under the heated pad whilst in use. Heat can damage the controller.
- Only operate the heated pad with the detachable controller supplied
- Do not immerse the heated pad in water.
- Not intended for medical use in hospital.
- Do not use this heated

pad on any parts of the body that are swollen, injured, or inflamed.

- If you have any medical conditions or concerns, consult your doctor before using this heated pad.
- If you experience persistent muscle or joints pain, please inform your doctor. Persistent pain can be a symptom of an illness.
- If you experience any pain or discomfort while using the heated pad, stop using it immediately.
- The heated pad can be used by children older than 3 years and younger than 8 years under supervision and with the control always set to minimum temperature value.
- The cable and control unit of the heated pad can give rise to risks of entanglement, strangulation, tripping or treading if not correctly arranged. The user shall make sure that excess ties and electric cords shall be arranged in a safe way.
- The heated pad is not to be used by persons (including children)

with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction.

- Children of less than 3 years and pets should be kept away unless continuously supervised.
- Children of less than 3 years must NOT use this product due to their inability to react to overheating
- Do not leave children alone or unattended in the area where the heated pad is in use.
- Cleaning and user maintenance shall not be made by children unless continuously supervised.
- The heated pad is not to be used by young children, unless the controls have been pre-set by a parent or guardian, or unless the child has been adequately instructed on how to operate the controls safely
- Avoid placing the heated pad in areas where it can become a tripping hazard.
- Do not wipe or wash

the heated pad with chemical solvents. Prior to cleaning ensure the heated pad is off and that the batteries have been removed.

- Do not use the heated pad when it is wet or damp
- Do not operate the heated pad with wet hands.
- Do not use this heated pad near water or in the immediate surroundings of a bath, shower or a swimming pool. Never place the heated pad where it may fall into a bathtub or other water container.
- Never leave the heated pad unattended while in use. Always turn the heated pad off and disconnected from the power supply when not in use.
- Do not use the heated pad while sleeping or place it near a person who is sleeping.
- Do not run cord under carpeting. Do not cover cord with throw rugs, runners, or similar coverings. Do not route cord under furniture or appliances. Arrange cord away from traffic area

and where it will not be tripped over.

- Avoid placing the power cord in areas where it can become a tripping hazard. Do not place the power cord under carpeting or cover it with rugs, runners or similar.
- To avoid a circuit overload when using this heated pad, do not operate another high-wattage appliance on the same electrical circuit.
- Always plug the heated pad directly into a wall socket. Avoid using extension leads where possible as they may overheat and cause a risk of fire.
- Do not use this heated pad with an external programmer, timer switch or any other device which would switch the heated pad on automatically unless this has already been pre-built into the heated pad by the manufacturer.
- Do not pull on the power cord. Never move, carry, or hang the heated pad by the power cord.
- Do not kink or wrap the power cord and plug around the heated pad,

as this may cause the insulation to weaken or split, particularly where it enters the heated pad.

- Ensure the mains cable is fully extended before use and arrange it so that it is not in contact with any part of the heated pad.
- Do not operate the heated pad if there are signs of damage to the heated pad, power cord, or any of the accessories supplied. If the supply cord is damaged, it must be replaced by the manufacturer, its service agent or a similarly qualified person in order to avoid a hazard.
- Do not attempt to repair or adjust any electrical or mechanical elements in the heated pad. Doing so may be unsafe and will void your warranty.
- Ensure that the mains supply (including voltage, frequency and power) complies with the rating label of the heated pad.
- Keep away and do not insert pins or sharp objects into the heated pad
- Examine the heated pad frequently and return to the supplier if has been

misused or damaged.

- Always unplug the appliance when not in use.
- Children must be supervised not to play with the appliance.
- Keep the appliance out of reach of children and pets.
- This appliance is intended for indoor household use and similar applications. Do not use outdoors.
- Do not use the appliance near petrol, paints, flammable gases, ovens, or other heating sources.
- The mains cable should not hang down from the surface on which the appliance is placed in order to prevent the appliance from being pulled down.

SPECIFICATION

- Model: CH-HP03G
- Voltage: 220-240V~
- Frequency: 50 Hz
- Power: 100 W
- Product Dimensions: 65x90 cm
- Net Weight: 0.8 kg
- Gross Weight: 1 kg
- Maximum Temperature: 65 °C

SYMBOLS EXPLAINED

-  Gentle cycle at a maximum of 40°C
-  Do not dry clean
-  Do not iron
-  Do not tumble dry
-  Do not bleach
-  Read instruction manual
-  Not to be used by children ages 3 or younger
-  Do not insert pins

 Equipment complies with the applicable EU safety requirements and directives.


 Do not use folded

 Equipment complies with the applicable UK safety requirements and regulations.

 Class II Double Insulated Equipment


PARTS

1. Heated pad
2. Connector one
3. Connector two
4. Controller
5. Power cable

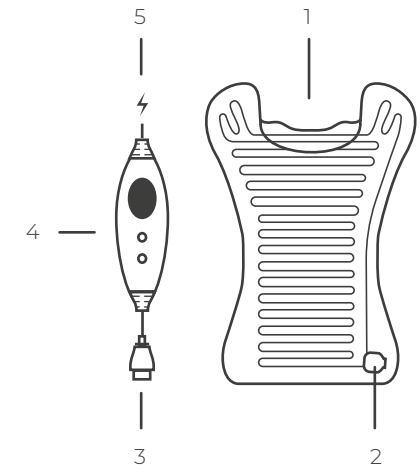CONTROLLER

1. LCD Display
2. Temperature setting indicator
3. On/Off button
4. Temperature setting button

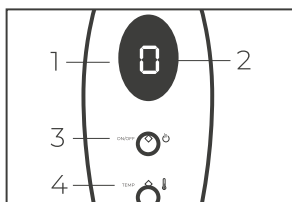These images are for reference only, please refer to your product.

If you find you're missing any parts from your order, please contact our customer care team by emailing help@cosihome.com for fast and friendly help.

OPERATING INSTRUCTIONS

Turning on

1. Insert connector two into connector one
2. Plug the power cable into a power socket.
3. Ensure the power socket is turned on.
4. Press and hold the on/off button for three seconds.
5. The heated pad will turn on and the temperature setting indicator number will illuminate on the controller's LCD display.
6. The heated pad will turn on in heat level 3.

Preheating

For best results, preheat your heated pad for at least 10 minutes before use.

1. Repeatedly press the temperature button setting until you reach heat level 6.
2. Leave the heated pad to preheat for at least 10 minutes.

Setting the Temperature

1. Press the temperature setting button repeatedly to cycle through the six available temperature settings.
2. The heat level will display on the controller's LCD display.

Automatic Turn Off

For your safety, your heated pad is fitted with a 90-minute automatic off timer that begins when you turn the heated pad on. After 90 minutes, the heated pad will turn off.

If you would like to continue using your heated pad after the 90-minute timer has elapsed, turn it back on and use as normal. If more than 10 minutes have elapsed, preheat your heated pad again.

Manual Turn Off

1. Press and hold the on/off button for three seconds.
2. Remove the power cable from the power socket.

Error Codes

Error Code	Cause	Solution
P	The controller is not connected properly.	Unplug the heated pad. Ensure the controller and or connector are connected properly. Plug it back in to power and use as normal.
E	The heated pad has been used while folded or crumpled and has overheated. As a safety measure, it has turned off.	Turn the controller off, unplug the heated pad, lay it down completely flat, and leave it to cool down.
C	The controller has been placed under the pad and has overheated. As a safety measure, it has turned off.	Turn the controller off, unplug the heated pad, lay it down completely flat, and leave it to cool down.

WASHING

⚠ WARNING: ALWAYS TURN THE HEATED PAD OFF, UNPLUG THE POWER CORD FROM THE ELECTRICAL OUTLET AND ALLOW THE HEATED PAD TO COMPLETELY COOL DOWN BEFORE HANDLING OR CLEANING IT.

⚠ CAUTION: ALWAYS DISCONNECT THE CONTROLLER AND THE POWER LEAD FROM THE HEATED PAD BEFORE WASHING. DO NOT ALLOW WATER OR OTHER LIQUIDS TO COME INTO CONTACT WITH THE

REMOTE CONTROL OR THE MAINS POWER PLUG.

⚠ CAUTION: DO NOT USE ANY ABRASIVE CLEANERS OR SOLVENTS, AS THIS CAN DAMAGE THE HEATED PAD.

⚠ CAUTION: WHEN SEPARATED FROM THE CONTROLLER, THE HEATED PAD CAN BE WASHED BUT NEVER DRY CLEANED.

- Do not attempt to clean the heated pad in any other way than that detailed by the manufacturer.

- Do not expose the controller or power lead to any moisture.
- For hand washing, use a mild washing powder or detergent and warm water. Soak for five minutes and then gently wash and rinse thoroughly in clean water.
- For machine washing, wash as wool at 40°C.

DRYING

⚠ CAUTION: DO NOT CONNECT THE HEATED PAD TO THE POWER SUPPLY IN ORDER TO DRY IT. DO NOT DRY THE HEATED PAD IN THE SUN OR USE AN IRON TO DRY IT AS THIS WILL DAMAGE THE INSULATION LAYER OF THE HEATED PAD.

- Never wring the heated pad after it has been washed.
- Do not place pegs on the wired areas of the heated pad as this will damage the internal circuit.
- Always air dry the heated pad.
- Always ensure the heated pad is 100% dry and has returned to its

original size before use.

STORAGE

- When storing the heated pad, allow it to cool down before folding.
- Do not crease the heated pad by placing items on top of it during storage.
- The heated pad is equipped with a detachable controller. Ensure the controller is disconnected from both the mains power and the heated pad and left in a safe, dry place, out of reach of children.
- To store your product, we recommend using the original or a similar sized box.
- Store the heated pad in a safe, clean and dry place, out of reach of children when not in use.
- Examine the heated pad frequently and prior to use after storage.

RECYCLING & DISPOSAL

Waste electrical and electronic products (WEEE), batteries, accumulators and packaging should not be disposed of with general household waste. This is to prevent possible harm to the environment and human well-being, and to help conserve our natural resources. Most electrical products with a plug, battery, or cable can be recycled. Contact your local council, household waste disposal service, or the store where you purchased the product for information regarding the collection schemes in your area.

Waste Electrical and Electronic Equipment (WEEE)

This symbol on the product and/or packaging indicates that Waste Electrical products (WEEE) under Directive 2012/19/EU which governs used electrical and

electronic appliances, should not be disposed of with general household waste.

The Mobius Loop

This symbol indicates that the product or packaging is capable of being recycled.

Take-Back Scheme

Waste electrical products (WEEE) can be sent to One Retail Group for recycling and disposal when you purchase a new product from One Retail Group, on a like-for-like basis. For more information on how to recycle this product please visit www.cosihome.com. The end-user of electrical and electronic products (WEEE) is responsible for separating old batteries and lamp bulbs before they are handed over to a designated disposal and recycling service. Waste batteries can also be disposed of in collection bins, which can be found in most local supermarkets and household waste recycling centres.

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© Copyright 2023
EU Authorised Representative: Brandrep Limited,