


**vive®**

# **BARIATRIC COMMODE**


User Manual  
LVA1059

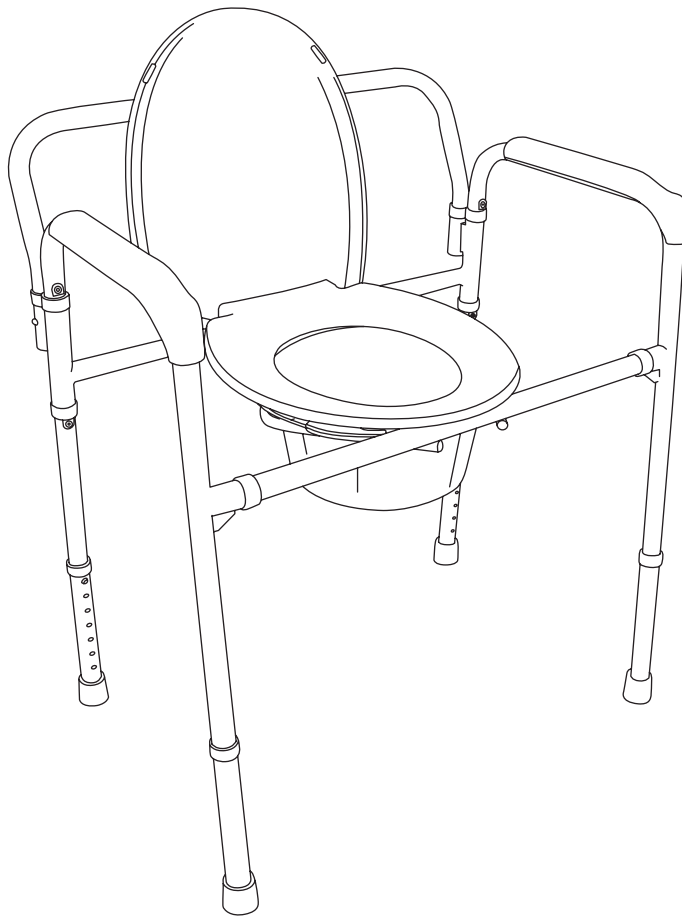

**vivehealth.com**

## PARTS LIST

- A) 2x Commode Frame
- B) 1x Bucket Frame
- C) 4x Legs
- D) 1x Backrest
- E) 1x Toilet Seat
- F) 1x Cover
- G) 1x Bucket & Lid

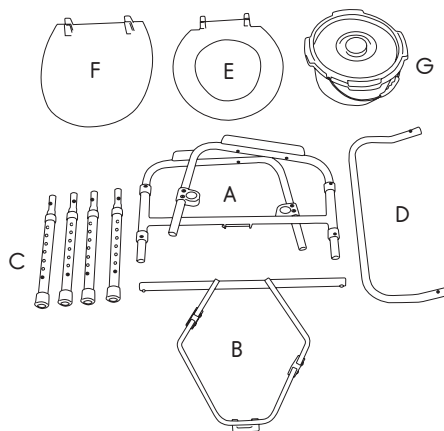


## SPECIFICATIONS


- Width: 21"
- Depth: 19"
- Adjustable Height (From Seat): 16.5" - 22.5"
- Adjustable Height (From Armrest): 24" - 30"
- Bucket Diameter: 11"
- Max User Weight: 350 LBS

## ASSEMBLY INSTRUCTIONS


1. Open the commode frame (A) and place on floor.

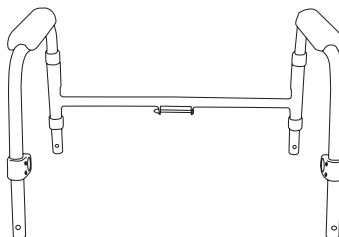


2. Slide the hook on the bucket frame (B) over the receptacle located on the back of the commode frame.

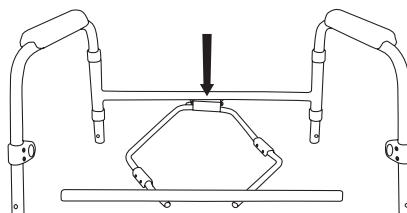

- Secure bucket frame (B) in place by inserting the front cross brace into the fastener receptacles located on the front of the commode frame. The push buttons will "click" in place when they are secure.

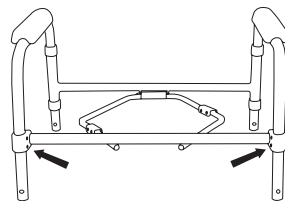


- Turn the frame upside-down. Insert the legs (C) into each receptacle by depressing the push pins and sliding into place. The brass push pins will have an audible click when they are in place and secure.

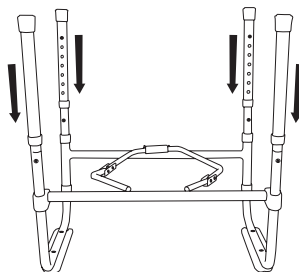


- Turn the commode right-side up. Insert the backrest (D) into the proper receptacles located on the back of the commode frame by depressing the push pins and sliding into place. You will hear an audible click when secure.


**Note:** Backrest should be facing outwards, away from the seat. If installed backwards, the toilet seat will not open completely.

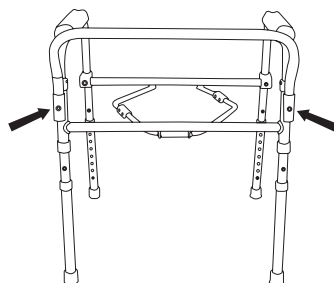

- Attach the toilet seat (E) and cover (F) to the back of the frame. The toilet seat hooks should fit inside the hooks of the cover.

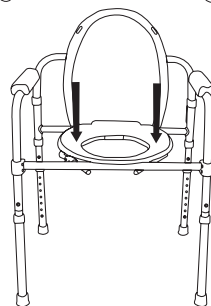

- Lower commode bucket into frame and return seat to down position. Now your seat is fully assembled

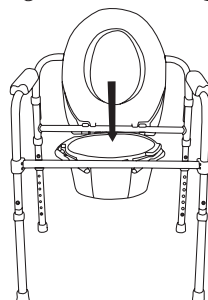

## WARNING

- Check periodically to make sure brass pins and fasteners are secure.
- Toilet seat must be in down position before use. Failure to do so may result in bodily injury or commode breaking.
- Do not install commode without first reading and understanding this user manual. Injury may occur if not doing so.
- Inspect all parts for shipping damage before assembly. If there is damage, please contact Vive Health Customer Support.
- Check rubber tips on the leg extensions for rips, wear or if they are missing. Immediately replace any if necessary.
- All four legs must be touching the floor simultaneously at all times to ensure safety.
- Users with limited physical strength should be supervised or assisted while using commode.

## ADJUSTABLE LEG HEIGHT

Before use, adjust height of each leg to your specific needs, making sure all legs are at the same height. When desired height is reached, lock legs into place by aligning buttons into proper holes.

## GUARANTEE REGISTRATION

You are protected by Vive Health's industry leading guarantees and customer service:



If you did not purchase through [vivehealth.com](https://vivehealth.com), please register at [vivehealth.com/register](https://vivehealth.com/register) to validate your guarantee.  
**Product Code: LVA1059**





We sincerely appreciate your business.  
We strive to provide you with the best  
quality products at great value.

If you have any questions please contact us.



[service@vivehealth.com](mailto:service@vivehealth.com)



1-800-487-3808

Monday - Friday 9am - 5pm EST