

The products in the warranty period(one year)because of the following circumstances the failure to provide free repairs.Or discretion to compensate the costs of repairs or replacement parts.


- 1、 Have found damage or purchase bad use.
- 2、 Arising from the use of operating in accordance with the specification of damaged natural persons(but misuse,or failure to use designated methods to natural disasters, such as earthquakes and fires damaged,the company not free repair responsibility).
- 3、 Non-expendable parts damage repairs.

Warning

- ① Not inconsistent with this machine in standard voltage of use.
- ① Surface water will not signal controllers to avoid controller failure.
- ① Not in the bathroom"swimming pools"hot springs and other moist places special use this machine.
- ① Please do not leave heaters placed in close traffic.
- ① In this product are advised not to use damaged shell or jacket to avoid the risk of injury should be timely replacement shell or jacket.
- ① Connecting power plugs in the sockets,you must check the power switches on the location.
- ① Please direct pulling plugs with wet hands Chaba power wiring or plugs.
- ① Stubbs power switches must be used here ,in the position and the correct weight plugs from outlets.
- ① Plugs or other components such as wires or damaged ,please do not arbitrarily packed this product.To contact local dealers for repairs or maintenance or repair point.
- ① The products have to work long hours may partially heat for heating sensitive time for the use of one-time too long.
- ① Among children under the age of 12 and prohibit the use of this product unclear.
- ① Pregnant women and people who do use physical period.

HEALTHY WITH YOU EVERYDAY

ROBOTIC CUSHION massage

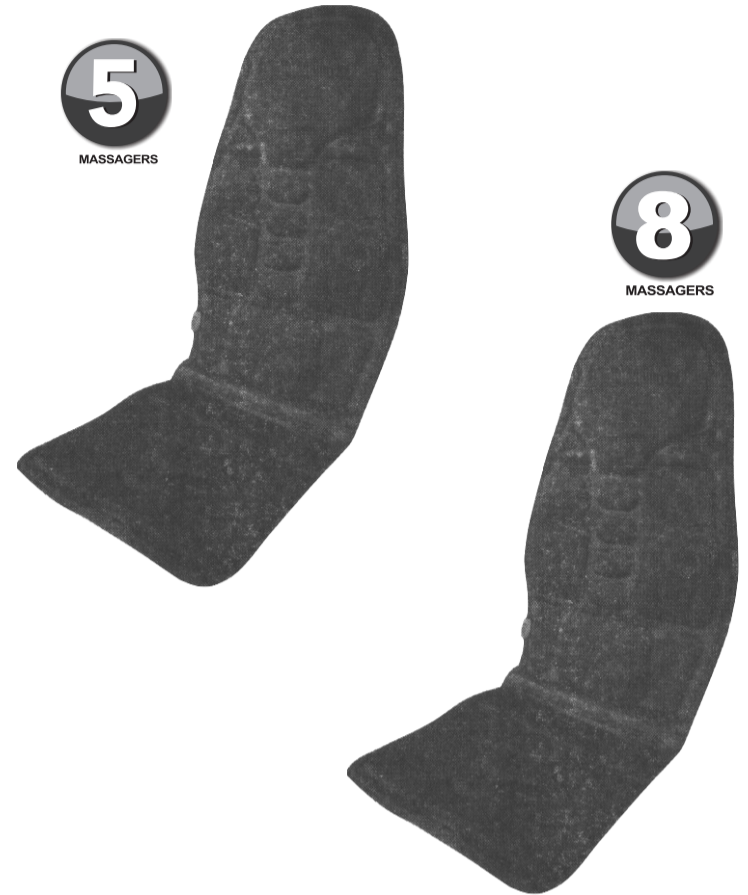

MASSAGE BODY ACUPUNCTURE POINT
PROMOTE BLOOD CIRCULATION
ELIMINATE STRESS

FEATURES

The products massage neck,back,waist,baokbone, leg,etc.can relieve fatigue,muscle tension,pressure, and give you vigor quickly.Having every parts of body massage can make you energetic.

1. Advanced comptrter chips and digital control technology are adopted.Make a break-through the traditional idea.with strong trapping, oscillates and massage two combine into one.
2. trapping,kneading,vibration,and magnetic curative effect, Initiation the strong trapping with massage control mode.
3. Imitate the speciality ma ssa gem an the trick in division.
4. Move-jing's magnetic field is combined.and more remarkable effects.
5. Health protection the entity of recreation appears.
6. It's safer and more reliable to operate the massager due to its automatically timing device and overheat protector
7. In accordance with physical engineering,advanced l uxurious fabrics.
8. Advanced intelligent control,Really accomplish of "the labour power bes difficult to do" enjoyment.
9. more kinds of massage mode,It is with facility to regulate.
10. more comvenient to operate,carry for the convenience. Furnish with car cigarette lighter plug.

OPERATION INSTRUCTIONS

1. Switch button: when pressing the button for the first time the massager starts, pressing the button again it turns off.
2. Mode button: when the massager starts,pressing the mode button one time to change mode different. the first mode to the ninth' s turning.and the light bes on.
3. Hand switch button: on unautomatic mode, pressing the button for the first time the massager turns off. pressing the button again it starts. and the light bes on.
4. Fixed time button: when the massager starts, pressing for the first time the strong change, ten minutes, twenty minutes,thirty minutes turn on change.and adoption digital to show.
5. Strength button: when the massager starts,pressing for the first time the strong change.one to eight. kinds of massage strength turn on change.
6. Fix button: on unautomatic mode,pressing the button for the first time the parts to change,back, waist,leg turn on change,and the light bes on.

Common Phenomenon processing

Phenomenon	Possible causes	Approach
Connected power, the products do not work	Switches are open	Will be open switches
	Whether access to undesirable power plugs and docketts	Inspection power sockets, plugs and sockets and identify good contacts
Massage movements, the sudden change in the voice	Whether the use of time in a row	Cut off the power and let rest for half an hour afer work products
Machine suddenly stopped	Whether the use of time in a row , Preventing automatic protection switching machines	Cut off the power and let rest for half an hour after work products

Note: the above checks have yet to find fault, and the products will be immediately sent to the company the maintenance services center maintenance.