

ARRIVAL
New

KNEE CARE DEVICE

KNEE-JOINT CARE SOLUTION

Life Can Feel Good Again



MANUAL

Please carefully read this manual before using the device

KNEE CARE DEVICE

Thank you for purchasing the Knee Care Device.
Please take a moment to read this guide and store it for future reference.

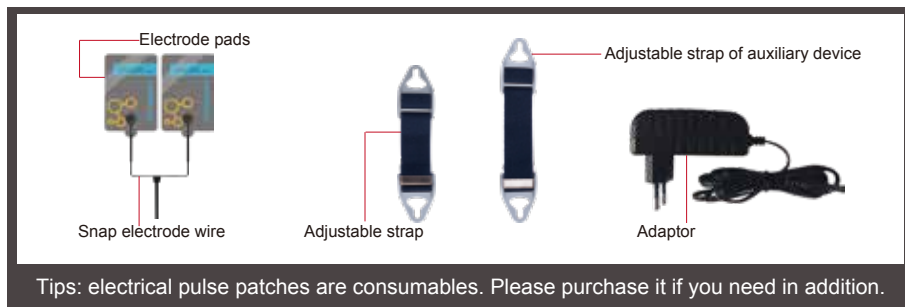
NOTE: This device is not intended for medical use. It is intended only to provide a relaxing massage at home.

STRUCTURE INTRODUCTION



- ① Main device ② Auxiliary device ③ Control display panel ④ Control line
- ⑤ Power light & LOGO ⑥ Power inlet ⑦ Auxiliary device control port
- ⑧ Electrical pulses output port ⑨ Heating pad ⑩ Vibrators

ACCESSORIES



Tips: electrical pulse patches are consumables. Please purchase it if you need in addition.

FEATURES

- Adjustable temperature 40°C-70°C (104°F-158°F)
- Vibration Massage
- Detachable and adjustable straps for knee
- Comfortable, ergonomically designed padded body
- Plug into an AC outlet


NOTE





- The default time for heat therapy is 30 minutes. Press the on/off button on the control panel again if you need more cycles.
- In different mode, some buttons for different functions will be locked. Please read HOW TO OPERATE to get more details.

CONTROL PANEL



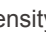

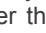




HOW TO OPERATE







1. Press  to turn on the device and Auto Mode(50°C, mixed strength vibration, level I EMS). Press again to turn off.

2. Press  to turn on the heating therapy. Press  to increase the temperature, Press  to decrease. Press  till the temperature is lower than 40°C, the heating therapy is turned off.

(PS: The temperature display flickers when it reaches the set temperature.)

3. Press  to turn on the ELECTRIC IMPULSE. Press  to increase the intensity level, press  to decrease the intensity level(press  till it is lower than the lowest level to turn off). Press  again to toggle between different electrotherapy modes. (PS: To avoid any uncomfortable, please increase the intensity slowly.)

4. Press  to toggle between low, high and mixed vibration strength. And the  Mode will be turned off.

5. Press  to turn off the 'Auto' mode. On the second pressing , the device will be restarted. In this condition, users will be allowed to press  or  to adjust the temperature, and press  or  to select different level of EMS.

Additional Description For EMS

1. To insert an external connector into the electrical pulse output, attach the pads to the location you'd like to massage (arms, legs etc.).

2. Reattach pulse pads to protective paper after using.

3. Please pay attention to distinguish the outlet.



USE OF SILICON ELECTRODE SLICES

- Insert one end of the 2-core connector wires into the electrode port and connect the other end to the electrode slices then use the massager as specified above.
- The electrode slices should be clean on their surfaces without any dust, oily or sticky substance. If their surfaces become dirty, flush with clean water and then air dry them for repeated use.
- The electrode slices can be used for about 50 times (20 minutes for each time). If their viscosity still cannot be recovered after several times of cleaning, consider purchasing new ones.
- Only the electrode slices provided with the massager should be used.
- Please dispose electrode slices according to your local environmental regulations rather than discarding them at will.

NOTE: It is recommended to use the EMS massager for one to two times a day. Please give your skin a break from the electrode pads every 20 minutes or so. Leaving pads in one place may cause skin irritation.

IMPORTANT INFORMATION

1. Only use the adapter that was supplied with the device.
2. Do not let the device fall from high places, and do not keep it around sharp objects.
3. Do not use this product when sleeping.
4. Do not use other medical electronic devices at the same time, such as pacemakers, heart-lung machines or electrocardiographs as this could be dangerous. If the patient is simultaneously using high frequency surgical equipment or other related equipment, the electrical pulse patches in the device may cause burns or damage the equipment. If the equipment is used near a microwave, the output of the device may become unstable.
5. Consult your physician before using this device, especially if you are pregnant, have skin sensitivities, heart disease, blood pressure abnormalities, malignant tumors, cerebrovascular diseases, diabetes, etc.
6. Please keep away from children.
7. When using the electric pulse function, do not move or change the position.

MAINTENANCE

- Store this device in a moisture-free, dry place when not in use.
- Do not wrap the power cord around the device.
- Do not use the power cord as a handle.
- Keep away from all solvents and harsh detergents.
- Unplug the device and allow it to cool before cleaning. Wipe with a soft, dry cloth as needed.
- Do not submerge the device in water or other liquids.
- Do not use abrasive cleaners, brushes, gasoline, kerosene, glass/furniture polish, or paint thinner to clean the massager.

PARAMETERS

NAME	KNEE CARE DEVICE
MATERIAL	ABS, TPU, EVA
SIZE	FREE
COLOUR	GRAY
VOLTAGE	12V DC
CYCLE	30 MINUTES
FUNCTION	EMS, MASSAGE, HEATING
APPLICABILITY	KNEE PAIN, MUSCLE RELEX, JOINT DISCOMFORT