
LIROMA

RED LIGHT THERAPY

Manual

INFRARED THERAPY LAMP (DOUBLE HEAD) - 507

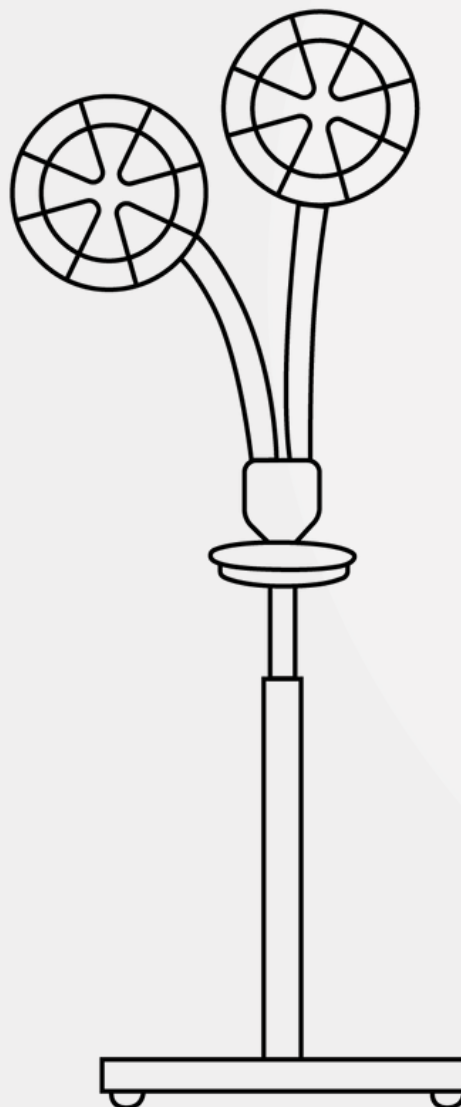Content

- Thank you for your purchase (2)
- Infrared therapy (2)
- Specifications (3)
- Infrared lamp double head - 507 (4)
- Benefits Liroma Infrared Therapy Lamps. (4)
- How to use it (5)
- Remote control (6)
- Warnings (7-8)
- Liroma (9)

1. Thank you for your purchase

Liroma wants its customers to feel healthy, fit and relaxed.

To achieve this, Liroma has a wide range of high-quality Infrared therapy products. Our products help with complaints such as: fatigue, poor circulation, pain in muscles and joints, pain due to underlying diseases, poor sleep, skin problems, stress and many other ailments. Our products target both mental and physical conditions.

2. Infrared therapy


Liroma's infrared radiator is for irradiating the human body only. Irradiation with infrared light transfers heat to the body. Blood flow to the irradiated skin increases and the metabolic process in the irradiated area is accelerated. Through the action of infrared light, healing processes in the body are stimulated and can be purposefully supported. Infrared lamps are great for muscle and joint pain, among other things. Infrared light can also be used as a complementary therapy in the treatment of throat, nose and ear disorders as well as to support facial and beauty care, especially for blemished skin.

First check with your doctor whether the use of an infrared lamp is medically useful and permissible in your case.

3. Specifications

Model	Infrared lamp double head - 507
Recommended user duration	5 - 12 minutes per session
Temperature emitter	40 C - 100 C
Dimensions	Maximum height 160cm
Mains connection	220-230V / 50Hz
Weight	11 kg
Recommended user distance	50 - 100 cm
Power	2x 300W
Warranty	2 years
Timer	0 - 60min
Operating modes	10
Infrared radiation type	IR-A

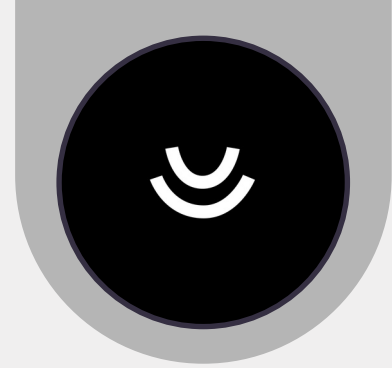4. Infrared lamp double head – 507

Liroma's infrared lamps are the only lamps that emit heat. The lamps have a strength that falls within the IR-A spectrum. This sets them apart from red light therapy lamps. Infrared lamps have a less profound effect than red light therapy lamps but still offer many health benefits.

IR-A

Infrared radiation has short wave, medium wave and long wave. Another name for short-wave infrared radiation is: IR-A. IR-A has the deepest effect in the skin compared to IR-B and IR-C. IR-A infrared lamps therefore also fall under therapeutic light. IR-A contains 700 to 1400 nm (nanometres). The wavelength is expressed in nanometres. Below 700nm we speak of visible light. Infrared light is not visible, but short-wave infrared light is closest to visible light and therefore has a very good therapeutic effect. Think of reducing muscle and joint pain. The closer the radiation is to visible light, the shorter the wavelengths and the more powerful the radiation is.

5. Benefits Liroma Infrared Therapy Lamps

- Pain relief: Infrared heat lamps can help relieve pain such as muscle pain, joint pain and chronic pain. The profound heat can help relax muscles and improve blood circulation, which can lead to pain relief.
- Improved blood circulation: Heat from infrared lamps can improve blood circulation by dilating blood vessels. This can help boost the supply of oxygen and nutrients to tissues, which can contribute to faster recovery.
- Relaxation and stress reduction: The heat from infrared lamps can help relax tense muscles and reduce stress. This can have a calming effect on the body and mind, leading to a sense of relaxation and well-being.
- Improved skin condition: Infrared therapy can improve skin condition by improving blood circulation to the skin and stimulating collagen production. This can help reduce wrinkles, improve skin texture and promote a healthy complexion.
- Relieving stiffness and flexibility: By heating the muscles and joints, infrared therapy can help reduce stiffness and improve flexibility. This can be especially beneficial for people with arthritis or other conditions that cause stiffness.
- Accelerated recovery: Infrared therapy can speed up the recovery process after injuries or intense exercise. The improved blood circulation and relaxing effect of heat can help reduce inflammation and promote tissue repair.
- Improved sleep quality: Some people find that infrared therapy can help them sleep better by promoting relaxation and reducing stress and tension in the body.

6. How to use it

The lamp has 10 Operating modes;

- ① Mode 1: 30 Watts
- ② Mode 2: 60 Watts
- ③ Mode 3: 90 Watts
- ④ Mode 4: 120 Watts
- ⑤ Mode 5: 150 Watts
- ⑥ Mode 6: 180 Watts
- ⑦ Mode 7: 210 Watts
- ⑧ Mode 8: 240 Watts
- ⑨ Mode 9: 270 Watts
- ⑩ Mode 10: 300 Watts

LIROMA


For the most optimal effect, we recommend using the lamp on the 300Watt stand at a distance of 60-80cm. This has been proven by research.

It is also recommended to use the lamp for a maximum of 12 minutes at 300Watt in one particular area. So the lamp may be used on several areas of the body in one day.

Setting Up and Using the Lamp

1) Remove all parts from the box.

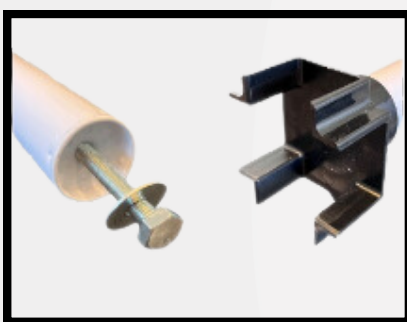2) Unscrew the screw from the bottom of the rod and attach the black part.


3) Assemble the cross base and attach the rod. Now tighten the screw securely.

4) Pull out the extension rod and tighten it securely.

5) Insert the lamp head into the extension rod.

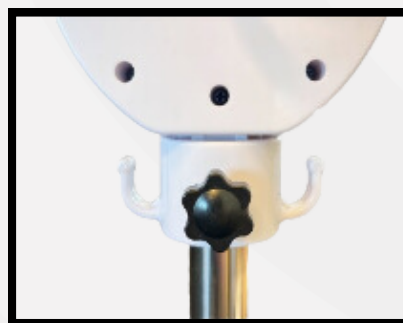
6) Tighten the screw securely

Would you rather see a video explaining how to assemble the lamp?


LIROMA
RED LIGHT THERAPY

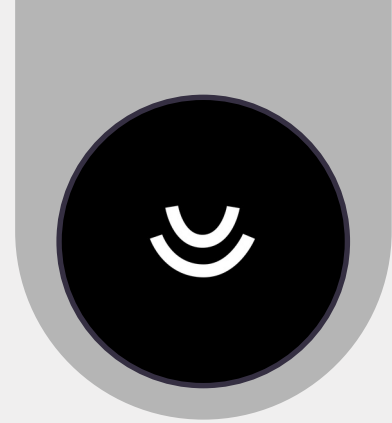7. Remote control

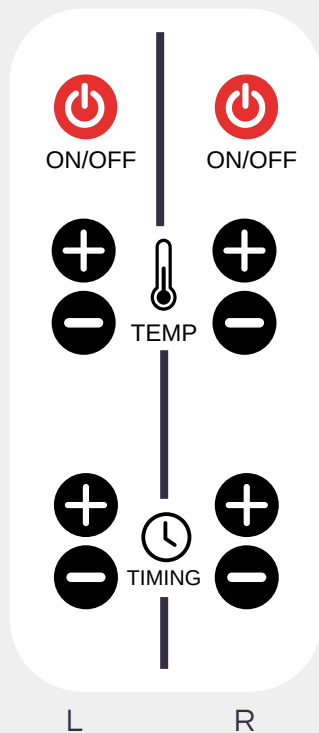
Press the button once, and the Infrared bulb will turn on. Press it again, and the light will turn off.

Adjust the temperature by clicking on the "+" sign. Each press of the button increases the lamp power by 30W.

Adjust the temperature by clicking on the "-" sign. Each press of the button decreases the lamp power by 30W.

Set the treatment duration by clicking on the "+". Each click increases the time by 5 minute.

Set the treatment duration by clicking on the "-". Each click decreases the time by 5 minute.

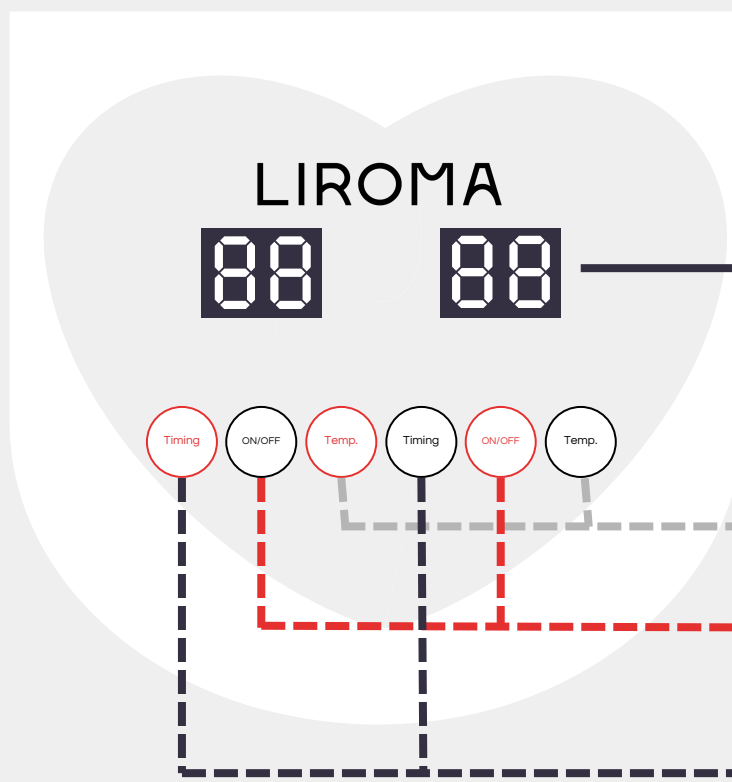Digital clock; Display of time duration (5-60min) & Modes (1-10).

Temperature adjustment; Each press of the button increases the setting by 1. Stages run from 1-10.

ON - OFF button; Press the button and the lamp turns on. Press the button again and the lamp goes off.

Time duration adjustment; Each press increases the time duration by 5 minutes. Clock runs from 5-60 minutes

9. Warnings

1. The housing of the lamp and the infrared emitter can become hot during use. Risk of burns if touched! Allow the lamp to cool down first.
2. Prolonged use can lead to skin burns.
3. Do not touch the infrared emitters and do not disassemble them while the plug is plugged into the outlet.
4. Children and adults with mental issues should use the lamp under supervision to prevent burns and electric shocks.
5. This device is not suitable for use by persons with physical, neurosensory, or mental limitations or lack of experience and/or understanding. Do not use the lamp in combination with electronic devices implanted in the body, such as artificial heart, artificial lungs, etc. Consult your doctor if in doubt.
6. Do not use the lamp on pregnant women, individuals with severe heart conditions, seriously ill patients, and other emergency patients. Consult your doctor if in doubt.
7. Do not use the lamp during menstruation, abnormal blood pressure, or special skin care. Consult your doctor if in doubt.
8. Do not use the device if you have an implant. Consult your doctor if in doubt.
9. Do not use the device on open or fresh wounds.
10. Do not use the device in combination with creams, ointments, or balms containing heat-producing ingredients. This may result in skin burns.
11. Medications, cosmetics, or food may, under certain circumstances, lead to hypersensitivity or allergic reactions of the skin. Stop the treatment if you experience this.
12. Are you taking medication? Always consult your doctor before use.
13. Never open the device under any circumstances. Never attempt to repair the device yourself. Serious injuries can occur. Your warranty will be void if you open or repair the device.
14. The experience of warmth may be limited or increased in the following cases:
 - In diabetics
 - In individuals suffering from insomnia, dementia, or concentration problems.
 - In individuals with allergies
 - In children and the elderly
 - After taking medications and alcohol

Infrared radiator

Prior to any cleaning and when changing the infrared heater, the appliance must be switched off, cooled down and the mains plug unplugged. Please note that only infrared heaters of the original Liroma type may be used during replacement.

Warning

- Do not touch the infrared emitters and do not unscrew them while the plug is plugged into the outlet.
- Disconnect the device from the mains to prevent electric shocks.
- For environmental reasons, the device should not be disposed of with regular household waste.


9. Warnings

This device is an electrical appliance. To prevent electric shocks and other electrical hazards, please follow these instructions:

1. Do not wash the product or parts of the product with water or other liquid agents.
2. Do not let your device fall into water or other liquids and do not place it in a location where it can fall into water or other liquids.
3. Do not expose it to running water.
4. If your device comes into contact with or is submerged in water or another liquid, immediately unplug it from the power outlet. If the device is used near water, we recommend using an outlet protected by a ground fault circuit interrupter.
5. Do not touch the device or plugs with wet hands.
6. Do not use or touch the device with wet or damp skin.
7. Do not pull or lift the device by the cord.
8. Do not use the device if the cord is damaged.
9. Do not use the device on or near heated surfaces.
10. Do not use the device in areas where it may be exposed to flammable or combustible products or fumes. A fire explosion may occur.
11. Do not store your device at temperatures higher than 140 degrees.
12. Do not use the device if it is damaged. Continuous use of a damaged device may result in electric shock or injury.
13. Do not use extension cords unless they are designed to handle the total power of the device being used.
14. Do not unplug the device by pulling directly on the cord. Grasp the plug itself, not the cord, when unplugging from the outlet.


If there is any damage or the device is not working properly, contact the seller as soon as possible or send an e-mail to info@liroma.nl.

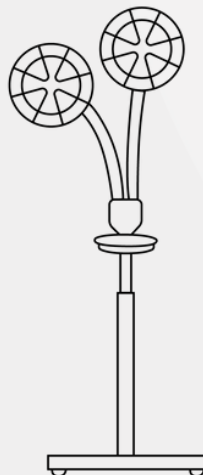LIROMA

Infrared Therapy

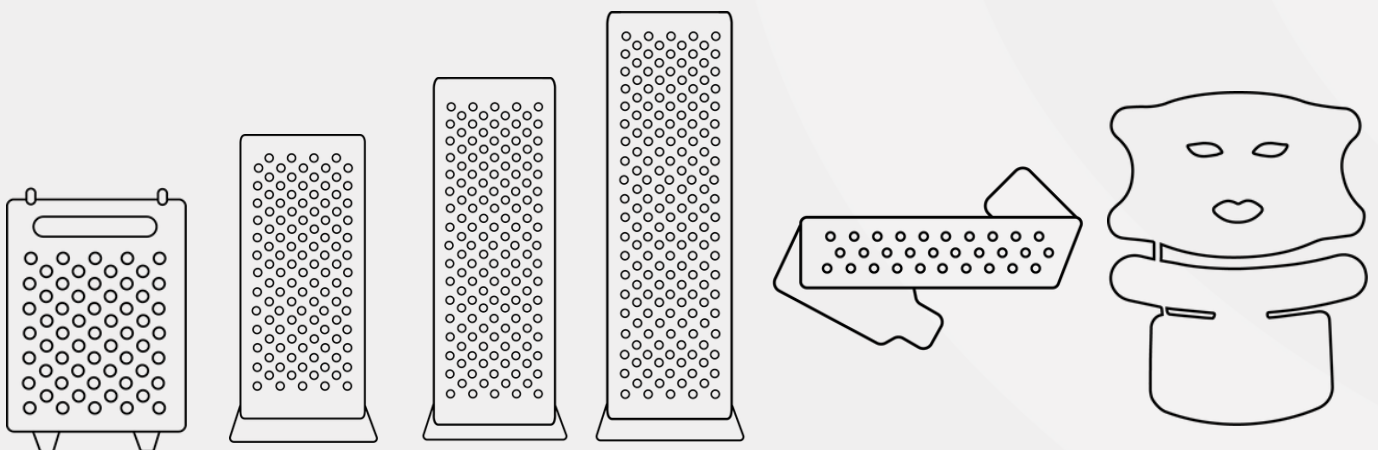Liroma has been a specialist in (infra)red light therapy for many years. We sell lamps that irradiate the entire body, as well as smaller lamps that you can target a body part. There are also (infra)red light devices specifically for the back or face. Check out our website for all our products.

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