

CONTROLLING SAMPLE DISPERSION IN MULTIDIMENSIONAL CAPILLARY LC-MS

R Henry¹, HJ Wirth¹, J Habsuda¹, Dan DiFeo², P Dawes¹ and E Dawes¹ ¹SGE International Pty Ltd, 7 Argent Place, Ringwood, Vic 3134, AUSTRALIA. ²SGE, Incorporated, 2007 Kramer Lane, Austin, Texas 78738, USA.

INTRODUCTION

A common technique in proteomics is to digest the whole protein fraction of a tissue sample and analyse the resulting mix of peptides. Samples generated this way tend to be very complex and due to a number of required sample preparation steps can be dilute. By using a multidimensional approach one can resolve a number of issues. The methods can be divided in two groups; sample preparation and improvement in resolution.

SAMPLE PREPARATION INCLUDES:

• PRECONCENTRATION/DESALTING

A column with a lower affinity (e.g. C4) can be positioned in front of an analytical C18 column. The column would be short with a relative large ID and is used to load a large, dilute sample fast onto the trap column while at the same time remove salts and other contaminations originating from the digest solution. The sample is then eluted in a concentrated and purified form onto the analytical column.

• DEPLETION/ENRICHMENT COLUMNS

These columns utilize a specific interaction to either remove or to capture certain parts of the sample. Any specific interaction or biorecognition can be utilized. Common examples are immuno affinity columns for the depletion of high abundance serum proteins (Figure 1), immobilized lectins for binding of carbohydrate structures in carboxylated peptides or immobilized metal ions (Fe(III) or Ga(III)) to bind phosphorylated peptides.

IMPROVEMENT OF RESOLUTION:

By employing a precolumn with opposite selectivity one can subdivide the sample in a number of fractions. The most commonly used combination is a strong cation exchanger followed by an analytical C18 column. The sample is loaded onto the SCX column and eluted with increasing salt concentrations in a step-wise fashion onto the analytical column (or a prior desalting column) (see Figure 2).

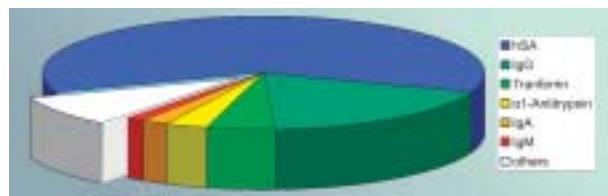


Figure 1. Abundance of proteins in human serum. Only a small portion of the "others" fraction is of any interest for proteomic research. Depletion of the 6 highest abundant proteins simplifies the analysis by reducing the masking effect of these proteins. Drawback is the possible co-precipitation of some proteins.

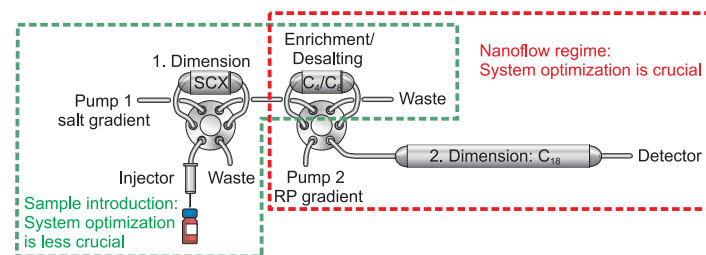


Figure 2. Possible setup for a multidimensional separation in proteomic science.

IMPORTANCE OF HIGH RESOLUTION

A mammalian cell expresses about 2000 proteins which when digested can form >10,000 peptides. Figure 3 shows the increase in complexity of the tryptic digests of Cytochrome C (a – Mw:12,200) and a monoclonal antibody (b – Mw: 150,000).

The task of separating the individual peptides becomes significantly more complex when mixtures of proteins are analyzed. When using an automated MS/MS detector, the software has to decide which masses (peptides) should be further investigated by the second dimension MS (sequence determination). To limit the amount of data the number of signals further analyzed are limited to, for example, 10. If there are more than ten peptides overlapping in one peak, the peptides from the lowest abundant protein are not analyzed and therefore lost. However, in a lot of cases the low abundant proteins are of higher scientific interest.

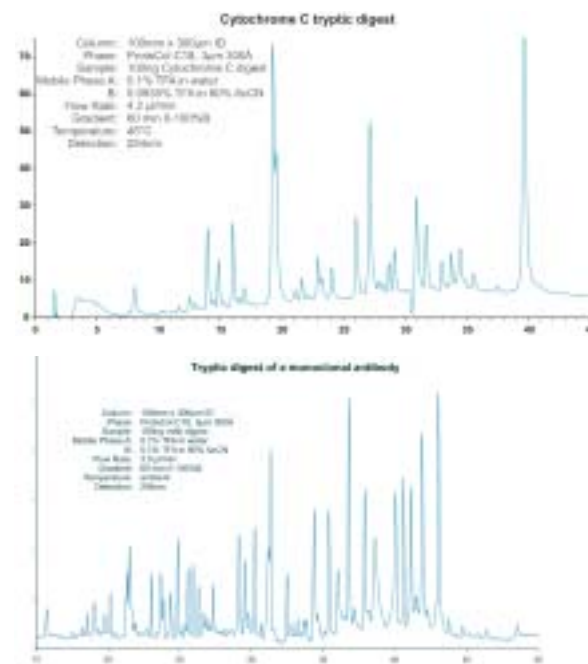


Figure 3. Chromatogram and sequence of Cytochrome C and IgG antibody

FACTORS CAUSING BAND BROADENING

There are a number of factors causing band broadening in Capillary LC:

• VOID VOLUMES

Void volumes can be introduced by non-ideal connections and excessive tubing length. The influence of tubing length on the chromatographic performance is shown in Figure 4.

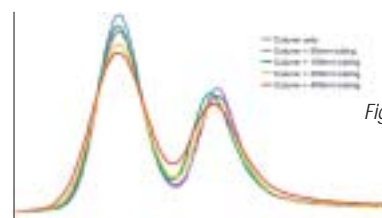


Figure 4. Effect of tubing length on resolution (tubing ID 50µm).

• DISTURBANCE IN THE LAMINAR FLOW

Under normal chromatographic conditions, the flow inside capillaries is laminar. Each connection can disturb the laminar flow and cause turbulence by having a gap between the connecting capillary ends or by misalignment of the bore caused by the tolerance on concentricity of capillary tubing (Figure 5).

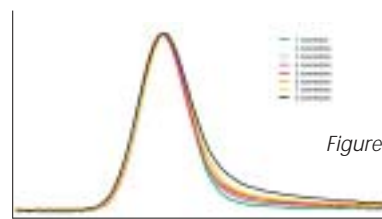


Figure 5. Effect of the number of connections on peak symmetry.

• NONSPECIFIC INTERACTIONS

A sample has the opportunity to interact with any surface within the flow path. If the surface chemistry differs from the selectivity of the stationary phase these interaction may lead to peak tailing and a reduction in chromatographic performance (see Table 1).

SOURCE OF NON-SPECIFIC INTERACTION	REMEDY
Exposed parts of the packing matrix	End-capping
Polymeric tubing material can be active for some molecules	Change tubing
Fused silica (Si-OH interacts with bases)	Chemical deactivation
Ferrous metal (tubing, frits) binds phosphate	Avoid Fe, coat element

Table 1. Origins on non-specific binding in a chromatographic system and method to minimize the effect.

THE PROTECOL™ APPROACH TO CAPILLARY LC

Unique features of the ProteCol range:

- Integrated connection tubing: reducing the amount of connections thus reducing the chance of user-introduced voids
- Zero volume butt connections: no additional volumes through connections
- Minimal volume filters and pre-columns
- All wetted surfaces are deactivated: minimization of nonspecific binding

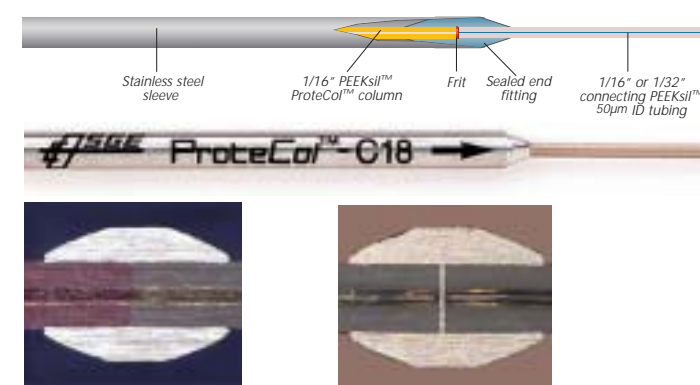


Figure 6. ProteCol column design and close up of cross-section of the column end.

