A close-up, high-angle shot of a woman's face, focusing on her eyes and nose. She has light-colored eyes and a slight smile. The background is a bright, ethereal blue with soft, wispy light rays emanating from the right side, creating a sense of depth and light. The overall tone is clean and clinical.

tria.®

**ACNE CLEARING BLUE LIGHT**

Instructions for Use

# WELCOME TO BLUE LIGHT TREATMENT WITH TRIA BEAUTY

## ENJOY HEALTHY, CLEAR SKIN WITH FEWER BREAKOUTS AT HOME

The Tria Acne Clearing Blue Light begins working immediately to eliminate existing blemishes and help prevent future breakouts with the same professional-level blue light technology once found only in leading dermatologist offices.

Unlike other acne treatments you may have tried, the clinically proven Tria Blue Light targets acne at the source, deep within the skin, to eliminate acne-causing bacteria safely and without irritation. Tria Blue Light is so gentle, you can use it every day to maintain consistently clear skin.

Before you use your Tria Acne Clearing Blue Light, please read these *Instructions for Use* to ensure best results. You can also find the *Quick Start Guide* on our website (see back cover).

# TABLE OF CONTENTS

<b>GLOSSARY</b> .....	2	<b>TAKING CARE OF THE DEVICE</b> .....	15
		Cleaning .....	15
<b>ABOUT ACNE</b> .....	3	Maintenance .....	15
		Storage .....	15
<b>GET TO KNOW YOUR TRIA ACNE CLEARING BLUE LIGHT</b> .....	4	Disposal .....	15
Display .....	5	Traveling with the Tria Acne Clearing Blue Light .....	15
		Life of the Device .....	15
<b>WHO SHOULD USE THE TRIA ACNE CLEARING BLUE LIGHT (INDICATIONS FOR USE)</b> .....	6	<b>TROUBLESHOOTING</b> .....	16
<b>GETTING STARTED</b> .....	8	<b>TECHNICAL SPECIFICATIONS</b> .....	17
<b>HOW TO USE THE TRIA ACNE CLEARING BLUE LIGHT</b> .....	10		
<b>WARNINGS</b> .....	12		
<b>RISKS</b> .....	13		

## GLOSSARY

**BACTERIA:** *Propionibacterium acnes* or *P. acnes* that causes acne breakouts. Lives on sebum secreted by the follicle.

**BLACKHEAD:** A small non-inflamed bump on the skin with a black center.

**CLARITY:** Freedom from blemishes and imperfections caused by more severe inflammatory reactions.

**COMEDONAL ACNE:** Blackheads and whiteheads.

**CYST:** A closed, liquid-filled structure under the skin. Acne cysts are present only in severe acne.

**INFLAMED:** Red and sensitive to the touch.

**NODULE:** A solid mass, like a knot, either raised or felt under the skin. Acne nodules are present only in severe acne.

**PIMPLE:** A red and sensitive, inflamed bump on the skin.

**PUSTULE:** A red and sensitive, inflamed bump on the skin filled with white or yellow pus.

**SEBUM:** Oil produced by sebaceous glands.

**TEXTURE:** Smoothness of your skin. When your complexion is free of bacteria, the result is a clearer, smoother texture.

**TONE:** Brightness and even skin colour. When your skin is healthy, it's brighter and the colour is consistently even.

**WHITEHEAD:** A small, non-inflamed bump on the skin with a white center.


## ABOUT ACNE

Acne begins deep under the skin with the buildup of oil in the sebaceous glands. When this oil and dead skin cells clog the hair follicle, *P. acnes* bacteria proliferate. The skin becomes inflamed and breakouts develop.

For years, dermatologists have treated acne with traditional topical and oral medications, which sometimes left the skin irritated or resulted in side effects that made long-term treatment difficult. Now, Tria Blue Light offers a true breakthrough in the treatment of acne by providing a safe, gentle and effective means of clearing and helping to prevent blemishes by eliminating the underlying cause at the source.

Using a gentle, non-UV high-energy light source, Tria Blue Light treats mild to moderate inflammatory acne by eliminating the bacteria that cause acne—deep within the skin—without affecting the surrounding tissue or irritating skin. Even those with sensitive skin types can treat blemishes long term, safely and effectively with the Tria Blue Light.

## GET TO KNOW YOUR ACNE CLEARING BLUE LIGHT
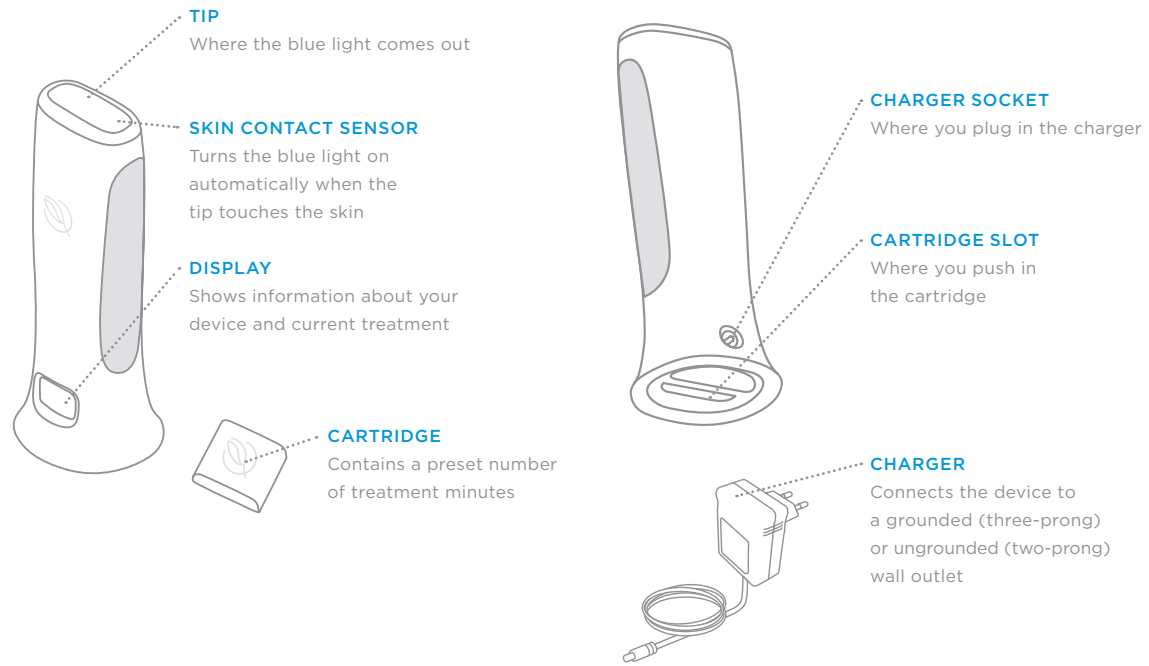
### DISPLAY

**BATTERY INDICATOR**  
Shows how charged the battery is: blinks when it is completely out of power and you need to recharge it.

**ELAPSED TIME**  
The large numbers show how long the current treatment has lasted (in minutes and seconds).

**MINUTES REMAINING**  
The small numbers show how many minutes of treatment time are left on the cartridge.

When there is no more time on the cartridge, or if there is no cartridge in the device, the treatment time disappears, the picture of the cartridge blinks, and the display says "REPLACE CARTRIDGE."

## WHO SHOULD USE THE TRIA ACNE CLEARING BLUE LIGHT (INDICATIONS FOR USE)

The Tria Acne Clearing Blue Light is indicated for the treatment of mild to moderate inflammatory acne. It can be used by individuals who want to eliminate acne-causing

bacteria, clear blemishes and acne breakouts, and improve the complexion without harshness. Use this chart to see if the Tria Acne Clearing Blue Light is right for you.

ACNE TYPE	USE THE TRIA ACNE CLEARING BLUE LIGHT?	DESCRIPTION	EXAMPLE
MILD ACNE	✓ YES	Several inflamed pimples (red and sensitive bumps)	
MODERATE ACNE	✓ YES	Many inflamed pimples and pustules (red and sensitive bumps filled with white or yellow pus)	
COMEDONAL ACNE	✗ NO	Blackheads and whiteheads (small bumps on the skin with a white or black center, not inflamed)	
SEVERE NODULAR OR CYSTIC ACNE	✗ NO	Many pimples and pustules with several nodules (solid masses of skin, like knots, either raised or felt under the skin) or cysts (closed, liquid-filled structures under the skin)	

### TEST ON SKIN BEFORE FIRST TREATMENT

Before your first treatment, test the Acne Clearing Blue Light on your forearm. This is important because your skin may be sensitive to light because of medications you may have taken or other medical conditions.

To perform the patch test, follow the step-by-step instructions beginning on page 8. It is normal for your skin to feel warm, but if you feel pain or your skin gets uncomfortably hot during treatment, stop using the Tria Acne Clearing Blue Light. The Tria Acne Clearing Blue Light is not right for you.

### WAIT 24 HOURS AND LOOK AGAIN AT THE TREATED SKIN

No redness or difference: If you don't notice any difference between the treated and untreated skin, you may begin the Tria Acne Clearing Blue Light treatments.

Redness or difference: If you notice any difference between the treated and untreated skin, such as redness, discoloration or discomfort, **SEE A HEALTHCARE PROFESSIONAL** and **DO NOT USE** the Tria Acne Clearing Blue Light again. The Tria Acne Clearing Blue Light is not right for you.

**NOTE:** Repeat this test if you think your skin sensitivity may have changed for any reason, for example as a result of changes in your medication or other medical conditions.


## GETTING STARTED

### THE FIRST TIME YOU TAKE THE TRIA ACNE CLEARING BLUE LIGHT OUT OF THE BOX, DO THESE 3 THINGS:

#### 1. CHARGE THE DEVICE

Plug the charger into the wall and into the back of the device. The display will light.

**NOTE:** You'll need to unplug the Tria Blue Light from the charger to perform your first treatment.

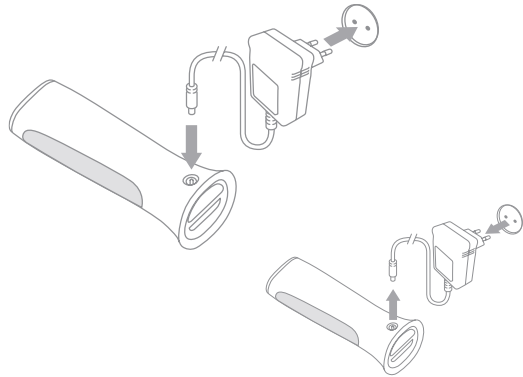
#### 2. INSERT A CARTRIDGE

To insert a cartridge, push it into the slot on the bottom of the device. The device will beep, and the display will light.

**NOTE:** The cartridge contains a pre-set number of treatment minutes. When the time on the cartridge runs out, you will need to replace it with a new cartridge.

To remove the cartridge, push it further into the slot until it pops out. You should take out the cartridge if you plan to store the device for a long time or carry it in a bag, so that the light does not accidentally come on and use up the cartridge time.
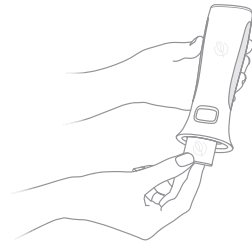
#### 3. TAP THE TIP OF THE DEVICE WITH THE PALM OF YOUR HAND

The light automatically comes on when the tip touches your skin and automatically turns off when the tip is not touching your skin. There is no On/Off switch.

During a treatment, when you move the tip back and forth across your skin, the light will stay on.

You may hear the internal cooling fan in the device, and the vents may feel cool. This is normal.

The display stays on for a while after the light turns off. This is so you can move the device somewhere else on your skin without losing track of your treatment time.

Now you're ready to start using the Tria Acne Clearing Blue Light.

## HOW TO USE THE TRIA ACNE CLEARING BLUE LIGHT

To use the Tria Acne Clearing Blue Light, follow these steps:


1. Clean and dry the skin you want to treat.
2. Place the Tria Acne Clearing Blue Light tip on your skin. This automatically turns on the blue light.
3. Move the Tria Acne Clearing Blue Light as if you are “painting” the skin with light. Listen for steady beeps every 90 seconds. Keep going for 6 minutes until the Tria Acne Clearing Blue Light beeps and turns off automatically.

**NOTE:** If you lift the device off your skin, the blue light will automatically turn off. Just put the tip back onto your skin to restart the treatment. Do not pause your treatment for more than 2 minutes, or the display will turn off and the treatment timer will be reset to zero.

**⚠ WARNING:** Do not use the Tria Acne Clearing Blue Light on or around your eyes or eyelids, and do not aim the Tria Acne Clearing Blue Light at your eyes or stare directly into the blue light. This could cause serious eye injury, including blindness or temporary loss of vision.

**⚠ WARNING:** Keep the tip moving. It is normal to feel your skin get warm, but **do not** hold the device in one place long enough for your skin to get uncomfortably hot, which could overheat the skin and causes skin injury. If your skin becomes uncomfortably hot, **stop** treating the area.

4. The Tria Acne Clearing Blue Light can also be used on the neck, chest or back. To treat on the neck chest or back, repeat steps 1-3.
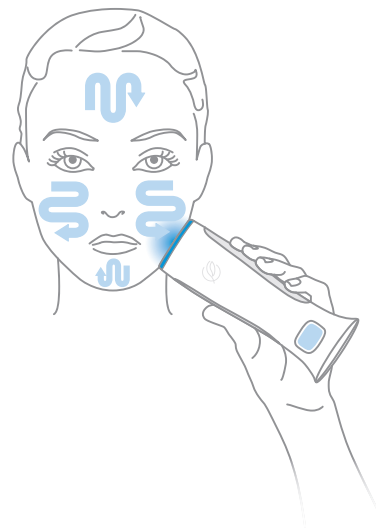
### HOW OFTEN TO USE THE ACNE CLEARING BLUE LIGHT

Treat your face with the Tria Acne Clearing Blue Light for 6 minutes once daily or 3 minutes twice daily (every 24 hours). For best results, do not skip treatments. If you do miss a treatment, do not perform an extra treatment. Continue to use the Tria Acne Clearing Blue Light daily.

### HOW LONG TO CONTINUE TREATMENTS

Continue to use the Tria Acne Clearing Blue Light for as long as you have visible acne and to help prevent future breakouts.

## WARNINGS

Caution – performance of the procedures other than those specified herein may result in hazardous radiation exposure.

### TO REDUCE THE RISK OF POTENTIALLY SERIOUS INJURY:

**⚠ DO NOT** use the Tria Acne Clearing Blue Light on or around your eyes or eyelids, and do not aim the Tria Acne Clearing Blue Light at your eyes or stare directly into the light. This could cause serious eye injury, including blindness or temporary loss of vision.

**⚠ STOP USING** the Tria Acne Clearing Blue Light and see a health care professional if you feel pain or your skin gets uncomfortably hot during treatment. Continuing to use Tria Acne Clearing Blue Light could result in skin injury such as temporary redness or blisters.

**⚠ STOP USING** the Tria Acne Clearing Blue Light and see a health care professional if you have skin redness that persists more than 24 hours after treatment. Failure to follow this warning may result in skin injury such as blisters or temporary pigment changes (hypo- and hyperpigmentation).

**⚠ DO NOT** use the Tria Acne Clearing Blue Light to treat severe nodular or cystic acne. Delaying proper treatment of severe nodular or cystic acne may be harmful. You should seek proper treatment without delay. See a healthcare professional.

**⚠ DO NOT** use the Tria Acne Clearing Blue Light to treat moles, warts, open sores, cancerous lesions or any skin condition other than acne. Delaying proper treatment of these conditions may

be harmful. You should seek proper treatment without delay. See a healthcare professional.

**⚠ DO NOT** treat the face more than twice a day (twice every 24 hours). Doing this could result in increased risk of skin injury such as temporary redness, blisters or pigment changes (hypo- and hyperpigmentation).

**⚠ DO NOT** let anyone under 13 use the Tria Acne Clearing Blue Light. The Tria Acne Clearing Blue Light has not been tested on children under 13 and may cause skin injury such as temporary redness or pigment changes (hypo- and hyperpigmentation) and may cause serious eye injury.

**⚠ DO NOT** share the Tria Acne Clearing Blue Light with other people. Sharing the device could transfer bacteria from one person to another and may cause harmful infection.

**⚠ DO NOT** use the Tria Acne Clearing Blue Light if the device is visibly damaged, and never attempt to open or repair the device. Doing this could result in a serious eye or skin injury or possibly fatal electric shock.

**⚠ DO NOT** use the Tria Acne Clearing Blue Light if the device is wet or if the device or power cord could fall into water while it's plugged in. Doing this could result in possibly fatal electric shock. These actions increase the risk of skin or eye injury. Please contact Customer Service for assistance.

## RISKS

The risks and benefits of the Tria Acne Clearing Blue Light treatments have been tested in a clinical study of 32 individuals with acne. If you use the Tria Acne Clearing Blue Light as directed, the risks are very low. In the clinical studies, none of the individuals treated (0 out of 32) experienced any serious side effects. However, the following mild side effects may occur with normal use:

**SKIN WARMTH:** The skin may feel warm during treatment.

**TEMPORARY SKIN REDNESS:** The skin may develop a slight redness that fades within 24 hours.

**TEMPORARY SKIN DISCOLOURATION:** The skin may develop a slight discoloration that lasts more than 24 hours. If this occurs, stop using the Tria Acne Clearing Blue Light.

**TEMPORARY SKIN DRYNESS:** The skin may develop slight dryness that goes away during the first 1-2 weeks of treatment.

**⚠** If the Tria Acne Clearing Blue Light is used improperly or it malfunctions, there are other risks. The following chart lists some of these risks.

ADVERSE EVENT	LIKELIHOOD OF ADVERSE EVENT	ADVERSE EFFECT	LIKELIHOOD OF ADVERSE EFFECT
General overuse of the device, such as performing many treatments per day	Unlikely	Skin injury, such as redness, blisters or pigment changes	Minor skin injury is unlikely; major skin injury is very unlikely.
Use of device on severe acne	Unlikely	User delays known effective treatments	Likely
Use of device on non-acne skin conditions	Unlikely	User delays known effective treatments	Likely

ADVERSE EVENT	LIKELIHOOD OF ADVERSE EVENT	ADVERSE EFFECT	LIKELIHOOD OF ADVERSE EFFECT
Use of device by children under 13	Unlikely	Same as for adults, skin injury, such as redness, blisters or pigment changes. Although, the severity of the effects may be different.	Unknown because use on children under 13 has not been tested
Use of device on or around eyes or eyelids	Very unlikely	Blindness or temporary loss of vision in the exposed eye	Very unlikely
Use of physically damaged or wet device or charger that exposes user to electrical hazard	Very unlikely	Shock or electrocution	Very unlikely (same as other common electrical devices)
Device malfunction with excessive blue light output	Very unlikely	Skin injury, such as pigment changes or blisters	Minor skin injury is unlikely; major skin injury is very unlikely.
Continued use despite side effects other than warmth, redness or dryness that quickly fades	Very unlikely	Skin injury, such as discolouration or blisters	Minor skin injury is unlikely; major skin injury is very unlikely.
Use of device despite skin sensitivity to light	Very unlikely	Skin injury, such as discolouration or blisters	Minor skin injury is unlikely; major skin injury is very unlikely.
Sharing a device with others who have common skin bacteria	Unlikely	Worse acne or a minor skin infection	Unlikely
Sharing a device with others who have serious skin infections	Very unlikely	Serious skin infection	Unlikely

## TAKING CARE OF THE DEVICE

### CLEANING

Clean the tip or the outside of the device whenever it appears dirty. Make sure the charger is not plugged in, then moisten a soft cloth with water and wipe down the tip or the outside of the device. **DO NOT** get water or liquid inside any part of the device.

### MAINTENANCE

The Tria Acne Clearing Blue Light does not need regular or special maintenance. If you think your Tria Acne Clearing Blue Light is not working properly or is damaged in any way, do not use the device. Contact Customer Service for assistance.

### STORAGE

Keep the Tria Acne Clearing Blue Light away from water, dampness or debris that may get inside the device. No special storage is needed.

### DISPOSAL

Please follow the regulations where you live for recycling or disposing of electrical devices.

### TRAVELING WITH THE TRIA ACNE CLEARING BLUE LIGHT

If traveling with the Tria Acne Clearing Blue Light, pack the device in its original packaging to prevent damage. The Tria Acne Clearing Blue Light will work properly at all common international voltages (any mains supply between 100 and 240 VAC at 50/60 Hz). You may need a plug adapter so that the charger that comes with the device can plug into the wall socket. Use only the accompanying charger unit to charge the Tria Acne Clearing Blue Light.

### LIFE OF THE DEVICE

The Tria Acne Clearing Blue Light expected service life is at least 3 years of normal use. There are no parts that need regular replacement. The Tria Acne Clearing Blue Light is working properly if light comes out of the tip when the tip touches your skin.

## TROUBLESHOOTING

Use this chart to help solve any problems. If you still have problems, please contact Customer Service for assistance.

PROBLEM	WHAT YOU SHOULD DO
I plugged in the charger and nothing happens. The display is off.	<p>Make sure the charger is plugged into the Tria Acne Clearing Blue Light and the wall socket.</p> <p>Check if the electrical outlet needs to be turned on with a wall switch or outlet switch.</p> <p>Check that the outlet works by plugging in another electrical device and confirming that it runs properly. If the other device doesn't work, plug the Tria Acne Clearing Blue Light into a different outlet.</p>
I touched the Tria Acne Clearing Blue Light tip to my skin, but I can't see any light.	<p>Look in a mirror while using the Tria Acne Clearing Blue Light on your face, or try using the Tria Acne Clearing Blue Light on your arm, and check if there is a halo of light around the tip of the device. If you see this halo, then the light is working. Never look directly into the light.</p> <p>Make sure that the sensor surrounding the Tria Acne Clearing Blue Light tip (see the picture on page 4) is making contact with your skin. Remember that the light comes on only when the sensor touches the skin.</p> <p>If you have treated your skin for a session of 3 minutes and the light went off automatically, then you are done treating.</p>

## TECHNICAL SPECIFICATIONS

**MODEL:** SPBL

### ELECTRICAL INPUT

**VOLTAGE:** 100-240 V - at 50/60 Hz

**CURRENT:** 0.30 A max

### LIGHT OUTPUT

**WAVELENGTH:** 412 nm

### PHYSICAL

**SIZE:** 65mm x 60mm x 170mm

**WEIGHT:** 270g


### CHARGER ELECTRICAL OUTPUT

**VOLTAGE:** 12 VDC

**CURRENT:** 0.5A

### SPECTRAL IRRADIANCE

**MAXIMUM OUTPUT:**  $6 \times 10^{-5}$  W/cm<sup>2</sup> (@412 nm)
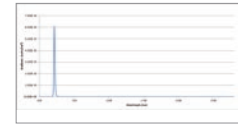
### DEVICE LABEL SAMPLE

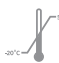
This device is:

1. A Class II ME Equipment
2. Internally powered in normal operation
3. IP22 rated to protect particulate >12.5mm entering the device and from dripping water.
4. Classified in Risk Group 1 in accordance with IEC 60601-2-57:2011. The charger provides isolation from supply mains.

Caution - UV emitted from this device.

## SYMBOL GLOSSARY

 Operating Temperature: 5°C to 30°C  
Operating Relative Humidity:  
Non-condensing 15-95%  
Operating Atmospheric Pressure:  
700 hPa to 1060 hPa

 Transport and Storage Temperatures:  
-20°C to 50°C  
Transport Relative Humidity:  
Non-condensing 0-95%

 This device is not suitable for use in the presence of flammable anesthetic mixture with air or with Oxygen or Nitrous Oxide.

Warning Symbol

Class II Symbol

WEEE Symbol

Type BF Symbol

Serial Number Symbol

CE Mark

Date of Manufacture Symbol

Manufacturer Symbol

Temperature Limitation Symbol

Accompanying Documentation  
Must be Consulted Symbol

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©2014 Tria Beauty, Inc., 4160 Dublin Blvd., Suite 200, Dublin, CA 94568. All rights reserved. Tria, Tria Beauty, and the Tria logo are trademarks of Tria Beauty, Inc.
MDSS GmbH